

ENGLISH



2020

SERVICES & ACCESSORIES
サービス & アクセサリー
SERVICE ET ACCESSOIRES
SERVICE UND ZUBEHÖR
服务或附件

MTS AFTERMARKET SOLUTIONS
be certain.



L'ENGAGEMENT MONDIAL, LA PRÉSENCE LOCALE 信守全球承诺, 专注本地服务 GLOBALE VERPFLICHTUNG,
REGIONALE BETREUUNG GLOBAL COMMITMENT, LOCAL DEDICATION グローバル コミットメント、地域レベルの貢献

Welcome to the MTS Services, Maintenance Parts and Accessories Catalog for 2020

Test and product development professionals rely on MTS to help increase lab productivity and accelerate time to market. Around the world, test engineers recognize that our solutions enable efficient test, measurement and simulation for a wide range of products and materials. Whether you need to evaluate advanced composites and alloys or test full-scale automobiles, airplanes, bridges or buildings, MTS testing solutions make possible many of the innovations that are changing how our world works.

Reliable Testing Solutions

Choosing the correct equipment, components and accessories to accomplish your test objectives is important because generating accurate test data requires reliable test equipment. MTS provides some of the most reliable and longest lasting test equipment on the market, and this catalog includes items to help you maximize equipment performance and longevity. Due to the wide range of solutions that MTS produces, including custom testing solutions, not every available item is listed within this catalog. If you cannot find a specific item, please contact your nearest MTS support office.

Responsive Global Service

With one of the largest, most experienced global service networks of any testing solutions provider, MTS can address your global testing needs quickly and efficiently. Local service and consulting representation throughout the world helps you maximize uptime and productivity. Additionally, MTS is adding connectivity to the lab to allow greater insight into equipment health and to streamline lab operations.

Table of Contents

INTRODUCTION 2-6

Welcome Letter	2
Table of Contents	3
How to Contact and Order from MTS	4
Contacting MTS	5
MTS Services	6

SERVICES 7-27

Training	7
Routine Maintenance	8-9
Fluid Assessment Services	9
Onsite Calibration	10-14
Factory Calibration Services	15
Metrology Calibration Services	16-17
Actuator Remanufacture Program	18-25
Software Support Plan	26-28
MTS Echo™ Intelligent Lab	29-32
Warranty Coverage	33

LOAD CELLS 34-41

Introduction	34
Servohydraulic Load Cells	35-36
Electromechanical Load Cells	37-38
Adapters Load Cells	39-41

EXTENSOMETERS 42-75

Introduction	42-43
How to Select an Extensometer	44
Application Index	45
Functional Index	46
Axial Extensometers	47-59
High Temperature Extensometers	60-62
Specialty Extensometers	63-67
MTS Advantage Video Extensometer	68-69
Laser Extensometer	70
Other Applications	71
Clip-On Displacement Gages	72
Displacement Gages	73
Extensometer Calibrator	74
DCPD	75

GRIPS & FIXTURES 76-121

Introduction	76-77
Polymer Matrix Composites Fixtures Index	78
Servohydraulic Grips & Fixtures	79-95
Composites Fixtures	96-102
Electromechanical Introduction	103
Advantage™ Grips	104-109
MTS Fundamental™ Grips & Fixtures	110-111
Bionix® Grips	112-117
Bionix EnviroBath	118-119
EM Extend® Kit	120-121

FURNACES & CHAMBERS 122-127

Introduction	122
Furnaces & Temperature Controller	123
Thermomechanical Fatigue (TMF)	124-125
Environmental Chambers	126-127

TESTING SYSTEM COMPONENTS 128-145

Introduction	128
201 Hydraulic Actuators	129-130
242 Hydraulic Actuators	131
244 Hydraulic Actuators	132
248 Hydraulic Actuators	133
249 Swivel Base & Rod Ends	134-135
252 Servovalves	136
215 & 216 Rotary Actuators	137
Rotary Actuator Options	138
515 Hydraulic Power Units	139
Multi-Pump Control Manager (MPCM)	140
293 Hydraulic Service Manifold	141
Isolation Service Manifold	142
Testline Modular Components	143
Additional Components	144

MTS MAINTENANCE PARTS 146-191

Introduction	146
Hydraulic Tools & Sample Bottles	147
Accumulator Parts	148
Actuator Parts	149-153
Load Frame Parts	154-155
Cables	156-159
Replacement Lamps	159
Adapter Cables for 3rd Part Equipment	160
Grip Parts	161
Hydraulic Hoses	162
Hydraulic Power Unit Parts	163-168
Hydraulic Service Manifold Parts	169-170
Filter Elements	171
Flushing & Shut-Off Valves	171
SWIFT® Parts	172-173
FlatTrac® & Rolling Road Systems Parts	174
Series 329 Road Simulator Parts & Durability Parts	175
329 System Parts	176
Model 329 4 DOF PC Road Simulator Parts	177-180
Model 329 6 DOF Road Simulator Actuator Parts	181
Model 329 6 DOF PC Road Simulator Parts	182-184
Model 329 4 DOF LT Road Simulator Parts	185-186
Model 329 6 DOF LT Road Simulator Parts	187-189
MAST Parts	190
320 Rebuild Kits	191

EXCHANGE & REPAIR PROGRAMS 192-207

Introduction	192
Just In Case™	193
Electronics	194-201
Onsite Exchange Programs	202-204
MTS Series 252 & 256 Servovalves	205-206
Remanufacture Programs - Actuators	207

COMPANY LOCATIONS 208-209

Company Locations	208-209
-------------------	---------

How to Contact and Order from MTS

With a global presence and local support, MTS has sales, service and training resources located near our customers. Within the next few pages of this catalog you will find a complete list of local offices and contact information. Simply call, e-mail or fax the nearest location to get in touch with a local representative.

MTS Teams Dedicated to Serve You

When you chose MTS for your test system solution needs, you gained access to a complete support network of teams who are ready to assist you. The teams on this page support the services and products offered within this catalog. There are many more MTS teams who support customers by consulting on, developing and installing custom or advanced testing systems.

ORDER SERVICES

This team researches repair and service parts and processes pick and ship orders that have no labor or custom content. Order Services also develops quotations for standard repair, maintenance, and service parts.

TECHNICAL SUPPORT

Support and assistance in answering technical questions that you may encounter when using your MTS testing system software and equipment is provided by Technical Support. See section below for information on how to contact this group.

SERVICE SALES

This group assists in arranging onsite services and consultation in developing service planning plus support for current service contracts. They will also provide quote support for products and services that require field service involvement.

SALES ENGINEERS

The Sales Engineers provide support and consultation when you have standard system or product expansion needs. They are readily available to provide you with information and assistance on new test systems, system upgrades, and new equipment for your full range of testing needs.

FIELD SERVICE

Your local source for onsite equipment installation is the Field Service team. They also perform onsite routine maintenance, calibrations and repair, or product exchange services for your equipment.

FIELD SERVICE COORDINATION AND SCHEDULING

This group will establish initial onsite service appointments, answer questions and help with scheduling priorities. They help accommodate schedule changes. Validate and clarify scope of work for onsite tasks to be performed while providing management of Field Service Resources.

CUSTOMER TRAINING

Standard and customized classroom training on MTS equipment and technology is available for those of all levels of technical ability. Training can be done at regional training centers or onsite. See the list of training centers in this section of the catalog to arrange for training.

Contacting MTS by...

PHONE

Call your local office listed on pages 208 and 209 of this catalog. Our North America Customer Care Center is also available at 1-800-328-2255 between the hours of 7:00 a.m. and 5:00 p.m. (USA Central time) Monday through Friday. In Europe, call toll-free number +800-81002-222.



FAX

The Regional Business Centers listed in this section have fax number access. When communicating via fax, please include information about the reason for your inquiry as well as how we should respond back to you.

INTERNET

Contact us online at www.mts.com. Use the website to learn more about MTS and the products and services you need to ensure the highest levels of testing performance. Simply select the Contact Us link from the menu in the top left of any page on the site to initiate an electronic information request. To contact MTS by e-mail, find the e-mail address for your local office in this section of the catalog, or send an e-mail to info@mts.com. If sending an e-mail, please include information about the reason for your inquiry as well as how we should respond back to you.

Contacting MTS

MTS Telephone Support

If you are a current MTS customer you may be familiar with our telephone support process. This process, which is described below, has been streamlined from what has been published in the past. If you are a new MTS service customer, please review these procedures to help us to meet your service needs.

CONTACTING TECHNICAL SUPPORT

MTS provides a full range of support services after your system is installed. If you have any questions about a system or product, contact Technical Support in one of the following ways.

TELEPHONE

- » Worldwide: 1-800-328-2255 - toll free in U.S.; +1-952-937-4000 - outside U.S.
- » Europe: +800-81002-222 - international toll free in Europe

INTERNET

- » www.mts.com > Contact Us (upper-right corner) > In the Subject field, choose "To escalate a problem; Problem Submittal Form"
- » E-mail, Worldwide: tech.support@mts.com
- » E-mail, Europe: techsupport.europe@mts.com

CALL PREPARATION

The Customer Care Center agent will ask you for some information to get you the help you need as quickly as possible. To help provide prompt support, please gather the following information prior to contacting MTS:

General information (essential)

- » Company name
- » Company address
- » MTS site number
- » Your name and phone number
- » All applicable support contract numbers

Describe the problem you are experiencing

- » How long has the problem been occurring?
- » Can you reproduce the problem?
- » Were any hardware or software changes made to the system before the problem occurred?

What to Expect When You Call

1. Your call will be registered by the Customer Care Center agent. The agent will ask for your site number.
2. The Customer Care Center agent may also ask you to verify information noted in the Call Preparation section.
3. What is the nature of your call today? The Customer Care Center agent will need to know if you are calling for technical support, to order parts, or to request a service call, to name a few.
4. If you have made a previous call regarding your issue, we can recall your file. You'll need to tell the Customer Care Center agent the following:
 - a. The MTS service request or the MTS quote number
5. The Customer Care Center agent, or another MTS employee, may ask you to perform certain tasks so we can identify the problem.
6. If you are calling to place an order, please have the following information ready:
 - a. MTS site number
 - b. Purchase order number, which you will need to complete your order.



MTS Services

Test professionals throughout the world rely on MTS Systems Corporation's innovative technologies, high-quality testing systems and applications expertise to optimize their testing programs. We have the experience to support your test equipment from pre-installation to de-commission and at every point in between. MTS has the service solutions to meet your needs for test schedule predictability, data integrity, system performance optimization and budget management.

ONSITE SERVICES

Although MTS builds some of the most rugged test solutions available, the constant motions and forces applied to test specimens can ultimately take their toll on the test systems as well. Our field service engineers have a worldwide reputation for applications expertise, and will respond to your request for support or repair quickly and efficiently. MTS can also assist with installation or movement of lab equipment – including disassembly, packing for transportation and installation at the new location. In addition, we offer consumables and spare parts for new-generation MTS equipment and most of our legacy systems.

ENGINEERING SERVICES

MTS offers a complete set of professional engineering services, including systems engineering, test consulting and facilities design services. MTS experts will listen to your test objectives, analyze your situation, and translate your desires to specific system requirements and an actionable plan. We can provide test designs, fixture engineering, control system evaluation, data collection and results analysis. By referencing the best practices of test labs worldwide, MTS can help you design test facilities, including hydraulic distribution systems. Let MTS help you develop long-range lab investment plans that support your business growth strategies or research plans.

CALIBRATION & ALIGNMENT

All test labs must calibrate their testing equipment to help ensure data accuracy. MTS provides top-quality calibration services accredited by A2LA to ISO/IEC 17025 standards. We can complete calibration at your location, or in our factory metrology labs. We also offer a range of services, including load frame alignment services, designed to help minimize data variance.

TRAINING

MTS training programs are designed to improve technician efficiency and maximize system performance. Expertly led and completely customizable, the courses provide hands-on learning to make sure you are thoroughly familiar with your test systems and know how to operate them effectively.

In addition to a broad selection of standard courses, MTS can customize courses to meet your specific lab needs and deliver the training at our Training Center or your workplace. MTS classroom training is available at Regional Training Centers located in the U.S., South Korea and Germany. All of our course offerings can be presented onsite at your facility.



MAINTENANCE & MONITORING

Making sure that equipment is ready for use when needed and test projects are completed on time without unexpected breakdowns are important aspects of test lab management. Based on decades of service experiences, MTS has a set of well-defined routine maintenance offerings tailored for specific systems and components, to help extend equipment life and provide you with confidence in your equipment. MTS has hydraulic fluid monitoring and maintenance programs that will help you optimize test system performance. Our condition based monitoring tools help you better understand equipment condition and anticipate potential issues before they become larger problems. And our real-time, remote monitoring software solution for the entire lab can be accessed anywhere by any Web-enabled device – smart phone, laptop or tablet, for increased visibility and productivity in the lab.

UPGRADE SOLUTIONS

As technology improves, an upgrade is often the most economical way of expanding your lab capabilities and extending the life of existing test equipment. MTS offers upgrades and replacements for all areas of your test system: mechanical components, controllers and software. Upgrade your old hydraulic power unit to a new, efficient one and save money. Our Software Support Plans make it easy to stay current with rapidly changing software technology. Within your contract period, you will automatically receive updates to all software covered in your contract.

Training

Course Selection

Factory Training Courses

Test System Operation & Application Theory Courses

	AMERICAS TRAINING CENTER	ASIA TRAINING CENTER	EUROPE TRAINING CENTER	SERVICE PRODUCT CODE
Damper Test System Operation	√			TRDTSOC
Durability Testing Technology	√	√		TRDTTC
Elastomer Testing on Controllers with MTS Series 793 Software	√	√	√	TRETC793C
Test Rig Design	√			TRTRDC

Software Operation Courses

	AMERICAS TRAINING CENTER	ASIA TRAINING CENTER	EUROPE TRAINING CENTER	SERVICE PRODUCT CODE
MTS Hardware Concepts & MTS Series 793 Software	√	√		TRCAHOMTWC
MTS Series 793 Advanced Software Operation			√	TR793ASOC
MTS Series 793 Software with MultiPurpose TestWare® (MPT) Test Design	√	√	√	TRCAMTSC
MTS Series 793 Software with MTS TestSuite™ Elite (mpe) Test Design	√	√	√	TRCAMTSMSC
MTS TestSuite Multipurpose Elite (mpe) Software	√	√	√	TRMTSMSC
MTS TestSuite TW Elite (TWE) Software	√	√	√	TRTSTWEC
MTS TestSuite mpe Advanced			√	TTRTSMPEADVC
MTS Acumen Operation with MTS TestSuite Multipurpose Elite Software	√			TTRACUMTSMPE
AeroPro Software Operation	√			TRAPOC
AeroPro Administrator	√		√	TRAPAC
Fatigue & Fracture with MTS TestSuite Software	√			TRFFSC
Component RPC Pro Software Operation	√	√	√	TRCFPSOC
RPC® Pro Software Operation	√	√	√	TRRPSOC
RPC Pro Advanced Software Operation	√	√	√	TRRPASOC
RPC Pro Fatigue Tools - Basic (add on course)	√	√	√	TRBRPPFATADDC
RPC Pro Fatigue Tools - Advanced (add on course)	√	√	√	TRARPPFATADDC

Training Alternatives*

MTS can provide onsite presentations of all course offerings. MTS also has the expert capability to develop custom courses on a broad range of test-related subjects not covered in the course schedule. (Please allow three months for onsite course scheduling).

ONSITE COURSE	SERVICE PRODUCT CODE
Per day onsite training - 1 to 8 students with 1 instructor	TTROD
Add per second instructor for class size above 8 - check availability for second instructor	TROAIC
Add per student above 8 (Class size limit of 15)	TROASC
REMOTE ONLINE COURSE	
Per hour online training - 1 to 8 students with 1 instructor	TTRRPH
Add per student above 8 (Class size limit of 15)	TTRRAS

Training/Consulting Package

DESCRIPTION	SERVICE PRODUCT CODE
MTS TestSuite MPE Training/Consulting Package	TTRTSMPETCP
MTS TestSuite TWE Training/Consulting Package (2 days software training plus 2 days consulting on custom test methods)	TTRTSTWETCP

Cancellation Fee

DESCRIPTION
If cancelled less than one week prior to course start

Price for a course held at MTS includes tuition, class material and lunch. MTS training center classroom size is limited, check for availability

*Price plus airfare and zone charges for instructor(s)

*For onsite courses with more than 8 students, contact MTS for second instructor availability

*Please allow 3 months for onsite course scheduling

We welcome your requests for custom classes or your language needs

Routine Maintenance

Hydraulic Power Unit

SERVICE PRODUCT CODE

HPU (500 hr)	MRMHPU1
HPU (1,000 hr)*	MRMHPU2
HPU (1,500 hr)	MRMHPU3
HPU (2,000 hr)*	MRMHPU4

Servohydraulic Load Frame

SERVICE PRODUCT CODE

SH Frame (500 hr)	MRMFRAME1
SH Frame (1,000 hr)	MRMFRAME2
SH Frame (1,500 hr)	MRMFRAME3
SH Frame (2,000 hr)*	MRMFRAME4

Static-Hydraulic Load Frame

SERVICE PRODUCT CODE

StH Criterion Frame (2,000 hr)*	MRMFRAMESTC
---------------------------------	-------------

Electromechanical Load Frame

SERVICE PRODUCT CODE

EM Frame (2,000 hr)	MRMFRAMEEM
EM Criterion Frame (2,000 hr)	MRMFRAMEEMC

Electromechanical Load Frame

SERVICE PRODUCT CODE

ED Acumen Frame (2,000 hr)	MRMACU1
----------------------------	---------

Controller

SERVICE PRODUCT CODE

Controller (2,000 hr)*	MRMCONTROL
------------------------	------------

Series 329 6DOF Road Simulator

SERVICE PRODUCT CODE

Per Corner (500 hr)	MRM3296D1
Per Corner (1,000 hr)	MRM3296D2
Per Corner (1,500 hr)	MRM3296D3
Per Corner (2,000 hr)*	MRM3296D4
Longitudinal Restraint Option 6DOF Per Corner (every 500 hrs)	MRM3296DLR
Static Support Option 6DOF Per Corner (every 500 hrs)	MRM3296DSS
XY Positioner Option 6DOF & 4DOF Per Corner (every 500 hrs)	MRM329XY
Brake Intensifier Option 6DOF & 4DOF Per System (every 2000 hrs)*	MRM329BI

Routine Maintenance Package pricing does not include materials. Each specific MTS equipment model may or may not require materials at each RM interval. Routine Maintenance Packages may not apply to all MTS equipment models. Please contact our service sales representative to verify equipment eligibility.

*Note: Filters and/or materials replaced during this service.

Routine Maintenance & Fluid Assessment Services

Tire Rolling Resistance - Hydraulic

SERVICE PRODUCT CODE

Tire Rolling Resistance - Hydraulic (500 hr)*	MRMRRH1
Tire Rolling Resistance - Hydraulic (1000 hr)*	MRMRRH2
Tire Rolling Resistance - Hydraulic (1500 hr)*	MRMRRH3
Tire Rolling Resistance - Hydraulic (2000 hr)*	MRMRRH4
Tire Rolling Resistance - Hydraulic Carriage Option (500 & 1500 hr)	MRMRRHO1
Tire Rolling Resistance - Hydraulic Carriage Option (1000 & 2000 hr)*	MRMRRHO2

Tire Rolling Resistance - Electric

SERVICE PRODUCT CODE

Tire Rolling Resistance - Electric (500 hr)*	MRMRRE1
Tire Rolling Resistance - Electric (1000 hr)*	MRMRRE2
Tire Rolling Resistance - Electric (1500 hr)*	MRMRRE3
Tire Rolling Resistance - Electric (2000 hr)*	MRMRRE4
Tire Rolling Resistance - Electric Carriage Option (500 & 1500 hr)	MRMRREO1
Tire Rolling Resistance - Electric Carriage Option (1000 & 2000 hr)*	MRMRREO2

Tire Tread Wear Simulation

SERVICE PRODUCT CODE

Tire Tread Wear - Hydraulic (500 hr)*	MRMTWH1
Tire Tread Wear - Hydraulic (1000 hr)*	MRMTWH2
Tire Tread Wear - Hydraulic (1500 hr)*	MRMTWH3
Tire Tread Wear - Hydraulic (2000 hr)*	MRMTWH4
Tire Tread Wear - Hydraulic Carriage Option (500 & 1500 hr)	MRMTWHO1
Tire Tread Wear - Hydraulic Carriage Option (1000 & 2000 hr)*	MRMTWHO2

Routine Maintenance Package pricing does not include materials. Each specific MTS equipment model may or may not require materials at each RM interval. Routine Maintenance Packages may not apply to all MTS equipment models. Please contact our service sales representative to verify equipment eligibility.
 *Note: Filters and/or materials replaced during this service.

MTS Fluid Care Program

SERVICE PRODUCT CODE

Program Management Per Sample Point	FFCPCONSULT
FSE Sampling Per Bottle (Use P/N 100-188-132)	
MTS Fluid Sample Tool (Use P/N 055-589-601)	

MTS Fluid Analysis

SERVICE PRODUCT CODE

D.I.Y. Sampling Per Bottle (Use P/N 100-030-731)	
FSE Sampling Per Bottle (Use P/N 100-188-132)	
MTS Sample Tool for D.I.Y. Sampling (Use P/N 055-589-601)	

Fluid Replacement Service*

SERVICE PRODUCT CODE

Replace hydraulic fluid* (systems <30 gpm)	MRMFLUIDCHG1
Replace hydraulic fluid* (systems 30 to 100 gpm)	MRMFLUIDCHG2

*Note: Fluid is not included and customer is responsible for disposal

Onsite Calibration Services

Calibration Standards Fee

SERVICE PRODUCT CODE

One charge per calibration incident, Load Cells 0-300 kN, LVDT, etc.
Load Cells 500-5000 kN and Torque*

CALSTDFEE

MTS Field Service is A2LA Accredited to perform on-site calibrations. Our Scope of Accreditation to ISO/IEC 17025-2005 covers many calibration types.

» Pricing is based on specimen and fixturing removed, leaving the system ready to be calibrated. Additional time required to prepare the system for calibration will be charged at current Field Service Engineer rates.

» Consult factory for models not listed.

» Calibrations are to current revision of applicable standard unless otherwise noted.

» ASTM E4 certification may not apply to electromechanical calibrations.

* This standards fee can vary by country.

Please contact your local MTS representative.

Load Cell to Relevant ASTM E4 or ISO 7500 Standards***

FIRST READOUT DEVICE

SERVICE PRODUCT CODE

Deadweight up to 100 lbs	CDWTO100
Axial 0 to 1 KIP Metric Equivalent 0 to 5 kN (Single range Acumen only/2 to 100% full-scale)	CT1KAA
Axial 0 to 5 KIP Metric Equivalent 0 to 25 kN** (Single range EM & Criterion/10 to 100% full-scale)	C5KAA
Axial 0 to 20 KIP Metric Equivalent 0 to 100 kN [up to 2 ranges (tension/compression)]	C20KAA
Axial 0 to 100 KIP Metric Equivalent 0 to 250 kN (tension/compression)]	C100KAA
Axial 0 to 240 KIP Metric Equivalent 0 to 500 kN [up to 2 ranges (tension/compression)]	C240KAA
Axial 0 to 500 KIP Metric Equivalent 0 to 1000 kN [up to 2 ranges (tension/compression)]	C500KAA
Axial 0 to 1000 KIP Metric Equivalent 0 to 5000 kN [up to 2 ranges (tension/compression)]	C1000KAA
Torque 0 to 60K inch-pound Metric Equivalent (0 to 6800 N.m) 1 range	C60KTA
Torque 0 to 200K inch-pound Metric Equivalent (0 to 22000 N.m) 1 range	C200KTA
Torque 0 to 750K inch-pound Metric Equivalent (0 to 85000 N.m) 1 range	C750KTA
ASTM E4 or ISO 7500*** readout verification, each additional readout device/additional bridges	CASTME4

** For multiple ranges, use the Axial 0 to 20 kip.

CRITERION AND EXCEED SYSTEMS (EM ONLY) - SINGLE RANGE

SERVICE PRODUCT CODE

Force Calibration, Deadweight up to 500 N EM only (Single range 0.5% / 10 to 100% full-scale)	CODW500NCEM
Axial Force Calibration 0 to 30 kN EM only (Single range 0.5% / 10 to 100% full-scale)	COA30KNCEM
Axial Force Calibration 0 to 100 kN EM only (Single range 0.5% / 10 to 100% full-scale)	COA100KNCEM
Axial Force Calibration 0 to 300 kN EM only (Single range 0.5% / 10 to 100% full-scale)	COA300KNCEM
Axial Force Calibration 0 to 600 kN EM only (Single range 0.5% / 10 to 100% full-scale)	COA600KNCEM

CRITERION AND EXCEED SYSTEMS (EM ONLY) - EXTEN RANGE

SERVICE PRODUCT CODE

Axial Force Calibration - Extended Range, 0 to 30 kN EM only (Extended range 0.5% / at points below 10% full-scale)	COA30KNCEMER
Axial Force Calibration - Extended Range, 0 to 100 kN EM only (Extended range 0.5% / at points below 10% full-scale)	COA100KNCEMER

*** Accreditation can vary by country. Please contact your local MTS representative to confirm.

MULTI-AXIS FORCE TRANSDUCER (LOAD CELL)

SERVICE PRODUCT CODE

Bi-axial (restraint) load cell	
Model 670.67B-03/04	CBIAXCAL1
Model 670.67B-10/11	CBIAXCAL2
Tri-axial 833 Elastomer load cell Utilizing hardware Matrix box (Matrix module and system controller required)	CTRIAXCAL1
Tri-axial 833 Elastomer load cell Utilizing SVW and calculated channels (System controller with conditioners required)	CTRIAXCAL2

Onsite Calibration Services

Displacement Transducers

TYPE	SERVICE PRODUCT CODE
Acumen Displacement Calibration	CDTACU
Servo-hydraulic displacement (LVDT)	CLVDT
Rotary position transducer	CRVDT
Delta P transducer	CDELTAPCAL
Electromechanical test system verification (includes both Speed and Displacement)	CSPDTMTS
Speed and Displacement Criterion and Exceed EM only (Single Range 0.5% / Displacement to ASTM E2309)	COSPDTCM

Reference AStM E2658 available upon request.

Clip-On Displacement Gages (COD)

Calibration of COD gages for use in compliance with ASTM** E399, E561, E1290, and E1820.

ASTM CLIP-ON DISPLACEMENT GAGE	SERVICE PRODUCT CODE
Calibration (Models Covered: 632.02, 632.03, 632.05, 632.06, & Equivalents)	CCODCAL

**ASTM E399 standard test method for plane-strain fracture toughness of metallic materials.

ASTM E561 standard practice for R-curve determination.

ASTM E1820 standard test method for measurement of fracture toughness.

Extensometer Calibration to Relevant ASTM E83/ISO 9513 Standards

	SERVICE PRODUCT CODE
Axial Extensometers (per extensometer) (Models Covered: 632.11, 632.12, 632.24, 632.25, 632.27, 632.29, 632.31, 632.90, 632.94, 634.12, 634.25, 634.28 & Equivalents)	CAXIALEXTCAL
Bi-Axial Extensometers (Models Covered: 632.8X & Equivalents)	CBIAXLEXTCAL
Diametral Extensometers (Models Covered: 632.18, 632.19, 632.20 & Equivalents)	CDIAEXTCAL
Averaging Axial Extensometers (Models Covered: 632.17 & Equivalents)	CAVGEXTCAL
Cross Sectional Strain Extensometers (Models Covered: 632.23 & Equivalents)	CCSSTREXTCAL
High Temperature Axial Extensometers (Models Covered: 632.4X, 632.5X, 632.6X & Equivalents)	CHITEMPEXTCL
Immersible Extensometers (Models Covered: 632.70 & Equivalents)	CIMMEXTCAL
Small Gage Length Extensometers (Models Covered: 632.13, 632.26 & Equivalents)	CSGLEXTCAL
Circumferential Extensometers (Models Covered: 632.92 & Equivalents)	CCIREXTCAL
High Elongation Extensometers (Models Covered: DXL-Ext & Equivalents)	CHIELEXTCAL
Biomedical Extensometers (Models Covered: 632.32 & Equivalents)	CBIMEDEXTCAL

Other Calibration Items

	SERVICE PRODUCT CODE
Non-MTS Calibration Custom Quote	
LDH Machine Calibrations 866.5X	CLDH866.5X
LX Cal Extensometer	CLXEXTCAL
AVX Cal Extensometer	CAVXEXTCAL
AVX Cal Extensometer Extra Working Distance	CAVXEXCALELEWD
FlexDAC Calibration ONSITE (1st chassis)	CFLEXDACCAL
FlexDAC Calibration ONSITE (each additional chassis)	CFLEXDACCALADD
FlexTest Calibration ONSITE 494.16 / 494.21 / 494.25 / 494.45 / 494.46 (Single AC/DC conditioner or A/D or D/A)	CFLXTEST494CAL
FlexTest Calibration ONSITE 494.26 (2 AC/DC conditioners)	CALFLXTST49426
FlexTest Calibration ONSITE 494.41 / 494.42 / 494.43 (All FT40)	CALFLXTST49440
FlexTest Calibration AC Option (per board or DUC, 494.26 = 2 DUCs)	CALFLXTST494AC

Static Alignment Verification to Relevant ASTM E1012/NASM 1312B

Dynamic Force Verification to Relevant ASTM E467/NASM 1312B

Onsite Calibration Services

Time and Material Calibrations

Some systems must be calibrated under a time and material method. An example would be a system without a load unit that requires the use of a fabricated reaction base. Standard load unit systems will always be calibrated under the fixed fee method. When a system is calibrated using the time and material method, these charges apply:

- » Time-and-material Standards fee for each standard used
- » Each additional load cell using the same load cell standard
- » Zone charges (based on applicable zone charge)
- » Labor (based on applicable labor rate)

Flat-Trac® Calibrations

Flat-Trac calibrations should be performed at least annually, or when any component in the chain of mechanical or electrical components comprising a selected software channel is replaced or repaired.

	New Matrix Generation Calibration	Current Matrix Calibration
	SERVICE PRODUCT CODE	SERVICE PRODUCT CODE
Free rolling calibration	CNMFTCFR	CCMFTCFR
Each additional free rolling matrix	CNMFTCFRA	CCMFTCFRA
Adder for spindle drive calibration	CNMFTCSD	CCMFTCSD
Each additional spindle drive matrix	CNMFTCSDA	CCMFTCSDA
Slip angle calibration (stand-alone)		CCMFTCSA
Camber angle calibration (stand-alone)		CCMFTCCA

New Matrix Generation pricing includes the following activities when applicable:

- » A/D and D/A calibrations
- » Camber angle calibration
- » Slip angle calibration
- » Calibrating the calibration standard (Flat-Trac II)
- » Force and moment calibration
- » Belt speed calibration
- » Generation of calibration matrix
- » Pre- and post-calibration data
- » Verify inertial compensation (Flat-Trac III)
- » Calibrate inflation pressure
- » Calibrate ambient and tire temperatures
- » Spindle speed calibration
- » Shunt calibration
- » Loaded radius calibration
- » Verification of matrix through calibration check

Onsite Calibration Services

Calibration for 329 Passenger Road Simulators

LOAD CELLS AND DISPLACEMENT TRANSDUCER CALIBRATION		Without brake SERVICE PRODUCT CODE	With brake SERVICE PRODUCT CODE
Corner		CCORNERA	CCORNERAB
1/2 Car		CHALFCARA	CHALFCARAB
Full Car		CFULLCARA	CFULLCARAB
LOAD CELL CALIBRATION (ONLY)			
Corner		CCORNERL	CCORNERLB
1/2 Car		CHALFCARL	CHALFCARLB
Full Car		CFULLCARL	CFULLCARLB
DISPLACEMENT TRANSDUCER CALIBRATION (ONLY)			
Corner		CCORNERD	CCORNERDB
1/2 Car		CHALFCARD	CHALFCARDB
Full Car		CFULLCARD	CFULLCARDB

Displacement Transducer Calibration (only)

LOAD CELLS AND DISPLACEMENT TRANSDUCER CALIBRATION		Without brake SERVICE PRODUCT CODE	With brake SERVICE PRODUCT CODE
Corner		CLTCORNERA	CLTCORNERAB
1/2 Truck		CLTHALFCARA	CLTHALFCARAB
Full Truck		CLTFULLCARA	CLTFULLCARAB
LOAD CELL CALIBRATION (ONLY)			
Corner		CLTCORNERL	CLTCORNERLB
1/2 Truck		CLTHALFCARL	CLTHALFCARLB
Full Truck		CLTFULLCARL	CLTFULLCARLB
DISPLACEMENT TRANSDUCER CALIBRATION (ONLY)			
Corner		CLTCORNERD	CLTCORNERDB
1/2 Truck		CLTHALFCARD	CLTHALFCARDB
Full Truck		CLTFULLCARD	CLTFULLCARDB

Note: zone charges are not included in pricing.

Onsite Calibration Services

Calibration of Multi-Axial Simulation Table (MAST™) Systems

6 CHANNEL MAST 323	END-TO-END-CALIBRATION	SERVICE PRODUCT CODE
Package 1	6 accelerometers	323C6MSTACC
Package 2	6 accelerometers and 6 LVDT	323C6MSTACCL
Package 3	6 accelerometers, 6 LVDT, and ASC D/A-A/D	323C6MSTACCLAD

6 CHANNEL MAST 353	END-TO-END-CALIBRATION	SERVICE PRODUCT CODE
Package 1	6 accelerometers	353C6MSTACC
Package 2	6 accelerometers and 6 LVDT	353C6MSTACCL
Package 3	6 accelerometers, 6 LVDT, and ASC D/A-A/D	353C6MSTACCLAD

Note:

- Prices include calibration standard fee but travel expenses (zone, airfare, etc.) are not included.
- 353 MAST Packages 2 & 3 require 2 FSEs - Prices above include labor for both FSEs, but not travel.

Onsite Accelerometer Calibration

TYPE	SERVICE PRODUCT CODE
Single Axis- Transducer only	CACC
Tri-Axis - Transducer only	CTRIAC
Single Axis - Transducer with conditioner	CACCCOND
Tri-Axis - Transducer with conditioner	CTRIACC

Static Alignment Verification

TYPE	SERVICE PRODUCT CODE
Static twelve gage alignment STATIC12	CSTATIC12

Perform to ASTM E1012

Dynamic Force Verification

TYPE	SERVICE PRODUCT CODE
Dynamic Force Verification	CDYNAMICCAL

Perform to ASTM E467

Factory Calibration Services

Factory Calibration

LOAD CELL

Calibration up to 220,000 lbf (As Found / As Left Data Provided)
Single Range / Single Readout Device*

	SERVICE PRODUCT CODE
Axial up to 22 kip (Metric Equivalent–up to 100 kN)	FFCFTA1
Additional Range / Readout Device	FFCADDR
Axial >22 kip to 240 kip(Metric Equivalent- >100 kN to 1000 kN)	FFCFTA2
Additional Range / Readout Device	FFCADDR

MULTI-AXIS LOAD CELL

Bi-axial (restraint) load cell	
Model 670.67B-03/04	FFCBARFT1
Model 670.67B-10/11	FFCBARFT2
Tri-axial 833 Elastomer load cell Utilizing hardware Matrix box (Matrix module and system controller required)	FFCTAFT1
Tri-axial 833 Elastomer load cell Utilizing SW and calculated channels (System controller with conditioners required)	FFCTAFT2

TORQUE TRANSDUCER

Calibration up to 12,000 lbf-in (As Found / As Left Data Provided)

Torque up to 12,000 lbf Metric Equivalent- up to 1300 Nm	
Single Range / Single Readout Device	FFCTCT1
Additional Range / Readout Device	FFCADDR

COD GAUGES

Model 632.02 & 632.03 Clip-on displacement gage	FFCCOD1
Model 632.06 1 Arm Bandit displacement gage	FFCCOD2

EXTENSOMETER**

Single Range / Single Readout Device* (As Found / As Left Data Provided)

Axial Extensometers (per extensometer) (Models Covered: 632.11, 632.12, 632.24, 632.25, 632.27, 632.29, 632.90, 632.94, 634.12, 634.25, 634.28, 634.31 & Equivalents)	FFCEA
Bi-Axial Extensometers (Models Covered: 632.8X & Equivalents)	FFCEBA
Diametral Extensometers (Models Covered: 632.18, 632.19, 632.20 & Equivalents)	FFCEDE
Averaging Axial Extensometers (Models Covered: 632.17 & Equivalents)	FFCEAA
Cross Sectional Strain Extensometers (Models Covered: 632.23 & Equivalents)	FFCECSS
High Temperature Axial Extensometers (Models Covered: 632.4X, 632.5X, 632.6X & Equivalents)	FFCEHTA
Immersible Extensometers (Models Covered: 632.70 & Equivalents)	FFCEI
Small Gage Length Extensometers (Models Covered: 632.13, 632.26 & Equivalents)	FFCESGL
Circumferential Extensometers (Models Covered: 632.92 & Equivalents)	FFCEC
High Elongation Extensometers (Models Covered: DXL-Ext & Equivalents)	FFCEHE
Biomedical Extensometers (Models Covered: 632.32 & Equivalents)	FFCEB
Additional Range/Readout Device	FFCADDR

* Single Range (Minimum of 9 data points from 2% to 100%. Applies to 493 and 494 product lines except 493.21 and 494.21)

** Not relevant to ASTM E83 or ISO9513 Standards

SWIFT and SWIFT EVO System Calibrations

SWIFT CALIBRATION PRICES ARE PER CALIBRATION*

MODEL	SERVICE PRODUCT CODE
SWIFT 10, 20, 30, & 40 transducer*	FFCSW1
SWIFT 45 & 50	FFCSW2

* PRICE is for Aluminum, Titanium or Stainless Steel units

* PRICE includes Spinning and Non-spinning accessories

* PRICES do not include rims and /or hub adapters

* PRICES are subject to change without notice

* SWIFT Loaner Program is no longer available



Metrology Calibration Services

Metrology Laboratory Services

LOAD CELLS

ASTM E74, EN 10002-3, AND ISO 376 CALIBRATIONS

SERVICE PRODUCT CODE

Load cell - unidirectional up to 25 kip	MMTFTU25
Load cell - bidirectional up to 25 kip	MMTFTB25
Load cell - unidirectional up to 100 kip	MMTFTU100
Load cell - bidirectional up to 100 kip	MMTFTB100
Load cell - unidirectional up to 240 kip	MMTFTU240
Load cell - bidirectional up to 240 kip	MMTFTB240
Load cell with dual bridge – add	MMTFTDBA
Various other Force related M&TE	MMTMTEF

TORQUE TRANSDUCER

Flange style transducers only. Additional costs may apply due to required fixturing

Torque transducers up to 12,000 lbf-in (Bidirectional, ASTM E2428)	MMTTQ1
Torque transducers up to 12,000 lbf-in (Unidirectional, ASTM E2428)	MMTTQ2
Torque transducers up to 12,000 lbf-in (Unidirectional, non-ASTM)	MMTTQ3
Torque transducers up to 100,000 lbf-in (Bidirectional, ASTM E2428)	MMTTQ4
Torque transducers up to 100,000 lbf-in (Unidirectional, ASTM E2428)	MMTTQ5
Various other Torque related M&TE	MMTMTEF

Notes:

Unidirectional Force: Compression and/or tension ascending only data points.

Bidirectional Force: Compression and/or tension ascending and descending data points.

Unidirectional Torque: Clockwise and/or counter-clockwise ascending only data points.

Bidirectional Torque: Clockwise and/or counter-clockwise ascending and descending data points.

Non-ASTM: Performance to manufacturers specifications.

VIBRATION

Accelerometer - single axis	MMTVIB1
Accelerometer - triaxial	MMTVIB2
Charge amplifiers/Signal conditioners	MMTVIB4
Various other Vibration related M&TE	MMTMTEV

ELECTRICAL INSTRUMENTS/DEVICES

Provide manufacturer and model number to confirm availability	MMTEID
---	--------

DIMENSIONAL AND MECHANICAL INSTRUMENTS/DEVICES

Provide manufacturer and model number to confirm availability for Rotary Encoders and dimensional related M&TE	MMTDMID
--	---------

TEMPERATURE AND HUMIDITY

Provide manufacturer and model confirm availability	MMTTH
---	-------

TIME AND FREQUENCY

Provide manufacturer and model number to confirm availability	MMTTF
---	-------

OTHER MTS PRODUCT CALIBRATIONS

MTS LX series laser extensometers (single standard calibration)	MMTLE1
709 alignment system (data acquisition and signal conditioning electronics)	MMTAS2
Extensometer Calibrators	MMTEC
Additional ranges or custom calibration	MMTADDRCC

Actuators Eligible for Remanufacture using *SureCoat Technology*

These may qualify for 10 day actuator quick-turn program, contact MTS for details

Model 244.xx

MODEL NUMBER	STROKE LENGTH	DIRECT REPLACEMENT NEW PRODUCTION ACTUATORS STILL AVAILABLE	SERVICE PRODUCT CODE
244.11	4-inch	Yes	CYLRM24411X04
244.11	6-inch	Yes	CYLRM24411X06
244.11	8-inch	Yes	CYLRM24411X08
244.11	10-inch	Yes	CYLRM24411X10
244.11	12-inch	Yes	CYLRM24411X12
244.11	14-inch	Yes	CYLRM24411X14
244.11	16-inch	Yes	CYLRM24411X16
244.11	18-inch	Yes	CYLRM24411X18
244.11	20-inch	Yes	CYLRM24411X20
244.12	4-inch	Yes	CYLRM24412X04
244.12	6-inch	Yes	CYLRM24412X06
244.12	8-inch	Yes	CYLRM24412X08
244.12	10-inch	Yes	CYLRM24412X10
244.12	12-inch	Yes	CYLRM24412X12
244.12	14-inch	Yes	CYLRM24412X14
244.12	16-inch	Yes	CYLRM24412X16
244.12	18-inch	Yes	CYLRM24412X18
244.12	20-inch	Yes	CYLRM24412X20
244.20	4-inch	Yes	CYLRM24420X04
244.20	6-inch	Yes	CYLRM24420X06
244.20	8-inch	Yes	CYLRM24420X08
244.20	10-inch	Yes	CYLRM24420X10
244.20	12-inch	Yes	CYLRM24420X12
244.20	14-inch	Yes	CYLRM24420X14
244.20	16-inch	Yes	CYLRM24420X16
244.20	18-inch	Yes	CYLRM24420X18
244.20	20-inch	Yes	CYLRM24420X20
244.21	4-inch	Yes	CYLRM24421X04
244.21	6-inch	Yes	CYLRM24421X06
244.21	8-inch	Yes	CYLRM24421X08
244.21	10-inch	Yes	CYLRM24421X10
244.21	12-inch	Yes	CYLRM24421X12
244.21	14-inch	Yes	CYLRM24421X14
244.21	16-inch	Yes	CYLRM24421X16
244.21	18-inch	Yes	CYLRM24421X18
244.21	20-inch	Yes	CYLRM24421X20
244.22	4-inch	Yes	CYLRM24422X04
244.22	6-inch	Yes	CYLRM24422X06
244.22	8-inch	Yes	CYLRM24422X08
244.22	10-inch	Yes	CYLRM24422X10
244.22	12-inch	Yes	CYLRM24422X12
244.22	14-inch	Yes	CYLRM24422X14
244.22	16-inch	Yes	CYLRM24422X16
244.22	18-inch	Yes	CYLRM24422X18
244.22	20-inch	Yes	CYLRM24422X20

(continued...)

Actuators Eligible for Remanufacture using *SureCoat Technology*

These may qualify for 10 day actuator quick-turn program, contact MTS for details

Model 244.xx (continued)

MODEL NUMBER	STROKE LENGTH	DIRECT REPLACEMENT NEW PRODUCTION ACTUATORS STILL AVAILABLE	SERVICE PRODUCT CODE
244.23	4-inch	Yes	CYLRM24423X04
244.23	6-inch	Yes	CYLRM24423X06
244.23	8-inch	Yes	CYLRM24423X08
244.23	10-inch	Yes	CYLRM24423X10
244.23	12-inch	Yes	CYLRM24423X12
244.23	14-inch	Yes	CYLRM24423X14
244.23	16-inch	Yes	CYLRM24423X16
244.23	18-inch	Yes	CYLRM24423X18
244.23	20-inch	Yes	CYLRM24423X20
244.31	4-inch	Yes	CYLRM24431X04
244.31	6-inch	Yes	CYLRM24431X06
244.31	8-inch	Yes	CYLRM24431X08
244.31	10-inch	Yes	CYLRM24431X10
244.31	12-inch	Yes	CYLRM24431X12
244.31	14-inch	Yes	CYLRM24431X14
244.31	16-inch	Yes	CYLRM24431X16
244.31	18-inch	Yes	CYLRM24431X18
244.31	20-inch	Yes	CYLRM24431X20
244.41	4-inch	Yes	CYLRM24441X04
244.41	6-inch	Yes	CYLRM24441X06
244.41	8-inch	Yes	CYLRM24441X08
244.41	10-inch	Yes	CYLRM24441X10
244.41	12-inch	Yes	CYLRM24441X12
244.41	14-inch	Yes	CYLRM24441X14
244.41	16-inch	Yes	CYLRM24441X16
244.41	18-inch	Yes	CYLRM24441X18
244.41	20-inch	Yes	CYLRM24441X20

For all other non-standard 244 actuator remanufactures, please contact MTS

Loadframes Eligible for Remanufacture using *SureCoat Technology*

Model 318

Standard, non-special models

MODEL NUMBER	CAPACITY	STROKE LENGTH	THREADS	PART KIT NUMBER	SERVICE PRODUCT CODE
318.10	100 kn	4-inch	English	100-484-438	CYR31810X4E100
318.10	100 kn	4-inch	Metric	100-484-439	CYR31810X4M100
318.10	100 kn	6-inch	English	100-484-440	CYR31810X6E100
318.10	100 kn	6-inch	Metric	100-484-441	CYR31810X6M100
318.10	100 kn	10-inch	English	100-484-442	CYR3181010E100
318.10	100 kn	10-inch	Metric	100-484-443	CYR3181010M100
318.25	100 kn	4-inch	English	100-484-240	CYR31825X4E100
318.25	100 kn	4-inch	Metric	100-484-241	CYR31825X4M100
318.25	100 kn	6-inch	English	100-484-242	CYR31825X6E100
318.25	100 kn	6-inch	Metric	100-484-243	CYR31825X6M100
318.25	100 kn	10-inch	English	100-484-264	CYR3182510E100
318.25	100 kn	10-inch	Metric	100-484-265	CYR3182510M100
318.25	250 kn	4-inch	English	100-484-266	CYR31825X4E250
318.25	250 kn	4-inch	Metric	100-484-267	CYR31825X4M250
318.25	250 kn	6-inch	English	100-484-268	CYR31825X6E250
318.25	250 kn	6-inch	Metric	100-484-269	CYR31825X6M250
318.25	250 kn	10-inch	English	100-484-270	CYR3182510E250
318.25	250 kn	10-inch	Metric	100-484-271	CYR3182510M250

For all crosshead mount and other non standard 318 remanufactures, please consult MTS.

Pricing is for standard parts only, labor and additional parts are extra.

Loadframes Eligible for Remanufacture using *SureCoat Technology*

Model 370

Standard, non-special models

MODEL NUMBER	CAPACITY	STROKE LENGTH	ACTUATOR RATING	BEARING TYPE	PART KIT NUMBER
370.10	100 kn	4 -inch	15 kn	Stepped	058-565-009
370.10	100 kn	4-inch	15 kn	Hydrostatic	058-565-010
370.10	100 kn	6-inch	15 kn	Stepped	058-565-011
370.10	100 kn	6-inch	15 kn	Hydrostatic	058-565-012
370.10	100 kn	10-inch	15 kn	Stepped	058-565-013
370.10	100 kn	10-inch	15 kn	Hydrostatic	058-565-014
370.10	100 kn	4-inch	25 kn	Stepped	058-565-015
370.10	100 kn	4-inch	25 kn	Hydrostatic	058-565-016
370.10	100 kn	6-inch	25 kn	Stepped	058-565-017
370.10	100 kn	6-inch	25 kn	Hydrostatic	058-565-018
370.10	100 kn	10-inch	25 kn	Stepped	058-565-019
370.10	100 kn	10-inch	25 kn	Hydrostatic	058-565-020
370.10	100 kn	4-inch	50 kn	Stepped	058-565-021
370.10	100 kn	4-inch	50 kn	Hydrostatic	058-565-022
370.10	100 kn	6-inch	50 kn	Stepped	058-565-023
370.10	100 kn	6-inch	50 kn	Hydrostatic	058-565-024
370.10	100 kn	10-inch	50 kn	Stepped	058-565-025
370.10	100 kn	10-inch	50 kn	Hydrostatic	058-565-026
370.10	100 kn	10-inch	67 kn	Stepped	058-565-027
370.10	100 kn	10-inch	67 kn	Hydrostatic	058-565-028
370.10	100 kn	4-inch	100 kn	Stepped	058-565-029
370.10	100 kn	4-inch	100 kn	Hydrostatic	058-565-030
370.10	100 kn	4-inch	100 kn	Stepped	058-565-031
370.10	100 kn	4-inch	100 kn	Hydrostatic	058-565-032
370.10/25	100/250 kn	6-inch	100 kn	Stepped	058-565-033
370.10/25	100/250 kn	6-inch	100 kn	Hydrostatic	058-565-034
370.10/25	100/250 kn	10-inch	100 kn	Stepped	058-565-035
370.10/25	100/250 kn	10-inch	100 kn	Hydrostatic	058-565-036
370.25	250 kn	6-inch	250 kn	Stepped	058-565-037
370.25	250 kn	6-inch	250 kn	Hydrostatic	058-565-038
370.25	250 kn	10-inch	250 kn	Stepped	058-565-039
370.25	250 kn	10-inch	250 kn	Hydrostatic	058-565-040
370.50	500 kn	6-inch	500 kn	Stepped	058-565-041
370.50	500 kn	6-inch	500 kn	Hydrostatic	058-565-042

370 Tabletop Load Frames

MODEL NUMBER	CAPACITY	STROKE LENGTH	ACTUATOR RATING	BEARING TYPE	PART KIT NUMBER
370.02	25 kn	4-inch	15 kn	Stepped	058-565-001
370.02	25 kn	4-inch	15 kn	Stepped	058-565-002
370.02	25 kn	6-inch	15 kn	Stepped	058-565-003
370.02	25 kn	6-inch	15 kn	Stepped	058-565-004
370.02	25 kn	6-inch	25 kn	Stepped	058-565-005
370.02	25 kn	6-inch	25 kn	Stepped	058-565-006
370.02	25 kn	4-inch	25 kn	Stepped	058-565-007
370.02	25 kn	4-inch	25 kn	Stepped	058-565-008

Pricing is for standard parts only, labor and additional parts are extra.



Standard Chrome Piston Actuator *Remanufacture Programs*

Model 201.xx

MODEL NUMBER	STROKE LENGTH	DIRECT REPLACEMENT NEW PRODUCTION ACTUATORS STILL AVAILABLE	SERVICE PRODUCT CODE
201.10	10-inch	Yes	CYLRM20110X10
201.10	20-inch	Yes	CYLRM20110X20
201.10	30-inch	Yes	CYLRM20110X30
201.10	40-inch	Yes	CYLRM20110X40
201.11	10-inch	Yes	CYLRM20111X10
201.11	20-inch	Yes	CYLRM20111X20
201.11	30-inch	Yes	CYLRM20111X30
201.11	40-inch	Yes	CYLRM20111X40
201.17	10-inch	Yes	CYLRM20117X10
201.17	20-inch	Yes	CYLRM20117X20
201.17	30-inch	Yes	CYLRM20117X30
201.17	40-inch	Yes	CYLRM20117X40
201.20	10-inch	Yes	CYLRM20120X10
201.20	20-inch	Yes	CYLRM20120X20
201.20	30-inch	Yes	CYLRM20120X30
201.20	40-inch	Yes	CYLRM20120X40
201.25	10-inch	Yes	CYLRM20125X10
201.25	20-inch	Yes	CYLRM20125X20
201.25	30-inch	Yes	CYLRM20125X30
201.25	40-inch	Yes	CYLRM20125X40
201.30	10-inch	Yes	CYLRM20130X10
201.30	20-inch	Yes	CYLRM20130X20
201.30	30-inch	Yes	CYLRM20130X30
201.30	40-inch	Yes	CYLRM20130X40
201.35	10-inch	Yes	CYLRM20135X10
201.35	20-inch	Yes	CYLRM20135X20
201.35	30-inch	Yes	CYLRM20135X30
201.35	40-inch	Yes	CYLRM20135X40
201.40	10-inch	Yes	CYLRM20140X10
201.40	20-inch	Yes	CYLRM20140X20
201.40	30-inch	Yes	CYLRM20140X30
201.40	40-inch	Yes	CYLRM20140X40
201.45	10-inch	Yes	CYLRM20145X10
201.45	20-inch	Yes	CYLRM20145X20
201.45	30-inch	Yes	CYLRM20145X30
201.45	40-inch	Yes	CYLRM20145X40
201.50	10-inch	Yes	CYLRM20150X10
201.50	20-inch	Yes	CYLRM20150X20
201.50	30-inch	Yes	CYLRM20150X30
201.50	40-inch	Yes	CYLRM20150X40
201.60	10-inch	Yes	CYLRM20160X10
201.60	20-inch	Yes	CYLRM20160X20
201.60	30-inch	Yes	CYLRM20160X30
201.60	40-inch	Yes	CYLRM20160X40
201.70	10-inch	Yes	CYLRM20170X10
201.70	20-inch	Yes	CYLRM20170X20
201.70	30-inch	Yes	CYLRM20170X30
201.70	40-inch	Yes	CYLRM20170X40

Standard Chrome Piston Actuator *Remanufacture Programs*

Model 204.xx

MODEL NUMBER	RATING	DIRECT REPLACEMENT NEW PRODUCTION ACTUATORS STILL AVAILABLE	SERVICE PRODUCT CODE
204.08	1.1 kip	No	RAA20408
204.09	2.2 kip	No	RAA20409
204.11	2.5 kip	No	RAA20411
204.12	4.0 kip	No	RAA20412
204.13	6.0 kip	No	RAA20413
204.21	6.0 kip	No	RAA20421
204.22	9.0 kip	No	RAA20422
204.23	12 kip	No	RAA20423
204.24	20 kip	No	RAA20424
204.25	24 kip	No	RAA20425
204.26	35 kip	No	RAA20426
204.31	50 kip	No	RAA20431
204.32	70 kip	No	RAA20432
204.41	100 kip	No	RAA20441
204.42	150 kip	No	RAA20442
204.51	3.3 kip	No	RAA20451
204.52	5.5 kip	No	RAA20452
204.61	11 kip	No	RAA20461
204.62	15 kip	No	RAA20462
204.63	22 kip	No	RAA20463
204.64	35 kip	No	RAA20464
204.71	55 kip	No	RAA20471
204.72	77 kip	No	RAA20472
204.81	110 kip	No	RAA20481
204.82	165 kip	No	RAA20482
204.91	220 kip	No	RAA20491

Note: These units are no longer in production; exact replacement of these units is not possible.
MTS does carry functional replacements.

Model 242.xx

MODEL NUMBER	STROKE LENGTH	DIRECT REPLACEMENT NEW PRODUCTION ACTUATORS STILL AVAILABLE	SERVICE PRODUCT CODE
242.00	2-inch	Yes	CYLRM24200X02
242.00	4-inch	Yes	CYLRM24200X04
242.00	6-inch	Yes	CYLRM24200X06
242.01	2-inch	Yes	CYLRM24201X02
242.01	4-inch	Yes	CYLRM24201X04
242.01	6-inch	Yes	CYLRM24201X06
242.02	2-inch	Yes	CYLRM24202X02
242.02	4-inch	Yes	CYLRM24202X04
242.02	6-inch	Yes	CYLRM24202X06
242.03	2-inch	Yes	CYLRM24203X02
242.03	4-inch	Yes	CYLRM24203X04
242.03	6-inch	Yes	CYLRM24203X06

Standard Chrome Piston Actuator *Remanufacture Programs*

Model 248.xx

MODEL NUMBER	STROKE LENGTH	DIRECT REPLACEMENT NEW PRODUCTION ACTUATORS STILL AVAILABLE	SERVICE PRODUCT CODE
248.01	2-inch	Yes	CYLRM24801X02
248.01	4-inch	Yes	CYLRM24801X04
248.01	6-inch	Yes	CYLRM24801X06
248.01	10-inch	Yes	CYLRM24801X10
248.02	2-inch	Yes	CYLRM24802X02
248.02	4-inch	Yes	CYLRM24802X04
248.02	6-inch	Yes	CYLRM24802X06
248.02	10-inch	Yes	CYLRM24802X10
248.03	2-inch	Yes	CYLRM24803X02
248.03	4-inch	Yes	CYLRM24803X04
248.03	6-inch	Yes	CYLRM24803X06
248.03	10-inch	Yes	CYLRM24803X10
248.04	2-inch	Yes	CYLRM24804X02
248.04	4-inch	Yes	CYLRM24804X04
248.04	6-inch	Yes	CYLRM24804X06
248.04	10-inch	Yes	CYLRM24804X10
248.05	2-inch	Yes	CYLRM24805X02
248.05	4-inch	Yes	CYLRM24805X04
248.05	6-inch	Yes	CYLRM24805X06
248.05	10-inch	Yes	CYLRM24805X10
248.11	2-inch	Yes	CYLRM24811X02
248.11	4-inch	Yes	CYLRM24811X04
248.11	6-inch	Yes	CYLRM24811X06
248.11	10-inch	Yes	CYLRM24811X10
248.12	2-inch	Yes	CYLRM24812X02
248.12	4-inch	Yes	CYLRM24812X04
248.12	6-inch	Yes	CYLRM24812X06
248.12	10-inch	Yes	CYLRM24812X10

Software Support Plans

The MTS Software Support Plan is a simple, cost-effective way to keep your organization's vital test systems running smoothly.

A Software Support Plan provides a variety of benefits for MTS software customers;

- » Resolve software issues quickly with priority access to MTS technical support via phone, email, or remote login.
- » Training provides greater knowledge and productivity for your test operations staff.*

- » Cloud-based Echo Test Monitoring; lab personnel can view real-time test system status messages from any web-enabled device, and set up alerts to notify the instant a test requires attention or reaches a defined status.
- » Keep your MTS software up-to-date and performing optimally, with updates downloadable through the MTS Echo Portal

**RPC/cRPC and AeroPro customers get tuition-free training for 1 seat in up to 2 courses per 12 month increment of contract period.*

TestSuite customers get access to online technical training videos, posted monthly and available on the Echo Portal.

When SSP is purchased, it must be purchased for all licenses and options purchased. Partial coverage is not permitted.

Application Software

SOFTWARE	MODULE	SERVICE PRODUCT CODE
TestSuite		
	Multipurpose Elite	MESTSME
	TW Elite Servohydraulic	MESTSTWESH
	TW Elite EM	MESTSTWECI
Elastomer		
	Elastomer Software	MES793EBE
Damper Test		
	Damper Software	MESDT
Wind/Civil Engineering Materials Test		
	Wind/Civil Materials Test Software (includes 793.61/62/63/64/66/67)	MES793GCEPM
AeroPro (793 controller system software included)		
	AeroPro Control	MESSAFTPC
	AeroPro Calculations in the Loop	MESSAFTAPCL
	Digital Data Interface (DDI)	MESAERODDICNCT
	AeroPro API (SDK)	MESSAFTAPAPI
	Data Reprocessor	MESSAFTDR
	C3 (Cross Coupling Compensation)	MESSAFTC3
	AeroPro Signal Based Command	MESSAFSBC
	AeroPro Mode Switching	MESSAFMS
	AeroPro Redundant Load Checker	MESSAFRLC
	793 DDI Access	MESSAF793DDI
	Position Display	MESSAFPD
	Trend Monitoring	MESSAFTTM
	Data Display	MESSAFDD
	Peak Valley Recorder	MESSAFPVR

Note: 793 controller system software SSP AND ALL 793 application software SSP must be purchased together for the SSP to be valid. Customer must have/ purchase the latest version of software in order to purchase SSP. An upgrade in software may require a FlexTest controller upgrade.

(continued...)

Software Support Plans

Application Software (continued)

SOFTWARE	MODULE	SERVICE PRODUCT CODE
cRPC Pro	GROUND VEHICLE DURABILITY TESTING & ROAD SIMULATION SOFTWARE	
	Pro Drive File Development and Test Package	MESCRPCPDFDTP
	Drive File Development Package	MESCRPCDFDP
	Simulation Package	MESCRPCSP
	Durability Test Package including one monitoring option)	MESCRPCDTP
	Data Validation Package	MESCRPCDVP
	Edit and Analysis Package	MESCRPCEAP
	Options	
	Time History Fatigue	MESCRPCTHF
	Fatigue Sensitive Editing	MESCRPCFSE
	Histogram Analysis	MESCRPCHA
	Component Test Generation	MESCRPCCTG
	Simultaneous Excitation	MESCRPCSE
	FRF Diagnostics	MESCRPCFRFD
	Vibration Option	MESCRPCVO
	Matrix Editing	MESCRPCME
	Advanced Editing	MESCRPCAE
	Advanced Analysis	MESCRPCAA
	Region Analysis	MESCRPCRA
	Batch Processing	MESCRPCBP
	Adaptive Inverse Modeling (Turbo RPC)	MESCRPCAIM
	Test Monitoring Option	MESCRPCTMO
	Additional User	MESCRPCAU
	Ride Comfort	MESCRPCRC
	Trigger Option	MESCRPCTO
	Matlab Interface	MESCRPCMI
	Modulation Analysis	MESCRPCMA
	Virtual Test Option	MESCRPCVTO
	Synthetic Road Generation	MESCRPCSRG
	Durability Test - w/o Monitor Option	MESCRPCDT
	Proving Ground Schedule Builder	MESCRPCPSB
	MDA Sequence Builder	MESCRPCMDA
First Road		
	First Road Software	MESLSAFR
STEX Pro		
	STEX PRO	MESSTEXPP

Software Support Plans

Application Software *(continued)*

SOFTWARE	MODULE	SERVICE PRODUCT CODE
RPC Pro		
	Project Manager	MESRPCPM
	Setup	MESRPCS
	Acquire	MESRPCA
	Analyze	MESRPCAN
	Model	MESRPCM
	Simulate	MESRPCSI
	Full Test	MESRPCFT
	Base Test <i>(one monitor included)</i>	MESRPCBT1
	Ride Comfort	MESRPCRC
	ERPC - Effective Road Profile	MESRPCERP
	Zero Crossing Peak Tool	MESRPCZCPT
	Edit	MESRPCE
	Virtual Test Option	MESRPCVTO
	Synthetic Road Generation	MESRPCSRG
	Point by point time history monitoring option <i>(Base Test is required)</i>	MESPTPMON
	Statistical trend monitoring option <i>(Base Test is required)</i>	MESTRDMON
	Spectral (frequency domain) monitoring option <i>(Base Test is required)</i>	MESSPECTMON
	Fatigue monitoring option <i>(Base Test is required)</i>	MESFTGMON

Note: For a customer with a FlexTest controller to have a supported RPC configuration they must have the latest 793 software. When selling an RPC SSP contract to FlexTest customers, a FlexTest upgrade may be required and a 793 system software SSP contract is required. See the 793 Controller System Software section of this document for pricing and requirements.

793 Controller Software

SOFTWARE	MODULE	SERVICE PRODUCT CODE
793 System Software		
	Supports FlexTest 40 and FlexTest SE Controllers without MPT	MESFTSENMP
	Supports FlexTest 60/100/200 and FlexTest GT Controllers without MPT	MESFTGTNMP
	Supports FlexTest 40 and FlexTest SE Controllers with MPT	MESFTSE
	Supports FlexTest 60/100/200 and FlexTest GT Controllers with MPT	MESFTGT

MTS Echo Intelligent Lab

MTS Echo® Test Monitoring

MTS Echo Test Monitoring

MTS Echo software is a real-time, remote monitoring solution for the entire lab that can be accessed anywhere by any Web-enabled device—smartphone, laptop or tablet. It allows customers to remotely view the status of their test systems, which can be used to increase productivity, reduce downtime and immediately share vital test information.

MTS Echo software is a subscription that is purchased per piece of equipment that is connected. Software subscriptions are billed annually in advance.

Subscription Software

SOFTWARE	SERVICE PRODUCT CODE
MTS Echo	ECHOEM

MTS Echo compatible equipment includes all test systems running on one of the following controllers with MTS Series 793 version 5.3 or later:

- » FlexTest® Models 40/60/100/200 Controller
- » FlexTest SE
- » FlexTest GT
- » FlexTest Ilm
- » MTS Criterion Systems when run with MTS TestSuite 3.6 or later
- » MTS SilentFlo™ G2 Hydraulic Power Units (HPUs) and G1 HPUs with upgraded HMI controls



MTS Echo Intelligent Lab

MTS Echo Health Monitoring - Hydraulic Power Units



All MTS Echo Health Monitoring packages provide:

- » Core sensors
- » Real-time monitoring of signals
- » User-defined alerts
- » Interactive trending charts
- » MTS 24/7 Proprietary Algorithm Protection
- » Health dashboard:
 - Actionable recommendations
 - Predictive trending for recommended maintenance

MTS Echo Test Tracking

MTS Echo Test Tracking offers a modern and easy way to keep customers informed of testing milestones. It reduces the time lab personnel spend on phone calls, email and manual reporting of test status. Echo Test Tracking provides a better experience for test requestors and makes it easier to manage and stay up to date on the progress of all the tests in the lab.

Echo Test Tracking works for the entire lab: MTS equipment, non-MTS equipment and analog equipment. And the test equipment does not need to be connected to the internet. The annual subscription includes test tracking for all test equipment at your site.

SOFTWARE

MTS Echo Test Tracking

SERVICE PRODUCT CODE

ECHOTT

Proactively Protect Your HPU From Common Failures

Health monitoring helps ensure the maximum operating performance and service life for your hydraulic power unit (HPU). MTS Echo Health Monitoring provides the tools to properly identify threats to hydraulic system health and keep your system in optimum condition.

With MTS Echo Health Monitoring, a modular set of options allows you to choose the level of monitoring that best meets your needs. With the Connectivity hardware as a base, you can add packages to monitor temperature, fluid condition, heat exchanger performance and more from any web-enabled device, including smartphones, laptops and tablets.

Each package provides the core sensors necessary for the real-time monitoring of signals. You will be able to set user-defined alerts to notify you of potential concerns. You will also have access to interactive trending data and a health dashboard that offers actionable recommendations for intervention and maintenance.

With 24/7 access to system performance measurements and trends over time, you can maintain your hydraulic systems in the most efficient manner. Condition-based monitoring allows for early detection of potential issues and prevents unnecessary downtime for unneeded maintenance. Convenient, secure connection to hydraulic system condition parameters will help you make the best decisions regarding HPU operation.

BENEFITS

- » Test schedule predictability
- » Test system lifecycle extension
- » Data integrity
- » Budget management

COMPATIBLE HPUS

- » MTS SilentFlo HPUs
- » MTS SilentFlo 505 G1 HPUs*
- » MTS Model 506 HPUs*

* With upgraded controls, HMI information can also be monitored

MTS Echo Intelligent Lab

MTS Echo Health Monitoring - Hydraulic Power Units

MTS Echo Health Monitoring

Standard Recommended Packages

Protect your HPU investment and ongoing lab operations with 24/7 Health Monitoring of your HPU. Start with the recommended packages and add-on for extended protection.

TEMPERATURE PACKAGE

Monitor HPU temperatures for early detection of an over-temperature alarm or fault. These fault notifications can prevent unplanned downtime and minimize the risk of downstream damage. Maintaining a stable temperature is critical to the performance of your test system. Elevated temperature is also indicative of potential system failure.

FLUID CONTAMINATION PACKAGE

Hydraulic power units and servovalves can be damaged by fluid contaminated with particles that are larger than the clearance between lubricated surfaces. These hard particles create more wear contaminant by continually scraping off softer metals, like copper, to further accelerate component failures. This package will help identify conditions that may lead to motor shelling, pressure instability and servovalve instability.

HEAT EXCHANGER FLOW PACKAGE

This package monitors several parameters including changes in the thermal transfer over time and water consumption. Trending the water consumption can identify failed water solenoids or indicate when someone has forgotten to turn on the water supply, giving the operator time to turn on the supply and cool the hydraulic fluid to prevent an over temp shutdown. When combined with the temperature package, Heat Exchanger efficiency also can be tracked.



HPU Health Monitoring Packages Aligned to Common Failure Mode

HEALTH MONITORING PACKAGES	OVER-TEMP FAULTS	HEAT EXCHANGER FAILURES	CONTAMINATION	MOTOR & PUMP FAILURES	CATASTROPHE PREVENTION	SERVICE PRODUCT CODE
Temperature	Yes	Yes, temperature	No	Wear & tear Fouling	Medium	ECHOHMT
Contamination	No	Yes, water saturation	Particulate contamination Water saturation	Wear & tear Fouling	High	ECHOHMF
Heat Flow Exchanger	Yes	Cooling System failures Heat Exchanger failures Cooling water flow	Water saturation	No	High	ECHOHMHXF

MTS Echo Intelligent Lab

MTS Echo Health Monitoring - Hydraulic Power Units

MTS Echo Health Monitoring

Additional Packages

Extend Fault protection with additional stand-alone packages.

POWER MONITORING PACKAGE

Expand your health protection by monitoring your HPU power, including 3- phase voltage, current and energy consumption. This package provides insight into HPU efficiency and enables early warning of pump or motor failures. This package can give you visibility to energy consumption and help improve energy management.




ACCUMULATOR PRE-CHARGE & FILTRATION PACKAGE

Manage hydraulic distribution by monitoring accumulators and filters. Remote monitoring will indicate which accumulators need charging so the hydraulic distribution system can readily manage flow demand. Monitoring filtration can help detect inadequacies that may lead to actuator scoring and potential component failure. The predictive trending and analytics provided in this package will optimize maintenance.



LEAK DETECTION PACKAGE

Protect your lab and the environment by monitoring for early detection of hose rupture, manifold failures and fitting leaks. Improve environmental health and safety while protecting against catastrophic hydraulic fluid loss.

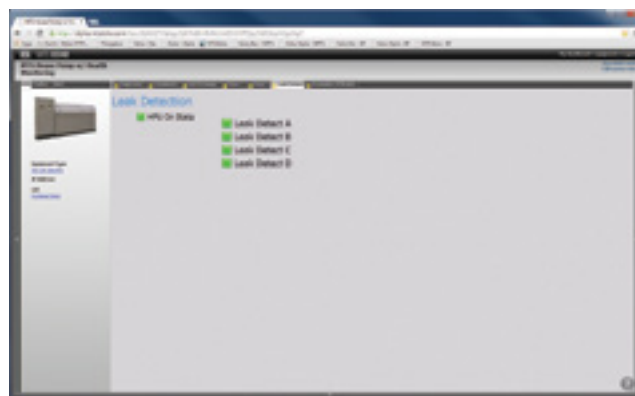

HPU Health Monitoring Packages Aligned to Common Failure Modes

HEALTH MONITORING PACKAGES	OVER-TEMP FAULTS	FAILURES MODES	CONTAMINATION	MOTOR & PUMP FAILURES	CATASTROPHE PREVENTION	SERVICE PRODUCT CODE
Power	No	Energy usage	No	Early fault detection possible	Low	ECHOHMP
Accumulator Pre-charge and Filtration	No	Output pressure accumulators Clogged filters	No	No	Medium	ECHOHMF
Leak Detection	No	Hose rupture Fittings Manifold failures	No	No	Low	ECHOMLD

Warranty Coverage

Overview

MTS builds some of the most rugged and reliable testing equipment on the market. We also stand behind our products and services with extensive warranty options. In addition to providing standard warranty coverage on new equipment and service parts, MTS also offers optional extended warranties for equipment and components to help ensure your equipment is ready for testing when you need it.

Standard Warranty – Equipment

MTS warrants its equipment to be free from defects in material and workmanship under normal use, given proper installation and maintenance, for 12 months from the date of shipment of your product. Defective equipment may be repaired or replaced at MTS' option.

Standard Warranty – Parts & Services

Warranty periods for service products start with the product shipment from MTS. Replacement, exchange, remanufacture, or repair of a system component does not extend the warranty coverage of the entire system. Coverage on system components is as follows:

PURCHASED PARTS

Parts with a model number	12 months
Parts without a model number	90 days

EXCHANGE, REMANUFACTURE OR REPAIR PROGRAMS

Exchange	12 months
Remanufacture	12 months
Repair	90 days

OTHER SERVICES

Factory calibrations	90 days
Field calibrations	90 days
Field service labor	90 days

Extended Warranty

The Extended Warranty Program is available as part of a service contract designed to protect you from incurring unplanned material and labor expenses for repairing MTS equipment. It includes the following:

- » Labor to troubleshoot and repair equipment
- » Parts to repair equipment
- » Exchange, remanufacture or repair as authorized by MTS

Extended Warranty Coverage

The Extended Warranty Program is available for purchase as part of an Assured Maintenance Plan (AMP). Under this plan option, MTS will troubleshoot and repair or replace affected equipment and/or components during a 12-month plan period.

Extended Warranty coverage is grouped into three classes:

- » Controllers (number of control channels)
- » Load frames and actuators
- » Hydraulic power units
- » It is easy to add coverage for your entire laboratory, or for only a portion of your MTS equipment

Extended Warranty Exclusions

- » MTS reserves the right to exclude any product from coverage. Please contact your sales engineer for complete details of product eligibility
- » Upgrades and updates required due to obsolete hardware and software cannot be part of an Extended Warranty Program
- » Hydraulic fluid is not covered under an Extended Warranty Program
- » MTS offers End-of-Life product communications as needed to assist in your equipment planning

MTS reserves the right to change products, services, and prices without notice.

If your system or product is more than 3 years old, a condition assessment may need to be performed by MTS to confirm availability. MTS serial number for each covered item must be recorded within the AMP Extended Warranty contract.

MTS reserves the right to reject claims for warranty where it is determined that failure is caused by Buyer – or third party made – modifications, improper maintenance, misuse, misapplication, improper or incomplete qualification, abuse of the Product, damage caused by connections, interfacing or use in unforeseen or unintended environment. These conditions will render warranties null and void.



LOAD CELLS

Highly accurate load cells deliver high stiffness, stability and linearity for a full range of both static and dynamic testing applications. They are constructed from aircraft quality military specification materials, which are heat treated to ensure uniform hardness and minimize distortion.



Introducing the *MTS Platinum Promise* for Load Cells

Our most popular load cells are in stock and ready to ship.

Look for this symbol to identify *Platinum Promise* load cells.



LOAD CELLS

Introduction	34
Load Cells for Servohydraulic Systems	35-36
Load Cells for Insight Systems	37
Load Cells for Electrodynamic Systems	38
Load Cell Adapters	39-41

Load Cells for Servohydraulic Systems

- » MTS axial, dynamic fatigue-rated load cells reduce errors caused by extraneous side loads or loading changes caused by geometry shifts in the specimen
- » Monolithic design incorporates high axial stiffness
- » Calibration, cables and attachment kits are sold separately
- » TEDs (Transducer Electronic Data Sheets) available – many models are available with TEDS, the self-identification capabilities that follow the IEEE 1451.4 standard
- » Dynamic performance – low deflection and high stiffness give you better dynamic performance
- » High output – provides you with excellent resolution and reading accuracy
- » Radially oriented beams – enable the unit to resist off-axis loads and moments for greater accuracy
- » High degree of component concentricity and parallelism – provide you with greater accuracy during your test setup
- » Unique proprietary wiring technique used on the bridge allows for minimal susceptibility to stray magnetic fields
- » Interchangeable – simple mounting makes installation and changing easy
- » Ability to stack load cells to test smaller/larger specimens in same load frame
- » Large choice of sizes – increases versatility
- » Dual bridge options available – contact MTS



These load cells are fatigue rated.



Axial, Dynamic Fatigue-Rated Load Cells Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.11B-01	250 N	M6 x 1	050-781-905
661.11B-02	500 N	M6 x 1	050-781-906
661.18F-01	1 kN	M12 x 1.25	050-459-902
661.18F-02	2.5 kN	M12 x 1.25	050-459-904
661.19F-01	5 kN	M12 x 1.25	045-438-102
661.19F-02	10 kN	M12 x 1.25	045-438-104
661.19F-03	15 kN	M12 x 1.25	045-438-106
661.19F-04	25 kN	M12 x 1.25	045-438-108
661.20F-01	25 kN	M27 x 2	045-430-004
661.20F-02	50 kN	M27 x 2	045-430-005
661.20F-03	100 kN	M27 x 2	045-430-006

Axial, Low Capacity Load Cells

Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.09B-20	10 N	M6 x 1	051-340-602
661.09B-20	10 N	M6 x 1	051-340-602

US Customary Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.09A-20	2 lbf	0.25 in - 28	051-340-601
661.09A-21	22 lbf	0.25 in - 28	051-340-603

US Customary Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.11A-01	50 lbf	0.25 in - 28	050-781-901
661.11A-02	100 lbf	0.25 in - 28	050-781-902
661.18E-01	220 lbf	0.50 in - 20	050-459-901
661.18E-02	550 lbf	0.50 in - 20	050-459-903
661.19E-01	1.1 kip	0.50 in - 20	045-438-101
661.19E-02	2 kip	0.50 in - 20	045-438-103
661.19E-03	3 kip	0.50 in - 20	045-438-105
661.19E-04	5 kip	0.50 in - 20	045-438-107
661.20E-01	5 kip	1.00 in - 14	045-430-001
661.20E-02	11 kip	1.00 in - 14	045-430-002
661.20E-03	22 kip	1.00 in - 14	045-430-003

Load Cells For Servohydraulic Systems



Axial, Dynamic Fatigue-Rated Load Cells

Single Bridge

Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.22D-01	250 kN	M36 x 2	042-361-802
661.23F-01	500 kN	M52 x 2	044-445-002
661.31F-01	1000 kN	M76 x 2	045-190-502
661.34F-01	1500 kN	M90 x 2	045-225-902
661.36D-03	2500 kN	M2 x 12	046-311-702

US Customary Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.22C-01	55 kip	1.50 in - 12	042-361-801
661.23E-01	110 kip	2.00 in - 12	044-445-001
661.31E-01	220 kip	3.00 in - 12	045-190-501
661.34E-01	330 kip	1.25 in - 12	045-225-901
661.36C-03	550 kip	1.50 in - 12	046-311-701

Axial, Dynamic Fatigue-Rated Load Cells

Single Bridge

Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
661.11H-01	250 N	M6 x 1	057-117-901
661.11H-02	500 N	M6 x 1	057-117-902
661.11H-03	1000 N	M6 x 1	057-117-903
661.18H-01	1000 N	M12 x 1.25	057-118-001
661.18H-02	2.5 kN	M12 x 1.25	057-118-002
661.19H-01	5 kN	M12 x 1.25	057-118-101
661.19H-02	10 kN	M12 x 1.25	057-118-102
661.19H-03	15 kN	M12 x 1.25	057-118-103
661.19H-04	25 kN	M12 x 1.25	057-118-104
661.20H-02	50 kN	M27 x 2	057-118-201
661.20H-03	100 kN	M27 x 2	057-118-202
661.22H-01	250 kN	M36 x 2	057-263-701
661.23H-01	500 kN	M52 x 2	057-263-801

Axial-Torsional Load Cells

Metric Units

MODEL	AXIAL FORCE	TORSIONAL FORCE	THREAD SIZE	PART NUMBER
662.20D-01	2.5 kN	25 N.m	M6 x 1	049-611-002
662.20D-03	10 kN	100 N.m	M8 x 1.25	049-381-502
662.20D-04	15 kN	150 N.m	M8 x 1.25	050-457-702
662.20D-05	25 kN	250 N.m	M10 x 1.5	050-027-102

US Customary Units

MODEL	AXIAL FORCE	TORSIONAL FORCE	THREAD SIZE	PART NUMBER
662.20C-01	550 lb	250 in.lb	1/4 in - 20	049-611-001
662.20C-03	2200 lb	1000 in.lb	5/16 in - 18	049-381-501
662.20C-04	3300 lb	1500 in.lb	5/16 in - 18	050-457-701
662.20C-05	5500 lb	2500 in.lb	3/8 in - 16	050-027-101

Axial-Torsional Load Cells with Integrated TEDS

Metric Units

MODEL	AXIAL FORCE	TORSIONAL FORCE	THREAD SIZE	PART NUMBER
662.20H-01	2.5 kN	25 N.m	M6 x 1	057-202-601
662.20H-03	10 kN	100 N.m	M8 x 1.25	057-202-701
662.20H-04	15 kN	150 N.m	M8 x 1.25	057-202-702
662.20H-05	25 kN	250 N.m	M10 x 1.5	057-202-801

Load Cells for Insight & Insight Renew Universal Test Systems

Refer www.mtsmonotonic.com for a complete list of all monotonic load cells for MTS Exceed and MTS Criterion Systems

- » Highly accurate MTS load cells for MTS Insight and Insight Renew are designed to offer world-class stiffness, stability, and linearity.
- » Provide overload and side load protection
- » Designed with built-in shunt resistors to facilitate regular verification of accuracy using calibration routines featured in MTS software.
- » Load cells feature TEDS (Transducer Electronic Data Sheets) self-identification capabilities that follow the IEEE 1451.4 standard



S-Beam Load Cells with TEDS for Insight Universal Test Systems

Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	PART NUMBER
S-Beam	5 N (1 lbf)	M6x1	056-932-601
S-Beam	10 N (2 lbf)	M6x1	056-932-602
S-Beam	25 N (5 lbf)	M6x1	056-932-603
S-Beam	50 N (10 lbf)	M6x1	056-932-604
S-Beam	100 N (20 lbf)	M6x1	056-932-605
S-Beam	250 N (50 lbf)	M6x1	056-932-606
S-Beam	500 N (110 lbf)	M6x1	056-932-701
S-Beam	1 kN (220 lbf)	M6x1	056-932-702
S-Beam	2 kN (450 lbf)	M6x1	056-932-703
Low Profile	125 N (25 lbf)	M6x1	056-932-801
Low Profile	250 N (50 lbf)	M6x1	056-932-802
Low Profile	500 N (110 lbf)	M6x1	056-932-803
Low Profile	1 kN (225 lbf)	M12x1.25	056-932-901
Low Profile	2.5 kN (550 lbf)	M12x1.25	056-932-902
Low Profile	5 kN (1,100 lbf)	M12x1.25	056-932-903
Low Profile	10 kN (2,200 lbf)	M12x1.25	056-932-904
Low Profile	30 kN (6,600 lbf)	M12x1.25	056-932-905
Low Profile	50 kN (11,000 lbf)	M16x1.5	056-933-201
Low Profile	100 kN (22,500 lbf)	M27x2	056-933-001
Low Profile	150 kN (33,750 lbf)	M27x2	056-933-002
Low Profile	200 kN (45,000 lbf)	M36x2	056-933-301
Low Profile	300 kN (67,500 lbf)	M36x2	056-933-101

Load Cells for MTS Acumen Electrodynamic Test Systems

Replacement MTS Acumen System Load Cell Kits with Built-in Accelerometer

- » Designed for the MTS Acumen Electrodynamic Test Systems, these axial, dynamic fatigue rated load cells can be mounted either to the actuator or to the T-slot table
- » Dynamic performance – low deflection and high stiffness give you better dynamic performance
- » High output – provides you with excellent resolution and reading accuracy
- » Accelerometer – Integrated acceleration sensor enables load cell to be mounted on moving actuator
- » Piloted load cell – simplifies load cell installation and minimizes the need for a separate alignment fixture
- » Radially oriented beams – enable the unit to resist off-axis loads and moments for greater accuracy
- » High degree of component concentricity and parallelism – provide you with greater accuracy during your test setup
- » Unique proprietary wiring technique used on the bridge allows for minimal susceptibility to stray magnetic fields
- » Interchangeable – simple mounting makes installation and changing easy
- » Ability to stack lower force transducers to test smaller specimens in same load frame
- » Integrated TEDS



Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	KIT PART NUMBER*	FACTORY CALIBRATED KIT PART NUMBER*
661.18SE-01	1.5 kN	M6 x 1	057-560-101	057-560-106
661.18SE-02	3 kN	M6 x 1	057-560-102	057-560-107

*Kits include load cell, mounting hardware, CD and protective case

Tandem Load Cell Kits with TEDS for MTS Acumen Systems

- » Compact axial, dynamic fatigue rated load cells mount easily to the MTS Acumen system load cells
- » Smaller force capacity load cell increases the resolution test data when testing at the lower end of the system load cell's capacity.
- » Dynamic performance – low deflection and high stiffness give you better dynamic performance
- » High output – provides you with excellent resolution and reading accuracy
- » Radially oriented beams – enable the unit to resist off-axis loads and moments for greater accuracy
- » High degree of component concentricity and parallelism – provide you with greater accuracy during your test setup
- » Unique proprietary wiring technique used on the bridge allows for minimal susceptibility to stray magnetic fields.
- » Interchangeable – simple mounting makes installation and changing easy
- » Easily mounted to electrodynamic system load cells.
- » Integrated TEDS



Metric Units

MODEL	FORCE CAPACITY	THREAD SIZE	KIT PART NUMBER*	FACTORY CALIBRATED KIT PART NUMBER*
661.11H	10 N	M6 x 1	057-560-111	057-560-112
661.11H	25 N	M6 x 1	057-560-113	057-560-114
661.11H	50 N	M6 x 1	057-560-115	057-560-116
661.11H-04	125 N	M6 x 1	057-560-105	057-560-110
661.11H-01	250 N	M6 x 1	057-560-104	057-560-109
661.11H-02	500 N	M6 x 1	057-560-103	057-560-108

*Kits include load cell, mounting hardware and spanner wrench

Load Cell Adapters

Adapter Introduction

MTS is pleased to offer a wide variety of load cell adapters to meet all of your special needs. Our adapters are made to the best possible specifications and adhere to the industry standards. In many cases, we offer more than one type of material so you may choose the lightweight aluminum or the heavy duty steel adapters. If there is an adapter you would like that you do not see listed, contact your sales engineer.

USE PIGGYBACK ADAPTERS FOR MOUNTING LOW CAPACITY LOAD CELLS

- » Adapters provide a convenient way to mount a lower capacity load cell to a large capacity load cell
- » Piggyback adapters eliminate the need to remove the larger load cell when a lower capacity test is needed for testing
- » Piggyback adapters are available for all load cell combinations



Piggyback

FROM THREAD SIZE	TO THREAD SIZE	PART NUMBER
0.5 in - 20	1 in - 14	100-082-366
0.25 in - 28	1 in - 14	100-082-361
0.25 in - 28	1.5 in - 12	100-082-608
0.25 in - 28	0.5 in - 20	100-082-319
M16 x 1.5	M24 x 1.5	100-105-984
1 in - 14	1.5 in - 12	100-082-598
M12 x 1.25	M6 x 1	100-093-064
1 in - 14	1 in - 14	100-179-743

AXIAL MONOTONIC AND DYNAMIC FATIGUE-RATED LOWER FORCE LOAD CELL KITS

- » Tandem Piggyback adapters have much of the same functionality as the standard piggyback adapters; however, the Tandem Piggyback adapters have locking rings
- » Locking rings allow you to preload the threads to eliminate the slack in the load train



Tandem Piggyback

FROM THREAD SIZE	TO THREAD SIZE	PART NUMBER
M6 x 1	M16 x 1.5	056-871-102
M12 x 1.25	M27 x 2	056-871-106
M12 x 1.25	M36 x 2	056-871-107
M12 x 1.25	M16 x 1.5	056-871-105
M16 x 1.5	M27 x 2	056-871-108
M16 x 1.5	M36 x 2	056-871-109
M27 x 2	M36 x 2	056-871-110
M6 x 1	M12 x 1.25	056-871-101
M6 x 1	M36 x 2	056-871-104
M6 x 1	M27 x 2	056-871-103

Load Cell Adapters

CLEVIS PIN

- » With a clevis pin, you can easily and quickly change out your grips and fixtures without decreasing the force of your test
- » MTS offers a wide variety of clevis pin devices that meet industry standards



Clevis Pin

THREAD SIZE	CONNECTION TYPE	PIN DIAMETER	CLEVIS DIAMETER	PART NUMBER
M27 x 2	Df	0.5 in	1.25 in	100-140-862
M12 x 1.25	Cf	0.25 in	0.625 in	100-140-788
M12 x 1.25	Df	0.5 in	1.25 in	100-140-546
1 in - 14	Df	0.5 in	1.25 in	100-088-645
M16 x 1.5	Df	0.5 in	1.25 in	100-092-928
0.25 in - 28	Df	0.5 in	1.25 in	100-092-552
M24 x 1.5	Df	0.5 in	1.25 in	100-092-905
M36 x 2	Df	0.5 in	1.25 in	100-092-937
0.5 in - 20	Df	0.5 in	1.25 in	100-088-650
M6 x 1	Df	0.5 in	1.25 in	100-092-919
M12 x 1.25	Bf	0.187 in	0.5 in	100-093-049
M6 x 1	Bf	0.187 in	0.5 in	100-092-923
0.25 in - 28	Bf	0.187 in	0.5 in	100-092-914
M6 x 1	Cf	0.25 in	0.625 in	100-092-910
0.25 in - 28	Cf	0.25 in	0.625 in	100-092-547

Load Cell Adapters

CONVERSION ADAPTER

- » MTS offers a wide variety of conversion adapters to help you quickly change from one type or style of clevis adapter to another without the costs associated with a full adapter



Conversion Adapter


THREAD SIZE	CLEVIS PIN CONNECTION	TYPE	PART NUMBER
30 mm	15f	Male to Female	100-093-030
0.5 in - 20	15f	Male to Female	100-093-001
1.25 in - 12	30f	Male to Female	100-092-991
0.625 in - 18	15f	Male to Female	100-092-996
1.25 in	Cf	Male to Female	100-092-880
0.25 in - 28	Df	Male to Female	100-081-395
1.50 in - 12	Df	Male to Female	100-082-558
0.5 in - 20	Cf	Male to Female	100-081-868
0.5 in - 20	Bf	Male to Female	100-081-862
1 in - 14	Cf	Male to Female	100-081-873

Insight Load Cell Mounting Kits

THREAD SIZE	WORKSTATIONS	PART NUMBER
M36 x 2	200 & 300	056-904-512
M27 x 2	200 & 300	056-904-511
M16 x 1.5	200 & 300	056-904-510
M27 x 2	100 & 150	056-904-509
M16 x 1.5	100 & 150	056-904-508
M12 x 1.25	100 & 150	056-904-507
M16 x 1.5	30 & 50	056-904-506
M12 x 1.25	30 & 50	056-904-505
M6 x 1	30 & 50	056-904-504
M12 x 1.25	5 & 10	056-904-503
M6 x 1	5 & 10	056-904-502
M6 x 1	1 & 2	056-904-501

EXTENSOMETERS

Extensometers are used for measuring and/or controlling strain in tension, compression and fatigue testing applications. Measuring on the specimen significantly increases the reliability of your data. Relying on the crosshead or actuator movement of your materials test system will measure too much – machine deflection, grip deflection, and possible slippage are all captured in the load train displacement measurement.

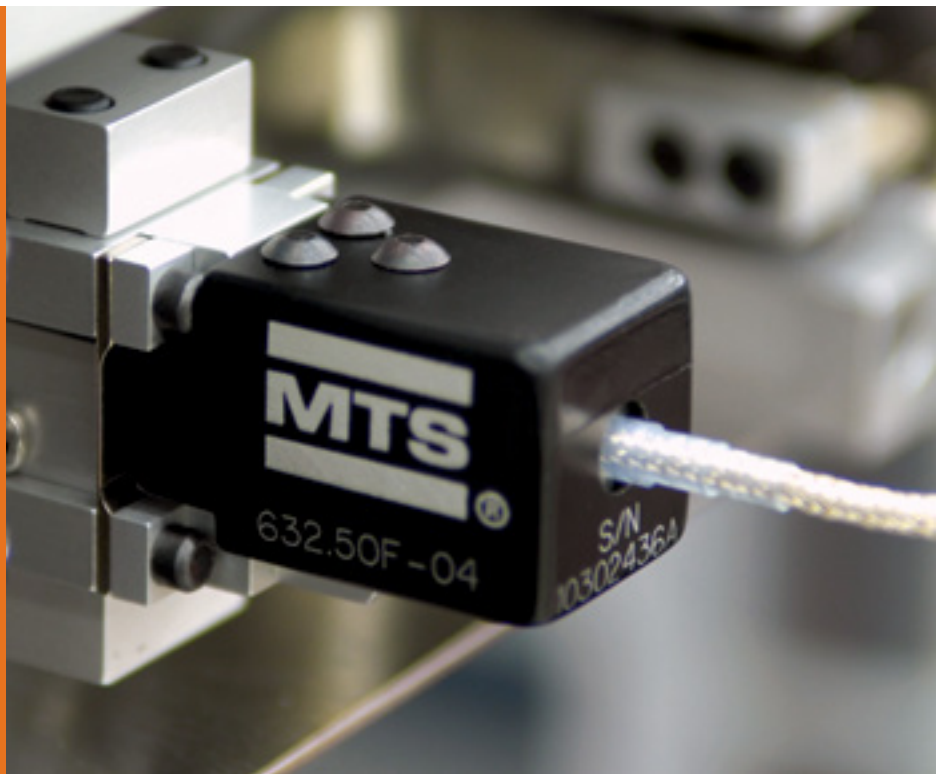

Extensometer Table of Contents

Introduction	42-43
How to Select an Extensometer	44
Extensometer Application Index	45
Extensometer Functional Index	46
Axial Extensometers	47-59
High Temperature Extensometers	60-62
Specialty Extensometers	63-67
MTS Advantage Video Extensometer	68-69
Laser Extensometer	70
Other Applications	71
Clip-On Displacement Gages	72
Displacement Gages	73
Extensometer Calibrator	74
MTS Direct Current Potential Drop (DCPD) Solution	75

Extensometers

Accuracy is Key

Pioneered to Perfection

MTS pioneered and perfected the cross-flexure design that ensures true center point bending of its extensometers. Precision, resistive-type foil strain gages are bonded to the MTS cross-flexure to form a 4-arm fully active Wheatstone bridge.

- » Bridges of 350 ohms and 1000 ohms are used
- » Output, depending on the unit is 2 mV/V to 6 mV/V
- » Nominal excitation range from 5 to 12 V DC for room temperature testing

HIGH PERFORMANCE, HIGH QUALITY

- » Proven stability over 30+ years of usage
- » Built-in over-travel protection provides repeatable accurate data from test initiation through specimen failure
- » Quickly and accurately set-up your extensometer using the zero set pin or zero stop
- » Cross-flexure design ensures the strain gages are subjected to true center point bending and the same force, time after time
- » Small, lightweight design – most are under 35 g mass


- » MTS extensometer models are available for both static and dynamic testing
- » MTS extensometers meet or exceed requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

CAN'T FIND WHAT YOU NEED?

- » Our extensometer portfolio contains hundreds and hundreds of additional dynamic rated extensometers. Contact your local sales representative or applications engineer to find the model that meets your exact needs
- » Custom models are also available upon request

MTS Extensometer Repair & Rework

Why go through the process of certifying a new extensometer when the one you have still does the job? Many of our extensometers survive tens of thousands of tests without a problem – so if the unit only needs a new zero set pin or cables let us repair it for you. Contact your local sales representative for details.

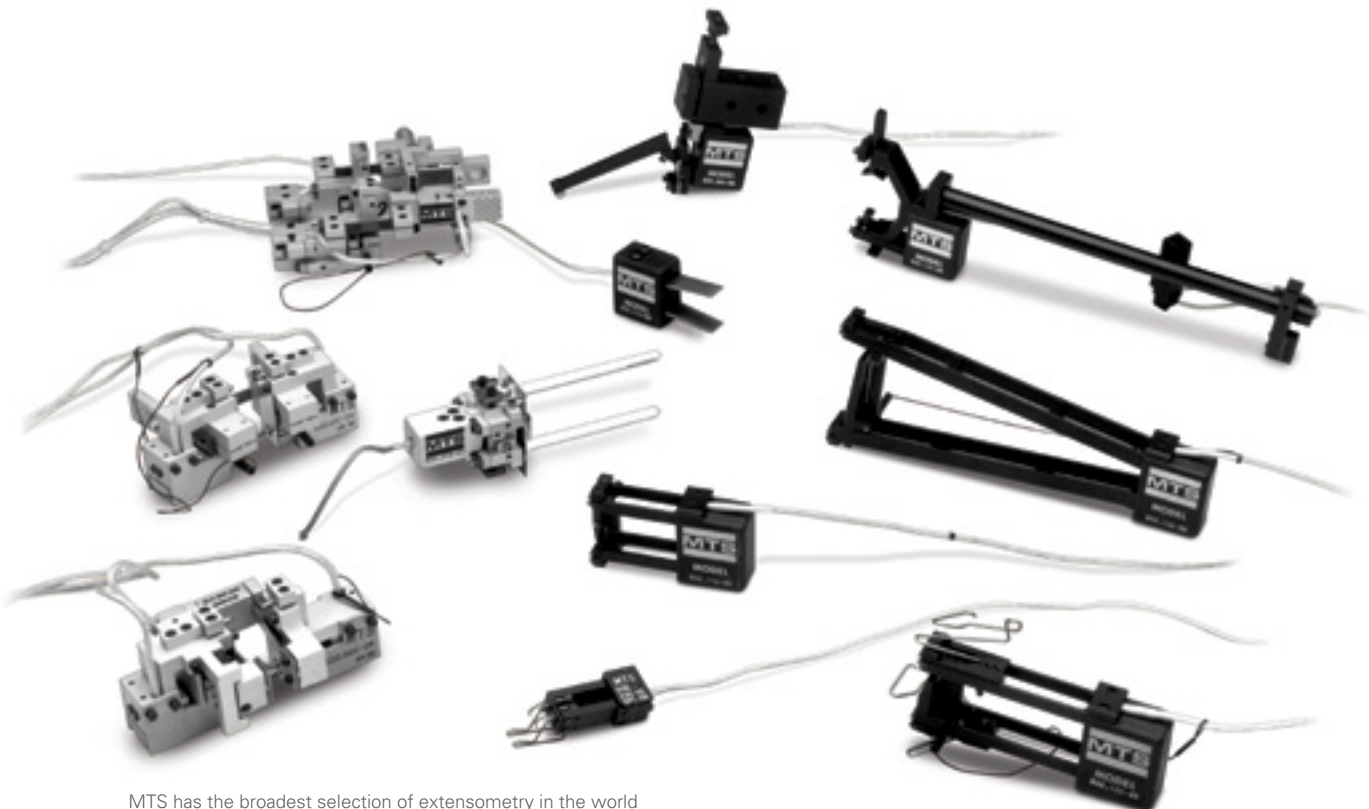

MTS has the broadest selection of extensometry in the world

How to Select an Axial Extensometer

Start With the Specimen

What shape is your specimen?

MTS offers contact extensometers with gage lengths between 3 mm (0.12 in) and 100 mm (3.93 in). Select an extensometer with a gage length slightly shorter than the gage section of your specimen.

Whether the specimen cross-section is round or flat determines which type of knife edges to use with the extensometer. Straight knife edges are standard in our extensometers and are recommended for testing round specimens. Knife edges with 3-point specimen contact are typically used when testing flat specimens and they include one single-point and one double-point contact knife edge.

What is your specimen made of?

Both the physical response of the test material and the test method itself, influence how far the specimen will travel. Ductile metals may see more than a 15% strain during a tensile test to failure but a fatigue test on the same metal may be limited to 1%. Plastic and rubber materials have higher elongation than metals and often require extensometers with a greater measurement range. Stiff materials, such as ceramics, should be tested with extensometers with smaller measurement ranges. Higher resolution measurements typically require extensometers with lower full scale strain ranges.

Next, Think About the Test Method

Are you performing a monotonic or dynamic test?

Monotonic Testing: MTS Fundamental Series 635 Extensometers are ideal for measuring strain in tension testing applications. Their design is optimized for monotonic testing. Series 635 Extensometers should not be used for testing in compression.



Dynamic Testing: Extensometers supported by and clipped directly onto a specimen exert a small bending stress on the specimen. Extensometers that are lighter weight and closer to the specimen create smaller bending moments. MTS Model 632.29 extensometers are ideal for measuring strain in wires and small organic specimens, because they have small gage lengths, low mass and mount very close to the specimen.

In dynamic testing, if the frequency of the test cycle matches the natural frequency of the extensometer, the extensometer can begin to vibrate increasing the noise in the strain measurement signal. Our dynamic extensometers are designed with high natural frequencies to ensure that this never happens. Dynamic extensometer MTS Model 632.27 has an optimized design for a higher natural frequency above other 634 and 632 models. The MTS Model 632.27 is designed for fatigue applications up to 150 Hz.



Will testing take place in a furnace, environmental chamber or liquid bath?

Check out MTS model 632.5x extensometers for high-temperature applications up to 1200°C (2200°F). These models are mounted outside a furnace or environmental chamber and are available with options for air or water cooling. Water cooling is the most effective for keeping the extensometer body at a constant temperature and maintaining low noise in the data measured.

Model 633.11 is designed for elevated temperature testing without external cooling. This model can be used up to 540°C (1000°F).

Most of our 632 and 634 models can be used in cryogenic applications as cold as -100°C (-150°F).

Are you following a standard test procedure?

Refer to the application index on the next page to find suggested model numbers for specific applications and ASTM standards.

When to Select a Non-Contact Extensometer

When your specimen response measurement requirements are more demanding than conventional contact measurement techniques, MTS non-contact extensometry is the preferred choice.

» **Perfect for Delicate Specimens:**

When testing a specimen that is fragile, brittle or irregularly shaped, MTS non-contact extensometry prevents damage caused by the contact forces of a clip on extensometer.

» **Take Multiple Measurements:**

The MTS Advantage Video Extensometer allows you to take up to 200 measurements on one specimen. One video extensometer can replace dozens of contact extensometers.

» **Protect Your Investment:**

Measurements through failure can be taken without any threat of damage to the video extensometer so tests can run without interruption.



Extensometer Application Index

Use this helpful index to determine which extensometers are best suited for your application. If the extensometer

listed is in blue, it is not shown in this catalog. Please contact us for more information.

Application	Extensometer Options	See Pages
METALS		
Tensile (ASTM E8)	632.13/.24 & 634.11/.12/.25/.28/.31; 635; AVX; FAX	48, 49, 51, 52, 53, 54, 57, 68
Compression (ASTM E9)	632.13/.26/.27 & 634.11/.12/.31	48, 49, 50, 53
KIC - Fracture Toughness (ASTM)	632.02 ASTM Standard; 632.03 MTS Enhanced; 632.65 (High Temperature)	72
Pre-Cracking (ASTM E399)	632.05 (High Frequency)	please call
JIC (ASTM E813)	632.03 with Long Travel Option	72
Crack Growth - Bend (ASTM E399)	632.03	72
Fatigue Crack Growth (ASTM E647)	632.03 with Long Travel Option; 632.05 (High Frequency)	72
High Cycle Fatigue	632.27 (Low Strain)	50
Low Cycle Fatigue (ASTM E606)	632.13/.18/.19/.20 & 634.11/.12/.31 & 632.50/.51/.53/.54	48, 49, 53, 60, 61, 66
Low Cycle Fatigue (High Temp)	632.50/.51/. 52	60, 61, (please call for .56)
Low Cycle Fatigue (High Temp)	632.53/.54	60
ADVANCED MATERIALS		
Bend (MIL SPEC 1942, JIS 1601R)	632.06; AVX	68, 73
Tensile/Compression	632.26; AVX; FAX	48, 57, 68
High Temperature (CMC's MMC's)	632.5X	60, 61
Ultra High Temperature Bend	632.70	please call
COMPOSITES		
Tensile (ASTM D3039)	632.17/.85 & 634.11/.31; AVX	49, 53, 63, 64, 68
Compression (ASTM D695)	632.17 & 634.11/.31; AVX	49, 53, 63, 68
±45° Shear Test (ASTM D3518)	632.17 & 634.11/.31	49, 53, 63
Flexure (ASTM D790)	632.06 & 634.11/.31	49, 53, 73
Tension/Tension Fatigue (ASTM D3479)	632.17/.85 & 634.11/.31	49, 53, 63, 64
PLASTICS		
Tensile (ASTM D638)	632.17/.24/.85 & 634.12/.25/.31, 635, LX 500/1500; AVX; FAX; AHX 850	49, 51, 52, 53, 54, 57, 58, 63, 64, 68, 69, 70
ROCK MECHANICS		
Triaxial Comp. Strength (ASTM D2664)	632.90/.92	please call
Triaxial Creep (ASTM D4406)	632.90/.92	please call
BIO-MATERIALS	632.29 & 634.11/.12/.25/.31, LX 500; AVX; FAX	47, 49, 51, 53, 57, 68, 70
AXIAL/TORSIONAL	632.68/.80	please call
HIGH TEMPERATURE		
Air Cooled	632.53/.54	60
Water Cooled	632.51/. 68	61, (please call for .73)
Uncooled	632.53/.54/. 70	60
Hostile Environments	632.10 & LX 500/1500	70
FRACTURE MECHANICS		
Fracture Toughness	632.02/.03; 632.65	72
Pre-Cracking	632.05 (High Frequency)	please call
BEND		
Displacement	632.06	73
Crack Growth	632.03	72
High Temperature	632.70	please call
FRACTURE		
Crack Growth	632.03 with Long Travel Option; 632.05 & 632.65	72
High Cycle	632.27 (Low Strain)	49
Low Cycle Mechanics	632.13/.18/.19/.20 & 634.11/.12/.31	48, 49, 53, 66
Low Cycle (High Temp)	632.50/.51/.53/.54	60, 61
TENSILE/COMPRESSION TESTING		
Tensile	632.13/.24 & 634.11/.12/.25/.28/.31; 635; AVX; FAX; AHX 850	48, 49, 51, 52, 53, 54, 57, 58, 68
Compression	632.13/.26/.27 & 634.12/.25/.31; AVX; FAX	48, 49, 50, 53, 57, 68
Poisson's Ratio	632.18/.19/.20/.23/.85 & 634.12/.25/.31; AVX; FAX	49, 51, 53, 57, 64, 65, 66, 68
Elastic Limit	632.26/.27/.29	47, 48, 50

Extensometer Functional Index

Axial Contact Extensometers

DYNAMIC MODELS (TENSION/COMPRESSION)

Metric Units

MODEL NUMBER	GAGE LENGTH	MEASURING RANGE: STRAIN (TRAVEL)
632.29F-20	3 mm	±8% (±0.24 mm)
632.29F-30	5 mm	-10%/+30% (-0.5 mm/+1.5 mm)
632.29F-30	6 mm	±4% (±0.24 mm)
632.26F-3X	8 mm	±6% (±0.48 mm)
632.26F-2X	8 mm	±15% (±1.2 mm)
632.13F-2X	10 mm	±15% (±1.5 mm)
632.26F-4X	12 mm	±9% (±1.1 mm)
632.27F-3X	25 mm	±2% (±0.5 mm)
632.27F-2X	25 mm	±4% (±1.0 mm)
634.11F-2X	25 mm	+20%/-10% (+5 mm/-2.5 mm)
634.12F-2X	25 mm	+50%/-10% (+12.5 mm/-2.5 mm)
634.25F-2X	50 mm	+50%/-10% (+25 mm/-5 mm)
634.31F-2X	10-50 mm	variable (+4 mm/-2 mm)

US Customary

MODEL NUMBER	GAGE LENGTH	MEASURING RANGE: STRAIN (TRAVEL)	SEE PAGE
632.29E-20	0.12 in	±8% (±0.0096 in)	47
632.29E-30	0.2 in	-10%/+30% (-0.02 in/+0.06 in)	47
632.29E-3X	0.24 in	±4% (±0.0096 in)	47
632.26E-3X	0.3 in	±6% (±0.018 in)	48
632.26E-2X	0.3 in	±15% (±0.045 in)	48
632.13E-2X	0.5 in	±15% (±0.075 in)	48
632.26E-4X	0.5 in	±9% (±0.045 in)	48
632.27E-3X	1.0 in	±2% (±0.02 in)	50
632.27E-2X	1.0 in	±4% (±0.04 in)	50
634.11E-2X	1.0 in	+20%/-10% (+0.2 in/-0.1 in)	49
634.12E-2X	1.0 in	+50%/-10% (+0.5 in/-0.1 in)	49
634.25E-2X	2.0 in	+50%/-10% (+1.0 in/-0.2 in)	51
634.31E-2X	0.5-2.0 in	variable (+0.2 in/-0.1 in)	53

MONOTONIC MODELS (TENSION ONLY)

Metric Units

MODEL NUMBER	GAGE LENGTH	MEASURING RANGE: STRAIN (TRAVEL)
635.25F-05	25 mm	+20% (+5 mm)
634.11F-5x	25 mm	+20% (+5 mm)
634.12F-5x	25 mm	+50% (+12.5 mm)
632.24F-50	25 mm	+100% (+25 mm)
635.50F-05	50 mm	+10% (+5 mm)
635.50F-10	50 mm	+20% (+10 mm)
635.50F-25	50 mm	+50% (+25 mm)
634.25F-5x	50 mm	+50% (+25 mm)
634.28F-24	50 mm	+100% (+50 mm)
635.100F-10	100 mm	+10% (+10 mm)


US Customary

MODEL NUMBER	GAGE LENGTH	MEASURING RANGE: STRAIN (TRAVEL)	SEE PAGE
			54
634.11E-5X	1.0 in	+20% (+0.2 in)	49
634.12E-5X	1.0 in	+50% (+0.5 in)	49
632.24E-50	1.0 in	+100% (+1.0 in)	52
634.25E-5X	2.0 in	+50% (+1.0 in)	51, 54
634.28E-24	2.0 in	+100% (+2.0 in)	52, 54
			54
			51
			52
			54

Axial Extensometers

Miniature Extensometers: <6 mm (<0.23 in) Gage Length

- » MTS Model 632.29 Extensometers are really small, because sometimes there just is not enough room. Most models are only 15.3 mm (0.60 in) tall and weigh roughly 6 g
- » Typical applications include testing of small wires, brittle or fragile materials, small organic specimens, or expensive materials
- » Includes a gage fixture for proper gage length setting and a 1000 Ohm fully-active Wheatstone bridge for high sensitivity
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

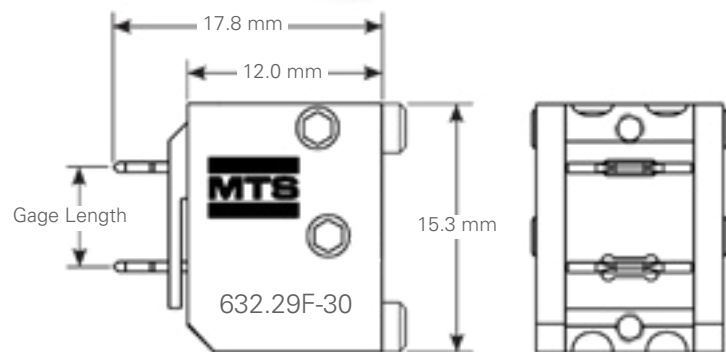

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
632.29F-20	3 mm	±8%	15.3 mm	17.8 mm	-100°C / 150°C	30 g		047-402-102
632.29F-30	5 mm	-10% / 30%	15.3 mm	36.8 mm	-100°C / 150°C	30 g	x	047-402-108
632.29F-30	6 mm	±4%	15.3 mm	17.8 mm	-100°C / 150°C	30 g		047-402-104

US Customary Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
632.29E-20	0.12 in	±8%	0.60 in	0.70 in	-150°F / 300°F	30 g		047-402-101

*These units may be used at 25°C (50°F) higher than the listed temperatures, but only for durations less than 24 hours.

Quick Attachment Kits


Improve your productivity by selecting the part number that includes a quick attachment kit. Includes adapters for flat and round specimens.



Axial Extensometers

Axial Extensometers with <25 mm (<1 in) Gage Length

- » Models 632.13 and 632.26 extensometers are ideal for testing small specimens
- » Measuring range is equally divided into tensile and compressive range
- » Includes knife edges and manual attachment fixtures for round specimens
- » Select the -21 or -31 model for cryogenic testing down to -269°C (-452°F)
- » Models include quick attachment kits for round and flat specimens where indicated. Additional quick attachment kits are available upon request
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

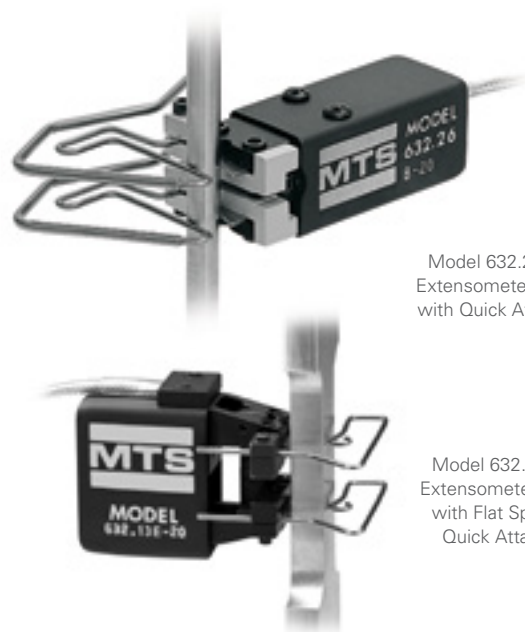

Model 632.26 Axial Extensometer (shown with Quick Attach Kit)

Model 632.13 Axial Extensometer (shown with Flat Specimen Quick Attach Kit)

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT	PART NUMBER
632.26F-3x	8 mm	±6%	17.5mm	44.5 mm				
632.26F-30					-100°C / 150°C	40 g		047-401-810
632.26F-30					-100°C / 150°C	40 g	x	047-401-821
632.26F-31					-269°C / 65°C	40 g		047-401-811
632.26F-2x	8 mm	±15%	17.5 mm	44.5 mm				
632.26F-20					-100°C / 150°C	40 g		047-401-804
632.26F-20					-100°C / 150°C	40 g	x	047-401-849
632.26F-23					-100°C / 175°C	60 g		047-401-806
632.13F-2x	10 mm	±15%	32.8 mm	41.1 mm				
632.13F-20					-100°C / 150°C	45 g		047-401-004
632.13F-20					-100°C / 150°C	45 g	x	047-401-010
632.13F-21					-269°C / 65°C	45 g		047-401-005
632.13F-23					-100°C / 175°C	65 g		047-401-006
632.26F-4x	12 mm	±9%	18.8 mm	44.5 mm				
632.26F-40					-100°C / 150°C	40 g		047-401-816
632.26F-40					-100°C / 150°C	40 g	x	047-401-822

US Customary Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT	PART NUMBER
632.26E-3x	0.300 in	±6%	0.69 in	0.175 in				
632.26E-30					-150°F / 300°F	40 g		047-401-807
632.26E-2x	0.300 in	±15%	0.69 in	0.175 in				
632.26E-20					-150°F / 300°F	40 g		047-401-801
632.13E-2x	0.500 in	±15%	1.29 in	1.62 in				
632.13E-20					-150°F / 300°F	45 g		047-401-001
632.13E-20					-150°F / 300°F	45 g	x	047-401-007
632.13E-21					-452°F / 150°F	45 g		047-401-002
632.26E-4x	0.500 in	±9%	0.74 in	1.75 in				
632.26E-40					-150°F / 300°F	40 g		047-401-813
632.26E-41					-452°F / 150°F	40 g		047-401-814
632.26E-43					-150°F / 350°F	60 g		047-401-815

*These units may be used at 25°C (50°F) higher than the listed temperatures, but only for durations less than 24 hours.

Axial Extensometers

Axial Extensometers with 25 mm (1 in) Gage Length

- » Ideal for a variety of applications including metals, plastics, rubbers, and composites
- » Low activation force without slipping on your specimen
- » High natural frequency and sensitivity make them useful in determining properties during failure
- » Select the -21 model for cryogenic testing down to -269°C (-452°F)
- » All part numbers include knife edges, quick attachment fixtures for flat and round specimens, and standard elastic attachment kit
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards



Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT	PART NUMBER
634.11F-2x	25 mm	+20% / -10%	33.0 mm	73.4 mm				
634.11F-21					-269°C / 65°C*	35 g	x	052-251-807
634.11F-24					-85°C / 120°C*	35 g	x	052-251-808
634.11F-25					-100°C / 175°C	45 g	x	052-251-809
634.11F-5x	25 mm	+20%	33.0 mm	73.4 mm				
634.11F-54					-85°C / 120°C*	35 g	x	052-251-811
634.11F-55					-100°C / 175°C	45 g	x	052-251-812
634.12F-2x	25 mm	+50% / -10%	33.0 mm	82.8 mm				
634.12F-21					-269°C / 65°C*	45 g	x	051-191-507
634.12F-24					-85°C / 120°C*	45 g	x	051-191-508
634.12F-25					-100°C / 175°C	55 g	x	051-191-509
634.12F-5x	25 mm	+50%	33.0 mm	82.8 mm				
634.12F-51					-269°C / 65°C*	45 g	x	051-191-510
634.12F-54					-85°C / 120°C*	45 g	x	051-191-511
634.12F-55					-100°C / 175°C	55 g	x	051-191-512

US Customary Units


MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT	PART NUMBER
634.11E-2X	1.00 in	+20% / -10%	1.30 in	2.89 in				
634.11E-21					-452°F / 150°F*	35 g	x	052-251-801
634.11E-24					-120°F / 250°F*	35 g	x	052-251-802
634.11E-25					-150°F / 350°F	45 g	x	052-251-803
634.11E-5X	1.00 in	+20%	1.30 in	2.82 in				
634.11E-54					-120°F / 250°F *	35 g	x	052-251-805
634.11E-55					-150°F / 350°F	45 g	x	052-251-806
634.12E-2X	1.00 in	+50% / -10%	1.30 in	3.26 in				
634.12E-21					-452°F / 150°F*	45 g	x	051-191-501
634.12E-24					-120°F / 250°F*	45 g	x	051-191-502
634.12E-25					-150°F / 350°F	55 g	x	051-191-503
634.12E-5X	1.00 in	+50%	1.30 in	3.26 in				
634.12E-54					-120°F / 250°F*	45 g	x	051-191-505

*These units may be used at 25°C (50°F) higher than the listed temperatures, but only for durations less than 24 hours.
See the MTS Fundamental Series 635 Extensometers on page 54 for an affordable monotonic only option.

Axial Extensometers

Low Strain, High Frequency Axial Extensometers with 25 mm (1 in) Gage Lengths

- » Small measuring ranges make Model 632.27 extensometers highly sensitive to detecting small displacement changes in ceramics and ceramic matrix composites
- » Extremely short arms and high natural frequencies ensure that this extensometer measures accurate strains at up to 150 Hz
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

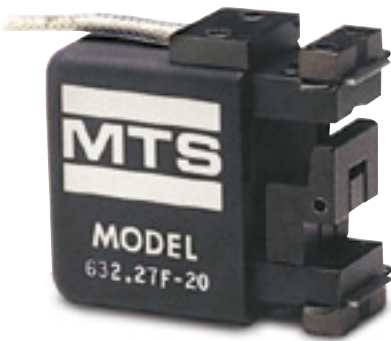

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
632.27F-2x	25 mm	±4%	33.0 mm	44.1 mm				
632.27F-20					-100°C / 150°C	60 g		047-401-904
632.27F-21					-269°C / 65 °C	60 g		047-401-905
632.27F-23					-100°C / 175°C	90 g		047-401-906
632.27F-3x	25 mm	±4%	33.0 mm	44.1 mm				
632.27F-23					-100°C / 150°C	60 g		047-402-004

US Customary Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
632.27E-2x	1 in	±4%	1.30 in	1.62 in				
632.27E-20					-150°F / 350°F	60 g		047-401-901
632.27E-3x	1 in	±2%	1.30 in	1.62 in				
632.27E-30					-150°F / 350°F	60 g		047-402-001

*These units may be used at 25°C (50 °F) higher than the listed temperatures, but only for durations less than 24 hours.

Testing at High Temperatures?

MTS offers a full range of high temperature extensometers that work with our chambers and furnaces.

See pages 60-62 for details.



Axial Extensometers

Axial Extensometers with 50 mm (2 in) Gage Length

- » Ideally suited for tensile and fatigue testing of plastics, rubber, and elastomers
- » Long travel and mechanical stops enable you to leave it on the specimen through failure without damaging the unit
- » Activation forces as low as 30 g prevent specimen damage without slipping
- » Select the -21 or -51 model for cryogenic testing down to -269°C (-452°F)
- » All part numbers include hardened knife edges, quick attachment fixtures for flat and round specimens, and standard elastic attachment kits
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards



Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
634.25F-2x	50 mm	+50% / -10%	58.4 mm	151.4 mm				
634.25F-21					-269°C / 65°C*	30 g	x	051-191-607
634.25F-24					-85°C / 120°C*	30 g	x	051-191-608
634.25F-25					-100°C / 175°C	45 g	x	051-191-609
634.25F-5x	50 mm	+50%	58.4 mm	151.4 mm				
634.25F-54					-85°C / 120°C*	30 g	x	051-191-611
634.25F-55					-100°C / 175°C	45 g	x	051-191-612

See the MTS Fundamental Series 635 Extensometers on page 56 for an affordable monotonic only option.

US Customary Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
634.25E-2x	2.00 in	+50% / -10%	2.3 in	5.96 in				
634.25E-21					-452°F / 150°F*	30 g	x	051-191-601
634.25E-24					-120°F / 250°F*	30 g	x	051-191-602
634.25E-25					-150°F / 350°F	45 g	x	051-191-603
634.25E-5x	2.00 in	+50%	2.3 in	5.96 in				
634.25E-51					-452°F / 150°F*	30 g	x	051-191-604
634.25E-55					-150°F / 350°F	45 g	x	051-191-606

*These units may be used at 25°C (50°F) higher than the listed temperatures, but only for durations less than 24 hours.

See the MTS Fundamental Series 635 Extensometers on page 56 for an affordable monotonic only option.

Axial Extensometers

Enhanced Travel Axial Extensometers

- » Designed for plastics, rubber, and elastomer testing, these extensometers can measure up to 100% strain
- » Utilizes a zero-stop for setting the initial gage length and simplifying mounting. Simply hold the arms together and attach the extensometer to the specimen
- » All part numbers include our quick attachment fixtures



Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
632.24F-5x	25 mm	+100%	33 mm	123.2 mm				
632.24F-50					-100°C / 150°C	30 g	x	047-401-602
634.28F-2x	50 mm	+100%	73.7 mm	155.6 mm				
634.28F-24					-85°C / 120°C	55 g	x	056-355-802

US Customary Units


MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
632.24E-5x	1 in	+100%	1.3 in	4.85 in				
632.24E-50					-150°F / 300°F	30 g	x	047-401-601
634.28E-2x	1 in	+100%	2.9 in	6.13 in				
634.28E-24					-120°F / 250°F	55 g	x	056-355-801

*These units may be used at 25°C (50°F) higher than the listed temperatures, but only for durations less than 24 hours.

Axial Extensometers

Axial Extensometers with Multiple Gage Lengths

- » Measure and control strain on a variety of different specimen types using this configurable extensometer
- » Features multiple gage lengths to provide flexibility to run a variety of tests without recalibration
 - 9 gage lengths available between 10 mm and 50 mm in the metric models
 - 6 gage lengths available in the US customary models
- » Easy gage length changeover reduces set-up time between tests
- » Quick attachment fixtures are included and can accommodate:
 - Round specimens with diameters from 2.5 to 14.0 mm (0.10 to 0.56 in)
 - Flat specimens ranging in thickness from 10.5 to 26.0 mm (0.42 to 1.02 in)
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B2 and ISO 9513 Class 0.5 standards

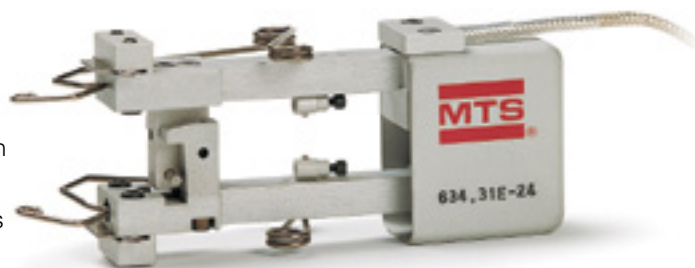

Metric Units

Available Gage Length	Measuring Range (% Strain, rounded to nearest 1 %)
10 mm	+40% / -20%
15 mm	+27% / -13%
20 mm	+20% / -10%
25 mm	+16% / -8%
30 mm	+13% / -7%
35 mm	+11% / -6%
40 mm	+10% / -5%
45 mm	+9% / -4 %
50 mm	+8% / -4%

US Customary Units

Available Gage Length	Measuring Range (% Strain, rounded to nearest 1 %)
0.500 in	+40% / -20%
0.750 in	+27% / -13%
1.000 in	+20% / -10%
1.500 in	+13% / -4%
1.750 in	+11% / -6%
2.000 in	+10% / -5%

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT (MIN/MAX)	LENGTH	TEMPERATURE RANGE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
634.31F-2x	See Above		32.5 mm / 59.4 mm	73.8 mm				
634.31F-21					-269°C / 65°C*	32 g	x	052-251-904
634.31F-24					-85°C / 120°C*	32 g	x	052-251-905
634.31F-25					-100°C / 175°C	40 g	x	052-251-906

Each 634.31E-2x extensometer includes components to set it to the six gage lengths listed.
The extensometer weighs between 36 and 42 g, depending on the mass of the components in use.

US Customary Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT (MIN/MAX)	LENGTH	TEMPERATURE RANGE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT FOR FLAT AND ROUND SPECIMENS	PART NUMBER
634.31E-2x	See Above		1.28 in / 2.34 in	2.9 in				
634.31E-21					-452°F / 150°F*	40 g	x	052-251-901
634.31E-24					-120°F / 250°F*	40 g	x	052-251-902
634.31E-25					-150°F / 350°F	50 g	x	052-251-903

Each 634.31E-2x extensometer includes components to set it to the six gage lengths listed.
The extensometer weighs between 36 and 42 g, depending on the mass of the components in use.

MTS Fundamental Series 635 Extensometers

Monotonic Extensometers for Tensile Testing

- » Ideal for measuring strain in tension applications
- » Optimized for monotonic testing with reliability, durability, and affordability in mind
- » May be left in place through specimen failure
- » All part numbers include hardened knife edges and quick attachment kits for flat and round specimens
- » Only available in metric gage lengths
- » Designed to meet requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards



Axial Extensometers for Monotonic Testing Only, MTS Fundamental Series 635

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE MIN/MAX	ACTIVATION FORCE	INCLUDES QUICK ATTACHMENT KIT	PART NUMBER
635.25F-05	25 mm	+20%	39.6 mm	77.5 mm	4°C / 50°C		x	057-863-506
635.50F-05	50 mm	+10%	59.2 mm	77.5 mm	4°C / 50°C		x	057-863-505
635.50F-10	50 mm	+20%	61.5 mm	77.5 mm	4°C / 50°C		x	057-863-504
635.50F-25	50 mm	+50%	69.1 mm	153.7 mm	4°C / 50°C		x	057-863-503
635.100F-10	100 mm	+10%	111.8 mm	77.5 mm	4°C / 50°C		x	057-863-502
635.100F-25	100 mm	+25%	111.8 mm	153.7 mm	4°C / 50°C		x	057-863-501

Accessories for Axial Extensometers

Gage Length Extenders

Expand the Range of Use of Your Existing Extensometer

- » Gage length extender kits modify the gage length of existing extensometers to perform strain measurements on specimens which have non-standard gage lengths
- » Travel of the extensometers is not affected by the addition of extenders, but the strain measuring range will be reduced
- » Each extender/extensometer combination requires separate calibration



Metric Units

EXTENSOMETER MODEL	GAGE LENGTH WITH EXTENSION	PART NUMBER
634.11/12 (F)		
	50 mm	051-191-901
	100 mm	051-191-902
	150 mm	051-191-903
	200 mm	051-191-904
	Kit of 50, 100, 150, & 200 mm extenders	051-191-905
634.25 (E)		
	100	051-191-908
	150	051-191-909
	200	051-191-910
	Kit of 100, 150, & 200 mm extenders	051-191-911

US Customary Units

EXTENSOMETER MODEL	GAGE LENGTH WITH EXTENSION	PART NUMBER
634.11/12 (E)		
	2 in	051-191-801
	4 in	051-191-802
	Kit of 2, 4, 6, and 8 in extenders	051-191-805
634.25 (E)		
	4 in	051-191-808
	Kit of 4, 6, and 8 in extenders	051-191-911

Extenders for additional models and gage lengths available upon request.

Quick Attachment Kits & Replacement Parts

- » Improve your productivity by buying a Quick Attachment Kit for your extensometer
- » Includes adapters for flat and round specimens

DESCRIPTION	PART NUMBER
Quick Attach Spring Set (includes 2 springs)	
For all 634.xx models except 634.31F	047-201-911
For Model 634.31F	047-201-920
Spring Retainer (holds 1 spring in place)	055-467-601

Extenders for additional models and gage lengths available upon request.

Accessories for Axial Extensometers

Specimen Attachment Kits


Large Diameter Specimen Attachment Kit

- » Recommended for specimens larger than 32 mm (1.25 in) in diameter
- » Provides a more effective attachment angle and increases the normal force pressing the extensometer onto the specimen
- » Increases stability of the extensometer on the specimen
- » Includes two remote spring attachment bracket assemblies that mount on the extensometer arms and an assortment of 16 tension springs
- » For use with models 632.11/.12/.25 and 634.11/.12/.25




Replacement Knife Edges for Extensometers

- » Replacements for those provided with extensometers
- » Stainless steel with a Rockwell hardness of Rc 58 to 60
- » Round specimen knife edges have a 90° included angle
- » Flat specimen knife edges provide three point specimen contact by means of one single point and one double point contact knife edge

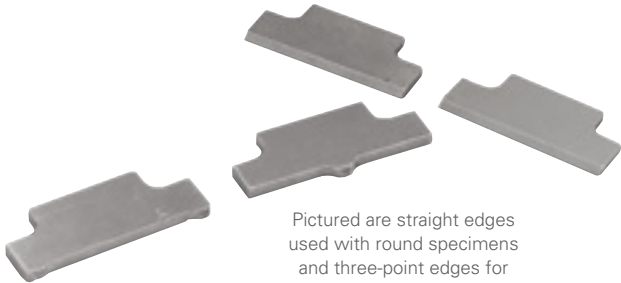

Flat Specimen Attachment Kit – Wire Form Method

- » Mounts a heavy duty, remote spring bracket assembly to the arms of the extensometer
- » Formed wire reacts with the force of the springs directly across from the knife edges

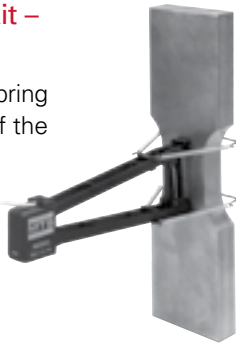

Flat Specimen Attachment Kit – Disc Method

- » Primarily used for specimens less than 12 mm (0.5 in) in width
- » Effectively makes specimen appear as a round specimen for extensometer mounting purposes



DESCRIPTION	634.XX PART NUMBER	632.XX PART NUMBER
Large Diameter Kit for 63x.11/.12/.25	053-063-001	038-863-801
Flat Specimen Attachment Kit – Disc Method	100-020-108	034-505-201
Flat Specimen Attachment Kit– Wire Form Method, for specimens up to 25 mm (1 in) wide	100-058-679	

Extenders for additional models and gage lengths available upon request.

Fundamental™ Automatic Extensometer (FAX)

Automatic Extensometer for Criterion® and Exceed® Universal Test Systems

- » Axial strain measurement for tensile testing on Universal Test Systems
- » High resolution of up to 0.05 μm
- » Axial with optional transverse deformation measurement
- » Designed for longevity and high-volume testing
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards



Two Mounting Options for the FAX

- » Fixed Mount attaches the FAX and moves simultaneously with the crossbeam of the Universal Test System
- » Rotational Mount allows the operator to quickly rotate the FAX out of the test area for easy access to the specimen



Fixed Mount



Rotational Mounts


Axial Model Specifications

							Transverse		
MODEL	ARM LENGTH	GAGE LENGTH	RELATIVE ERROR	MEASURING RANGE	AXIAL RESOLUTION	THICKNESS OR DIAMETER RANGE	SPECIMEN WIDTH	MEASURING RANGE	RESOLUTION
FAX1352									
Axial Only	350 mm	100-200 mm	±0.5%	0-80 mm	≤0.2 μm	Flat: 0.2-40 mm Round: 0.2-40 mm			
Axial & Transverse (FAX1352-T)	350 mm	100-200 mm	±0.5%	0-80 mm	≤0.2 μm	Flat: 0.2-5 mm	10-25 mm	Up to 25% of specimen width	0.05 μm
FAX1452									
Axial Only	450 mm	100-200 mm	±0.5%	0-100 mm	≤0.5 μm	Flat: 0.2-40 mm Round: 0.2-40 mm			
Axial & Transverse (FAX1452-T)	450 mm	100-200 mm	±0.5%	0-100 mm	≤0.5 μm	Flat: 0.2-5 mm	10-25 mm	Up to 25% of specimen width	0.05 μm

Advantage™ High Elongation Extensometer (AHX850)

High Elongation Extensometer for Criterion® Universal Test Systems

- » Accurately measure strain in specimens prone to large displacement
- » Displacement up to 850 mm (33.5 in)
- » Six preset gage lengths: 10 mm (0.4 in), 20 mm (0.8 in), 25 mm (1.0 in), 50 mm (2.0 in), 75 mm (2.95 in), 100 mm (3.9 in)
- » High resolution optical digital encoder eliminates noise, signal drift, and output changes
- » Balanced design with optimized arm and head weights
- » Swings away from test area when not in use
- » Optional fixed or rotational mount
- » Compatible with MTS Criterion® Universal Test Systems
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

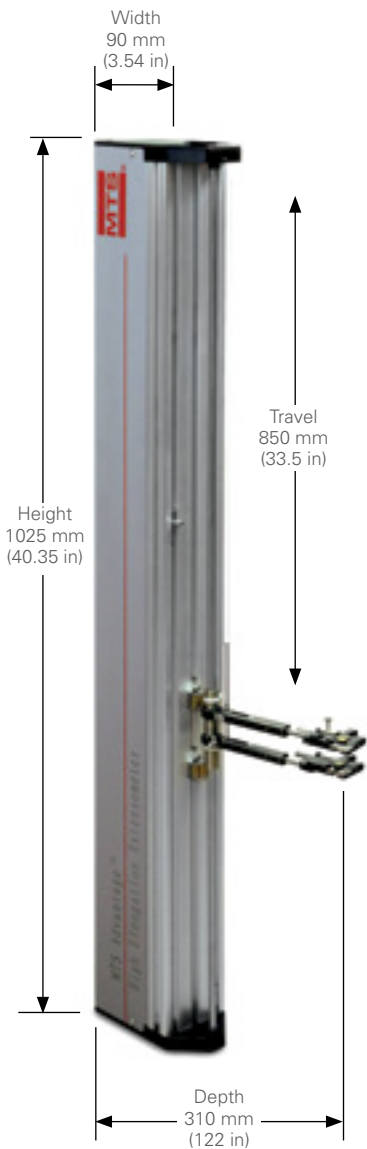


MODEL	GAGE LENGTH	MEASUREMENT RANGE	HEIGHT	LENGTH	WEIGHT	TEMPERATURE RANGE	GRIPPING FORCE	ARM OPENING
AHX850	10 mm (0.4 in) to 100 mm (3.9 in)	850 mm (33.5 in)	1025 mm (40.35 in)	310 mm (12.2 in)	6 kg (13.2 lbs)	5°C to 50°C (41°F to 122°F)	Adjustable spring	30 mm (1.18 in)

Long Travel Extensometer (LTX850)

High Elongation Extensometer for Exceed® Universal Test Systems

- » Used to measure displacement of materials, such as polymers and elastomers, which exhibit high elongation while in tension
- » Durable high strength aluminum structure
- » Dual independent digital input channels for upper and lower arms ensure accurate and reliable measurement
- » Changeable knife edges, adjustable gripping forces, and balanced head and arm weight allow smooth following of material strain change with minimal stickiness
- » Six preset gage lengths: 10 mm (0.4 in), 20 mm (0.8 in), 25 mm (1.0 in), 50 mm (2.0 in), 75 mm (2.95 in), 100 mm (3.9 in)
- » Choose fixed or rotational mounts
- » Compatible with Exceed® Universal Test Systems
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

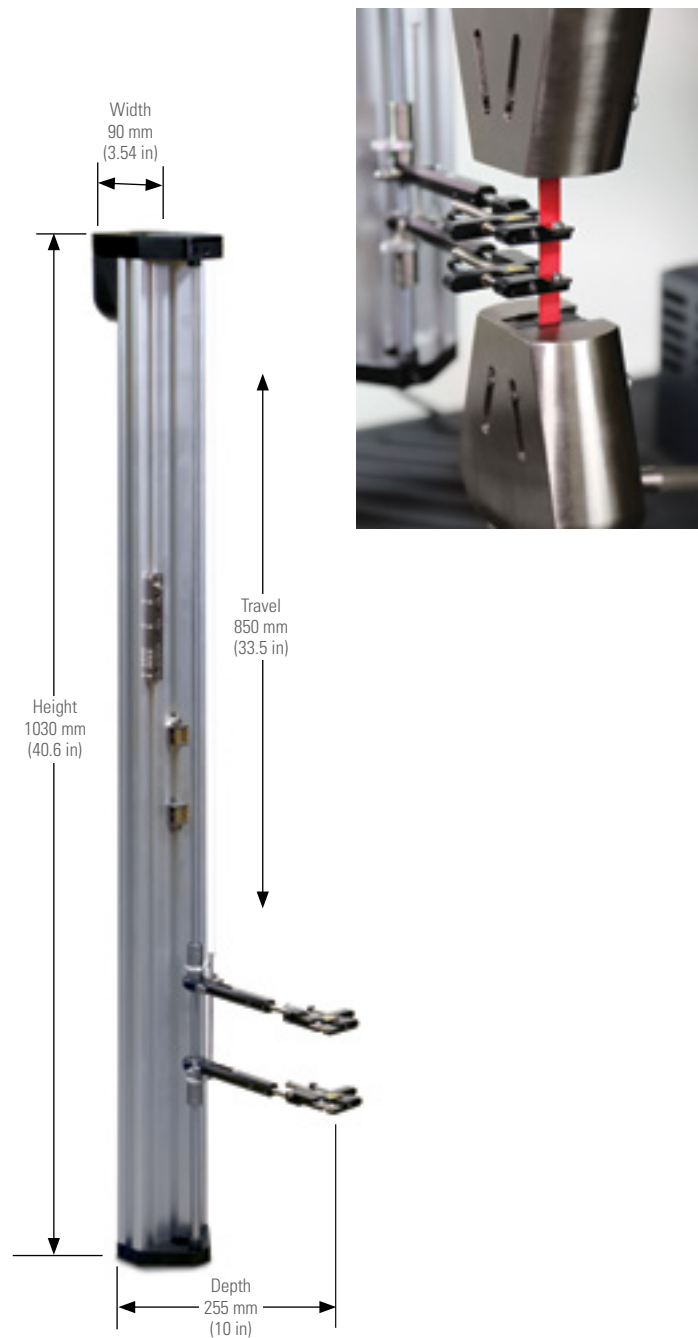


MODEL	GAGE LENGTH	MEASUREMENT RANGE	HEIGHT	LENGTH	WEIGHT	TEMPERATURE RANGE	GRIPPING FORCE	ALLOWED SPECIMEN SIZE
LTX850	10 mm (0.4 in) to 100 mm (3.9 in)	850 mm (33.5 in)	1030 mm (40.6 in)	90 mm (3.54 in)	5.3 Kg (11.7 lb)	5°C to 50°C (41°F to 122°F)	Adjustable spring	Width: 20 mm (0.8 in) Thickness: 30 mm (1.2 in)

High Temperature Extensometers

Air Cooled High Temperature Extensometers

- » Designed for measuring axial strain in high temperature tension, compression, and through-zero fatigue testing applications up to 1200°C (2200°F)
- » Includes ceramic extension rods that extend into the furnace and minimize rod creep, thermal expansion errors, and conduction losses from the specimen
- » Contact force varies by model. Select from models with 100 g, 300 g, 400 g, or 780 g contact forces
- » Air cooling required
- » Additional models and rod geometries available upon request
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

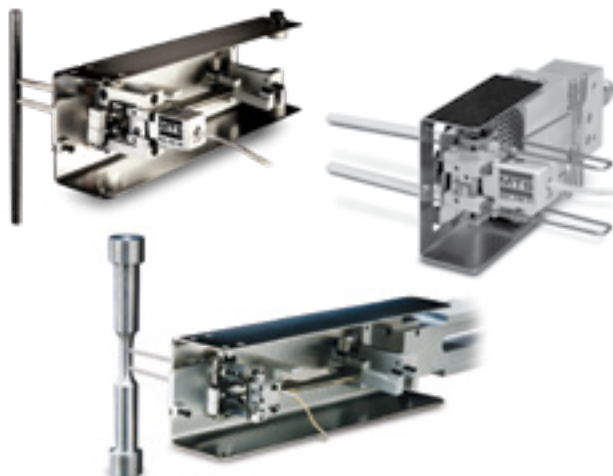

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	MAXIMUM TEMPERATURE	CONTACT FORCE	CERAMIC EXTENSION RODS INCLUDED	PART NUMBER
632.50F-08	10 mm	+20%/-10%	1200°C	400 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-538 ¹
632.54F-11	12 mm	+20%/-10%	1200°C	100 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-304 ¹
632.53F-14	12 mm	+20%/-10%	1200°C	300 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-105 ¹
				400 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-106 ¹
632.50F-04	12 mm	+20%/-10%	1200°C	780 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-402-539 ¹
632.54F-11	25 mm	+10% / - 5%	1200°C	100 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-303 ¹
632.50F-01	25 mm	+15% / -15%	1200°C	400 g	Qty 2, V-chisel Edge Rods, 94 mm (3.7 in) long	047-402-548 ¹
632.53F-11	25 mm	+10%/-5%	1200°C	300 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-111 ¹
				400 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-112 ¹

US Customary


MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	MAXIMUM TEMPERATURE	CONTACT FORCE	CERAMIC EXTENSION RODS INCLUDED	PART NUMBER
634.54E-14	0.50 in	+20% / -10%	2200°F	100 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-302 ¹
632.53E-14	0.50 in	+20%/-10%	2200°F	300 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-102 ¹
				400 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-103 ¹
632.50E-04	0.50 in	+20%/-10%	2200°F	780 g	Qty 2, V-chisel Edge Rods, 84 mm (3.35 in) long	047-402-525 ¹
632.53E-14	0.50 in	+20%/-10%	2200°F	300 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-103 ¹
				400 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	
632.54E-14	1.00 in	+10%/-5%	2200°F	100 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-301 ¹
632.53E-11	1.00 in	+10%/-5%	2200°F	300 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-108 ¹
				400 g	Qty 4, V-chisel Edge Rods, 85 mm (3.5 in) long	047-404-109 ¹

¹ Mount is not required when using with 653.xx.B furnace, for other installations contact MTS

High Temperature Extensometers

Water Cooled High Temperature Extensometers

- » Designed for measuring axial strain in high temperature tension, compression, and through-zero fatigue testing applications up to 1200°C (2200°F)
- » Water cooling is more effective at keeping the extensometer body at a constant temperature and reducing noise in the data measured
- » Includes ceramic extension rods that extend into the furnace and minimize rod creep, thermal expansion errors, and conduction losses from the specimen
- » Contact force varies by model. Select from models with 400 g or 780 g contact forces.
- » Water cooling required
- » Additional models and rod geometries available upon request
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

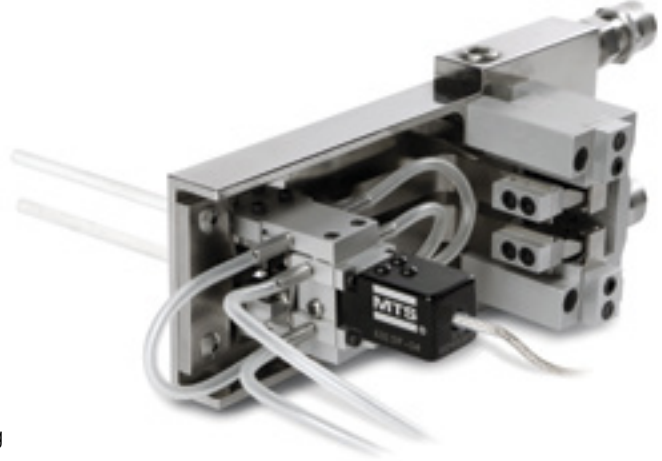

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	MAXIMUM TEMPERATURE	CONTACT FORCE	CERAMIC EXTENSION RODS INCLUDED	PART NUMBER
632.51F-04	12 mm	+20% / -10%	1200°C	780 g	Qty 2, V-chisel Edge Rods, 84 mm (3.35 in) long	047-402-651 ¹
632.51F-01	25 mm	+15% / -10%	1200°C	400 g	Qty 2, V-chisel Edge Rods, 84 mm (3.35 in) long	047-402-652 ¹

US Customary

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	MAXIMUM TEMPERATURE	CONTACT FORCE	CERAMIC EXTENSION RODS INCLUDED	PART NUMBER
632.51E-04	0.5 in	+20% / -10%	2200°F	780 g	Qty 2, V-chisel Edge Rods, 84 mm (3.35 in) long	047-402-650 ¹

¹ Mount is not required when using with 653.xx.B furnace, for other installations contact MTS.

High Temperature Extensometers

Elevated Temperature Axial Extensometer for Use in Chambers

- » Designed for measuring axial strain in elevated temperature testing environments without external cooling
- » Provides stable output, with low creep, for both tension and fatigue testing applications at temperatures from 40°C to 540°C (100° to 1000°F) and up to 85% relative humidity.
- » Fatigue rated up to 30 Hz
- » Compatible with controllers that use 494 or 493 series hardware. Not available for MTS electromechanical systems
- » Packed in a rugged storage case that contains the instrument, converter module, spare parts, springs, attachment devices, and tools



Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX ³	ACTIVATION FORCE	INCLUDES ATTACHMENT KIT
633.11M-15	25.00 mm	+16 / -8%	34.7 mm	35.8 mm	40°C / 540 °C	60 g	x

US Customary

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX ³	ACTIVATION FORCE	INCLUDES ATTACHMENT KIT
633.11L-15	1.000 in	+16 / -8%	1.4 in	2.2 in	100°F to 1000°F	60 g	x

¹ Relative humidity in excess of 85% may cause erroneous output. This erroneous output will disappear when chamber is elevated above ambient.


² Maximum operating frequency stated for sinusoidal displacements of 0.05 mm (0.002 in) amplitude or less.

³ May be used 25°C (50°F) higher than specified for short durations less than 24 hours.

Specialty Extensometers

Averaging Axial Extensometer

- » Simultaneously measures axial deflection on the opposite sides of the specimen and then sums those measurements to provide a single average strain output
- » Measuring strain on both sides of the specimen reduces errors in strain measurements caused by bending strains, resulting in more accurate elastic modulus values
- » Units can be used to perform testing on solid or tubular composite metal specimens of almost any configuration
- » Easily adjustable from 3 mm to 32 mm (0.12 in to 1.26 in) in diameter for round specimens
- » Easily adjustable for up to 51 mm (2 in) width for flat specimens
- » Low strain range detects small displacement changes in ceramics and composites
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

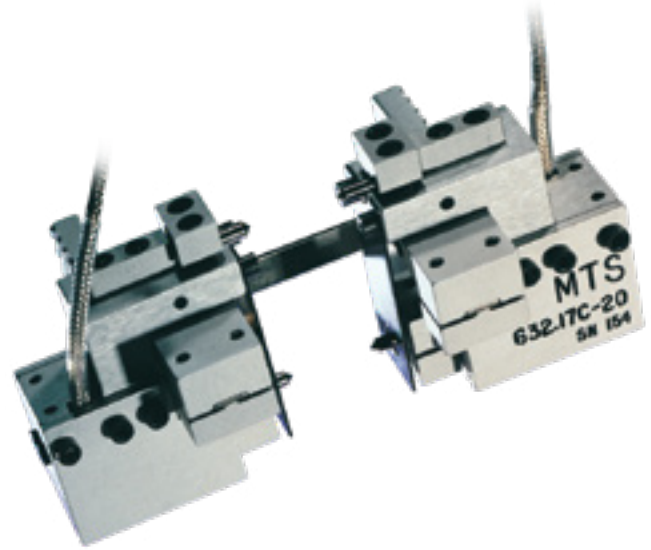

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	DEPTH	TEMPERATURE RANGE MIN/MAX*	PART NUMBER
632.17F-x0							
632.17F-40	10 ±0.05 mm	+12% / -5%	63.5 mm	101.6 mm	38.1 mm	-100°C / 150°C	047-404-802
632.17F-20	25 ±0.05 mm	+5% / -2%	63.5 mm	101.6 mm	38.1 mm	-100°C / 150°C	047-401-102

US Customary Units


MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	DEPTH	TEMPERATURE RANGE MIN/MAX*	PART NUMBER
632.17E-x0							
632.17E-30	0.5 ±0.002 in	+10% / -4%	2.5 in	4.0 in	1.5 in	-150°F / 300°F	047-404-801
632.17E-20	1.0 ±0.002 in	+5% / -2%	2.5 in	4.0 in	1.5 in	-150°F / 300° F	047-401-101

*These units may be used at 25°C (50 °F) higher than the listed temperatures, but only for durations less than 24 hours.

Specialty Extensometers

Biaxial Extensometer

- » Measure the cross-sectional, diametral, or average axial strain with a single extensometer
- » Readily attach to many shapes and sizes of specimens
- » If you are worried about mechanical crosstalk between the axial and transverse channels, select Model 632.85x-05. This model operates only with conical point contacts.
- » If you have very thin specimens (i.e. sheet metal), select Model 632.85x-14. This model is designed for use with either conical point or v-chisel specimen contacts. Two output signals require separate conditioning for the averaged axial and the cross-sectional signals.
- » Meets or exceeds requirements for calibration according to ASTM E83 Class B1 and ISO 9513 Class 0.5 standards

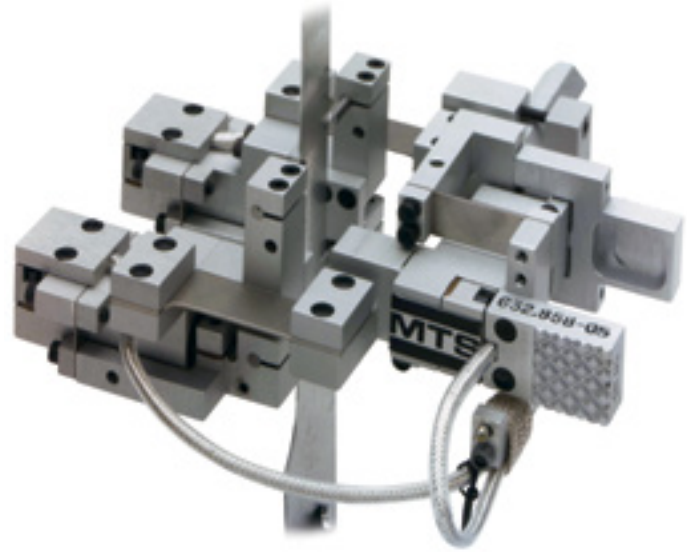

Metric Units

MODEL	GAGE LENGTH	AXIAL TRAVEL	TRANSVERSE TRAVEL	HEIGHT	LENGTH	DEPTH	TEMPERATURE RANGE MIN/MAX*	PART NUMBER
632.85F-xx								
632.85F-05	25±0.05 mm	+1.25 to -0.50 mm	±0.5 mm	101.6 mm	38.1 mm	101.6 mm	-100°C / 150°C	047-403-502
632.85F-14	25±0.05 mm	+1.25 to -0.50 mm	0.5 mm	101.6 mm	38.1 mm	101.6 mm	-100°C / 150°C	047-403-506

US Customary Units

MODEL	GAGE LENGTH	AXIAL TRAVEL	TRANSVERSE TRAVEL	HEIGHT	LENGTH	DEPTH	TEMPERATURE RANGE MIN/MAX*	PART NUMBER
632.85E-xx								
632.85E-05	1.000±0.002 in	+0.050 to -0.020 in	±0.02 in	4.0 in	1.5 in	4.0 in	-150°F / 300°F	047-403-501
632.85E-14	1.000±0.002 in	+0.050 to -0.020 in	±0.02 in	4.0 in	1.5 in	4.0 in	-150°F / 300°F	047-403-505

*These units may be used at 25°C (50 °F) higher than the listed temperatures, but only for durations less than 24 hours.

Specialty Extensometers

Cross-Sectional Strain Extensometer

- » Measures cross-sectional strain and when used with a long gage length or enhanced travel extensometers (pages 51 and 52), it is ideal for determining the “R” value for plastics and sheet metal
- » Free floating feature enables it to travel with the specimen as it is elongated during axial loading



Metric Units

MODEL	GAGE WIDTH	TRAVEL	SPECIMEN THICKNESS	TEMPERATURE RANGE MIN/MAX*	CLAMP FORCE	PART NUMBER
632.23F-x0						
632.23F-20	12.7 ± 0.5 mm	-4 to 0 mm	0.5 to 5 mm	-100°C / 150°C		047-401-502
632.23F-30	20 ± 1 mm	-4 to 0 mm	0.5 to 5 mm	-100°C / 150°C		047-401-504
632.23F-30	25 ± 1 mm	-4 to 0 mm	0.5 to 5 mm	-100°C / 150°C		047-401-506

US Customary Units

MODEL	GAGE WIDTH	TRAVEL	SPECIMEN THICKNESS	TEMPERATURE RANGE MIN/MAX*	CLAMP FORCE	PART NUMBER
632.23E-x0						
632.23E-20	0.5 ± 0.02 in	-0.16 to 0 in	0.02 to 0.2 in	-150°F / 300°F		047-401-501

*These units may be used at 25°C (50 °F) higher than the listed temperatures, but only for durations less than 24 hours.

Specialty Extensometers

Diametral Extensometer

- » Ideal for measuring diametral strain in tension or compression
- » Designed for measuring cross sectional area changes on round specimens and helping to measure Poisson's ratio
- » Models available for testing at ambient, cryogenic, and elevated temperatures
- » Contact MTS for options for transverse measurements of specimens with rectangular cross sections



Metric Units

MODEL	GAGE DIAMETER	TRAVEL	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	CLAMP FORCE	PART NUMBER
632.18F-2x	6.1 mm to 26 mm	±2 mm	12.7 mm	191 mm			
632.18F-20					-100°C / 150 °C	400 g	047-401-204
632.18F-21					-265°C / 65° C	400 g	047-401-205
632.19F-2x	3.6 mm to 13 mm	±1 mm	12.7 mm	191 mm			
632.19F-20					-100°C / 150 °C	400 g	047-401-304
632.20F-2x	2 mm to 7.9 mm	±0.5 mm	12.7 mm	191 mm			
632.20F-20					-100°C / 150 °C	400 g	047-401-404

US Customary Units


MODEL	GAGE DIAMETER	TRAVEL	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	CLAMP FORCE	PART NUMBER
632.18E-2x	0.24 in to 1.02 in	±0.08 in	0.5 in	7.5 in			
632.18E-20					-150°F / 300 °F	400 g	047-401-201
632.19E-2x	0.14 in to 0.52 in	±0.04 in	0.5 in	7.5 in			
632.19E-20					-150°F / 300 °F	400 g	047-401-301
632.20E-2x	0.14 in to 0.31 in	±0.02 in	0.5 in	7.5 in			
632.20E-23					-150°F / 300 °F	400 g	047-401-403

*These units may be used at 25°C (50 °F) higher than the listed temperatures, but only for durations less than 24 hours.

Specialty Extensometers

Immersible Extensometer

- » Provide accurate strain measurement while totally immersed in water or saline solution
- » Designed for both tension and fatigue testing up to 30 Hz
- » Accurately translates specimen displacement to the hermetically sealed LVDT using a patented parallel flexure system
- » Gage Length extenders to extend gage length to 200 mm (8.0 in)
- » AC conditioning required

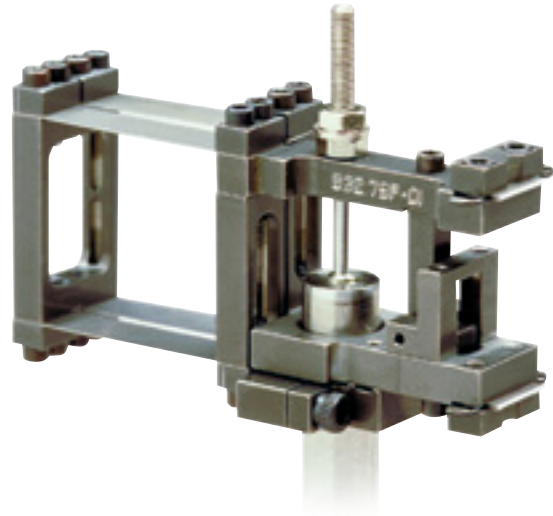

Metric Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	PART NUMBER
632.79F-01	25 mm	±20%	76 mm	69.1 mm	-15°C to 85°C	60 g	050-855-602

US Customary Units

MODEL	GAGE LENGTH	MEASURING RANGE (% STRAIN)	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX*	ACTIVATION FORCE	PART NUMBER
632.79E-01	1 in	±20%	3.0 in	2.7 in	4°F to 185°F	60 g	050-855-601

MTS Advantage™ Video Extensometer

The MTS Advantage Video Extensometer delivers the flexibility, accuracy and repeatability needed in applications where it is difficult to measure strain using conventional means. This extensometer combines several features that make it an outstanding choice to measure strain, rotation and displacement of plastics, metals, composites and high-elongation materials.


- » Magnetic-return support arm for easy specimen access and precise positioning
- » Sixteen quick-attach measurement heads supporting a wide array of gage lengths and strain ranges
- » Pass/fail calibration verification blocks and software calibration wizard to confirm setup to desired ASTM E83 or ISO 9513 class
- » Mounting options on left, right, front or back of load frame

Advantages of Video Extensometry

- » Measures samples through failure without the risk of damage to your extensometer
- » Perfect for fragile, brittle or irregularly shaped specimens
- » Allows up to 200 measurements in real-time
- » Allows repeatable measurements to be made after the physical test is complete
- » Ideal option for measuring in high-temperature or submerged


User-friendly Software Enables Quick Measurements

- » Drag and drop a full selection of measurement devices onto your specimen
- » Select from 10 different measurement devices to determine position, displacement, distance or rotation
- » Virtual strain measurement devices convert motion into axial strain, dual average strain, Poisson's ratio, shear strain and 2D strain maps
- » Understand results quickly with intuitive graphical displays
- » Ethernet communication between systems sends ± 10 V signals to the MTS TestSuite TW software

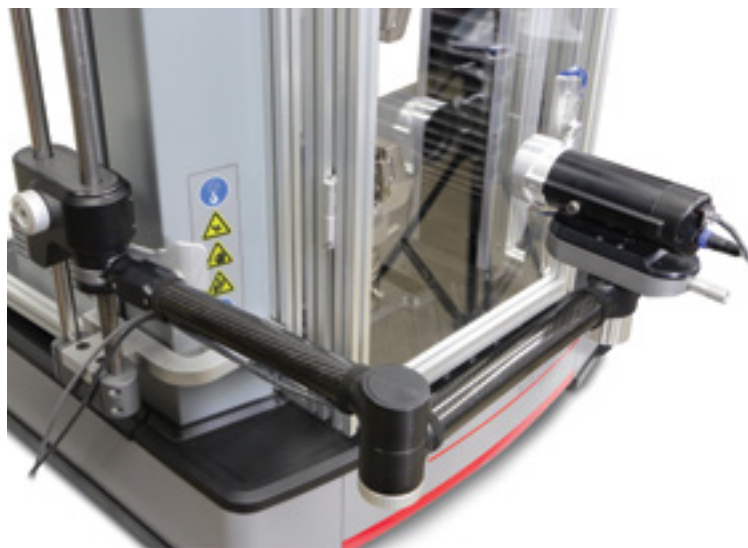

Upgrade to 2D Digital Image Correlation (DIC)

- » Create planar strain maps with real-time measurements
- » Output 2D contour plots where color gradients show levels of strain
- » Define the x-axis and show E_{xx}, E_{yy} or E_{xy}
- » Identify areas of high stress, crack option or other discontinuities
- » Toggle between all measurement options without reprocessing a test

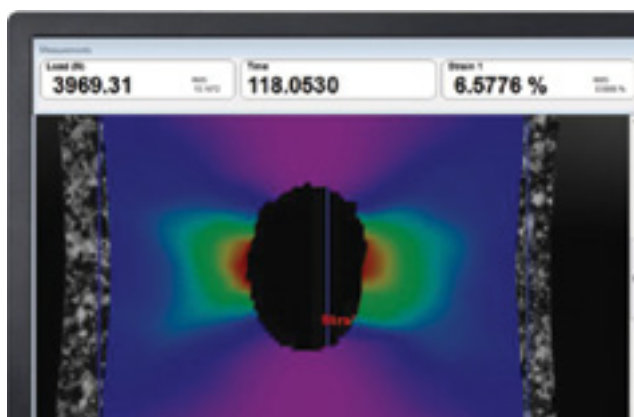

2D Strain Map

Measurement Head Specifications

Each AVX Measurement Head includes an inspected, traceable validation block, which can be used to confirm calibration at the beginning of a series of tests. Validation blocks are serialized and include an ISO 17025 Accreditation Calibration Report.

XT-100 Series Measurement Heads

Ideal for determining higher strain (>10%) materials properties such as yield point & elongation, and for long gage lengths. All models are capable of meeting Class B-2 (ASTM E-83) & Class 0.5 (ISO 9513) at the specified gage lengths and strain ranges (>10%). They are also capable of meeting Class B-1 (ASTM E-83) over most of their operating range (gage lengths where maximum axial tensile strain less than 600%).

The XT-100 series operate at measurement rates from 0.1 - 500 Hz.



AVX Measurement Head	Maximum Axial Tensile Strain Range (%) at Specified Gage Length (mm) ¹					Maximum Axial Compressive Strain Range (%) at Specified Gage Length (mm) ¹					Maximum Transverse Gage Length (mm) ²	Typical Extension Resolution (μm) ³	Minimum Specimen Width for Measurements (mm) ⁴		Maximum Tracking Speed (mm/min) ²
	10	25	50	100	200	10	25	50	100	200			Axial	Transverse	
XT-101	350	100	-	-	-	80	80	50	-	-	14	0.2	1.5	4.0	1500
XT-102	600	200	60	-	-	80	80	80	-	-	22	0.3	2.3	6.4	2500
XT-103	900	320	120	30	-	-	80	80	80	-	32	0.4	3.4	9.4	4000
XT-104	1000	490	210	70	-	-	80	80	80	-	46	0.6	4.8	13	5500
XT-105	-	800	370	150	40	-	-	80	80	80	74	1.0	7.6	21	9000
XT-106	-	1000	510	220	70	-	-	80	80	80	97	1.3	10	28	11500
XT-107	-	-	800	380	160	-	-	-	80	80	149	1.9	15	43	18500

1. Strain ranges assume a distance between the grips of twice the GL, with a centrally positioned gage. Actual strain ranges may be greater or less than these values, depending on gage positioning, grip separation and specimen behavior.

2. Maximum transverse GL and tracking speed is quoted at a measurement rate of 100 Hz.

3. Resolution is based on typical lab performance.

4. For Class B-1, minimum specimen width for transverse gage lengths should be increased by 65%.

XT-200 Series Measurement Heads

Ideal for determining low strain materials properties (from 0.01%), such as tensile & compressive modulus, Poisson's ratio & R-value. All models are capable of meeting Class B-1 (ASTM E-83) & Class 0.5 (ISO 9513) at the specified gage lengths and strain ranges.

The XT-200 series are our highest accuracy measurement heads. These models operate at measurement rates from 0.1 - 30 Hz. The XT-250 series are suitable for many high accuracy dynamic applications, and operate at measurement rates from 0.1 - 500 Hz.



AVX Measurement Head	Maximum Axial Tensile Strain Range (%) at Specified Gage Length (mm) ¹			Maximum Axial Compressive Strain Range (%) at Specified Gage Length (mm) ¹			Maximum Transverse Gage Length (mm) ²	Typical Extension Resolution (μm) ³	Minimum Specimen Width for Measurements (mm)		Maximum Tracking Speed (mm/min) ²
	10	25	50	10	25	50			Axial	Transverse	
XT-201	75	-	-	80	-	-	17	0.05	0.4	1.7	50
XT-202	120	10	-	80	30	-	19	0.12	0.9	4.3	150
XT-203	190	35	-	80	80	-	31	0.09	0.7	3.2	100
XT-204	300	80	5	80	80	25	45	0.13	1.0	4.6	150
XT-205	460	145	40	80	80	80	63	0.18	1.4	6.4	250
XT-251	100	5	-	80	10	-	6	0.08	0.6	2.8	700
XT-253	250	60	-	80	80	-	10	0.14	1.1	5.0	1300
XT-254	390	120	25	80	80	80	15	0.20	1.6	7.3	1900
XT-255	580	190	65	80	80	80	21	0.28	2.2	10.3	2600

1. Specified strain ranges assume a distance between the grips of twice the GL, with a centrally positioned gage. Actual strain ranges may be greater or less than these values, depending on gage positioning, grip separation and specimen behavior.

2. Maximum transverse GL and tracking speed is quoted at a measurement rate of 15Hz (XT-20x series) and 100 Hz (XT-25x series).

3. Resolution based on typical lab performance.

Laser Extensometers

MTS offers the LX 500 and LX 1500 Extensometers. These high performance devices are completely self-contained, easily transportable, and easy to use with their digital display and four-button keypad. They can be operated from the rear control panel or remotely via the RS-232 interface. The unique scanning laser beam technique, used to measure elongation entirely from one side of the specimen, eliminates the need for a separate receiver mounted on the opposite side of the specimen.

FEATURES

- » Selectable averaging of 2 to 512 scans
- » Measures the initial gage length of the specimen directly
- » Can measure strain simultaneously over three segments within the specimen's gage length
- » Meets ASTM B1 extensometer classification (25 mm/1 in and longer gage length)
- » Measurements are updated 100 times per second.
- » Analog output port for closed loop strain control or for input to a data acquisition board, chart recorder, and so forth
- » Parallel beam design accommodates operation through chamber glass, a clear bath, water, and so forth
- » Certified for compliance with the Center for Devices and Radiological Health as a Class II product




MODEL	MAXIMUM TRAVEL	RESOLUTION	NON-LINEARITY MAXIMUM	REPEATABILITY	ANALOG OUTPUT	PART NUMBERS* 115V±10%	PART NUMBERS* 230V±10%
LX 500	8-127 mm (0.3-5.0 in)	0.001 mm (0.0001 in)	±0.025 mm (±0.001 in)	±0.003 mm (±0.0001 in)	16 bit selectable range, ±10 V full scale	055-283-101	055-283-102
LX 1500	8-381 mm (0.3-15 in)	0.01 mm (0.001 in)	±0.15 mm (±0.006 in)	±0.03 mm (±0.001 in)	16 bit selectable range, ±10 V full scale	056-103-001	056-103-002

*Includes tripod, mounting plates, instruction manual, line cord, and tripod head.

ACCESSORIES	MEASUREMENT	PART NUMBERS
Retroreflective Tape	3.2 mm x 7.6 m (1/8 in x 25 ft)	100-039-188
Retroreflective Tape	6.4 mm x 7.6 m (1/4 in x 25 ft)	100-039-189
Retroreflective Tape, Submersible (water)	3.2 mm x 7.6 m (1/8 in x 25 ft)	100-214-873
Retroreflective Target, High Temp. (900°F), Non-adhesive	2.5 mm x 76.2 mm (0.1 in x 3 in)	100-039-190
Laser Cart, Vertical Adjustment	913 mm x 1556 mm (36 in x 61 in)	100-165-138

Other Applications

With over 101 different models, MTS has the largest offering of extensometers in the world. In fact, here's a list of some of our other models that you won't find in this catalog. Contact us if you're interested in any of these capabilities!

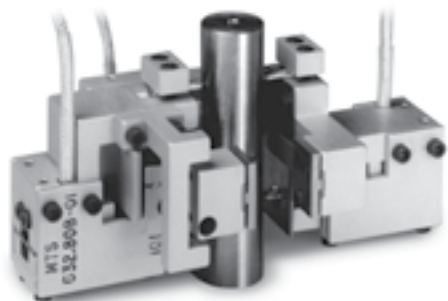

HOSTILE ENVIRONMENT DISPLACEMENT GAGE (MODEL 632.10)

Use for salt fog or boiling water reaction simulations. Standard temperature range to 290°C (550°F), with options to 400°C (750°F).

BIOMEDICAL SOFT TISSUE OR SIMILAR MATERIAL (MODEL 632.32)

Ultra-low contact force with serrated knife edges.

HIGH-TEMPERATURE AXIAL APPLICATIONS (MODELS 632.50/51/59)

Many options available. Air-cooled, water-cooled, vacuum, and temperatures to 2300°C (4100°F).

HIGH-TEMPERATURE DIAMETRAL/TRANSVERSE APPLICATIONS (MODELS 632.60/61)

Standard models to 1000°C (2200°F), with options available to 1500°C (2700°F).



HIGH-TEMPERATURE AXIAL TORSIONAL APPLICATIONS (MODELS 632.68 AND 632.80)



Simultaneous measurement of axial and torsional strain. The Model 632.80 has a temperature range of -100°C to 150°C (-150°F to 300°F), while the high-temperature version Model 632.68 is rated up to 1200°C (2200°F).

AXIAL ROCK MECHANICS APPLICATIONS (MODELS 632.90/94)

Models available for use in air, or in high pressure, single measurement, and averaged axial applications.

CIRCUMFERENTIAL APPLICATIONS (MODEL 632.21/92)

Unique roller chain design provides low hysteresis.

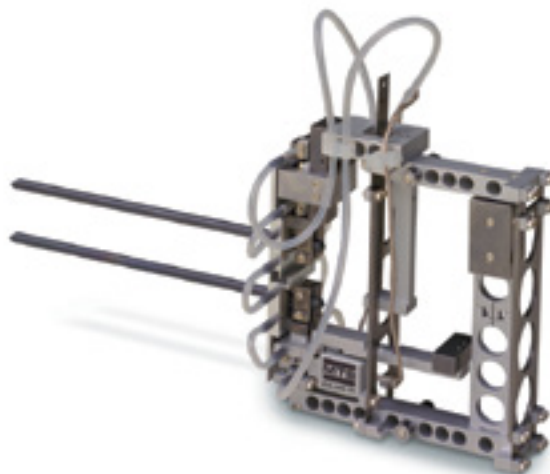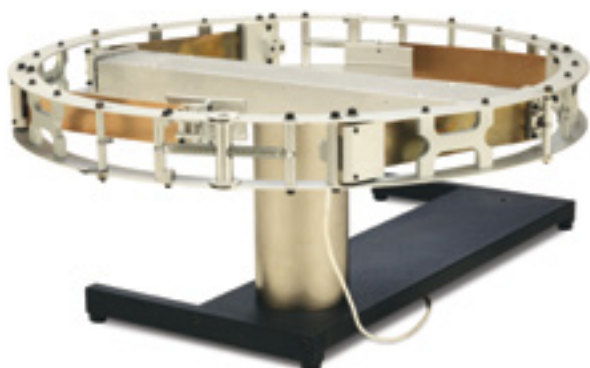

Clip-On Displacement Gages

- » Double cantilever displacement-type gages are primarily used in fracture mechanics testing, K_{IC} , J_{IC} , crack growth and R-curve determination
- » Useful for measuring the displacements of components and specimens with odd geometries, performing elongation studies, and testing pipes, bars and rods
- » Compatible with a variety of specimens including compact tension, WOL, round compact tension, C-shaped, bend, and other common specimens or panels
- » Can be immersed in a wide variety of non-conductive fluids without special precautions
- » Model 632.02 Clip-On gages when testing to ASTM E399
- » Model 632.03 Clip-On Gages are our most accurate featuring improved slot geometry and three-point contact knife edges. Options are available with longer travel for J-integral and crack growth tests
- » Replacement knife edges available upon request
- » Additional models available, contact your sales representative for more options



Standard equipment includes the gage and storage case, two reusable hardened knife edges for attachment to the specimen, and an instruction manual.

Pictured are 632.02 and 632.03/05 Clip-On Gage knife edges.

Metric Units

MODEL/OPT	GAGE LENGTH	TRAVEL	MAXIMUM FREQUENCY	COMPRESSED FORCE	HEIGHT/LENGTH	TEMPERATURE MIN/MAX	PART NUMBER
632.02F-2x	5 mm	+3 / -1 mm	50 Hz	2500 g	30.5 / 48.3 mm		
632.02F-20						-100°C / 150°C	047-400-104
632.03-2x	12 mm	4 mm	100 Hz	2500 g			
632.03F-20						-100°C / 150°C	047-400-204
632.03-3X							
Option 001	2 mm	4 mm				-100°C / 150°C	047-400-310
Option 002	2.5 mm	4 mm				-100°C / 150°C	047-400-311

US Customary Units

MODEL/OPT	GAGE LENGTH	TRAVEL	MAXIMUM FREQUENCY	COMPRESSED FORCE	HEIGHT/LENGTH	TEMPERATURE MIN/MAX	PART NUMBER
632.02E-2x / -	0.2000 in	+0.100 / -0.05 in	50 Hz	2500 g	1.2 / 1.9 in		
632.02E-20						-150°F / 300°F	047-400-101
632.03E-20	0.475 in	0.15 in	100 Hz	2500 g			
						-150°F / 300°F	047-400-201
632.03E-3X							
Option 001	0.075 in	0.15 in				-150°F / 300°F	047-400-301
Option 002	0.10 in	0.15 in				-150°F / 300°F	047-400-302


Displacement Gages

- » Versatile design measures specimens, fixtures, shafts, or the movement of any other component during a test
- » Designed for use where small deformations must be measured such as bend tests or tests with unusual geometries
- » Over-travel protection included:
 - Gage arm releases if it travels too far or when a side load is applied.
 - Gage arm can be reset without changing the body position
- » Comes with an adjustable mounting block
- » Displacement measurements can be made with the gage in contact against a specimen or an active component in the force train



MODEL/OPT	TRAVEL	ZERO ADJUST	FORCE AT ZERO	HEIGHT	LENGTH	TEMPERATURE RANGE MIN/MAX	PART NUMBER
632.06H-2x / --	±4 mm (±0.16 in)	8 mm (0.30 in)	150 g	101 mm (4 in)	101 mm (4 in)		
632.06H-20						-100°C / 150°C (-150°F / 300°F)	047-400-501
632.06H-23						-100°C / 175°C (-150°F / 350°F)	047-400-503
632.06H-3x/ 001	±2 mm (±0.08 in)	8 mm (0.30 in)	150 g	101 mm (4 in)	101 mm (4 in)		
632.06H-30						-100°C / 150°C (-150°F / 300°F)	047-400-504
632.06H-33						-100°C / 175°C (-150°F / 350°F)	047-400-506
632.06H-3x/ 002	±1 mm (±0.04 in)	8 mm (0.30 in)	250 g	101 mm (4 in)	69 mm (2.7 in)		
632.06H-30						-100°C / 150°C (-150°F / 300°F)	047-400-507
632.06H-3x/ 003	±6 mm (±0.25 in)	8 mm (0.30 in)	150 g	101 mm (4 in)	101 mm (4 in)		
632.06H-30						-100°C / 150°C (-150°F / 300°F)	047-400-510
632.06H-33						-100°C / 175°C (-150°F / 350°F)	047-400-512
632.06H-3x/ 004	±8 mm (±0.32 in)	8 mm (0.30 in)	100 g	101 mm (4 in)	140 mm (5.5 in)		
632.06H-30						-100°C / 150°C (-150°F / 300°F)	047-400-513
632.06H-3x/ 005	±12.5 mm (±0.50 in)	8 mm (0.30 in)	80 g	101 mm (4 in)	190 mm (7.5 in)		
632.06H-30						-100°C / 150°C (-150°F / 300°F)	047-400-516
632.06H-31						-270°C / 65°C (-452°F / 150°F)	047-400-517
632.06H-33						-100°C / 175°C (-150°F / 350°F)	047-400-518
632.06H-3x/ 006	±25 mm (±0.50 in)	8 mm (0.30 in)	35 g	101 mm (4 in)	343 mm (13.5 in)		
632.06H-30						-100°C / 150°C (-150°F / 300°F)	047-400-519
632.06H-31						-270°C / 65°C (-452°F / 150°F)	047-400-520

Extensometer Calibrator

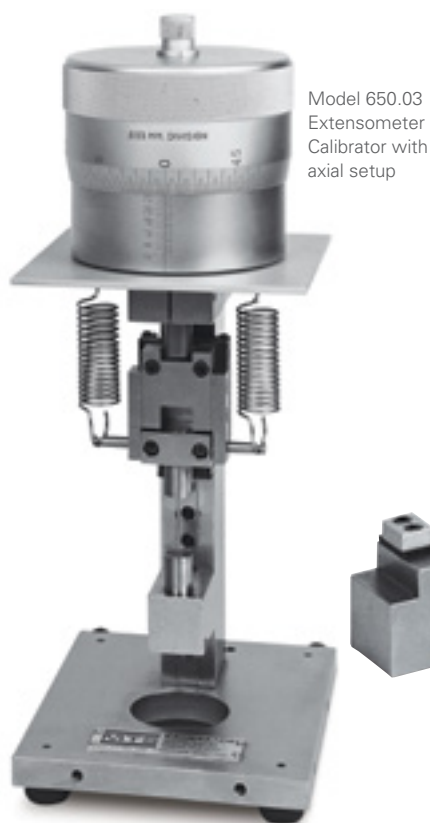

Model 650.03
Extensometer
Calibrator with
axial setup

Model 650.03

The MTS Model 650.03 Extensometer Calibrator is a precise instrument designed and built for rugged, long-lasting durability. It gives you freedom from wear; vital for the high accuracy required for precise calibrator verification of high resolution instruments.

The micrometer head has a spindle with a large diameter lead screw to reduce wear and resulting backlash and to provide higher strength and stiffness. Larger mounting diameters reduce the unit stress required in

mounting the micrometer head to the frame. The micrometer head is large for easy readability, accurate adjustment, and improved repeatability of gage calibration.

1. Standard Axial includes Extensometer models 632.11/12/13/24/25/26/27/31, 634.11/12/25/31 set-up with straight knife edges and attachment for round specimens. Some specific knife edges or attachment kits may require other sizes of specimen adapters. The installation drawing for the extensometer will identify the specimen size range.

2. Standard High Temperature Axial includes 632.53/54 with standard length rods, 632.50/51 with length options 075, 076, 077, 078, 079 and 080. Calibration require rods with either conical point (with Dimpled Specimen) or V-chisel rod end geometry (with optional 9.5 mm round).

EXTENSOMETER MODEL	632.02	632.03/05	STANDARD AXIAL (1)	EXTENDED GL AXIAL	632.06	632.06 OPTION 006	632.18/19/20	632.29	632.17/85 AXIAL	632.85 TRANSVERSE	STANDARD HI TEMP AXIAL (2)	632.92, 632.11 OR .12 W/ 632.21 KIT
--------------------	--------	-----------	--------------------	-------------------	--------	-------------------	--------------	--------	-----------------	-------------------	----------------------------	-------------------------------------

650.03 Calibrator frame in fitted case	031-508-302	Required for all
Micrometer head, U.S. customary (0.00005 in graduation)	031-508-205	Required for all, choose one
Micrometer head, Metric (0.001 mm graduation)	031-508-206	Required for all, choose one

CALIBRATOR KITS

Calibrator Kit, Basic Specimen Adapters	100-049-249	•	•	opt					•	•		•
Calibrator Kit, Model 632.03/05 Knife Edges	100-049-250		•									
Calibrator Kit, Model 632.02 Knife Edges	100-049-251	•										
Calibrator Kit, Standard, 12.7 mm (0.50 in), Round Spec.	100-049-270			•	•							
Calibrator Kit, Frame Invert	100-049-252				•	•	•					•
Calibrator Kit, Extended, 12.7 mm (0.50 in), Round Spec.	100-049-253				•							
Calibrator Kit, Dimpled, 9.5 mm (0.375 in) Round Spec.	100-054-156											•
Calibrator Kit, 9.5 mm (0.375 in) Round Specimen	100-049-255			opt						opt		
Calibrator Kit, 3.2 mm (0.125 in) Round Specimen	100-049-256			opt				•				
Calibrator Kit, Mounting Bracket for Model 632.53/54	100-054-133											•
Calibrator Kit, Model 632.18/19/20 Diametral Extensometers	100-049-258						•					
Calibrator Kit, Displacement Gages	100-049-259				•	•						
Calibrator Kit, Model 632.06 Option 006 Adapter	100-054-796					•						
Calibrator Kit, Model 632.17/85 Axial	100-049-260								•			
Calibrator Kit, Model 632.85 Transverse	100-049-261									•		
Calibrator Kit, Model 632.92/21 Circumferential	100-049-262											

Find the model of the extensometer, then follow the column down to determine the kit(s) required to calibrate the extensometer.

MTS Direct Current Potential Drop (DCPD) Solution

Integrated Software & Hardware Package Allowing Accurate and Efficient Measurement of Fatigue Crack Growth

Direct Current Potential Drop (DCPD) measures the change in resistance of a specimen, which correlates directly with crack growth. It is an effective way to collect fatigue crack growth data in high-temperature tests and other challenging environments. DCPD is used to overcome the challenges that make contact extensometry and compliance calculations impractical for crack length measurement. Unfortunately, to get the best results, some DCPD solutions require excessive set-up time while others limit flexibility.

The MTS DCPD Solution features integrated software and hardware for test control, data analysis, and results reporting. It greatly simplifies DCPD test configuration using the powerful automation tools built into MTS TestSuite software. Integrated into MTS FlexTest 40 Controllers and MTS TestSuite Software, the MTS DCPD solution eliminates the challenges that can arise with non-integrated components from multiple suppliers. This provides a convenient and comprehensive way to incorporate DCPD into your material testing programs.

This solution uses the reversing DCPD method. By monitoring load and potential (voltage) data from a loaded specimen and a reference specimen, it determines crack length by calculating ratios of potential on both. Applying a calibration equation to the ratios produces the crack length. The current is switched on and off (reversed) to correct for thermocouple effects.

Overall

Current:	Adjustable 0 - 20 Amps DC Maximum
Voltage:	5 Volts DC Maximum
Gain:	x 5,000 or x 50,000
Gain Accuracy:	0.5% Typical, 2% Maximum
Filter:	50 Hz and 300 Hz
Bandwidth:	300 Hz Maximum
Output:	±10 Volts DC

Main Chassis

Gain:	1x or 10x
-------	-----------

Preamp

Gain:	5,000x
-------	--------



Comprehensive DCPD Solution

This solution offers a comprehensive approach for adding DCPD to your materials testing program. It includes the software and hardware required for integrated DCPD test development and execution. The MTS DCPD Testing Solution includes:

- » High-fidelity MTS DCPD Measurement System (2- or 7-channel configuration)
- » Modular Preamplifier(s)
- » MTS TestSuite DCPD Module (includes all test templates, report templates and calculations for ASTM-compliant fatigue crack growth assessments)
- » MTS TestSuite Fracture Analyzer software
- » Setup specimen kit for verifying hardware setup and test configuration
- » Setup, installation and equipment verification by MTS
- » Additional training, consulting and template modifications are available.

KEY FEATURES

- » Compatible with MTS Series 318 and MTS Landmark® test systems for electrical isolation
- » Integrates with MTS FlexTest Controller (Model 40) and MTS TestSuite Software
- » Designed for low noise
- » Modular Preamplifier maintains signal integrity
- » Compliant to CE, CSA and NRTL
- » Available with advanced consultation, template changes and training.

GRIPS & FIXTURES

For every material test – from basic tensile or fracture mechanics to ultra high temperature – MTS offers the best specimen test interface option. From fixtures designed to accommodate extensometers to compression platens for rock mechanics, our grips and fixtures will meet your testing needs.



Grips & Fixtures

Introduction	76-77
Polymer Matrix Composites Fixtures Index	78
Servohydraulic Grips & Fixtures	79-95
Composite Fixtures	96-102
Electromechanical Introduction	103
Advantage™ Grips	104-109
MTS Fundamental™ Grips & Fixtures	110-111
Bionix® Grips	112-117
Bionix EnviroBath	118-119
EM Extend® Kit	120-121

MTS Grips

Durable and Versatile

When it comes to grips, MTS offers the best in the business. Our Series 647 Hydraulic Wedge Grips, for example, set the standard for ease-of-use, accuracy, and repeatability. We offer grips designed for specialized applications such as high-temperature, low-cycle fatigue and fracture testing.

For electrodynamic systems, MTS has an expanding selection of grips designed to meet typical application requirements, supplementing the existing low force offering.

Three categories of electromechanical accessories are available to help you match your application needs, whether it be demanding R&D testing of advanced composites and alloys; standard testing of plastics, textiles, and other materials; or complex testing for the medical markets. For a broader selection of electromechanical grips, and to see the offering for static hydraulic applications, please refer to the "Accessories for MTS Criterion System" or "Accessories for MTS Exceed" catalogs.

Whether you need a 10 N manual grip or a 2500 kN hydraulic grip, we've got it all.

Available options for most MTS grips include:

- » Choice of gripping surfaces
- » Mounting hardware
- » Water cooling kits (hydraulic grips)
- » Alignment fixtures, software, and strain-gaged specimens to assure proper load train alignment.

A small sampling of specimens that MTS grips can accommodate. MTS grips are precise and repeatable, ensuring that your specimen is held in the same manner each time.



Fixtures

MTS fixtures often go beyond the industry standard design to make them easier to use, more durable, and to provide you with more precise test data. Some can be modified to accommodate special measurement devices such as extensometers. Others may be constructed out of special materials to withstand higher temperatures, corrosive environments, or simply to provide longer life.

Can't find what you need?

We offer many more grips and fixtures. Contact your local sales representative or applications engineer to find the model that meets your exact needs.

Polymer Matrix Composites Fixtures Application Index

A Comprehensive Array of Polymer Matrix Composites / Fibre Reinforced Plastics Accessories

MTS complements its electromechanical and servohydraulic testing lines with a comprehensive array of accessories to fulfill a full spectrum of polymer matrix composite material testing – from basic quality control, to demanding research and development applications.

CAN'T FIND WHAT YOU NEED?

We offer many more grips and fixtures. We can offer higher temperature version of many of the fixtures. Contact your local sales representative to find the model that meets your exact needs.

	Application	Test Standard	Fixture Option	See Page
Laminae & Laminate	Tensile	ISO 527-4 & 5, ASTM D3039, EN 2561, EN 2597	Model 647 Side-Loading Hydraulic Wedge Grip	79
			MTS Advantage Wedge Action Grips	108
	Compression	ISO 14126 Method 1A	Modified Celanese Compression Loading Fixture	96
		ISO 14126 Method 1B	IITRI Compression Loading Fixture	96
		ASTM D3410		
		ASTM D6641	Combined Loading Compression (CLC) Test Fixture	97
	Flexure	ISO 14125	Model 642.01 3- & 4-Point Bend Fixture with Roller Assembly Size 10 mm Diameter	94
			Model WA104A; WA204A & ZWA304 3-Point Bend Fixture with Loading Edge R5 Supporting R2 or R5	95
		ASTM D7264	Model 642.01 or 642.10 3- & 4-Point Bend Fixture with Roller Assembly Size 10 mm Diameter	94
			Model WA104A; WA204A & ZWA304 3-Point Bend Fixture with Loading Edge & Supporting R5	95
		EN 2562	Model 642.10 3- & 4-Point Bend Fixture with Roller Assembly Size 25 mm (loading) & 10 mm (support) Diameter	94
		EN 2746	Model 642.01 3- & 4-Point Bend Fixture with Roller Assembly Size 10 mm (loading) & 4 mm (support) Diameter	94
			Model WA104A; WA204A & ZWA304 3-Point Bend Fixture with Loading Edge R5 Supporting R2	95
	Shear	ISO 14129, ASTM D3518	Model 647 Side-Loading Hydraulic Wedge Grips	79
			MTS Advantage Wedge Action Grips	108
		ASTM D5379	V-Notched Beam (Iosipescu) Shear Fixture	98
		ASTM D7078	V-Notched Rail Shear Test Fixture	97
	Interlaminar Shear	ISO 14130	Short Beam Shear and Three-Point Flexure Fixture	98
			Model WA104A; WA204A & ZWA304 with Loading Edge R5 Supporting R2	95
		ASTM D2344	Short Beam Shear and Three-Point Flexure Fixture	98
		EN 2377	Model WA104A with Loading Edge R3 or R5 Supporting R2; WA204A & ZWA304 with Loading Edge R5 Supporting R2	95
	Fracture Mechanics	ASTM D6671	Mixed Mode Bending Fixture	99
	Fatigue (tension / tension)	ISO 13003, ASTM D3479	Model 647 Side-Loading Hydraulic Wedge Grip	79
Structural	Tension (open & filled hole)	ASTM D5766, ASTM D6742, ASTM D7615	Model 647 Side-Loading Hydraulic Wedge Grip	79
	Compression (open & filled hole)	ASTM D6484, ASTM 6742, BS 07260, ASTM D7615	Open / Filled Hole Compression Fixture	99
	Compression After Impact	ASTM D7137	Compression After Impact Test Fixture	100
Sandwich Structures	Tension	ASTM C273, ASTM C394	Flatwise Plane Shear Fixture, Tensile Mode	100
	Compression	ASTM C273, ASTM C394	Flatwise Plane Shear Fixture, Compression Mode	101
		ASTM C364	Adjustable Edgewise Compression Fixture	101
	Flexure / Shear	ASTM D5467, ASTM C393, ASTM D7249, ASTM D7250	Three- & Four-Point Sandwich Beam Flexure / Shear Fixture	102
Adhesives	Peel	ASTM D1781	Climbing Drum Peel Fixture	102

Multipurpose, Side Loading Hydraulic Wedge Grips

The MTS Series 647 Hydraulic Wedge Grips are versatile, easy-to-load grips for a wide range of tensile and fatigue applications. The symmetrical housing design ensures an even specimen loading across the entire face of the wedge. The lateral movement of the wedges won't change the gripping position on the specimen once the grips are activated.

Features

- » These grips clamp onto your specimen in the same position, test after test, to minimize the bending strains that can invalidate your test results
- » Tension and fatigue capability
- » Adjustable pressure allows these grips to be used for testing a variety of materials
- » A wide variety of wedges are available to meet your requirements
- » Side loading capability for easy specimen insertion.

Four Versions of the 647 Family

- » Axial – This model has all of the standard features described above
- » Aluminum – The 647.02 (Axial model only) grip housings are made of aluminum. This lightweight design provides excellent response in high frequency applications and makes them ideal for tabletop use. It weighs only 7 kg (15 lb)
- » Axial-Torsional – This capability is available in three capacity ranges
- » All Temperature – These models allow for temperatures to 540°C (1000°F). Contact MTS for additional information.

For use on non-hydraulic load frames or for high temperature applications below -7°C (20°F) or above +66°C (+150°F), see our hydraulic grip supplies on pages 90-91.

All grips are sold as pairs.

All wedges and attachment kits are sold separately.



TECH NOTE

Attachment Kits

The term "attachment kit" can refer to any of the hardware required to connect your grips to your system. This is an important item to remember since the attachment kits are generally sold separately. With a "matched" system (the load frame, load cell, and actuator all have the same capacity rating, like 100 kN), a basic attachment kit is required. This can consist of either spiral washers and a mounting stud, a cap screw, or pin and clevis depending on the specific grip and load frame. Inserts or loading washers may be required if there are variances from a matched system. If you have questions, please contact MTS for assistance with defining the correct attachment kit.

Multipurpose, Side Loading Hydraulic Wedge Grips

Axial Model 647 Grip Specifications

MODEL	DYNAMIC FORCE	STATIC FORCE	PRESSURE	TEMPERATURE* MIN/MAX	OVERALL HEIGHT	DIAMETER
647.02B	25 kN (5.5 kip)	31 kN (7 kip)	21 MPa (3000 psi)	-40°C/177°C (-40°F/350°F)	131 mm (5.2 in)	150 mm (6.0 in)
647.10A	100 kN (22 kip)	120 kN (27 kip)	21 MPa (3000 psi)	-40°C/177°C (-40°F/350°F)	188 mm (7.4 in)	203 mm (8.0 in)
647.25A	250 kN (55 kip)	333 kN (75 kip)	69 MPa (10,000 psi)	-40°C/177°C (-40°F/350°F)	249 mm (9.8 in)	266 mm (10.5 in)
647.50A	500 kN (110 kip)	550 kN (120 kip)	69 MPa (10,000 psi)	-18°C to 65°C (0°F to 150°F)	291 mm (11.5 in)	330 mm (13.0 in)
647.100A	1000 kN (220 kip)	1200 kN (264 kip)	69 MPa (10,000 psi)	-18°C to 65°C (0°F to 150°F)	414 mm (16.2 in)	444 mm (17.5 in)
647.250	2500 kN (550 kip)	2750 kN (610 kip)	69 MPa (10,000 psi)	-18°C to 65°C (0°F to 150°F)	819 mm (32.3 in)	737 mm (29.0 in)

*Temperatures above 77°C/150°F require a stand-alone grip supply and extension rods.

MODEL	WEIGHT	METRIC/US CUSTOMARY STUD SIZE	PART NUMBER
647.02B	7 kg (15 lb)	M12 x 1.25 (1/2"-20)	056-078-605
647.10A	30 kg (67 lb)	M27 x 2 (1"-14)	047-080-605
647.25A	77 kg (170 lb)	M36 x 2 (1 1/2"-12)	047-080-905
647.50A	148 kg (325 lb)	M52 x 2 (2"-12)	047-595-505
647.100A	386 kg (850 lb)	M76 x 2	053-137-201
647.250	1153 kg (3335 lb)	N/A	Contact MTS

Wedges and attachment kits sold separately.



Axial-Torsional Model 647 Grip Specifications

MODEL	DYNAMIC FORCE	STATIC FORCE	PRESSURE	TEMPERATURE* MIN/MAX	OVERALL HEIGHT	DIAMETER
647.02B-22	25 kN (5.5 kip)	220 N.m (2000 in.lb)	21 MPa (3000 psi)	-40°C/150°C (-40°F/300°F)	135 mm (5.3 in)	150 mm (6.0 in)
647.10A-05	100 kN (22 kip)	550 N.m (5000 in.lb)	21 MPa (3000 psi)	-40°C/121°C (-40°F/250°F)	188 mm (7.4 in)	203 mm (8.0 in)
647.10A-11	100 kN (22 kip)	1,100 N.m (10,000 in.lb)	21 MPa (3000 psi)	-40°C/121°C (-40°F/250°F)	224 mm (8.8 in)	203 mm (8.0 in)
647.25A-22	250 kN (55 kip)	2,200 N.m (20,000 in.lb)	69 MPa (10,000 psi)	-18°C/65°C (0°F/150°F)	343 mm (13.5 in)	266 mm (10.5 in)

*Temperatures above 77°C/150°F require a stand-alone grip supply and extension rods.

MODEL	WEIGHT	MOUNTING	PART NUMBER
647.02B-22	8 kg (16 lb)	40 mm dia.	100-026-042
647.10A-05	27 kg (60 lb)	M68 x 2 (LH)	049-157-301
647.10A-11	34 kg (75 lb)	M68 x 2 (LH)	049-817-001
647.25A-22	95 kg (210 lb)	M92 x 3 (LH)	056-124-001

Wedges for Series 647 Wedge Grips Offer a Variety of Surfaces

MTS Wedges Come in a Variety of Surfaces to Meet Your Requirements

- » Diamond tip – aggressive surface for gripping soft materials (steel, plastic)
- » Sawtooth – designed for increased holding capacity in tension tests
- » Surfalloy – grit incorporated onto the wedge surface for testing brittle samples
- » Serrated – universal surface for soft metals on vee wedges

MTS employs a unique wedge design that significantly reduces the stress concentration on the specimen, enabling even very brittle composites to be gripped securely without grip-induced failure.

MTS 647 Wedges can be ordered with a water-cooled option. These wedges are used with a Water Cooling Kit that is purchased separately.

Model 647 All-Temperature Wedges are available for the all-temperature grips. Contact MTS for additional information.

Water Cooling Kit

There are two versions of the Water Cooling Kit. Both versions include hoses, fittings, metering valves and a manifold that attaches magnetically to your load frame or other flat, steel surface. One version also includes water flow switches. Some MTS products, such as the Model 653 furnace, support integration with the flow switches, enabling protection of the equipment in the event of a water supply failure. Please check with an MTS Applications Engineer to determine if your equipment is compatible with the flow switch protection circuit.

Wedges are sold in matched sets of four.


Wedges for Model 647.02B Axial, Axial-Torsional Grips and 10 kN Advantage Pneumatic Grips

SPECIFICATIONS

Temperature Range: -40°C to 177°C (-40°F to 350°F)



The diamond impressions in the specimen provide good gripping without causing failure.

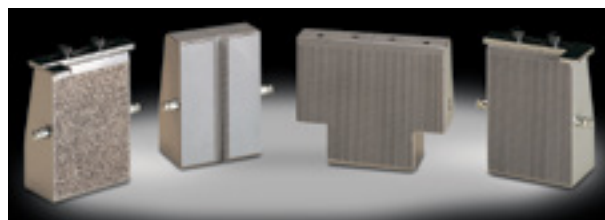

Model 647.02B Flat Specimen Wedges

SURFACE	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Sawtooth	0–7.2 mm (0–0.28 in)	25 mm (1 in)	050-507-906
Sawtooth	7.1–14.4 mm (0.28–0.57 in)	25 mm (1 in)	050-507-907
Sawtooth	13.4–20.5 mm (0.53–0.81 in)	25 mm (1 in)	050-507-901
Sawtooth	18.8–25.9 mm (0.74–1.02 in)	25 mm (1 in)	050-507-905
Surfalloy	0–7.2 mm (0–0.28 in)	25 mm (1 in)	050-507-917
Surfalloy	7.2–14.4 mm (0.28–0.57 in)	25 mm (1 in)	050-507-918
Surfalloy	13.4–20.5 mm (0.53–0.81 in)	25 mm (1 in)	050-507-904

Insertion depth 38 mm (1.5 in)

Model 647.02B Round Specimen Wedges

SURFACE	SPECIMEN DIAMETER	PART NUMBER
Surfalloy	10 mm	050-507-912
Surfalloy	15 mm	050-507-913
Surfalloy	25 mm	050-507-914
Surfalloy	0.5 in	050-507-915

Insertion depth 38 mm (1.5 in)

Model 647.02B Vee Specimen Wedges

SURFACE	FOR SPECIMEN DIAMETERS WHEN SIDE/TOP LOADING	PART NUMBER
Diamond tip	3.0–8.1/9.4 mm (0.12–0.32/0.37 in)	050-507-908
Diamond tip	8.9–10.9/15.2 mm (0.35–0.43/0.60 in)	050-507-909
Diamond tip	14.0–18.0/20.8 mm (0.55–0.71/0.82 in)	050-507-910
Diamond tip	20.1–22.9/27.2 mm (0.79–0.90/1.07 in)	050-507-911

Insertion depth 38 mm (1.5 in)

Wedges for Model 647.10 for Axial and Axial-Torsional Grips

Model 647.10 Flat Specimen Wedges

SURFACE	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Diamond tip	0–7.6 mm (0–0.30 in)	44 mm (1.75 in)	041-842-101
Diamond tip	7.1–14.2 mm (0.28–0.56 in)	44 mm (1.75 in)	041-842-102
Diamond tip	11.7–19.1 mm (0.46–0.75 in)	44 mm (1.75 in)	041-842-109
Wide diamond tip	0–7.6 mm (0–0.30 in)	76 mm (3 in)	046-198-604
Wide diamond tip	7.1–14.2 mm (0.28–0.56 in)	76 mm (3 in)	046-198-603
Surfallo	0–7.6 mm (0–0.30 in)	44 mm (1.75 in)	041-842-108
Surfallo	7.1–14.2 mm (0.28–0.56 in)	44 mm (1.75 in)	041-842-111
Surfallo	11.7–19 mm (0.46–0.75 in)	44 mm (1.75 in)	041-842-121
Wide surfallo	0–7.6 mm (0–0.30 in)	76 mm (3 in)	046-198-602
Wide surfallo	7.1–14.2 mm (0.28–0.56 in)	76 mm (3 in)	046-198-601

Insertion depth 63.5 mm (2.5 in)

WATER-COOLED FLATS	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Diamond tip	0–7.6 mm (0–0.30 in)	44 mm (1.75 in)	046-838-702
Diamond tip	7.1–14.2 mm (0.28–0.56 in)	44 mm (1.75 in)	046-838-701
Diamond tip	11.7–19.1 mm (0.46–0.75 in)	44 mm (1.75 in)	046-838-713
Surfallo	0–7.6 mm (0–0.30 in)	44 mm (1.75 in)	046-838-705
Surfallo	7.1–14.2 mm (0.28–0.56 in)	44 mm (1.75 in)	046-838-706

Insertion depth 63.5 mm (2.5 in)

Model 647.10 Round Specimen Wedges

SURFACE	SPECIMEN DIAMETER	PART NUMBER
Surfallo	12 mm	041-842-134
Surfallo	15 mm	041-842-135
Surfallo	20 mm	041-842-136
Surfallo	0.5 in	041-842-132
Surfallo	0.75 in	041-842-133
Surfallo	1.0 in	041-842-149

Insertion depth 63.5 mm (2.5 in)

WATER-COOLED ROUNDS	SPECIMEN DIAMETER	PART NUMBER
Surfallo	12 mm	046-838-716
Surfallo	15 mm	046-838-717
Surfallo	20 mm	046-838-718
Surfallo	0.5 in	046-838-714

Insertion depth 63.5 mm (2.5 in)

Model 647.10 Vee Specimen Wedges

SURFACE	FOR SPECIMEN DIAMETER WHEN SIDE/TOP LOADING	PART NUMBER
90° Serrated	5.8–10.2/11.9 mm (0.23–0.40/0.47 in)	041-842-103
90° Serrated	10.9–12.7/16.5 mm (0.43–0.50/0.65 in)	041-842-104
90° Serrated	17–17/22.9 mm (0.67–0.67/0.9 in)	041-842-110

Insertion depth 63.5 mm (2.5 in)

WATER-COOLED VEES	FOR SPECIMEN DIAMETER WHEN SIDE/TOP LOADING	PART NUMBER
90° Serrated	5.8–10.2/11.9 mm (0.23–0.40/0.47 in)	046-838-703
90° Serrated	10.9–12.5/16.5 mm (0.43–0.50/0.65 in)	046-838-704

Insertion depth 63.5 mm (2.5 in)

WATER COOLING KIT	PART NUMBER
Water Cooling Kit, with flow switches	057-697-502
Water Cooling Kit, without flow switches	057-697-505

Required for Series 647 Water Cooled Wedges

Includes 3.6 m (12 ft) of hose for connection to water source.

Attaches magnetically to load frame base or other flat steel surface.

Wedges for Model 647.25 for Axial and Axial-Torsional Grips

Model 647.25 Flat Specimen Wedges

SURFACE	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Diamond tip	1.0–11.9 mm (0.04–0.47 in)	50 mm (2 in)	041-842-201
Diamond tip	6.1–17.0 mm (0.24–0.67 in)	50 mm (2 in)	041-842-202
Diamond tip	15.0–25.9 mm (0.59–1.02 in)	50 mm (2 in)	041-842-203
Wide diamond tip	1.0–11.9 mm (0.04–0.47 in)	100 mm (4 in)	046-198-804
Wide diamond tip	6.1–17.0 mm (0.24–0.67 in)	100 mm (4 in)	046-198-806
Wide diamond tip	15.0–25.9 mm (0.59–1.02 in)	100 mm (4 in)	046-198-805
Surfallo	1.0–11.9 mm (0.04–0.47 in)	50 mm (2 in)	041-842-207
Surfallo	6.1–17.0 mm (0.24–0.67 in)	50 mm (2 in)	041-842-208
Surfallo	15.0–25.9 mm (0.59–1.02 in)	50 mm (2 in)	041-842-209
Wide surfallo	1.0–11.9 mm (0.04–0.47 in)	100 mm (4 in)	046-198-817
Wide surfallo	6.1–17.0 mm (0.24–0.67 in)	100 mm (4 in)	046-198-802
Wide surfallo	15.0–25.9 mm (0.59–1.02 in)	100 mm (4 in)	046-198-803

Insertion depth 89 mm (3.5 in)

WATER-COOLED FLATS	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Diamond tip	1.0–11.9 mm (0.04–0.47 in)	50 mm (2 in)	045-966-201
Diamond tip	6.1–17.0 mm (0.24–0.67 in)	50 mm (2 in)	045-966-202
Diamond tip	14.9–25.9 mm (0.59–1.02 in)	50 mm (2 in)	045-966-204
Surfallo	1.0–11.9 mm (0.04–0.47 in)	50 mm (2 in)	045-966-205
90° Serrations	1.0–11.9 mm (0.04–0.47 in)	50 mm (2 in)	045-966-208
90° Serrations	6.1–17.0 mm (0.24–0.67 in)	50 mm (2 in)	045-966-209
90° Serrations	15.0–25.9 mm (0.59–1.02 in)	50 mm (2 in)	045-966-210

Insertion depth 89 mm (3.5 in)

Model 647.25 Round Specimen Wedges

SURFACE	SPECIMEN DIAMETER	PART NUMBER
Surfallo	15 mm	041-842-231
Surfallo	20 mm	041-842-232
Surfallo	30 mm	041-842-233
Surfallo	0.5 in	041-842-234
Surfallo	1 in	041-842-235

Insertion depth 89 mm (3.5 in)

Model 647.25 Vee Specimen Wedges

SURFACE	FOR SPECIMEN DIAMETER WHEN SIDE/TOP LOADING	PART NUMBER
Serrated	6.4–10.2/13.5 mm (0.25–0.40/0.53 in)	041-842-206
Serrated	10.7–16.8/19.9 mm (0.42–0.66/0.78 in)	041-842-204
Serrated	16.8–20.0/26.2 mm (0.66–0.79/1.03 in)	041-842-205

Insertion depth 89 mm (3.5 in)

WATER-COOLED VEES	FOR SPECIMEN DIAMETER WHEN SIDE/TOP LOADING	PART NUMBER
Serrated	16.8–20.0/26.2 mm (0.66–0.79/1.03 in)	045-966-203
Serrated	6.4–10.2/13.5 mm (0.25–0.40/0.53 in)	045-966-206
Serrated	10.7–16.8/19.9 mm (0.42–0.66/0.78 in)	045-966-207

Insertion depth 89 mm (3.5 in)

WATER COOLING KIT	PART NUMBER
Water Cooling Kit, with flow switches	057-697-502
Water Cooling Kit, without flow switches	057-697-505

Required for Series 647 Water Cooled Wedges
Includes 3.6 m (12 ft) of hose for connection to water source.
Attaches magnetically to load frame base or other flat steel surface.

Wedges for Model 647.50, 647.100, and 647.250 Axial Grips

Model 647.50 Flat Specimen Wedges

SURFACE	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Diamond tip	0–10.9 mm (0.00–0.43 in)	100 mm (4 in)	047-641-606
Diamond tip	10.2–21.1 mm (0.40–0.83 in)	100 mm (4 in)	047-641-607
Diamond tip	20.3–31.2 mm (0.80–1.23 in)	100 mm (4 in)	047-641-608
Diamond tip	30.5–41.4 mm (1.20–1.63 in)	100 mm (4 in)	047-641-609
Diamond tip	40.6–51.6 mm (1.60–2.03 in)	100 mm (4 in)	047-641-610
Surfallo	0–10.9 mm (0.00–0.43 in)	100 mm (4 in)	047-641-611
Surfallo	10.2–21.1 mm (0.40–0.83 in)	100 mm (4 in)	047-641-612
Surfallo	20.3–31.2 mm (0.80–1.23 in)	100 mm (4 in)	047-641-613
Wide diamond tip	0.0–10.9 mm (0.00–0.43 in)	200 mm (8 in)	048-966-302
Wide diamond tip	10.1–21.0 mm (0.40–0.83 in)	200 mm (8 in)	048-966-303
Wide diamond tip	20.3–31.2 mm (0.80–1.23 in)	200 mm (8 in)	048-966-304
Wide diamond tip	30.4–41.4 mm (1.20–1.63 in)	200 mm (8 in)	048-966-305
Wide diamond tip	40.6–51.5 mm (1.6–2.03 in)	200 mm (8 in)	048-966-306
Wide diamond tip	0.0–10.9 mm (0.00–0.43 in)	200 mm (8 in)	048-966-307
Wide diamond tip	10.1–21 mm (0.40–0.83 in)	200 mm (8 in)	048-966-301

Insertion depth 89 mm (3.5 in)

Model 647.50 Vee Specimen Wedges

SURFACE	FOR SPECIMEN DIAMETERS WHEN SIDE/TOP LOADING	PART NUMBER
90° Serrated	6.4–12.7/15.5 mm (0.25–0.50/0.61 in)	047-641-601
90° Serrated	15.2–18.3/24.4 mm (0.60–0.72/0.96 in)	047-641-602
90° Serrated	24.1–25.1/33.5 mm (0.95–0.99/1.32 in)	047-641-603
90° Serrated	33.0–33.0/42.4 mm (1.30–1.30/1.67 in)	047-641-604
90° Serrated	41.9–41.9/51.3 mm (1.65–1.65/2.02 in)	047-641-605

Insertion depth 89 mm (3.5 in)

Model 647.100 Flat Specimen Wedges

SURFACE	SPECIMEN THICKNESS	USABLE WIDTH	PART NUMBER
Diamond tip	0–23 mm (0–0.89 in)	100 mm (4 in)	053-137-403
Diamond tip	23–45 mm (0.89–1.77 in)	100 mm (4 in)	053-137-402
Diamond tip	45–67 mm (1.77–2.64 in)	100 mm (4 in)	053-137-401

Insertion depth 89 mm (3.5 in)

Model 647.100 Vee Specimen Wedges

SURFACE	FOR SPECIMEN DIAMETERS WHEN SIDE/TOP LOADING	PART NUMBER
90° Serrated	6–10/10 mm (0.24–0.39/0.39 in)	053-137-407
90° Serrated	9.9–20.8/25.1 mm (0.39–0.82/0.99 in)	053-137-406
90° Serrated	25.1–29.5/44.4 mm (0.99–1.16/1.75 in)	053-137-405
90° Serrated	44.4–50.8/63.0 mm (1.75–2.0/2.5 in)	053-137-404

Insertion depth 89 mm (3.5 in)

Model 647.250 Wedge and Liner Sets

WEDGE TYPE	SPECIMEN SIZE RANGE	PART NUMBER
Flat*	0–29.5 mm (0–1.16 in) 26.8–54 mm (1.06–2.16 in) 51.8–79 mm (2.04–3.11 in) 76.8–104 mm (3.02–4.09 in)	057-367-101
Vee*	30–44.0/57.2 mm (1.18–1.73/2.25 in) 55.8–68.9/80.5 mm (2.20–2.71/3.17 in) 78.6–93.7/104 mm (3.10–3.69/4.09 in)	057-367-102
Flat & Vee	Includes the flat and vee wedges described above and single set of liners	057-367-103

* Flat wedge size ranges refer to specimen thickness

** Vee wedge diameter ranges refer to top/side loading

Hydraulic Collet Grips Are Ideal for Fatigue Testing

The MTS Series 646 Collet Grips offer the best value of any fatigue testing grips available. The design of the Series 646 Grips allows for excellent repeatability from one specimen to the next, minimizing bending strains. Loading is simply a matter of sliding the specimen into the smooth collet. There are no expensive threads or button ends to machine and no additional bolts or collars to install. The Series 646 Grips were designed primarily for round specimens. A special high temperature collet package is available for the 646.10 that utilizes the extension rods supplied with our Model 680 High Temperature Grips. This package will provide you with the same high-temperature capabilities as the 680 grips, and requires the same specimen adapters (see page 89).

See pages 90-91 for our grip supplies.
Pressure 45 MPa (6500 psi)



Axial Model 646 Specifications

MODEL	FORCE	TEMPERATURE	HEIGHT	DIAMETER	MOUNTING THREAD SIZE	PART NUMBER
646.10B	100 kN (22 kip)	-40°C to 65°C (-40°F to 150°F)	114 mm (4.50 in)	171 mm (6.75 in)	M27x2 (1"-14)	042-481-807
646.25B	250 kN (55 kip)	-40°C to 65°C (-40°F to 150°F)	171 mm (6.75 in)	254 mm (10.0 in)	M36x2 (1 1/2"-12)	042-481-705

Grips are sold in pairs.

Axial-Torsional Model 646 Specifications

MODEL	AXIAL FORCE	TORSIONAL FORCE	TEMPERATURE	HEIGHT	DIAMETER	MOUNTING THREAD SIZE	PART NUMBER
646.10B	100 kN (22 kip)	1100 N.m (10,000 in.lb)	-40°C to 65°C (-40°F to 150°F)	147 mm (5.8 in)	171 mm (6.75 in)	M68x2	043-816-502
646.25B	250 kN (55 kip)	2200 N.m (20,000 in.lb)	-40°C to 65°C (-40°F to 150°F)	231 mm (9.1 in)	254 mm (10.0 in)	M92x3	043-816-301

Grips are sold in pairs.

Metric Collet Sizes

GRIP MODEL	SPECIMEN TYPE	PART NUMBER
646.10-XX	Round 10 mm diameter	038-058-917
646.10-XX	Round 12 mm diameter	038-058-903
646.10-XX	Round 15 mm diameter	038-058-902
646.10-XX	Round 30 mm diameter	038-058-901
646.25-XX	Round 20 mm diameter	038-059-203
646.25-XX	Round 25 mm diameter	038-059-207
646.25-XX	Round 40 mm diameter	038-059-223

US Customary Collet Sizes

GRIP MODEL	SPECIMEN TYPE	PART NUMBER
646.10-XX	Round 0.25 in diameter	038-058-905
646.10-XX	Round 0.50 in diameter	038-058-907
646.10-XX	Round 0.75 in diameter	038-058-909
646.10-XX	Round 1.0 in diameter	038-058-908
646.25-XX	Round 1.0 in diameter	038-059-202

WATER COOLING KIT	PART NUMBER
Water Cooling Kit, with flow switches	057-697-502
Water Cooling Kit, without flow switches	057-697-505

Required for Series 647 Water Cooled Wedges
Includes 3.6 m (12 ft) of hose for connection to water source.
Attaches magnetically to load frame base or other flat steel surface.

HIGH TEMPERATURE PACKAGE COMPONENTS	PART NUMBER
Furnace Extension Kit for 646.10	045-337-601
Hand Pump	100-254-831
Specimen Adapters	See page 89

Fatigue-rated Side Loading Pneumatic Wedge Grips

The MTS Series 645 Fatigue-rated Pneumatic Wedge Grips provide the versatility and reliability you won't find in other pneumatic grips. The symmetrical housing design ensures even specimen loading across the entire face of the wedge. Based on proven MTS designs, the lateral movement of the wedges won't change the gripping position on the specimen when the grips are energized.

FEATURES AND BENEFITS

- » Ideally suited to MTS Acumen electrodynamic load frames, as well as other materials test systems
- » These grips clamp onto your specimen in the same position, test after test, to minimize the bending strains that can invalidate your test results
- » Tension and fatigue capability
- » Adjustable pressure allows these grips to be used for testing a variety of materials
- » A wide variety of wedges are available to meet your requirements
- » Side loading capability for easy specimen insertion.

All grips are sold as pairs.

All wedges and attachment kits are sold separately.



Axial Model 645 Specifications

MODEL	FORCE	TEMPERATURE	PRESSURE	OVERALL HEIGHT	DIAMETER	PART NUMBER
645.002	2 kN (0.44 kip)	-40°C to 200°C (-40°F to 400°F)	0.6 MPa (80 psi)	131 mm (5.15 in)	104 mm (4.1 in)	100-242-422
645.005	5 kN (1.1 kip)	-40°C to 200°C (-40°F to 400°F)	0.6 MPa (80 psi)	133 mm (5.25 in)	147 mm (5.8 in)	100-242-417
645.012	12 kN (2.7 kip)	-40°C to 200°C (-40°F to 400°F)	1.1 MPa (160 psi)	131 mm (5.16 in)	168 mm (6.6 in)	100-415-143

Grips are sold in pairs.

Model 645 Wedges

DESCRIPTION	FOR SPECIMEN DIAMETERS WHEN SIDE/TOP LOADING	USABLE WIDTH
Flat Specimen Sawtooth Wedge Set	0-7.1 mm (0-0.28 in)	050-507-938
Flat Specimen Sawtooth Wedge Set	3.3-10.7 mm (0.13-0.42 in)	050-507-939
Flat Specimen Sawtooth Wedge Set	5.3-12.4 mm (0.21-0.49 in)	050-507-940
Flat Specimen Surfalloy Wedge Set	0-7.1 mm (0-0.28 in)	050-507-941
Flat Specimen Surfalloy Wedge Set	3.3-10.7 mm (0.13-0.42 in)	050-507-942
Flat Specimen Surfalloy Wedge Set	5.3-12.4 mm (0.21-0.49 in)	050-507-943
Vee Wedge Set Diamond tip	3-7.9 mm (0.12-0.31 in)	050-507-944
Vee Wedge Set Diamond tip	7.1-7.9/12.7 mm (0.28-0.31/0.5 in)	050-507-945
Vee Wedge Set Diamond tip	10.9-13.2/16.5 mm (0.43-0.52/0.65 in)	050-507-946
Vee Wedge Set Diamond tip	14.5-15.7/20.3 mm (0.57-0.62/0.8 in)	050-507-947

Insertion depth 89 mm (3.5 in)

Accessories

DESCRIPTION	PART NUMBER
Fundamental Pneumatic Grip Supply	100-231-865
Mounting & Cables, Fundamental Grip Supply to Acumen	100-263-782

Alignment Fixture, Software, and Strain-Gaged Specimens

Save Time and Get Better Data

Alignment Software for up to 12-Gage Specimens

MTS Alignment hardware and software acquires, analyzes, and displays bending strains for alignment and bending strain verification purposes. Its graphical interface quickly allows you to align your system or to verify how much bending strain you currently have. One of the key attributes of this software is that it can separate the bending strain of the specimen from that of the load train.

The best method of determining bending strain is to use a strain-gaged specimen. This product interfaces directly with a strain-gaged specimen and displays bending strain graphically on the screen of your PC. Continual scanning of bending strains allows you to align your system or verify the bending strain while cycling load. When finished, an Excel report can be generated to certify your system.

Features

- » Supports 12-gage specimens for best results. Also accommodates 4, 6, 8, and 9 gage specimens
- » Data acquisition and conditioning is external to the PC. This allows for much more PC flexibility, including laptops
- » This alignment system ensures compliance with industry standards: ASTM E1012, GES400 (NADCAP), GE450 and ISOTC 164SC5WG11
- » The "Alignment Wizard" verbally instructs the user how to perform concentric and angular adjustments
- » Accepts a calibrated load signal for a bending percent and bending strain versus load validation graph
- » Separates the specimen eccentricities from the system alignment.

Alignment Fixture

The MTS Model 609 Alignment Fixture for material testing systems saves time compared to manual alignment methods.


The main feature of the Model 609 Alignment Fixture is its capability to perform alignment adjustments while the load train is fully preloaded. This eliminates any inaccuracies involved in trying to account for the small changes in alignment that frequently occur during the preloading process.



In addition, because the fixture remains preloaded at all times, previous alignment adjustments are not lost when small changes in alignment are required.

The Model 609 Alignment Fixtures are compatible with MTS load frames, excluding MTS Exceed. The alignment fixtures are also readily adaptable to other load frames.


Testing Results Hitting the Target?

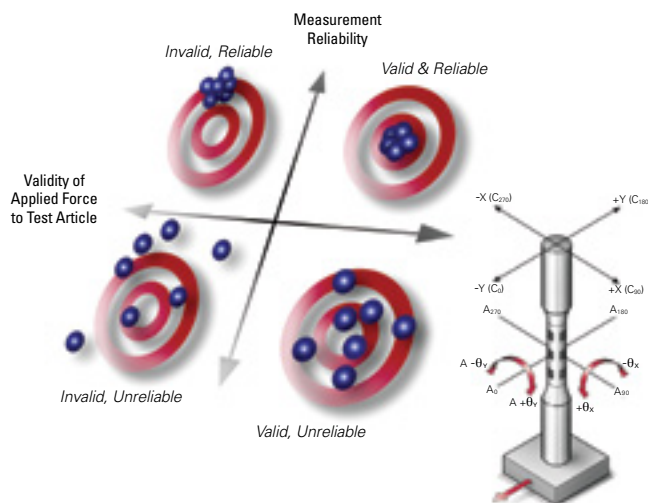

Specimens for Alignment Use

MTS offers round and flat standardized strain-gaged specimens for alignment use. Each specimen has 12 gages so you can correct both concentric and angular misalignments.

The gages are bonded with oven-cured epoxy for long life and are covered to prevent damage from handling. The specimen material is 4340 steel and is heat treated to safely allow stresses to 413 MPa (60,000 psi).

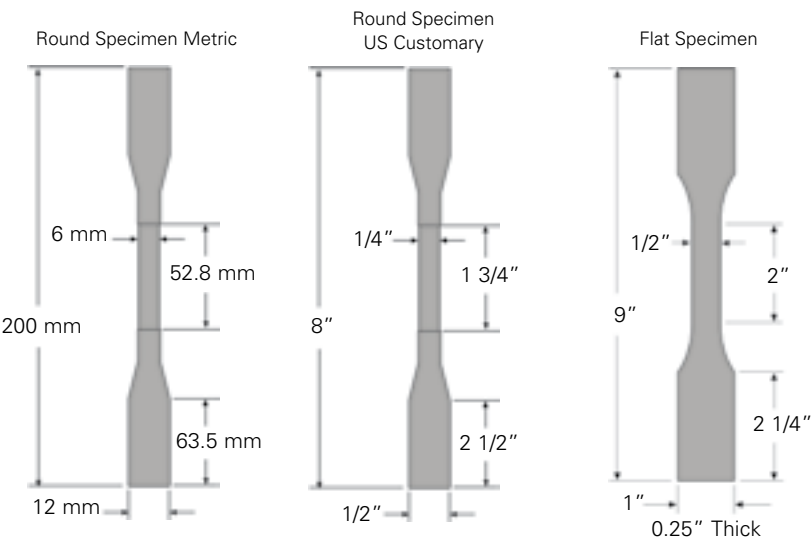


Alignment Fixture for Landmark, Model 318, and Model 312

MODEL	LOAD FRAME FORCE CAPACITY	HEIGHT	DIAMETER	STUD SIZE/LENGTH	METRIC PART NUMBER	US CUSTOMARY PART NUMBER
609.02A-01	25 kN (5.5 kip)	72 mm (2.81 in)	121 mm (4.75 in)	M12 x 1.25/355 mm (1/2"-20/14.0 in)	049-083-402	049-083-401
609.10A-01	100 kN (22 kip)	72 mm (2.81 in)	121 mm (4.75 in)	M27 x 2/343 mm (1"-14/13.5 in)	049-083-502	049-083-501
609.25A-01	250 kN (55 kip)	80 mm (3.12 in)	162 mm (6.38 in)	M36 x 2/462 mm (1 1/2"-12/18.2 in)	049-083-602	049-083-601
609.50A-01	500 kN (110 kip)	181 mm (7.11 in)	254 mm (10 in)	M52 x 2/686 mm (2"-12/27 in)	051-499-202	051-499-201

Insertion depth 89 mm (3.5 in)
Contact MTS for Acumen, Exceed, or Criterion Alignment fixture options.

REQUIRED ADAPTER KIT ORDERING INFORMATION
FOR THE MODEL 609 ALIGNMENT FIXTURE (FOR LOAD
FRAMES OTHER THAN MTS' MODELS 370, 318, AND 312)

- A Load cell thread size
- B Load cell thread depth
- C Crosshead height
- D Crosshead through-hole diameter
- E Dimensional information of any counter-bores in either the top or bottom of the crosshead

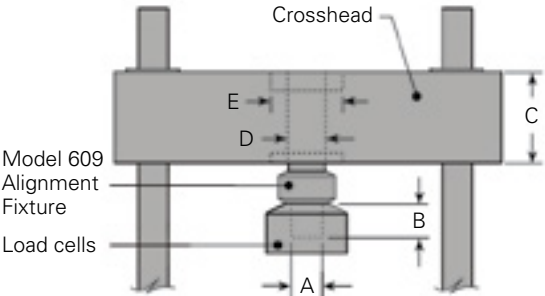

SPECIMENS	PART NUMBER
Round (Metric)	056-651-703
Round (US Customary)	056-651-701
Flat	056-651-702

709.20E Alignment Data Acquisition and Conditioning Units

709.20E Alignment Data Acquisition and Conditioning Units

VOLTAGE/PLUG	PART NUMBER	CERTIFIED PART NUMBER
115 V AC, CE	057-677-601	057-677-606
230 V AC, US Plug, CE	057-677-602	057-677-607
230 V AC, Euro Plug, CE	057-677-603	057-677-608
230 V AC, US Plug, CE	057-677-605	057-677-610

For other plug options, contact MTS

ALIGNMENT SOFTWARE*	PART NUMBER
709.20E-03 Alignment Software Perpetual License	056-866-201
709.20E-04 Alignment Software Annual License	056-866-202
709.20E-05 Alignment Software 5 year Term Annual Renewable License	056-866-203

*Requires XP/Vista/Windows 7.

High-Temperature LCF Grips

MTS Model 680 High-Temperature Grips were expressly designed for high-temperature, low-cycle fatigue testing, but can also be used for other types of tension-compression or tension-only testing. The grip heads and extension rods are constructed of high-temperature, super alloys selected to extend into the heated zone of the furnace, minimizing the thermal gradients within the specimen. Accessories for these grips include a hand pump (required), your choice of buttonhead or threaded specimen adapters, and a water cooling kit.

Specifications

MODEL NUMBER	FORCE RATING*		WEIGHT (each)	DIMENSIONS (see photo)			THREAD SIZE	PART NUMBER
	@700°C (1300°F)	@1000° (1832°F)		A	B	C		
680.01B	15 kip	15 kip	20 lb	1.5 in	4.6 in	4.9 in	14-Jan	044-057-103
680.01B	68 kN	8.9 kN	9 kg	38 mm	117 mm	124 mm	M27x2	044-057-103

**Note: Force rating decreases linearly with temperature between the ratings shown*

Accessories

TYPE	PART NUMBER
Hand Pump (Required)	100-254-832
Water Cooling Kit, with flow switches	057-697-502
Replacement Seal Kit	040-985-901

SPECIMEN ADAPTERS	DIMENSIONS	PART NUMBER
Buttonhead	0.25 in dia. shank, 0.75-0.85 in dia. button	041-901-902
Buttonhead	0.40 in dia. shank, 0.75-0.85 in dia. button	041-901-904
Buttonhead	0.50 in dia. shank, 0.75-0.85 in dia. button	041-901-901
Buttonhead	12 mm dia. shank, 21.1-21.6 mm dia. button	041-901-903
Threaded	1/2 in-13 UNC-2B	041-901-801
Threaded	M12 x1.75 mm	041-901-802



Grip Supplies and Intensifiers

MTS Model 685 self-contained, hydraulic Grip Supplies have been engineered for both performance and ease of use. There are two grip supplies and one grip intensifier to choose from.

STANDARD FEATURES INCLUDE:

- » Directional control valve for each grip
- » Center valve detent, allowing unparalleled control of over gripping
- » Continuous positive pressure design, providing high pressure stability over the entire operating range
- » Separate flow control valve for control of grip engagement speed
- » Independent grip circuits eliminate crosstalk.
- » Easy to maintain and service
- » Accommodate a wide range of electrical connections

MODEL 685.22 AND MODEL 685.10 STAND ALONE HYDRAULIC GRIP SUPPLIES

The 685.22 and 685.10 units feature a self-contained hydraulic pump, a 0.75 kW (1 hp) electric motor, a 11.3 l (3 gal) reservoir, a 10-micron absolute return line filter, and hoses for connection to grips. These units are furnished with individual directional control valves for upper and lower grips. The grip supplies use a special hydraulic fluid which allows the grips to be used in environmental chambers at elevated temperatures. They are designed to run continuously, which results in good pressure stability and easy adjustment of the output pressure. Grip closure rate is also adjustable. Since the grip supplies are self-contained systems, they allow the use of hydraulic grips on non-hydraulic test systems.



MODEL 685.60 HYDRAULIC GRIP INTENSIFIER

The 685.60 Hydraulic Grip Intensifier, which utilizes an innovative fluid-to-fluid intensification system, is used with grips that require a pressure higher than the normal system hydraulic pressure. Two versions are available with factory adjusted output pressures of 45 MPa (6,500 psi) and 69 MPa (10,000 psi). The output pressure is adjustable from 10 MPa (1,500 psi) up to the grip supply output rating.

To minimize the set-up time of larger grips, a high volume upgrade kit is available.



MODEL**	DESCRIPTION	PRESSURE RANGE	SETTING REPEATABILITY	OPERATING TEMPERATURE	WEIGHT	ELECTRICAL POWER REQUIREMENTS	PART NUMBER
685.10E-05	Stand-alone grip supply	10-70 MPa (1,500-10,000 psi)	±0.7 MPa (±100 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	115 V (60 Hz)	057-509-601
685.10E-06	Stand-alone grip supply	10-70 MPa (1,500-10,000 psi)	±0.7 MPa (±100 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	100-115 V (50 Hz)	057-509-602
685.10E-07	Stand-alone grip supply	10-70 MPa (1,500-10,000 psi)	±0.7 MPa (±100 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	208-230 V (60 Hz)	057-509-603
685.10E-08	Stand-alone grip supply	10-70 MPa (1,500-10,000 psi)	±0.7 MPa (±100 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	200-240 V (50 Hz)	057-509-604
685.22D-05	Stand-alone grip supply	0.7-21 MPa (100-3,000 psi)	±0.07 MPa (±10 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	115 V (60 Hz)	057-598-001
685.22D-06	Stand-alone grip supply	0.7-21 MPa (100-3,000 psi)	±0.07 MPa (±10 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	100-115 V (50 Hz)	057-598-002
685.22D-07	Stand-alone grip supply	0.7-21 MPa (100-3,000 psi)	±0.07 MPa (±10 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	208-230 V (60 Hz)	057-598-003
685.22D-08	Stand-alone grip supply	0.7-21 MPa (100-3,000 psi)	±0.07 MPa (±10 psi)	-40 to 177°C (-40 to 350°F)	76 kg (170 lb)	200-240 V (50 Hz)	057-598-004
685.60C-02*	Hydraulic system intensifier	10-45 MPa (1,500-6,500 psi)	±0.7 MPa (±100 psi)	-18 to 66°C (0 to 150°F)	39 kg (87 lb)	—	057-597-901
685.60C-03*	Hydraulic system intensifier	10-70 MPa (1,500-10,000 psi)	±0.7 MPa (±100 psi)	-18 to 66°C (0 to 150°F)	39 kg (87 lb)	—	057-597-902

*Optional kit available to minimize set-up time of larger grips.

**Pressure Stability: ±100 psi (1% of operating range). Continuous pump operation.

FILTRATION

- » 10 micron

HYDRAULIC FLUID

- » Mobile 525 SHC (685.10 & 685.22)

MAXIMUM AMBIENT TEMPERATURE

- » 32°C (90°F)

DIMENSIONS (SEE PHOTO)

- » Depth: 43.18 cm (17 in)
- » Width: 44.45 cm (17.5 in)
- » Height: 91.44 cm (36 in)



Hydraulic Grip Control Kit

Load Frame Mounted Grip Control Unit

The MTS Model 685.53 Grip Control Unit accurately controls pressure to hydraulic grips. Each grip can be actuated independently. The kit includes the grip control unit, hydraulic hoses to the hydraulic power supply, hydraulic hoses to the grips, and a universal mounting bracket. The control unit requires 21 MPa (3,000 psi) input pressure and can control the pressure to the grips from 0.7 MPa to 20 MPa (100 psi to 3,000) psi with a resolution of ± 0.07 MPa (± 10 psi). This range and resolution allow the user to select the proper grip pressure so it holds the specimen sufficiently without damaging it. The stability of the controlled pressure is ± 0.5 MPa (± 25 psi) yielding no extraneous load fluctuations within your test. The clamping speed of the grips is also adjustable.



MODEL	LOAD FRAME (COLOR)	PRESSURE CAPACITY	USED FOR GRIP MODEL NUMBERS	PART NUMBER
685.53A-01 *	Load Frame Mounted Grip Control (brown)	21 MPa /3,000 psi	647.02/.10	044-667-001
685.53C-02	Load Frame Mtd. Grip Control 318.10/.25 (gray)	21 MPa /3,000 psi	647.02/.10	050-195-402

* JIC 6 drain hose is required if 318 load frame does not have lift and locks.

Fracture Mechanics Clevis Grips

The MTS Model 640.20B Fracture Mechanics Grips are manufactured per ASTM E-399 specifications. They are the most popular grips used for fracture toughness and fatigue crack growth rate tests. They are constructed from high-strength, aircraft-quality 4340 steel and are machined to the exacting tolerances required by ASTM. The package includes two grips, two .4W pins, and four spring retainers.



MODEL	SPECIMEN WIDTH	STATIC FORCE	PIN DIAMETER	MOUNTING THREAD
640.20B-X1	12.7 mm (0.50 in)	60 kN (13.6 kip)	12.2 mm (0.48 in)	M27 x 2 (1"-14)
640.20C-X3	25.4 mm (1.00 in)	60 kN (13.6 kip)	12.2 mm (0.48 in)	M27 x 2 (1"-14)

Temperature Rating

-129°C to 177°C

MODEL	METRIC PART NUMBER
640.20B-21	007-172-745
640.20C-03	007-172-743

Temperature Rating

-200°F to 350°F

MODEL	METRIC PART NUMBER
640.20B-01	007-172-741
640.20C-03	007-172-743

Compression Platens

Model 643 Compression Platens

MTS Model 643 Compression Platens are made from case-hardened alloy steel with hard chrome plating. The platens have a smooth face with etched concentric rings, enabling the specimen to be centered visually for better test results. The upper platen is available with a spherical seat for improved alignment and ensuring even pressure across the entire surface of the specimen. Model 643 compression platens are offered in three configurations: Fixed & Spherical, Two Fixed, and One Fixed.

Units range in diameter from 60 mm (2.4 in) to 300 mm (12 in) and are designed for a unit stress of 689 MPa (100,000 psi) static (275 MPa (40,000 psi) dynamic) centered on the bearing surfaces. The usable temperature range is -129°C to 177°C (-200°F to +350°F).

Specifications

MODEL	MAXIMUM SPECIMEN DIAMETER	MOUNTING THREAD INSERT SIZES
643.06B	60 mm (2.4 in)	M12 x 1.25 /1/2"-20
643.10B	100 mm (4 in)	M27 x 2 /1"-14
643.15B	150 mm (6 in)	M27 x 2 /1"-14
643.20B	200 mm (8 in)	M27 x 2 /1"-14
643.30B	300 mm (12 in)	M52 x 2 /2"-12

MODEL	DESCRIPTION	PART NUMBER
643.06B-01	Fixed & Spherical	050-292-603
643.06B-03	Two Fixed	050-292-503
643.06B-05	One Fixed	050-612-103
643.10B-01	Fixed & Spherical	050-292-803
643.10B-03	Two Fixed	050-292-703
643.10B-05	One Fixed	050-612-203
643.15B-01	Fixed & Spherical	050-293-003
643.15B-03	Two Fixed	050-292-903
643.15B-05	One Fixed	050-612-303
643.20B-01	Fixed & Spherical	053-042-503
643.20B-03	Two Fixed	053-042-803
643.20B-05	One Fixed	053-042-703
643.30B-01	Fixed & Spherical	050-293-203
643.30B-03	Two Fixed	050-293-103
643.30B-05	One Fixed	050-612-403



Bend Fixtures


MTS Model 642 Bend Fixtures are configured to meet a variety of testing requirements. The fixtures have adjustable spans with easy-to-use, permanently attached scales for equal positioning of the rollers. The hardened

rollers ensure test result accuracy by reducing undesirable loading and frictional forces on the specimen. All models can be used for both 3- and 4-point tests.



MODEL 642.001

- » 1, 2, 3, 4, and 5 mm rollers are included
- » Adjustable spans
- » Metric and US Customary scales.

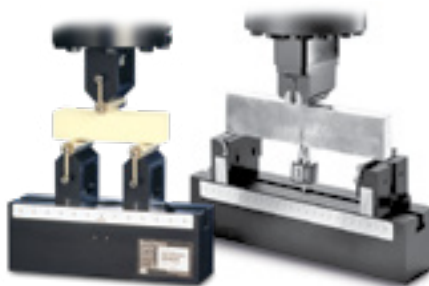

MODELS 642.01 AND 642.10:

- » Meet ASTM D790 and ASTM D7246 flex test
- » Model 642.01 also meets ISO 14125 and EN 2746
- » Model 642.10 also meets ASTM E399 and EN 2562
- » Flexible configurations provide either a region of constant stress or a line of maximum stress
- » Adjustable spans
- » Metric and US Customary scales.



MODEL 642.25:

- » Well-suited for ASTM E-399 fracture toughness tests
- » Adjustable spans
- » Metric and US Customary scales.

MODEL	DESCRIPTION	UPPER FIXTURE SPAN	LOWER FIXTURE SPAN	DYNAMIC FORCE*	TEMPERATURE RANGE	COMBINED HEIGHT**	PART NUMBER
642.001	3 & 4 point bend fixture		14-60 mm (0.6-2.4 in)	0.9 kN (200 lbf)	-129°C to 149°C (-200°F to 300°F)	61 mm (2.4 in)	100-027-125
642.01A-01	3 point bend fixture	NA	24-152 mm (0.94-6.0 in)***	10 kN (2.2 kip)	-129°C to 149°C (-200°F to 300°F)	172 mm (6.8 in)	051-427-701
642.01A-02	3 & 4 point bend fixture	24-76 mm (0.94-3.0 in)***	24-152 mm (0.94-6.0 in)***	10 kN (2.2 kip)	-129°C to 149°C (-200°F to 300°F)	243 mm (9.6 in)	051-427-801
642.10B-01	3 point bend fixture	NA	38-305 mm (1.5-12.0 in)****	100 kN (22 kip)	-129°C to 177°C (-200°F to 350°F)	273 mm (10.75 in)	050-032-601
642.10B-02	3 & 4 point bend fixture	38-152 mm (1.73-6.0 in)****	38-305 mm (1.5-12.0 in)****	100 kN (22 kip)	-129°C to 177°C (-200°F to 350°F)	356 mm (14.00 in)	050-032-701
642.25B-01	3 point bend fixture	NA	79-610 mm (3.12-24.0 in)*****	250 kN (55 kip)	-129°C to 177°C (-200°F to 350°F)	470 mm (18.50 in)	050-876-201
642.25B-02	3 & 4 point bend fixture	50.8-203 mm (2-8 in)*****	79-610 mm (3.12-24.0 in)*****	250 kN (55 kip)	-129°C to 177°C (-200°F to 350°F)	660 mm (26.00 in)	050-876-301

* Static and dynamic force rating depends upon roller diameter.

** Dimension depends upon roller diameter. Largest roller diameter shown.

*** Dimension depends upon roller diameter. 6.35 mm (0.25 in) roller diameter shown.

**** Dimension depends upon roller diameter. 25.4 mm (1 in) roller diameter shown.

***** Dimension depends upon roller diameter. 50.8 mm (2 in) roller diameter shown.

Model 642.01

Roller Assemblies*

DIAMETER	PART NUMBER
4 mm	051-284-607
5 mm	051-284-601
10 mm	051-284-603
0.25 in	051-284-602
0.50 in	051-284-604

Model 642.10 Roller Assemblies*

DIAMETER	PART NUMBER	DIAMETER	PART NUMBER
5 mm	049-578-501	0.25 in	049-578-502
10 mm	049-578-503	0.375 in	049-578-510
15 mm	049-578-505	0.50 in	049-578-504
20 mm	049-578-507	0.75 in	049-578-506
25 mm	049-578-509	1.00 in	049-578-508

Model 642.25 Roller Assemblies*

DIAMETER	PART NUMBER	DIAMETER	PART NUMBER
20 mm	050-875-202	0.75 in	050-875-201
30 mm	050-875-204	1.00 in	050-875-203
40 mm	050-875-207		
50 mm	050-875-209		



*Includes one roller and attachment springs.

Order quantity 3 for 3-point bend and 4 for 4-point bend configurations.

Roller assemblies listed above are not included with bend fixtures and must be purchased separately.

MTS Exceed® 3-Point Bend Fixtures

- » Designed to use with MTS Exceed load frames; can also be used on MTS Criterion® load frames with appropriate conversion adapters
- » Loading edge and supports can be changed to optional parts or customized designs
- » Fast and accurate specimen positioning with centering device
- » Adjustable stepless lower span on the support beam

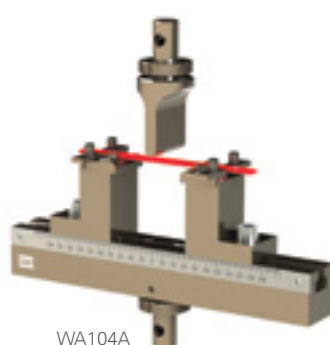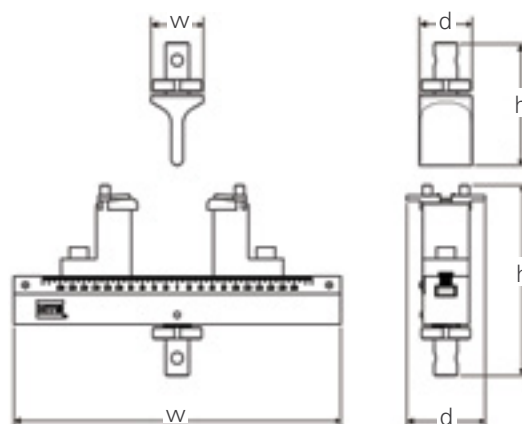

WA104A



WA204A



ZWA304

METHOD		STANDARD	FIXTURE OPTIONS
Laminae & Laminate	Flexure	ISO 14125 (3P)	Model WA104A; WA204A & ZWA304 with Loading Edge R5 Supporting R2 or R5
		ISO 14130	Model WA104A; WA204A & ZWA304 with Loading Edge R5 Supporting R2
		ASTM D7264	Model WA104A; WA204A & ZWA304 with Loading Edge R5 Supporting R5
		EN 2377	Model WA104A with Loading Edge R3 or R5 Supporting R2; WA204A & ZWA304 with Loading Edge R5 Supporting R2
		EN 2746	Model WA104A; WA204A & ZWA304 with Loading Edge R5 Supporting R2

Specifications

MODEL	WA104A	WA204A	ZWA304
Description	10 kN Bend fixture, plastics	20 kN Bend fixture, plastics	30 kN Bend fixture, SST
Rated Force	10 kN	20 kN	30 kN
Temperature Range	Room temperature	Room temperature	-70°C to 350°C
Weight (upper part)	500 g	670 g	510 g
Weight (lower part)	4.95 kg	9.22 kg	4.7 kg
Dimensions (h*w*d) (upper part)	106 mm × 42 mm × 42 mm	108 mm × 42 mm × 42 mm	108 mm × 42 mm × 42 mm
Dimensions (h*w*d) (lower part)	151 mm × 280 mm × 77 mm	180 mm × 340 mm × 88 mm	180 mm × 190 mm × 88 mm
Loading Edge	R5	R5	R5
Supporting	R2	R2	R2
Maximum Span	160 mm	200 mm	80 mm
Maximum Specimen Width	40 mm	45 mm	45 mm
Part Number	100-302-794	100-302-795	100-302-798

Optional Loading Edge

MODEL	EDGE RADIUS	WIDTH	COMPATIBLE FIXTURE
WA104A-06Ab	R2	40 mm	WA104A
WA104A-06Bb	R3	40 mm	WA104A
WA104A-06Cb	R7.6	40 mm	WA104A
WA204A-10Ac	R2	45 mm	WA204A
WA204A-10Bb	R10	45 mm	WA204A

Optional Supporting

MODEL	EDGE RADIUS	WIDTH	COMPATIBLE FIXTURE
WA104A-08Ab	R5	40 mm	WA104A
WA204A-06Ab	R5	45 mm	WA204A
ZWA304-04A	R5	45 mm	ZWA304

Composite Fixtures

Modified Celanese Compression Loading Fixture


- » Recommended to test in accordance with ISO 14126 Method 1A
- » Constructed out of high quality stainless steel
- » Design based on the University of Wyoming Modified Celanese Compression Test Fixture
- » Supported specimen dimensions:
 - Maximum width: 12.7 mm (0.5 in)
 - Thickness (with tabs): 3.8 - 6.35 mm (0.15 - 0.25 in)
 - Length: 114.3 mm (4.5 in)
- » Includes wedges with flame sprayed high friction surface
- » Requires compression platens for mounting (*purchased separately*)



STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	PART NUMBER
88 kN (20 kip)	-85 to 122°C (-120 to 250°F)	≈ 7.3 kg (16 lbs)	Ø 89 mm (3.5 in) x 191 mm (7.5 in)	100-351-817

IITRI Compression Loading Fixture

- » Recommended to test in accordance with ASTM D3410 and ISO 14126 Method 1B Constructed out of high quality stainless steel
- » Supported specimen dimensions:
 - Maximum width: 25.4 mm (1 in)
 - Maximum thickness (with tabs): 15.2 mm (0.6 in)
 - Length: 140 mm (5.5 in)
- » Includes sets of wedges to accommodate specimen thicknesses from 5.1 - 10.2 mm (0.2 - 0.4 in)
Wedges that support other specimen thicknesses are available on request.
- » Requires threaded adapters or compression platens for mounting (*purchased separately*)

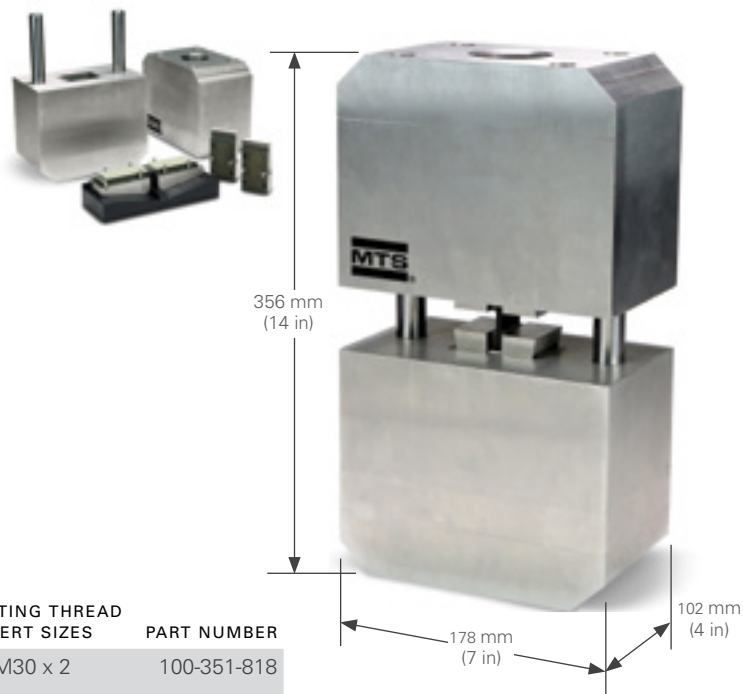


STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
250 kN (55 kip)	-85 to 122°C (-120 to 250°F)	≈ 36 kg (80 lbs)	178 mm (7 in) x 102 mm (4 in) x 356 mm (14 in)	M30 x 2	100-351-818

Composite Fixtures

Combined Loading Compression (CLC) Test Fixture

- » Recommended to test in accordance with ASTM D6641
- » Constructed out of high quality stainless steel
- » Supported specimen dimensions:
 - Maximum width: 25.4 mm (1 in)
 - Maximum thickness (with tabs): 12.7 mm (0.5 in)
 - Length: 140 mm (5.5 in)
- » Requires compression platens for mounting (*purchased separately*)

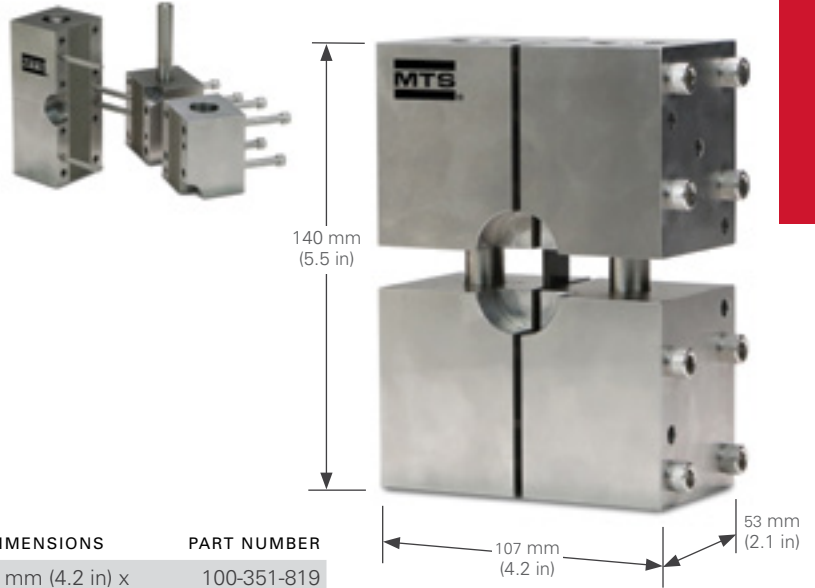

STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	PART NUMBER
89 kN (20 kip)	-85 to 122°C (-120 to 250°F)	≈ 6.8 kg (15 lbs)	107 mm (4.2 in) x 53 mm (2.1 in) x 140 mm (5.5 in)	100-351-819

V-Notched Rail Shear Test Fixture

- » Recommended to test in accordance with ASTM D7078
- » Constructed out of high quality stainless steel
- » Supported specimen dimensions:
 - Width: 55.6 mm (2.2 in)
 - Maximum thickness: 12.7 mm (0.5 in)
 - Maximum length: 76 mm (3.0 in)
- » Requires threaded adapter for top and bottom mounting (*purchased separately*)




STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
44 kN (10 kip)	-152 to 318°C (-240 to 600°F)	≈ 7.7 kg (17 lbs)	102 mm (4 in) x 64 mm (2.5 in) x 165 mm (6.5 in)	M30 x 2	100-351-818

Composite Fixtures

V-Notched Beam (Iosipescu) Shear Fixture


- » Recommended to test in accordance with ASTM D5379
- » Constructed out of high quality stainless steel
- » Supported Specimen Dimensions:
 - Width: 19 mm (0.75 in)
 - Thickness: 0.76 - 12.7 mm (0.03 -0.5 in)
 - Length: 76 mm (3.0 in)
 - Notch: 90 degree with 1.27 mm (0.05 in) radius minimum
- » Includes adjustable wedges
- » Requires threaded adapter for top and compression platen for bottom mounting (*purchased separately*)

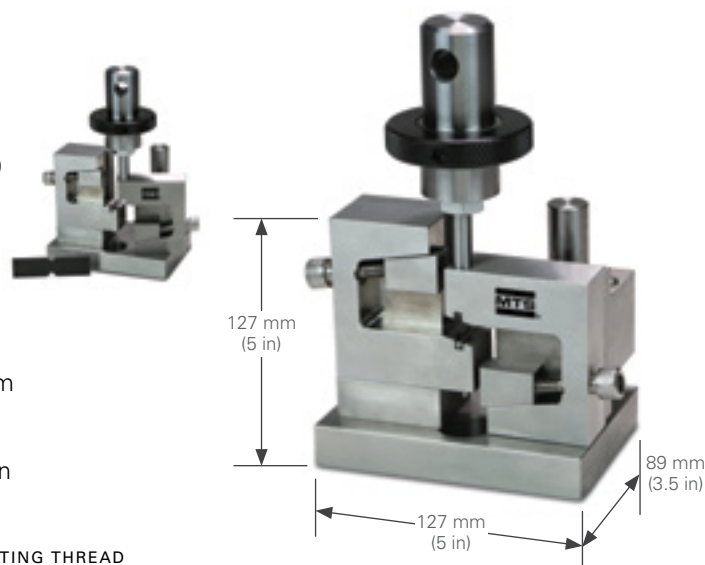

STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
44 kN (10 kip)	-85 to 122°C (-120 to 250°F)	≈ 6.8 kg (15 lbs)	127 mm (5 in) x 89 mm (3.5 in) x 127 mm (5 in)	0.5 - 20 in	100-087-239

Short Beam Shear and Three-Point Flexure Fixture

- » Recommended to test in accordance with ASTM D2344
- » Constructed out of high quality stainless steel
- » Supported Specimen Dimensions:
 - Maximum Width: 38 mm (1.5 in)
 - Maximum Thickness: 50 mm (2 in)
 - Maximum Length: 152 mm (6 in)
- » Adjustable support span
- » Requires male clevis adapter for top and threaded adapters for bottom mounting (*purchased separately*)

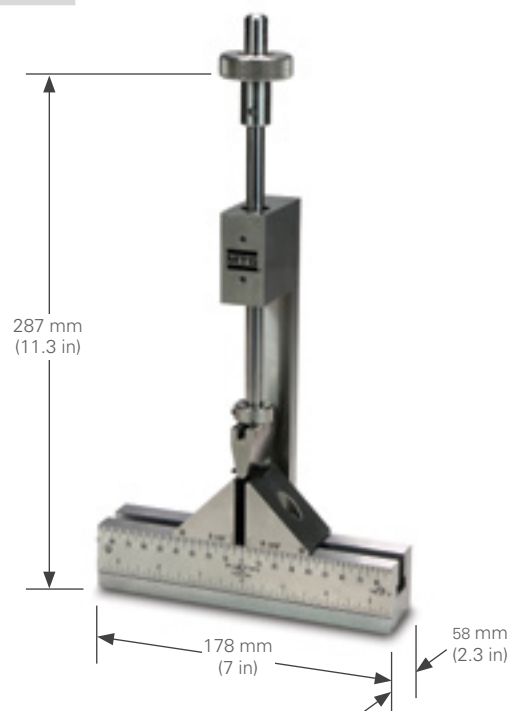

LOWER FIXTURE SPAN	LOADING NOSE DIAMETER	SUPPORTS DIAMETER	LOADING NOSE / SUPPORTS WIDTH
3.2 - 152 mm (0.125 - 6 in)	6.35 mm (0.25 in)	3.175 mm (0.125 in)	38 mm (1.5 in)

STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	TOP MOUNTING MALE CLEVIS	BOTTOM MOUNTING THREAD INSERT SIZES	PART NUMBER
8.9 kN (2 kip)	-85 to 122°C (-120 to 250°F)	≈ 6.8 kg (15 lbs)	178 mm (7 in) x 58 mm (2.3 in) x 287 mm (11.3 in)*	12 mm (Type O)	0.5 - 20 in	100-351-821

* Plus any specimen up to 51 mm (2 in)


Additional Loading Nose and Supports

STANDARD	MATERIAL	LOWER FIXTURE SPAN	LOADING NOSE DIAMETER	SET OF SUPPORTS DIAMETER	LOADING NOSE / SUPPORTS WIDTH	PART NUMBER
ISO 14130	High quality stainless steel	4 - 152 mm (0.157 - 6 in)	10 mm (0.394 in)	–	38 mm (1.5 in)	100-352-347
			–	4 mm (0.157 in)		100-352-348

Composite Fixtures

Mixed Mode Bending Fixture

- » Recommended to test in accordance with ASTM D6671
- » Constructed out of high quality stainless steel and aluminum
- » Supported specimen dimensions:
 - Maximum width: 38 mm (1.5 in)
 - Maximum thickness: 6.35 mm (0.25 in)
 - Maximum length: 228 mm (9.0 in)
- » Includes 5 sets of specimen hinges
- » Requires threaded adapter for top and compression platen for bottom mounting (*purchased separately*)

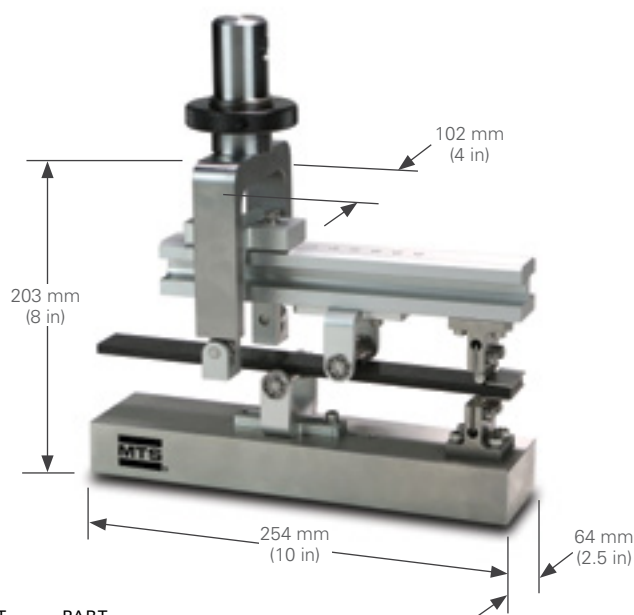


STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
4.4 kN (1 kip)	-85 to 122°C (-121 to 250°F)	≈ 7.3 kg (16 lbs)	254 mm (10 in) x 102 mm (4 in) x 203 mm (8 in)	0.25 - 28 in	100-351-822

Open / Filled Hole Compression Fixture

- » Recommended to test in accordance with ASTM D6484, ASTM D6742 and BS 07260
- » Constructed out of high quality stainless steel
- » Supported specimen dimensions:
 - Width: 38 mm (1.5 in)
 - Maximum thickness: 12.7 mm (0.5 in)
 - Maximum length: 305 mm (12 in)
- » Requires compression platens or hydraulic grips for mounting (*purchased separately*)

Note: Fixture thickness for gripping = 30 mm (1.18 in) + specimen thickness

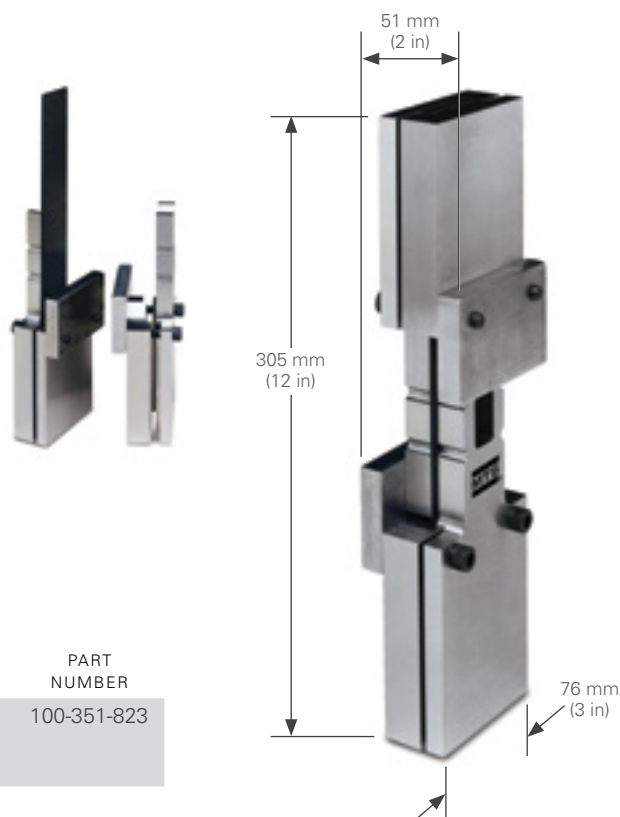


STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	PART NUMBER
222 kN (50 kip)	-152 to 318°C (-240 to 600°F)	≈ 6.8 kg (15 lbs)	76 mm (3 in) x 51 mm (2 in) x 305 mm (12 in)	100-351-823

Composite Fixtures

Compression After Impact Test Fixture


- » Recommended to test in accordance with ASTM D7137
- » Constructed out of high quality stainless steel
- » Supported specimen dimensions:
 - Width: 102 mm (4 in)
 - Thickness: 3.175 - 12.7 mm (0.125 - 0.500 in)
 - Length: 152 mm (6 in)
- » Requires threaded adapter for top and compression platen for bottom mounting (*purchased separately*)

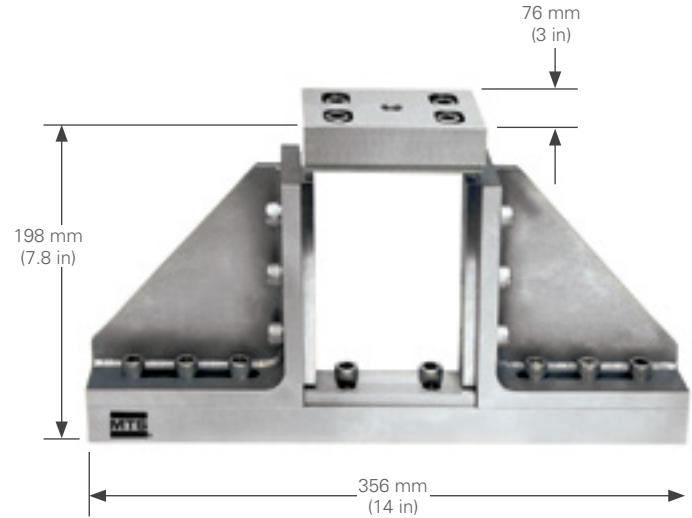

STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
222 kN (50,000 lbs)	-152 to 318°C (-240 to 600°F)	≈ 16 kg (35 lbs)	356 mm (14 in) x 76 mm (3 in) x 198 mm (7.8 in)	0.5 - 13 in	100-351-824

Flatwise Plane Shear Fixture, Tensile Mode

- » Recommended to test in accordance with ASTM C273 and ASTM C394 (Fatigue)
- » Constructed out of high quality stainless steel
- » Includes three sets of aluminum bonding plates
- » Supported specimen dimensions:
 - Maximum width: 76 mm (3 in)
 - Thickness: 6.3 - 19.1 mm (0.25 - 0.75 in)
(optional plates for thicker samples on request)
 - Maximum length: 229 mm (9 in)
- » Includes 5 sets of specimen hinges
- » Requires threaded adapter for top and bottom mounting (*purchased separately*)

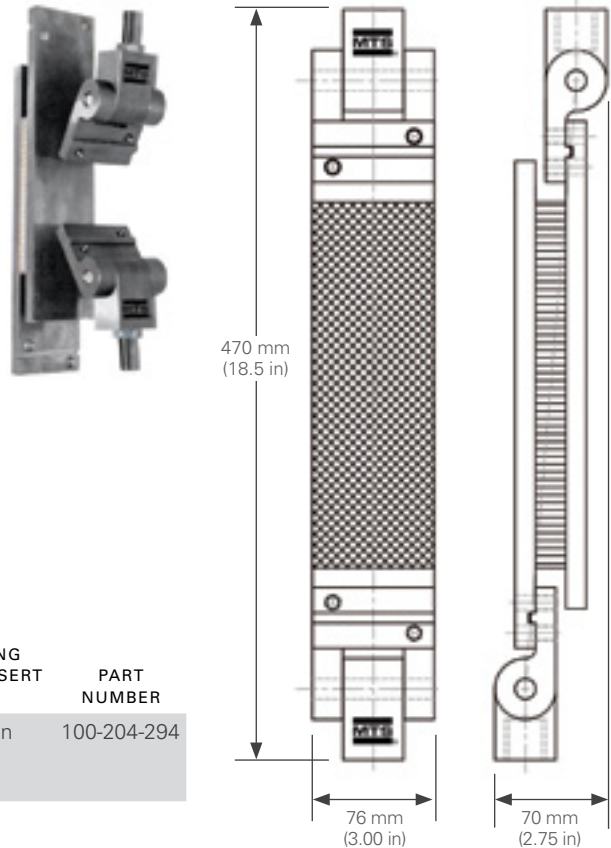


STATIC FORCE	TEMPERATURE RATING*	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
89 kN (20 kip)	-152 to 318°C (-152 to 600°F)	≈ 14.5 kg (32 lbs)	76 mm (3 in) x 70 mm (2.75 in) x 470 mm (18.5 in)	1 - 14 in	100-204-294

* Temperature Range of Aluminum Bonding Plates -29 to 49°C (-20 to 120°F)

Composite Fixtures

Flatwise Plane Shear Fixture, Compression Mode

- » Recommended to test in accordance with ASTM C273 and ASTM C394 (Fatigue)
- » Constructed out of high quality stainless steel
- » Includes three sets of aluminum bonding plates
- » Supported specimen dimensions:
 - Maximum width: 76 mm (3 in)
 - Thickness: 6.3 - 19.1 mm (0.25 - 0.75 in) (optional plates for thicker samples on request)
 - Maximum length: 229 mm (9 in)
- » Requires threaded adapter for top and bottom mounting (*purchased separately*)

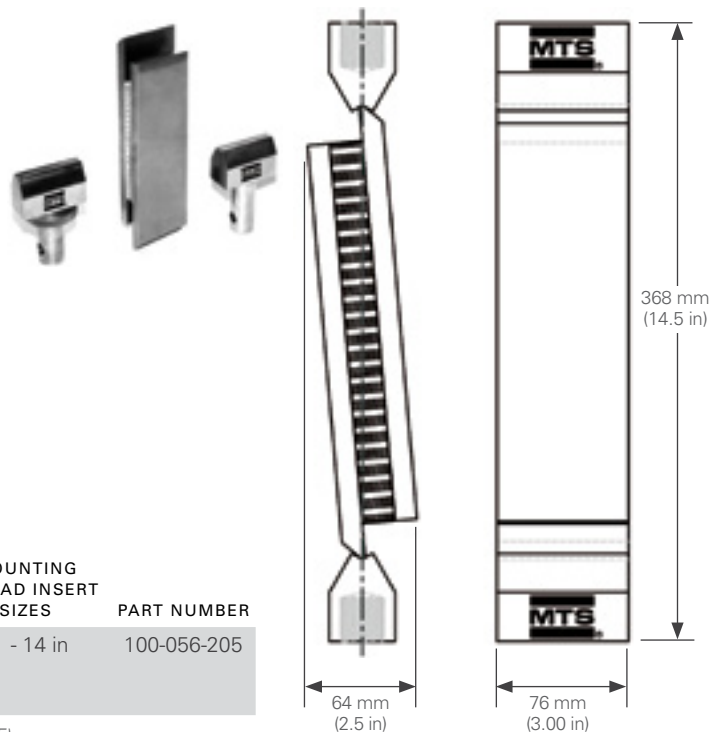


STATIC FORCE	TEMPERATURE RATING*	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
89 kN (20 kip)	-152 to 318°C (-152 to 600°F)	≈ 14.5 kg (32 lbs)	76 mm (3 in) x 64 mm (2.5 in) x 368 mm (14.5 in)	1 - 14 in	100-056-205

* Temperature Range of Aluminum Bonding Plates -29 to 49°C (-20 to 120°F)

Adjustable Edgewise Compression Fixture

- » Recommended to test in accordance with ASTM C364
- » Supported specimen dimensions:
 - Maximum width: 65 mm (2.5 in)
 - Maximum thickness: 22 mm (0.86 in)
- » Faces:
 - Maximum opening: 22 mm (0.86 mm)
 - Width: 65 mm (2.5 in)
 - Material: Rubber
- » Minimize system errors with self-aligning mechanism that secures the specimen
- » Male Clevis Adapters (40 mm) for top and bottom mounting (included)

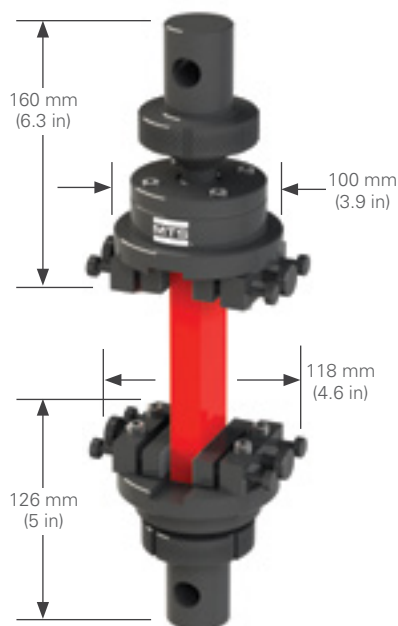


STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS (W*D*H)	MOUNTING MALE CLEVIS	MODEL NUMBER	PART NUMBER
50 kN (11.2 kip)	Room Temperature	≈ 6.8 kg (15 lbs)	118 mm (4.6 in) x 100 mm (3.9 in) x 286 mm (11.3 in)	(Upper part) 40 mm/ (Lower part) 40 mm	DKF1005089.03	100-302-784

Composite Fixtures

Three & Four Point Sandwich Beam Flexure / Shear Fixture

- » Recommended to test in accordance with ASTM C393, ASTM D5467, ASTM D7249 and ASTM D7250
- » Constructed out of high strength steel with a durable black oxide finish (except for rollers and pads)
- » Supported specimen dimensions:
 - Maximum width: 100 mm (4 in)
 - Maximum length: 610 mm (24 in)
- » Adjustable loading and support spans
- » Loading and support bars are supplied with loading pins and flat steel loading blocks held in alignment with springs (rubber pads not included)
- » Requires threaded adapter for top and bottom mounting (*purchased separately*)

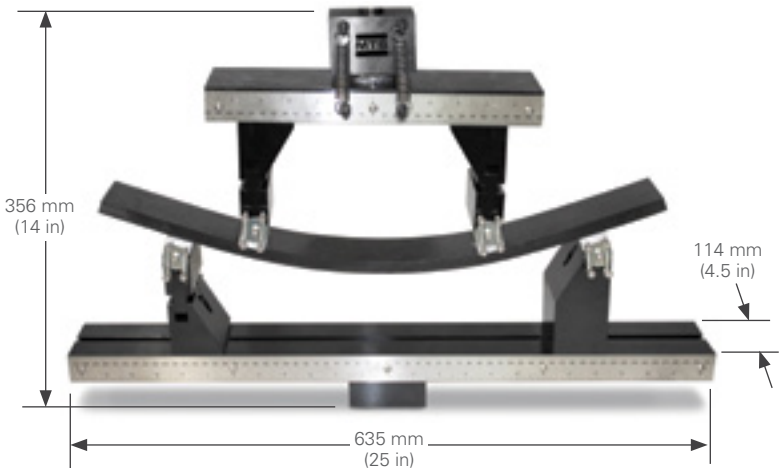


UPPER FIXTURE SPAN	LOWER FIXTURE SPAN	LOADING PINS DIAMETER	SUPPORT PINS DIAMETER	LOADING & SUPPORT PINS WIDTH
51- 305 mm (2 - 12 in))	152 - 610 mm (6 - 24 in)	25.4 mm (1 in)	25.4 mm (1 in)	100 mm (4 in)

STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING THREAD INSERT SIZES	PART NUMBER
11 kN (2.5 kip)	-85 to 122°C (-120 to 250°F)	≈ 52 kg (114 lbs)	635 mm (25 in) x 114 mm (4.5 in) x 356 mm (14 in)	1 - 14 in	100-351-826

Climbing Drum Peel Fixture with Roller Type Grips

- » Recommended to test in accordance with ASTM D1781
- » Constructed out of high quality stainless steel with an aluminum drum
- » Supported specimen dimensions:
 - Width: 25.4 - 102 mm (1 - 4 in)
 - Thickness: 0.762 - 25.4 mm (0.03 - 1 in)
 - Length: 254 mm (10 in)
- » Includes Type D Male Clevis Adapters for top and bottom mounting

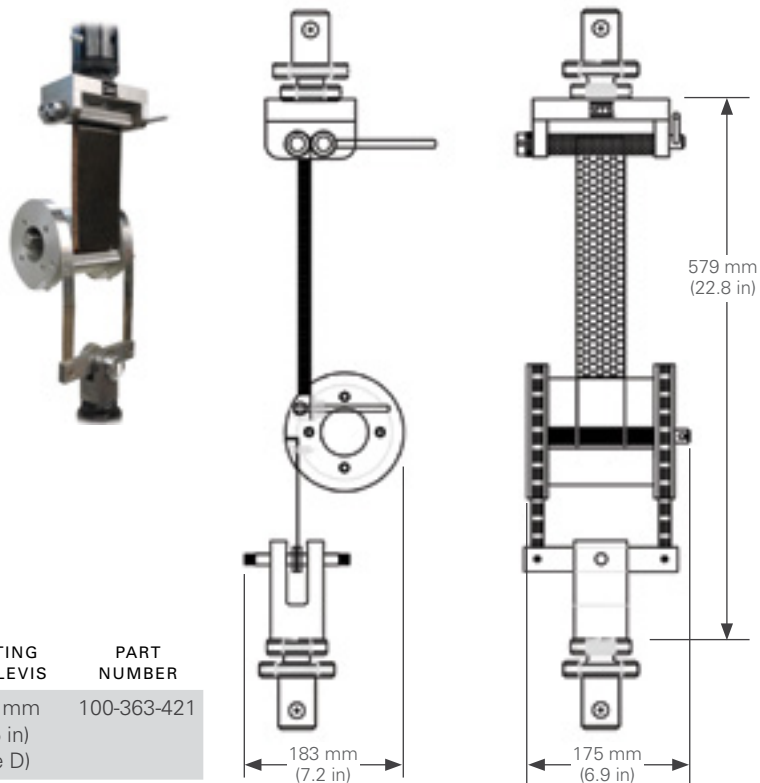

STATIC FORCE	TEMPERATURE RATING	WEIGHT	DIMENSIONS	MOUNTING MALE CLEVIS	PART NUMBER
2.2 kN (0.5 kip)	-29 to 49°C (-20 to 120°F)	≈ 18 kg (40 lbs)	183 mm (7.2 in) x 175 mm (6.9 in) x 579 mm (22.8 in)	31.75 mm (1.25 in) (Type D)	100-363-421

A Comprehensive Array of Electromechanical Accessories

Grips and Fixtures for Electromechanical Systems

MTS complements its electromechanical testing lines with a comprehensive array of accessories to fulfill a full spectrum of material and small component testing – from basic quality control, to complex biomedical simulations,

to demanding research and development applications. This catalog includes three distinct accessory families to accommodate your specific and evolving testing needs:



MTS Advantage Accessories

Highly versatile and full-featured wedge, pneumatic, and screw action grips for demanding R&D testing of advanced composites and alloys. Ideal for the specific needs of the high-end researcher, this accessory family accommodates a very broad range of clamping force and temperature requirements and features numerous control and grip face options. MTS stands behind the MTS Advantage family line with a three-year warranty – the best in the industry!

MTS Fundamental Accessories

Basic affordable grips and fixtures for standard testing of plastics, textiles, rubber, wire, rope, and more. These accessories feature a universal adapter design and optional threaded frame adapters to facilitate easy installation onto both MTS electromechanical and servohydraulic load frame systems, as well as other electromechanical test systems. Anti-rotation features and integrated alignment pins enhance test accuracy and repeatability.



Bionix Accessories

Affordable and extremely durable grips, fixtures, and platens for accurately replicating biomaterial and medical device service environments and extending the utility of Bionix electromechanical and servohydraulic test solutions. These accessories feature a universal adapter design and optional threaded frame adapters. This facilitates easy installation onto both MTS electromechanical and servohydraulic load frame systems as well as other electromechanical test systems. Anti-rotation features and integrated alignment pins enhance test accuracy and repeatability.

Can't find what you need?

We offer many more grips and fixtures. Contact your local sales representative or applications engineer to find the model that meets your exact needs.

10 N MTS Advantage Pneumatic Grips

Grips Ideal for Thin Strips, Wires and Films with Low Breaking Strength

These 10 N capacity (2 lbf) pneumatic grips are designed for low tensile applications and work well on thin specimens with low breaking strengths like the gold wire in electronics, tungsten wire, thermocouples, gels, and plastic films.

A big advantage these grips have over other low capacity grips is that they have been designed and machined to precise tolerances, thereby eliminating side loads.

Other features include a lightweight design to let you test with your low-capacity load cells without the grips using up too much of its capacity. It also has a self-aligning feature which lets these grips check alignment even while actuated.

These grips feature dual acting, swivel faces. With dual action, both faces move simultaneously to the centerline of the grip and eliminate the bending strains on your specimen that can invalidate your test results. The ability to swivel lets these faces conform to any variation in your specimen's geometry. The faces shown below are sold separately.

These pneumatic grips require an air supply to provide filtered, dry air and an air grip controller. They are compatible with servohydraulic and electromechanical test systems.



FEATURES

- » Designed with plenty of room to manipulate your specimens
- » Swivel faces to accommodate varying specimen widths
- » No side loads to influence your specimen.

10 N (2 lbf) Pneumatic Grips

TEMPERATURE RATING	AIR PRESSURE	FORCE	CLAMPING THICKNESS	MAX SPECIMEN WEIGHT	PART NUMBER
-40°C to 200°C (-40°F to 400°F)	0.6 MPa (80 psi)	31 N (7 lbf)	5 mm (0.20 in)	270 g (9.6 oz)	100-032-017

Type Bm upper grip mounting, type Dm lower.

Faces** for 10 N (2 lbf) Pneumatic Grips


FACE SURFACE	WIDTH	HEIGHT	PART NUMBER
Smooth	15 mm (0.60 in)	8 mm (0.30 in)	056-163-701
Serrated	15 mm (0.60 in)	8 mm (0.30 in)	056-163-702
Rubber*	15 mm (0.60 in)	8 mm (0.30 in)	056-163-703

* Room temperature only.

** Refer to the TechNote on page 105 for guidance.

Electromechanical Attachment Scheme

CLEVIS PIN CONNECTION TYPE	TYPE B	TYPE C	TYPE D	15	30
Load	Up to 10 N	Up to 200 N	Up to 150 kN	Up to 200 N	Up to 150 kN
Clevis Diameter (X)	0.50 in	0.625 in	1.25 in	15 mm	30 mm
Pin Diameter (Y)	0.187 in	0.25 in	0.50 in	6 mm	14 mm
Upper Grip Connection (male)	Bm	Cm	Dm	15 m	30 m
Lower Grip Connection (male)	Dm	Dm	Dm	30 m	30 m
Upper Load Frame Connection (female)	Bf	Cf	Df	15 f	30 f
Lower Load Frame Connection (female)	Df	Df	Df	30 f	30 f

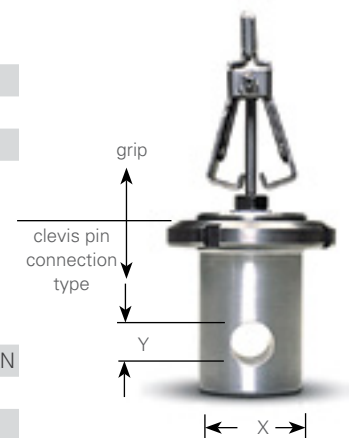

100/200 N and 2 kN MTS Advantage Pneumatic Grips

Designed for testing thin sheets, films, and tapes, these 100/200 N and 2 kN (20/40 lbf and 400 lbf) pneumatic grips will let you perform tests with your low capacity load cells without taking up too much of your load cells capacity.

Both models feature dual acting grip faces that simultaneously move to the centerline of the grip. This assures correct specimen alignment and removes the bending strains that can invalidate your test results.

The design of both grips lets you quickly and easily load and align your specimens. This lets you spend less time loading specimens and more time testing. They also feature an area above the grip faces that is large enough for you to load the specimens while holding them with your fingers.

Both grips use standard pin connectors for the connection to the baseplate and load cells and can be used with servohydraulic or electromechanical test systems.

FEATURES

- » Lightweight design so you can perform low force tests without using up too much of your low-capacity load cells capacity
- » Dual acting grip faces which move to the centerline of the grip for better specimen alignment
- » Large area above the grip faces so you can load your specimens while holding them with your fingers
- » Integrated clamp/release switch located close to the specimen for intuitive inspection.



These 100 N (20 lbf) (top) and 2 kN (400 lbf) (bottom) capacity pneumatic grips can be used to test a variety of thin films, sheets, and tapes. Grip faces are shown on the next page.

Standard Specifications

AIR PRESSURE

0.6 MPa (80 psi)

100/200 N and 2 kN Pneumatic Grips

DESCRIPTION	CLAMPING FORCE	MAXIMUM SPECIMEN THICKNESS	TEMPERATURE RANGE	DIMENSIONS		UPPER GRIP WEIGHT	PART NUMBER
				HEIGHT	WIDTH		
100/200 N (20/40 lbf)	620 N (140 lbf)	10 mm (0.39 in)	-40°C to 100°C (-40°F to 212°F)	180 mm (7.1 in)	114 mm (4.5 in)	1.0 kg (2.0 lb)	100-036-576
2 kN (400 lbf)	7100 N (1600 lbf)	12 mm (0.47 in)	-40°C to 200°C (-40°F to 400°F)	226 mm (8.9 in)	131 mm (5.2 in)	3.2 kg (7 lb)	100-280-342

* Type Cm upper grip mounting, type Dm lower. (See page 104)

** Type Dm upper grip mounting, type Dm lower. (See page 104)

100/200 N and 2 kN MTS Advantage Grips Faces

Faces for the 100/200 N and 2 kN (20/40 lbf and 400 lbf) capacity pneumatic grips are available with a variety of surface coatings so you can use your grips to test a variety of materials.

To find out more about which surface works best on the material you're testing, read the TechNote on this page.

The 100/200 N and 2 kN (20/40 lbf and 400 lbf) capacity grips share the same faces and are sold in sets of four.

Faces for 100/200 N and 2 kN (20/40 lbf and 400 lbf) Pneumatic Grips

FACE SURFACE	HEIGHT	WIDTH	PART NUMBER
Rubber (smooth)*	12.5 mm (0.5)	25 mm (1 in)	056-163-825
Smooth	25 mm (1 in)	25 mm (1 in)	056-163-801
Corrugated	25 mm (1 in)	25 mm (1 in)	056-163-802
Serrated	25 mm (1 in)	25 mm (1 in)	056-163-803
Diamond Tip	25 mm (1 in)	25 mm (1 in)	056-163-804
Rubber (matte)*	25 mm (1 in)	25 mm (1 in)	056-163-805
Rubber (corrugated)*	25 mm (1 in)	25 mm (1 in)	056-163-806
Rubber (smooth)*	25 mm (1 in)	25 mm (1 in)	056-163-807
Line Contact (r=4.75 mm)	25 mm (1 in)	25 mm (1 in)	056-163-808
Line Contact w/Rubber Opposite* (r=4.75 mm)	25 mm (1 in)	25 mm (1 in)	056-163-826
Line Contact w/Rubber Opposite* (r=2.5 mm)	25 mm (1 in)	25 mm (1 in)	056-163-829
Smooth	25 mm (1 in)	75 mm (3 in)	056-163-809
Corrugated	25 mm (1 in)	75 mm (3 in)	056-163-810
Serrated	25 mm (1 in)	75 mm (3 in)	056-163-811
Diamond Tip	25 mm (1 in)	75 mm (3 in)	056-163-812
Rubber (matte)*	25 mm (1 in)	75 mm (3 in)	056-163-813
Rubber (corrugated)*	25 mm (1 in)	75 mm (3 in)	056-163-814
Rubber (smooth)*	25 mm (1 in)	75 mm (3 in)	056-163-815
Diamond	25 mm (1 in)	150 mm (5.9 in)	056-163-828
Rubber*	25 mm (1 in)	150 mm (5.9 in)	056-163-830
Smooth	38 mm (1.5 in)	58 mm (2.3 in)	056-163-816
Corrugated	38 mm (1.5 in)	58 mm (2.3 in)	056-163-817
Serrated	38 mm (1.5 in)	58 mm (2.3 in)	056-163-818
Diamond Tip	38 mm (1.5 in)	58 mm (2.3 in)	056-163-819
Rubber (matte)*	38 mm (1.5 in)	58 mm (2.3 in)	056-163-820
Rubber (corrugated)*	38 mm (1.5 in)	58 mm (2.3 in)	056-163-821
Rubber (smooth-EPDM)*	38 mm (1.5 in)	58 mm (2.3 in)	056-163-822
Line Contact (r=4.75 mm)	38 mm (1.5 in)	58 mm (2.3 in)	056-163-823
Line Contact w/Rubber Opposite* (r=4.75 mm)	38 mm (1.5 in)	58 mm (2.3 in)	056-163-827
Grab Test	38 mm (1.5 in)	58 mm (2.3 in)	056-163-824

* Rubber grip faces are for use only at room temperature.
Rubber faces can be reconditioned. Contact Order Services.

TECH NOTE


Choosing the Right Surface

IMPROVE YOUR TEST RESULTS

Grip-induced failure and slippage can cost you time and money, so it's important to pick out the right faces for the material you're testing. MTS manufactures faces for our pneumatic grips with a variety of surfaces to accommodate your specimen. Here's a brief description of the surfaces offered:

- » **Corrugated** (top left) for smooth specimens (textile, fabrics, tissue, leather).
- » **Serrated** (top right) for gripping soft materials without causing failure (paper, board, aluminum, copper and steel wire, soft steels).
- » **Rubber Coated** (black colored) to offer extra protection for your thin specimens (films, fibers, aluminum).
- » **Diamond Tipped** (bottom right) to provide you with an aggressive surface for gripping soft materials (soft steels, rigid plastic, wood).
- » **Grab Test** (center) according to international standards, including ASTM D2208 (leather), D5034 (fabrics), D1683 (woven fabrics) and ISO 5082 (fabrics).
- » **Line Contact** (bottom left) faces are for rubber, latex, and cable sheath.

Faces are also available in a wide width for some grip models.
All of the faces for the grips in this catalog are sold in sets of four.

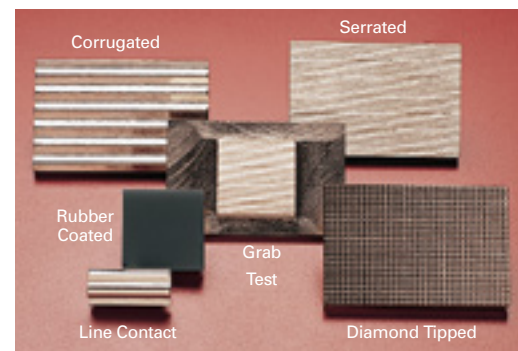

10 kN MTS Advantage Pneumatic Grip

The 10 kN (2.2 kip) capacity pneumatic grips are designed for performing tensile tests on a variety of materials including elastomers, plastics, rigid and semi-rigid films, and sheets.

These grips feature dual acting grip faces that simultaneously move to the centerline of the grip. This feature assures correct specimen alignment and removes the bending strains on your specimens that can invalidate test results.

In addition to these two unique features, these grips save you time by offering quick and easy specimen insertion. This lets you spend less time inserting and aligning specimens and more time testing. Both models also feature a constant clamping load to protect your specimens from damage due to slippage.

A pneumatic controller (see below), available with foot-control to free up both hands for more accurate specimen alignment, and an air supply to provide dry, filtered air are required to operate these grips. The 10 kN MTS Advantage Pneumatic grips share the same faces as the Model 647.02B, which can be found on page 81. They are available for flat, round, or vee specimens with a variety of surfaces. They can be used with servohydraulic or electromechanical test systems.

FEATURES

- » Dual acting grip faces which move to the centerline of the grip for better specimen alignment
- » Easy specimen insertion to save time and a constant clamping load to protect your specimen from damage due to slippage.

Specifications

Temperature range:	-40°C to 200°C (-40°F to 400°F)
Air pressure:	0.6 MPa (80 psi)
Upper grip mounting:	Type Dm (See page 104)
Lower grip mounting:	Type Dm (See page 104)



DESCRIPTION	CLAMPING FORCE	MAXIMUM SPECIMEN THICKNESS	OVERALL HEIGHT	WIDTH	UPPER GRIP WEIGHT	PART NUMBER
10 kN (2.2 kip) Pneumatic Grips	17 kN (4 kip)	10-20 mm (0.4-0.8 in)	183 mm (7.2 in)	200 mm (7.9 in)	17.3 kg (16 lb)	100-034-623

MTS Advantage Pneumatic Grip Controllers

for use with Single Acting Pneumatic Grips

- » Provides Open/close function for MTS Advantage Pneumatic Grips
- » Regulates air pressure and controls air flow to upper and lower grips
- » Choose either a hand switch or foot switch. Magnetic mounts on the hand switch enable convenient mounting on the load frame column or base. Rugged cover on the foot switch prevents accidentally opening or closing the grips. Control sequence of the hand and foot switches is user programmable.
- » Also enables grips to be opened or closed using Test Suite Software*
- » Compatible 6 mm (1/4 in) diameter air tubing. Includes adapter to work with 4 mm (5/32 in) diameter tubing.
- » Connects to load frame E-stop
- » Operates on worldwide voltages (90-264 V AC, 47-63 Hz).



DESCRIPTION	MAXIMUM DRY AIR INTAKE	MAXIMUM OUTPUT	PART NUMBER
Grip Controller	0.85MPa (125psi)	0.7MPa (100psi)	100-393-631
Hand Switch	Not Applicable	Not Applicable	100-393-629
Foot Switch	Not Applicable	Not Applicable	100-393-630

*May require TestSuite software update to operate. Not compatible with Test Works 4.

MTS Advantage Wedge Action Grips

MTS Advantage Wedge Action Grips are versatile general-purpose grips in which the faces remain stationary during loading. This makes it especially useful for applications where screw or pneumatic grips do not provide sufficient clamping force, or where compressive or buckling forces are not desirable during specimen insertion. It works with servohydraulic and electromechanical machines and even accommodates the side insertion of specimens.

FEATURES

- » Quick and easy interchangeable faces for round and flat specimens
- » Self-tightening during test reduces slipping
- » Specimen positioning aids
- » Side loading design

- » Standard pinned adapter for easy installation and removal
- » Suitable for use in environmental chambers
- » Improved serrations secure specimen with minimal clamping force.

APPLICATIONS

- » Composites
- » Metals
- » Plastics
- » Polymers
- » Wood

FUNCTION

Wedges

- » Spring and mechanical retraction
- » Easy access to wedges for quick changeover

Preload

- » Uses right-hand/left-hand thread mechanism for reducing effort



Grip Interface

- » Type D upper and lower mounting (except for 300 kN).
- » 300 kN mounting is M36x2 thread

MTS Advantage Wedge Action Grips

TENSILE CAPACITY	WEIGHT	TEMPERATURE RATING	PART NUMBER
10 kN (2,200 lbf)	4.5 kg (10 lb)	-130°C (-200°F) up to 315°C (600°F) at 7.5 kN (1.6 kip)	056-079-501
30 kN (6,700 lbf)	5.5 kg (12 lb)	-130°C (-200°F) up to 315°C (600°F) at 22 kN (5 kip)	052-862-001
50 kN (11,000 lbf)	7 kg (15 lb)	-130°C (-200°F) up to 315°C (600°F) at 37 kN (8 kip)	054-951-001
100 kN (22,000 lbf)	15 kg (33 lb)	-130°C (-200°F) up to 315°C (600°F) at 75 kN (16 kip)	056-079-801
150 kN (33,000 lbf)	19.6 kg (43 lb)	-130°C (-200°F) up to 315°C (600°F) at 112 kN (24 kip)	053-536-901
300 kN (67,000 lbf)	54 kg (118 lb)	-130°C (-200°F) up to 315°C (600°F) at 213 kN (48 kip)	056-144-702

Flat Wedges

SPECIMEN THICKNESS	CAPACITY	PROFILE	DIMENSIONS	PART NUMBER
0-7.9 mm (0-0.312 in)	10 kN (2.2 kip), 30 kN (6.7 kip), 50 kN (11 kip)	Serrated Steel	50 mm x 25 mm	053-140-801
36 mm-13.2 mm (0.236 in-0.520 in)	10 kN (2.2 kip), 30 kN (6.7 kip), 50 kN (11 kip)	Serrated Steel	50 mm x 25 mm	053-140-802
0-9 mm (0-0.355 in)	100 kN (22 kip), 150 kN (33 kip), 300 kN (67 kip)	Serrated Steel	50 mm x 50 mm	053-537-401
6.4 mm-16 mm (0.25 in-0.63 in)	100 kN (22 kip), 150 kN (33 kip), 300 kN (67 kip)	Serrated Steel	50 mm x 50 mm	053-537-402

Vee Wedges

SPECIMEN DIAMETER MINIMUM-MAXIMUM SIDE/TOP LOADING	CAPACITY	PROFILE	DIMENSIONS	PART NUMBER
3 mm-7.9 mm/7.9 mm (0.118 in-0.312 in/0.312 in)	10 kN (2.2 kip), 30 kN (6.7 kip), 50 kN (11 kip)	Serrated Steel	50 mm x 25 mm	053-140-803
7 mm-9.5 mm/12.7 mm (0.276 in-0.375 in/0.500 in)	10 kN (2.2 kip), 30 kN (6.7 kip), 50 kN (11 kip)	Serrated Steel	50 mm x 25 mm	053-140-804
11.5 mm-12.7 mm/16 mm (0.453 in-0.500 in/0.630 in)	10 kN (2.2 kip), 30 kN (6.7 kip), 50 kN (11 kip)	Serrated Steel	50 mm x 25 mm	053-140-805
3.2 mm-5.8 mm/7.6 mm (0.125 in-0.230 in/0.300 in)	100 kN (22 kip), 150 kN (33 kip), 300 kN (67 kip)	Serrated Steel	50 mm x 50 mm	053-537-405
5 mm-10.4 mm/12.5 mm (0.196 in-0.410 in/0.492 in)	100 kN (22 kip), 150 kN (33 kip), 300 kN (67 kip)	Serrated Steel	50 mm x 50 mm	053-537-403
12.2 mm-15.5 mm/19.5 mm (0.482 in-0.610 in/0.768 in)	100 kN (22 kip), 150 kN (33 kip), 300 kN (67 kip)	Serrated Steel	50 mm x 50 mm	053-537-404

MTS Advantage Screw Action Grips

The MTS Advantage Screw Action Grip offers versatility and a high clamp force. This makes it especially good for applications where pneumatic grips do not provide sufficient clamping force. It works with servohydraulic and electromechanical machines and even accommodates the side insertion of specimens.

FEATURES

- » Quick and easy interchangeable faces
- » Side loading design
- » Suitable for use in environmental chambers
- » An alignment guide to ensure concentricity
- » Right or left hand operation
- » Ability to test lap shear specimens
- » Faces are compatible with MTS pneumatic grips
- » Flexible mounting allows adapters to be easily changed
- » Accommodates threaded configurations
- » Approximately twice the clamp force of comparably rated pneumatic grips
- » Interchangeable, resilient pucks to allow follow-through action, to compensate for specimen neck-down, minimize damage to delicate specimens, and increase clamp load for difficult specimens
- » Can be used in a fixed (one side follows through) and non-fixed (both sides follow through) configuration
- » Room to hold flexible specimens above and below the faces
- » Knurled screw for hand tightening of specimen during installation
- » Low profile grip body increases test space
- » Faces pivot for self alignment and reduced likelihood of breakage at the specimen face contact
- » Temperature range of -129°C to 200°C (-200°F to 400°F) when using the aluminum version of the pucks described above.

APPLICATIONS

- » Metals
- » Plastics*
- » Polymers
- » Wood

* Independently actuated clamps provide accurate modulus for plastic specimens.



The MTS Advantage Screw Action grips come with 2 socket head cap screws.

MTS Advantage Screw Action Grips

CAPACITY	SPECIMEN THICKNESS	WEIGHT	UPPER MOUNTING/ LOWER MOUNTING*	PART NUMBER
100 N	13 mm	0.36 kg (0.80 lbs)	Type Cm Type Dm	055-426-701
2 kN	19 mm	0.95 kg (2.1 lbs)	Type Dm Type Dm	055-426-801
5 kN	25 mm	2.1 kg (4.6 lbs)	Type Dm Type Dm	055-426-901
10 kN	25 mm	4.9 kg (10.8 lbs)	Type Dm Type Dm	100-030-185

*See page 105

5 kN and 10 kN Grip Faces

(see page 106 for the 100 N and 2 kN grip faces)

FACE SURFACE	HEIGHT	WIDTH	PART NUMBER
Smooth	50 mm (2 in)	75 mm (3 in)	056-163-901
Corrugated	50 mm (2 in)	75 mm (3 in)	056-163-902
Serrated	50 mm (2 in)	75 mm (3 in)	056-163-903
Diamond	50 mm (2 in)	75 mm (3 in)	056-163-904
Rubber (matte)*	50 mm (2 in)	75 mm (3 in)	056-163-905
Rubber (corrugated)*	50 mm (2 in)	75 mm (3 in)	056-163-906
Rubber (smooth)*	50 mm (2 in)	75 mm (3 in)	056-163-907
Line Contact	50 mm (2 in)	75 mm (3 in)	056-163-908
Grab Test	50 mm (2 in)	75 mm (3 in)	056-163-909
Line Contact w/rubber opposite	50 mm (2 in)	75 mm (3 in)	056-163-910

*Rubber grip faces are for use at room temperature.

To Order

- 1 Select the required grip capacity.
- 2 Choose the grip faces that meet your application (see the TechNote on page 106). For the 100 N and 2 kN grips, refer to the faces on page 106. For the 5 kN and 10 kN grips, use the faces on this page.
- 3 Determine if you need upper and lower adapters.

MTS Fundamental Pneumatic Bollard Grips

MTS Fundamental Pneumatic Bollard Grips are value-priced, “horn” style grips designed to reduce stress concentration on specimens and avoid grip-induced failures. They are suitable for testing cords, filaments, fibers, fine wire, and yarns at temperatures ranging from -10°C to 80°C (14°F to 176°F). Specimen slippage is minimized by a pneumatic gripping mechanism and clamping forces are easily controlled by adjusting air supply pressure. Operation of these grips requires a dry, filtered air supply and a pneumatic controller supplied by MTS (see MTS Advantage Pneumatic Grip Controllers, page 107).

FEATURES

- » “Horn” style design
- » Pneumatic gripping mechanism – requires controller
- » Temperature range: -10°C to 80°C (14°F to 176°F)
- » Applications: Cords, filaments, fibers, fine wire, yarns.



MTS Fundamental Pneumatic Bollard Grips

CAPACITY	MAX SPECIMEN DIAMETER	UPPER/LOWER MOUNTING	APPLICATIONS	PART NUMBER
200 N (40 lbf)	1 mm (0.04 in)	Cm/Dm	cotton threads, natural silks, spandex yarn, precious metals	100-121-778
2 kN (450 lbf)	4 mm (0.16 in)	Dm/Dm	technical yarns, thin copper cables, magnetic tape	100-140-565
5 kN (1 kip)	5 mm (0.20 in)	Dm/Dm	nylon tennis racquet strings, high strength composite	100-139-065
5 kN (1 kip)	16 mm (41.3 in)	Dm/Dm	cord, yarn, wire cord, yarn, wire	100-034-764 100-034-765

MTS Fundamental Coefficient of Friction Grips

MTS Fundamental Coefficient of Friction Grips employ a pulley and string mechanism to measure the force required to pull a friction sled over a material specimen. They are suitable for meeting a variety of material testing standards, including ASTM D1894 (plastic film), TAPPI T542 (paper and cardboard), TAPPI T549 (non-fibrous materials), and TAPPI T816 (corrugated). The upper grip is not included.

FEATURES

- » Pulley and string mechanism
- » Applications: Plastic film, paper and cardboard, non-fibrous materials, corrugated

FIXTURE	PART NUMBER
Coefficient of Friction	100-087-526



MTS Fundamental 90° Peel Fixture

MTS Fundamental 90° Peel Fixtures employ precision bearings to maintain perpendicular axial alignment while measuring the force required to peel adhered material from a clamped-down substrate. They are suitable for meeting a variety of material testing standards, including ASTM D1876, FINAT FTM1/2/3, and AFERA 4015 T4. The upper grip is not included.

FEATURES

- » Precision bearings
- » Applications: ASTM D1876, FINAT FTM1/2/3, and AFERA 4015 T4.

FIXTURE	PART NUMBER
90° Peel	100-170-762




Bionix Vise Action Grips

Bionix Vise Action Grips are miniature, medium, and heavy-duty general purpose tensile grips that use a clamping screw to tighten a vise; each order includes a set of rubber faces. Available in a variety of load capacities and sizes, they are suitable for testing bandages, bio-textiles, diapers, and plastic films at temperatures ranging from -10°C to 50°C. Comprised of durable stainless steel, they can be safely deployed within a saline environmental bath. These accessories feature a universal adapter design and an optional threaded frame adapter to facilitate easy installation onto both MTS electromechanical and servohydraulic load frame systems, as well as other electromechanical test systems. Anti-rotation features and integrated alignment pins enhance test accuracy and repeatability.

FEATURES

- » Stainless steel – suitable for saline environmental bath
- » Set of rubber faces standard with each order
- » Universal adapter design (Type D standard)
- » Optional threaded frame adapter
- » Anti-rotation features and integrated alignment pins
- » Temperature range: -10°C to 50°C
- » Applications: Bandages, bio-textiles, diapers, plastic films.

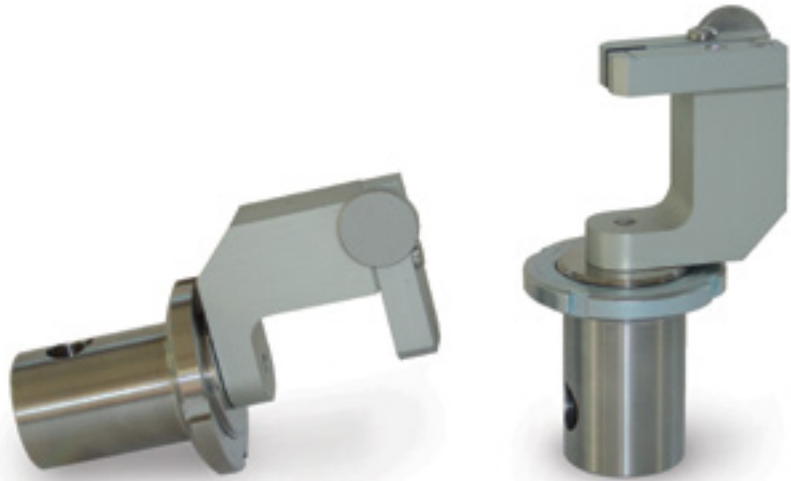

Model 100



Model 2000/5000

MTS Bionix Vise Action Grips

MODEL	MAXIMUM LOAD	FACE PROFILE	MAXIMUM SPECIMEN THICKNESS	MAXIMUM SPECIMEN WIDTH	WEIGHT (UPPER)	HEIGHT (FROM CLEVIS CENTER)	PART NUMBER
100	100 N	Rubber	2 mm	10 mm	560 g	100 mm	100-186-411
2000	2 kN	Rubber	3.5 mm	63 mm	1000 g	95 mm	100-174-783
5000	5 kN	Rubber	10 mm	80 mm	3100 g	140 mm	100-186-413

Bionix Scissors Action Grips

Bionix Scissors Action Grips are self-tightening, and self-aligning grips featuring diamond faces. They are suitable for testing irregular bio-materials, bone, cartilage, tendons, and replacement bio-medical components at temperatures ranging from -130°C to 250°C. Comprised of durable stainless steel, they can be safely deployed within a saline environmental bath. These accessories feature a universal adapter design and an optional threaded frame adapter to facilitate easy installation onto both MTS electromechanical and servohydraulic load frame systems, as well as other electromechanical test systems. Anti-rotation features and integrated alignment pins enhance test accuracy and repeatability.

FEATURES

- » Stainless steel – suitable for saline environmental bath
- » Self-tightening and self-aligning
- » Diamond grip faces
- » Universal adapter design (Type D standard)
- » Optional threaded frame adapter
- » Anti-rotation features and integrated alignment pins
- » Temperature range: -130°C to 250°C
- » Applications: Irregular bio-material, bone, cartilage, tendons, replacement bio-medical components.



Model 2000

Bionix Scissors Action Grips

MODEL	MAXIMUM LOAD	MAXIMUM SPECIMEN THICKNESS	MAXIMUM SPECIMEN WIDTH	WEIGHT (UPPER)	HEIGHT (EACH CLEVIS CENTER)	PART NUMBER
2000	1 kN	10 mm	25 mm	850 g	160 mm	100-181-625

Bionix Stainless Steel Compression Platens

Bionix Stainless Steel Compression Platens are compression platens with precision-ground, hardened surfaces. Available in a variety of sizes, they are suitable for testing bio-materials, bone, cartilage, tendons, and replacement bio-medical components at temperatures ranging from -130°C to 250°C. Comprised of durable stainless steel, they can be safely deployed within a saline environmental bath. These accessories feature a universal adapter design and an optional threaded frame adapter to facilitate easy installation onto both MTS electromechanical and servohydraulic load frame systems, as well as other electromechanical test systems. Specimen centering grooves, anti-rotation features, and integrated alignment pins enhance test accuracy and repeatability.

FEATURES

- » Stainless steel – suitable for saline environmental bath
- » Precision-ground, hardened surface
- » Universal adapter design (Type D standard)
- » Optional threaded frame adapter
- » Specimen centering grooves
- » Anti-rotation features and integrated alignment pins
- » Temperature range: -130°C to 250°C
- » Applications: Irregular bio-material, bone, cartilage, tendons, replacement bio-medical components.



Bionix Stainless Steel Compression Platens

MODEL	MAXIMUM LOAD	MAXIMUM SPECIMEN DIAMETER	WEIGHT (UPPER)	HEIGHT (FROM CLEVIS CENTER)	PART NUMBER
50	10 kN	50 mm	681 g	64 mm	100-182-229
100	10 kN	100 mm	1820 g	64 mm	100-182-227
150	10 kN	150 mm	3725 g	64 mm	100-182-226

Bionix Roller Action Grips

Versatile Bionix Roller Action Grips are quick-loading and self-tightening roller grips suitable for testing bandages, bio-textiles, diapers, synthetics, and flexible polymers at temperatures ranging from -10°C to 50°C for the 1 kN grips and from -130°C to 250°C for the 2kN and 5 kN grips. Available in a variety of load capacities and sizes, these accessories feature a universal adapter design and an optional threaded frame adapter to facilitate easy installation onto both MTS electromechanical and servohydraulic load frame systems, as well as other electromechanical test systems. Anti-rotation features and integrated alignment pins enhance test accuracy and repeatability.

FEATURES

- » Quick-loading, self-tightening
- » Universal adapter design (Type D standard)
- » Optional threaded frame adapter
- » Anti-rotation features and integrated alignment pins
- » Temperature range: -10°C to 50°C for 1 kN; -130°C to 250°C for 2 kN and 5 kN
- » Applications: Bandages, bio-textiles, diapers, synthetics, flexible polymers.



Model 1000



Model 2000/5000

Bionix Roller Action Grips


MODEL	MAXIMUM LOAD	MAXIMUM SPECIMEN THICKNESS	MAXIMUM SPECIMEN WIDTH	WEIGHT (UPPER)	HEIGHT (FROM CLEVIS CENTER)	TEMPERATURE RANGE	PART NUMBER
1000	1 kN	3 mm	50 mm	1300 g	123 mm	-10°C to 50°C	100-185-262
2000	2 kN	4 mm	50 mm	1425 g	111 mm	-130°C to 250°C	100-184-839
5000	5 kN	7 mm	75 mm	2260 g	121 mm	-130°C to 250°C	100-184-841

Bionix Manual Bollard Style Grips

Bionix Manual Bollard Grips are “horn” style grips designed to reduce stress concentration on specimens and avoid grip-induced failures. They are suitable for testing cords, filaments, fibers, fine wire, and yarns at temperatures ranging from -130°C to 250°C. These accessories feature a universal adapter design and an optional threaded frame adapter to facilitate easy installation onto both MTS electromechanical and servohydraulic load frame systems, as well as other electromechanical test systems. Anti-rotation features and integrated alignment pins enhance test accuracy and repeatability.

FEATURES

- » “Horn” style design
- » Universal adapter design (Type D standard)
- » Optional threaded frame adapter
- » Anti-rotation features and integrated alignment pins
- » Temperature range: -130°C to 250°C
- » Applications: Cords, filaments, fibers, fine wire, yarns.

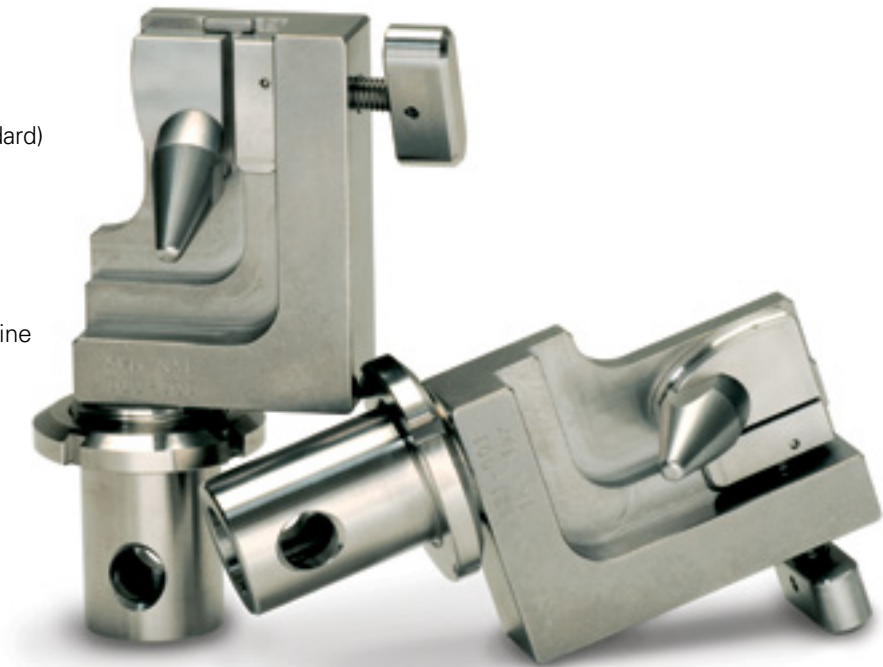

Model 1000

Bionix Manual Bollard Style Grips

MODEL	MAXIMUM LOAD	MAXIMUM SPECIMEN DIAMETER	WEIGHT (UPPER)	HEIGHT (EACH FROM CLEVIS CENTER)	PART NUMBER
1000	1 kN	3 mm	1010 g	124 mm	100-185-264


Bionix Spring Action Grips

Bionix Spring Action Grips are a reliable choice for low-force biomedical testing applications. Available with a variety of mounting threads and jaw surfaces, these lightweight grips are suitable for testing films, papers, and textiles at temperatures ranging from -75°C to 200°C (-103°F to 400°F). Comprised of durable stainless steel, they can be safely deployed within a saline environmental bath. All Bionix Spring Action Grip sets include two grip assemblies, two spare springs, two spare wires for link pin retainers, rubber jaw pad material, and a shipping case.



FEATURES


- » Stainless steel – suitable for saline environmental bath
- » Available with a variety of mounting threads and jaw surfaces
- » Spring action follow-through accommodates specimen neckdown
- » Faces pivot for self alignment and reduced likelihood of breakage at the specimen/face contact
- » Temperature range: -75°C to 200°C (-103°F to 400°F)
- » Applications: Films, papers, textiles
- » Force of jaw tips (jaws parallel): 32 N $\pm$ 4.4 N (7.2 lbf $\pm$ 1.1 lbf).

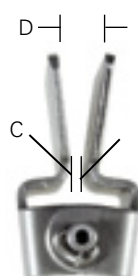

Model 10A,
fully open




Model 10A,
jaws parallel

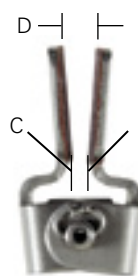

Model 10B,
(with rubber
face), fully open



Model 10B, (with
rubber face), jaws
parallel

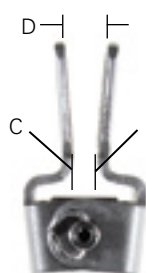

Model 10B,
fully open,
(with no rubber
face material)



Model 10B,
(with no rubber
face material),
jaws parallel

Bionix Spring Action Grips

MODEL	MOUNTING THREAD	SPAN "B"	SPAN "C"	SPAN "D"	JAW SURFACE	PART NUMBER
10A	M6 x 1 mm	0.0	1.0 mm (0.039 in)	2.6 mm (0.103 in)	Flat Smooth	056-077-501
10A	1/4 in -28UNF	0.0	1.0 mm (0.039 in)	2.6 mm (0.103 in)	Flat Smooth	056-077-502
10B	M6 x 1 mm	0.0	1.0 mm (0.039 in)	2.6 mm (0.103 in)	Rubber Face	056-077-503
10B	1/4 in -28UNF	0.0	1.0 mm (0.039 in)	2.6 mm (0.103 in)	Rubber Face	056-077-504
10B	M6 x 1 mm	1.6 mm (0.064 in)	2.6 mm (0.103 in)	4.2 mm (0.167 in)	Flat Smooth	056-077-503
10B	1/4 in -28UNF	1.6 mm (0.064 in)	2.6 mm (0.103 in)	4.2 mm (0.167 in)	Flat Smooth	056-077-504

Bionix EnviroBath

The versatile Bionix EnviroBath facilitates accurate and efficient mechanical testing of medical device and biomaterial specimens in fluids heated to body temperatures. Engineered for biomedical and general material test applications, the Bionix EnviroBath is easy to set up, operate, and maintain and is available in a range of volume configurations to accommodate a broad variety of test specimens. A universal adapter design ensures full compatibility with all MTS electromechanical, electrodynamic and servohydraulic load frame systems, as well as other electromechanical test systems.

- » Variety of volume configurations: 1 l, 6 l, 10 l
- » Supports saline and protein-based fluids
- » Utilizes a wide selection of Bionix grips and fixtures
- » Reliable fluid temperature control mechanism
- » Compatible with video or laser extensometer
- » Optional protein-based fluid system.



Easy to disassemble to clean.



General Specifications

TEMPERATURE

Range:	5°C above ambient to 40°C
Stability:	±2°C at 37°C

POWER (V AC, HZ, A)

US:	120 V AC, 50/60 Hz, 11 A
Europe:	240 V AC, 50 Hz, 10 A

Bionix EnviroBath

	VOLUME	AXIAL FORCE	TORSIONAL FORCE	INTERNAL DIMENSIONS WXHxD	EXTERNAL DIMENSIONS WXHxD
1	1 l	2.4 kN	28 N.m	100 x 200 x 56 mm 4 x 8 x 2.2 in	180 x 295 x 95 mm 7 x 11.5 x 3.65 in
6	6 l	2.4 kN	28 N.m	130 x 480 x 100 mm 5 x 19 x 3.9 in	205 x 575 x 140 mm 8 x 22.5 x 5.5 in
10	10 l	10 kN	100 N.m	215 x 305 x 150 mm 8.5 x 12 x 5.8 in	295 x 395 x 190 mm 11.5 x 15.5 x 7.4 in

Bionix EnviroBath

Options

	1 LITER	6 LITER	10 LITER
Sprayer Option		✓	✓
Protein Based Fluid	✓	✓	✓
Horizontal	✓	✓	
Digital Temperature Monitor	✓	✓	✓
Upper Pull Rod (SST)	✓	✓	✓

Compatibility Matrix*

GRIPS AND FIXTURES	1 LITER	6 LITER	10 LITER	PART NUMBER
Bionix Vise Action Grips, 100 N	✓	✓	✓	100-186-411
Model 642.001 Bend Fixture (SST)	✓	✓	✓	100-201-456 (S/H & E/D)
Bionix Manual Thumb Screw Grip, 100 N	✓	✓	✓	100-167-987
Bionix Compression Platens, 40 mm	✓	✓	✓	100-203-455
Bionix Spring Action Grips		✓	✓	056-077-50x
Bionix Compression Platens, 50 mm		✓	✓	100-182-229
Bionix Vise Action Grips, 2 kN		✓	✓	100-174-783
Bionix Roller Action Grips, 1 kN		✓	✓	100-185-262
Model 642.01 Bend Fixture (SST)			✓	100-203-453 (S/H & E/D)
Bionix Compression Platens, 100 mm			✓	100-182-227
Model 686.01A-31 Axial Torsional Grips (SST)			✓	051-882-201

Note: EnviroBath 10 is compatible with customer-supplied spinal fixture per ASTM F1717-01.

Some grips available in titanium. Contact MTS for additional information.

* Vertical Orientation

Electromechanical Load Frames

	1 LITER	6 LITER	10 LITER
MTS Single Column Standard Length	✓		✓
MTS Single Column Extended Length	✓	✓	✓
MTS Dual Column Standard and Extended Length	✓	✓	✓
Legacy MTS load frames and Instron E/M load frames			Consult MTS

Electrodynamic Load Frames

	1 LITER	6 LITER	10 LITER
MTS Acumen 1	✓		✓
MTS Acumen 3	✓		✓


Servohydraulic Load Frames

	1 LITER	6 LITER	10 LITER
Bionix Tabletop Axial or Axial/Torsional	✓ ¹	✓ ^{1, 2}	✓
MTS Landmark ³	✓	✓	✓
Legacy MTS load frames ³			Consult MTS

¹ Additional adapter required.

² Extended columns required.

³ Cross-head mounted actuator.

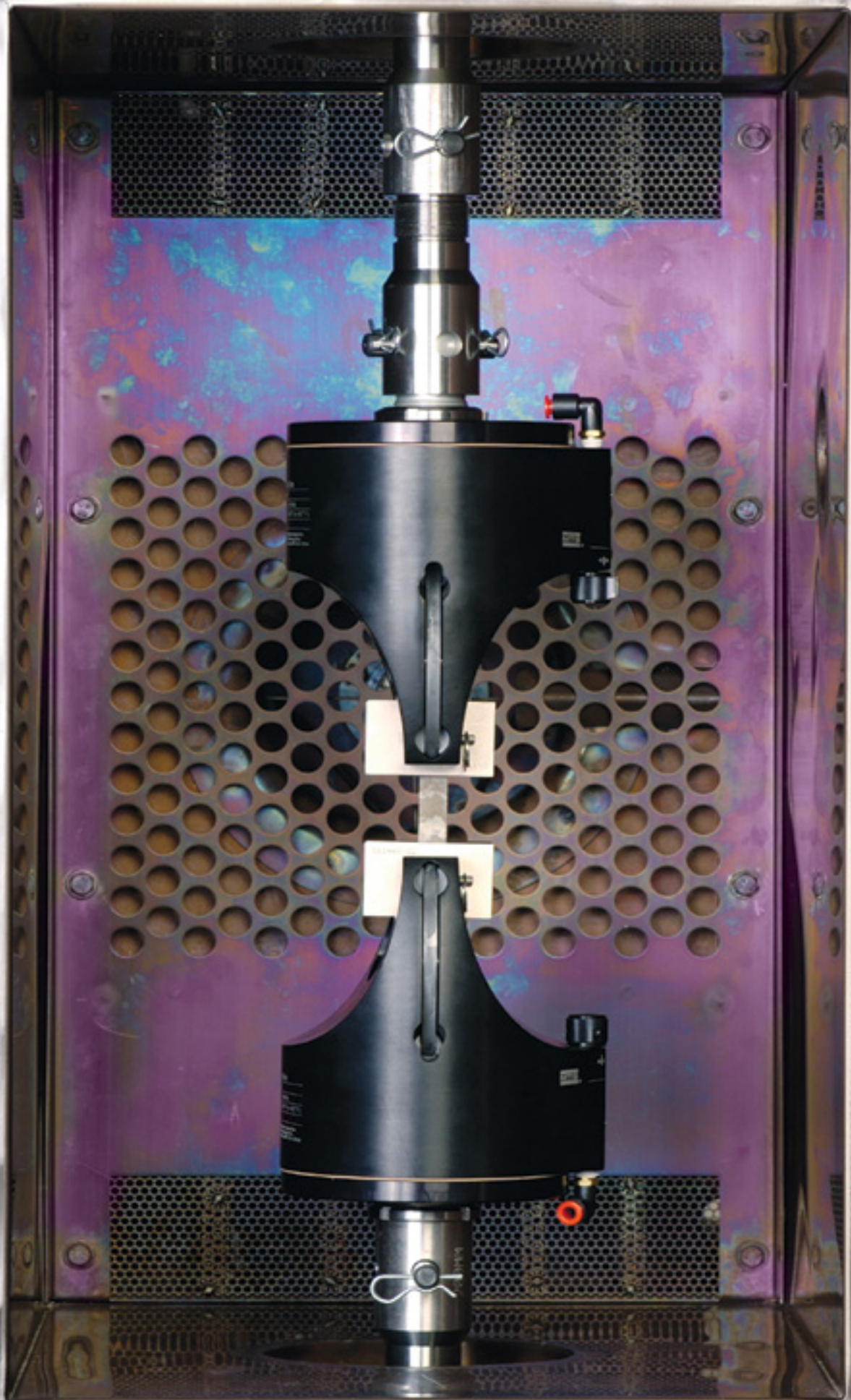


MTS EM Extend Kits

MTS load frame extension rod kits enhance your testing capabilities within environmental chambers. The rods are designed to apply tension and compression in an environmental chamber from temperatures of -130°C (-200°F) to 315°C (600°F). The Type C and Type D kits consist of five different extension lengths. These five different lengths can be configured into 23 different combination lengths. This configurability allows users to test low and high profile accessories or short and tall specimens. Compression platens are usually low profile and will need more extension lengths. Pneumatic grips take up more space and will require fewer lengths.

A load frame rod extension kit maximizes the internal test space of the environmental chamber for grips and fixtures. The upper extension rod can be adjusted in 25 mm increments. By combining the five different lengths, 23 combination heights are ultimately available providing pin-to-pin dimensions ranging from 100 mm (3.9 in) to 825 mm (32.4 in).

The Type D extension kit is compatible with MTS and Instron® clevis pin couplings. It is designed for tension or compression loads from 2 kN to 150 kN. All extension lengths are stainless steel and have ports for water cooling.

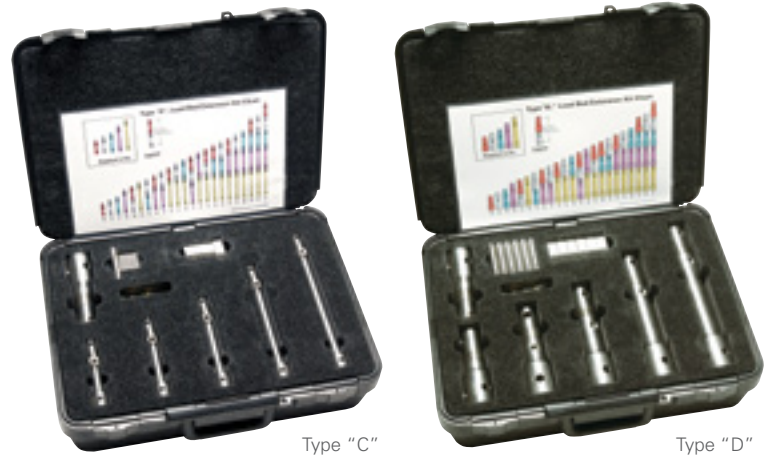

Load extension kits are packaged in a protective case. From left to right, the lengths are denoted as L1 through L5.

FEATURES

- » Water-cooling ports
- » Lightweight design minimizes tare weight
- » Tight tolerances ensure angular and concentric alignment are retained
- » Rated for tension and compression
- » Stainless steel construction minimizes heat transfer
- » Preloaded joints.

BENEFITS


- » Highly configurable to test various size specimens
- » Reduces test setup time.

MTS EM Extend Kit

MODEL	TYPE C	TYPE D
Part Number	100-151-425	100-150-815
Force Rating	200 N (40 lb)	150 kN (33,000 lb)
Temperature Rating	-130°C to 315°C (-200°F to 600°F)	-130°C to 315°C (-200°F to 600°F)
Length 1 (L1)*	100 mm (3.9 in)	100 mm (3.9 in)
Length 2 (L2)*	125 mm (4.9 in)	125 mm (4.9 in)
Length 3 (L3)*	150 mm (5.9 in)	150 mm (5.9 in)
Length 4 (L4)*	200 mm (7.8 in)	200 mm (7.8 in)
Length 5 (L5)*	250 mm (9.8 in)	250 mm (9.8 in)
Weight 1	0.18 kg (0.40 lb)	0.93 kg (2.06 lb)
Weight 2	0.20 kg (0.45 lb)	1.1 kg (2.47 lb)
Weight 3	0.23 kg (0.50 lb)	1.25 kg (2.77 lb)
Weight 4	0.29 kg (0.65 lb)	1.57 kg (3.47 lb)
Weight 5	0.38 kg (0.85 lb)	1.93 kg (4.27 lb)

FURNACES & CHAMBERS

Our high-temperature testing solutions leverage decades of collaboration with industry pioneers and innovative technologies proven to deliver accurate repeatable results. These solutions precisely replicate real-world loading of high-temperatures and other environments.

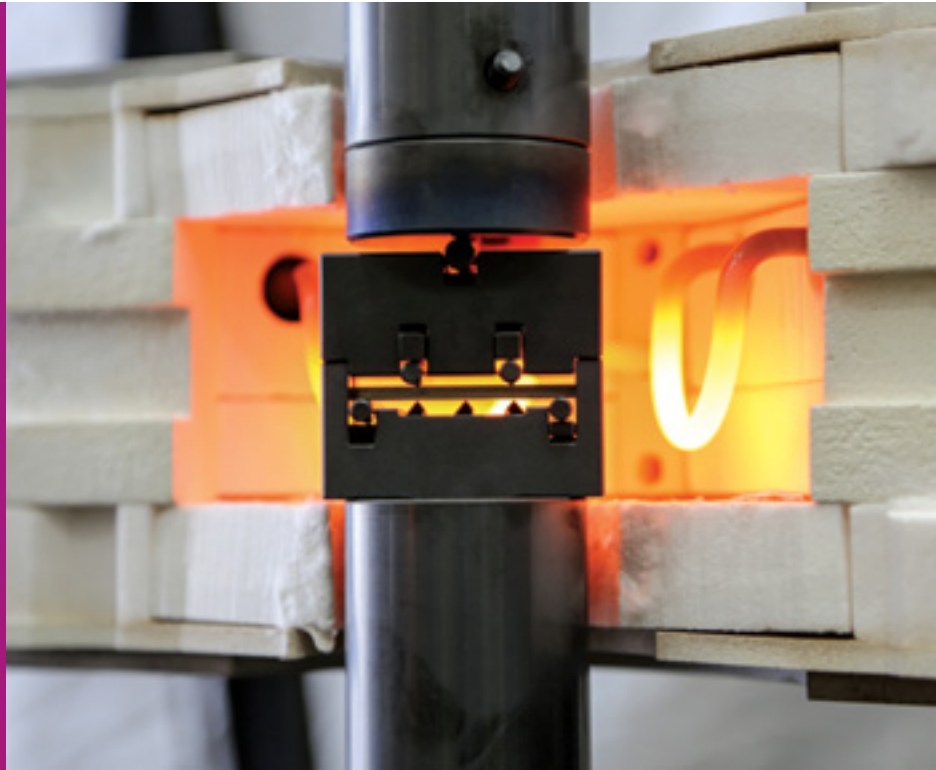

Furnaces & Chambers

Introduction	122
Furnaces & Temperature Controller	123
Thermomechanical Fatigue (TMF)	124-125
Environmental Chamber with Built-In Controller	126-127

Reliable Furnace and Temperature Controller

Combine Performance and Affordability

MTS Model 653 High-Temperature Furnaces are ideal for a wide variety of high-temperature tests, including tension, compression, bend and fatigue testing of metals, composites, ceramics and many other materials. The furnaces have a user-friendly design that helps operators perform setup tasks more quickly – without compromising test quality or accuracy.

Model 653 furnaces streamline test setup and specimen changeover with a clamshell design that allows operators to work from the front of the furnace at all times. This design also improves furnace alignment, which is critical to achieving proper gradients. The durable, low-maintenance insulation for Model 653.03 and 653.04 furnaces comes pre-cut, so there is no need to handle the material in the lab. It fits tighter around the pushrods, reducing heat loss and eliminating the need for wool. The flow switches of the grip water cooling circuit can be integrated with the Model 409 temperature controller, shutting down the furnace to protect the equipment and specimen in the event of a water supply failure.



MODEL 653.XX FURNACE FEATURES:

- » Capable of achieving temperatures up to 1400°C (2550°F) (non-testing environment)
- » Capable of achieving 1000°C (1800°F) in validated testing conditions with standard solution
- » Single or multiple zone heating
- » Clamshell design streamlines test setup, furnace alignment, and specimen changeover
- » Silicon carbide heating elements and alumina fiber insulation system for low heat loss and long life
- » Multiple furnace heights to accommodate diverse test requirements
- » Mounting bracket for a variety of MTS load frames is included
- » Designed to accommodate MTS high-temperature axial extensometers

MODEL 409.83 TEMPERATURE CONTROLLER FEATURES:

- » Multiple mounting options – on included furnace mounting bracket or optional stand for floor or table
- » Compact, ergonomic design
- » Multiple level, self-tuning PID control
- » SCR power relays included
- » Digital communications available.

Furnaces

MODEL	TEMPERATURE MAX/MIN*	OVERALL HEIGHT	HOT ZONE HEIGHT	HOT ZONE WIDTH & DEPTH	NUMBER OF ZONES
653.01	1400°C/100°C	55 mm	19 mm	50 x 50 mm	1
653.02	1400°C/100°C	85 mm	50 mm	50 x 50 mm	2
653.03	1400°C/100°C	126 mm	90 mm	62.5 x 62.5 mm	2
653.04**	1400°C/100°C	220 mm	185 mm	62.5 x 62.5 mm	3

* Nominal temperatures may vary depending on specimen geometry and material. 1400°C is achieved in a non-testing environment.

** Supports testing to ASTM E606-04e1, BSI 7270, JIS Z2279, AFNOR A03-403 or ISO 12106 requirements.

Note: When ordering, please indicate voltage requirements and provide necessary load frame dimensions in order to determine system integration requirements.

THERMAL GRADIENT VERIFICATION KIT COMPONENTS (OPTIONAL)	PART NUMBER
Thermocouple Data Acquisition Kit	100-255-716
Type K Thermocouple Connector Kit	100-255-714
Thermocouple Specimen	057-743-803
Thermocouple Welding Kit	100-256-698



Thermomechanical Fatigue (TMF)

Subsystem with Induction Heating

The MTS TMF Subsystem is ideal for inducing thermal cycling in a metallic (magnetically susceptible) specimen. A combination of induction heating and air cooling allows thermomechanical fatigue tests to be performed in or out of phase with a load force. The TMF subsystem has a user-friendly design that helps operators perform repeatable test setup tasks quickly. Air cooling jets and the induction coil have easy, repeatable adjustments, meaning that specimen changes do not require as much setup time to reestablish the initial conditions. A specimen centering tool provides quick initial position for the coil. The coil is mounted on an x, y, z positioner that has manual adjustment knobs with numerical indicators allowing 50 mm of travel in 3-axes for fine positioning relative to the specimen. A quick release slide allows the coil to be easily pulled from the test space with the ability to return to its original position. Air cooling jets are also easily folded back or removed to gain access to the specimen with the same ability to snap back into the original position.

The induction coil is coated in electrically insulating material to provide an additional safety factor against electric shock. The coil, generator and RF oscillator are water cooled and interlocked for water flow. A second interlock specimen thermocouple also provides indication in the event that the primary control thermocouple becomes detached. TestSuite software templates for TMF include test setup, to establish thermal gradients, and test run to perform and report results on the TMF test.

TMF SUBSYSTEM FEATURES:

- » 10 kW RF Generator with temperature controller and secondary temperature interlock controller
- » Air-cooling jets with detachable arms and adjustable flow profiles
- » 10°C/second heating and cooling with higher open loop heating rates possible
- » Up to 1200°C
- » X,Y,Z coil with position indication and quick release slide-back mount
- » Type K or Type R thermocouple compatibility
- » Water chiller for RF generator and coil
- » Tool kit for coil centering and adjustment
- » Standard 5-turn electrically insulated induction coil
- » TestSuite Templates for setup and test control
- » Instrumented specimen for equipment checkout and validation
- » 16 channel thermocouple data acquisition module
- » Thermocouple connector kit to interface with specimen
- » Designed to accommodate MTS high-temperature axial extensometers.

When ordering, please indicate power choice of water chiller 50 Hz or 60 Hz, thermocouple type (R or K).



Environmental Chambers

Engineered for Dependable Operation

Every aspect of Series 651 Environmental Chambers was designed using decades of MTS experience developing test systems for numerous industries. These environmental chambers may include:

- » Circulating fan helps to ensure small specimen temperature gradients
- » Removable “U-plug” sections allowing chambers to be moved into place after the specimen is mounted and instrumented
- » Doors that open to 180° or can be lifted away for convenient specimen access
- » Multi-panel, tempered glass windows are sealed to reduce fogging and frosting
- » Internal light to illuminate the test area
- » Mobile carts to prevent interference between chambers and the load frame

Model 651.05F-01 works in a variety of general materials testing applications up to 350°C (662°F) with floor-standing models of the MTS Acumen test system.

Model 651.05F-02 is ideal for general materials testing applications up to 350°C (662°F) with tabletop models of the MTS Acumen test system.

Model 651.06E-03 handles general materials testing up to 315°C (600°F) with floor-standing MTS Landmark test systems.

Model 651.10E-04 offers support for general materials testing up to 540°C (1000°F) with floor-standing MTS Landmark test systems.

Model 651.06E-04 is designed to accommodate large specimens or specimens that require large grips within the chamber. Works up to 315°C (600°F) with floor-standing MTS Landmark test systems.



Can't Find What You Need?

MTS also supplies chambers with mechanical refrigeration, humidity control, or salt fog testing. Contact your local applications engineer for more information.

Specifications

Environmental Chamber Specifications

MODEL NUMBER	651.05F-01	651.05F-02	651.06E-03	651.06E-04	651.10E-04
Load frame	Acumen 1 Acumen 3 Acumen 12	Landmark Bionix Acumen 3 Acumen 12	Landmark	Landmark	Landmark
Maximum temperature***	350°C (662°F)	350°C (662°F)	315°C (600°F)	340°C (1000°F)	315°C (600°F)
Chamber part number	058-052-704 (US Plug) 058-052-703 (WW Plug)	058-052-702 (Acumen, US Plug) 058-052-701 (Acumen, WW Plug) 058-052-602 (w/table top legs, US Plug) 058-052-601 (w/table top legs, US Plug)	056-097-206 (works with FlexTest SE Controller) 056-097-205 (works with FlexTest GT Controller) 056-097-207 (works with 494 FlexTest Controller)	056-097-005 (works with FlexTest SE Controller) 056-097-004 (works with FlexTest GT Controller) 056-097-008 (works with 494 FlexTest Controller)	056-097-105 (works with FlexTest SE Controller) 056-097-104 (works with FlexTest GT Controller) 056-097-106 (works with 494 FlexTest Controller)
Cart part number¶	100-294-319# 100-297-346	100-294-319# 100-297-346	100-160-646 100-213-156**	100-160-646 100-213-156**	100-160-646 100-213-156**
Test space (nominal)	216 mm W, 228 mm D, 228 mm H (8.5 in W, 9 in D, 9 in H)	286 mm W, 305 mm D, 453 mm H (11.25 in W, 12 in D, 17.8 in H)	356 mm W, 432 mm D, 559 mm H (14 in W, 17 in D, 22 in H)	356 mm W, 432 mm D, 813 mm H (14 in W, 17 in D, 32 in H)	356 mm W, 432 mm D, 660 mm H (14 in W, 17 in D, 26 in H)
Removable "U-Plug" sections	Not included	Included	Included	Included	Included
Typical application	Material testing, including fatigue and fracture, tension, compression, flex/bend, etc.				
Air temperature performance*					
Heating	Ambient to +177°C (+350°F) in 45 minutes	Ambient to +177°C (+350°F) in 45 minutes	Ambient to +315°C (+600°F) in 30 minutes	Ambient to +315°C (+600°F) in 30 minutes	Ambient to +540°C (+1000°F) in 70 minutes
LN2 cooling	Ambient to -101°C (-150°F) in 45 minutes	Ambient to -101°C (-150°F) in 45 minutes	Ambient to -129°C (-200°F) in 26 minutes	Ambient to -129°C (-200°F) in 26 minutes	Ambient to -129°C (-200°F) in 26 minutes
Temperature gradient †, ‡	±3°C or less	±3°C or less	±2°C (±5°F)	±2°C (±5°F)	±2°C (±5°F) from -129° to 260°C (-200° to 500°F) ±5°C (±10°F) from -260° to 540°C (-500° to 1000°F)
Temperature stability†	±2°C (±5°F)	±2°C (±5°F)	±1°C (±3°F)	±1°C (±3°F)	±1°C (±3°F)
Cooling requirements					
LN2 pressure	0.15 MPa (22 psi nominal) 0.21 MPa (30 psi max)	0.15 MPa (22 psi nominal) 0.21 MPa (30 psi max)	0.15 MPa (22 psi nominal) 0.34 MPa (50 psi max)	0.15 MPa (22 psi nominal) 0.34 MPa (50 psi max)	0.15 MPa (22 psi nominal) 0.34 MPa (50 psi max)
Recovery time to set point temperature					
Door open 1 minute			5 minutes	4 minutes	8 minutes
Door open 2 minute			6 minutes	5 minutes	10 minutes
Heater circuit power	208/230 V AC	208/230 V AC	208/230 V AC	208/230 V AC	208/230 V AC
Requirements	50/60 Hz, single-phase; 20 A circuit	50/60 Hz, single-phase; 20 A circuit	50/60 Hz, single-phase;	50/60 Hz, single-phase;	50/60 Hz, single-phase;
Weight	82 kg (180 lbs)	108 kg (250 lbs)	191 kg (4200 lbs)	200 kg (440 lbs)	195 kg (430 lbs)
Internal light	None	Internal 120 V bulb (40 watt recommended)	Internal 110 V bulb (220 V optional) (40 watt recommended)	Internal 110 V bulb (40 watt recommended)	Internal 110 V bulb (40 watt recommended)

* Performance data derived with chamber empty and access holes blocked. Consult MTS Systems Corporation for temperature performance for testing specimens having high thermal mass.

** For Landmark 370.50 only.

*** Typical maximum temperature for general materials testing applications. Many factors, including thermal mass of grips and fixtures, impact maximum operating temperature.

¶ Consult MTS for cart information specific to your configuration.


Cart stand for table top load units.

† After 30 minute stabilization time.

‡ Due to large variety of possible test setups an air gradient is given. Specimen gradient is much less.

TEST SYSTEM COMPONENTS

MTS test system components are designed to maximize the safety, reliability and energy-efficiency of hydraulic power generation, actuation and distribution. These components are specifically engineered to meet the unique and rigorous demands of mechanical testing, and since decisions related to hydraulic power and distribution have a profound impact on lab productivity, they play an increasingly important role in the lab.

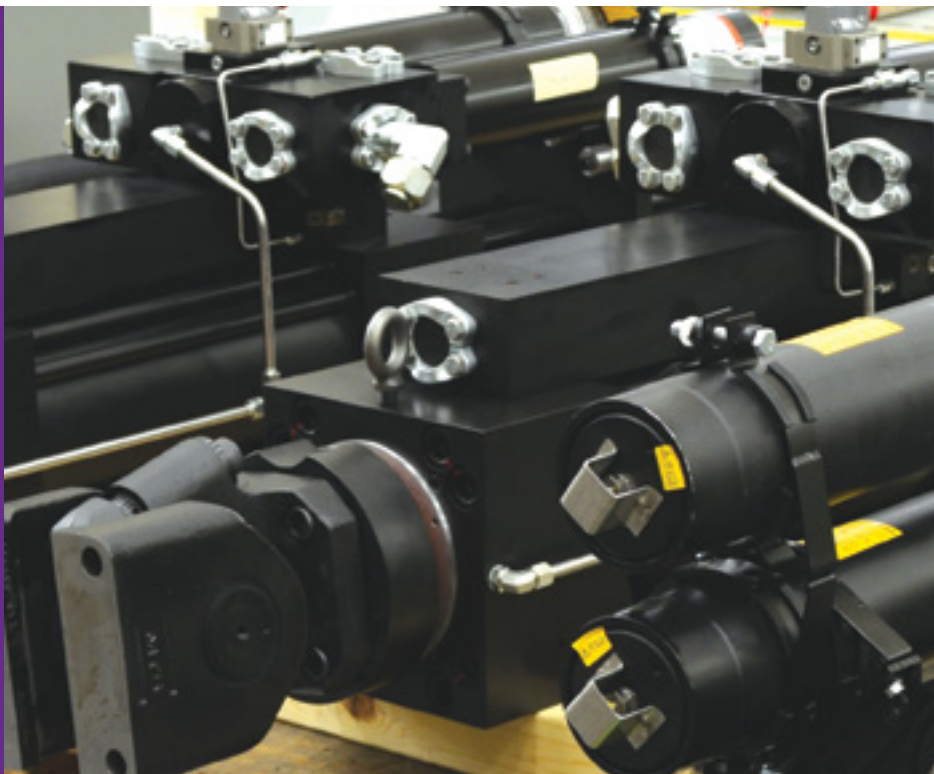


Testing System Components

Introduction	128
201 Hydraulic Actuators	129-130
242 Hydraulic Actuators	131
244 Hydraulic Actuators	132
248 Hydraulic Actuators	133
249 Swivel Base & Rod Ends	134-135
252 Servovalves	136
215 & 216 Rotary Actuators	137
Rotary Actuator Options	138
515 Hydraulic Power Units	139
Multi-Pump Control Manager (MPCM)	140
293 Hydraulic Service Manifold	141
Isolation Service Manifold	142
Testline Modular Components	143
Additional Components	144

MTS Series 201 Hydraulic Actuators

MTS Series 201 Actuators are single-ended, double-acting, fatigue-rated force generators targeted at long stroke and mid-level performance dynamic structural applications. This actuator is ideal for low frequency testing and simulation applications.

- » Proprietary seal and bearing designs that reduce friction
- » Available with or without displacement transducer to meet specific application needs
- » Low-friction seals comprise only 1% of rated load for energy-efficient operation
- » Optional hydraulic cushions to protect the actuator from the effects of high-speed or high-inertial impacts
- » Optional MTS SureCoat™ Rod Finishing Technology to maximize performance and increase actuator life
- » Ready integration with a wide variety of MTS servohydraulic components

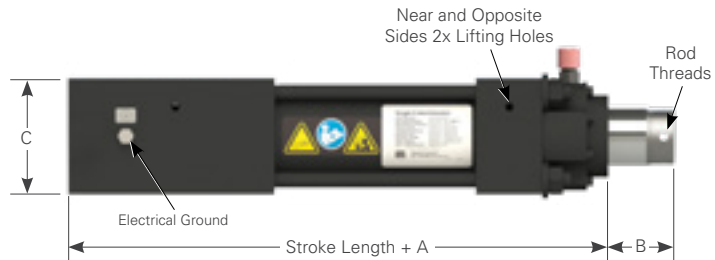

Specifications

201 Cylinder Dimensions

MODEL	Stroke Options*				Transducer Ready		Solid Rod		ROD EXTENSION FULLY RETRACTED (B)		REAR END CAP (C)		REAR END CAP MOUNTING PATTERN	
	MINIMUM		MAXIMUM		CYLINDER BODY LENGTH (A)		CYLINDER BODY LENGTH (A)						Primary	Secondary***
	mm	in	mm	in	mm	in	mm	in	mm	in	mm	in		
201.10	100	4	1000	40**	304.8	12.00	228.6	9.00	25.4	1.00	76.2	3.00	Type 1, 249XB.M25	
201.11	100	4	1500	60**	314.2	12.37	238.0	9.37	50.8	2.00	82.6	3.25	Type 1, 249XB.M25	
201.17	100	4	1500	60**	314.2	12.37	238.0	9.37	50.8	2.00	88.9	3.50	Type 1, 249XB.M25	
201.20	100	4	1500	60**	314.2	12.37	238.0	9.37	50.8	2.00	114.3	4.50	Type 1, 249XB.M70	
201.25	100	4	1500	60**	263.4	10.37	238.0	9.37	50.8	2.00	127.0	5.00	Type 1, 249XB.M160	
201.30	100	4	1500	60	286.3	11.27	286.3	11.27	50.8	2.00	165.1	6.50	Type 1, 249XB.M160	
201.35	100	4	1500	60	305.6	12.03	286.3	11.27	50.8	2.00	190.5	7.50	Type 1, 249XB.M340	Type 1, 249XB.M160
201.40	100	4	1500	60	330.2	13.00	–	–	50.8	2.00	225.0	8.86	Type 1, 249XB.M340	
201.45	100	4	1500	60	359.4	14.15	–	–	50.8	2.00	241.3	9.50	Type 2, 249XB.M500	Type 1, 249XB.M340
201.50	100	4	1500	60	391.2	15.40	–	–	50.8	2.00	292.1	11.50	Type 2, 249XB.M500	
201.60	100	4	1500	60	410.2	16.15	–	–	50.8	2.00	355.6	14.00	Type 2, 249XB.M1000	
201.70	100	4	1500	60	475.5	18.72	–	–	50.8	2.00	406.4	16.00	Type 2, 249XB.M1000	
201.80	100	4	1500	60	499.1	19.65	–	–	50.8	2.00	457.2	18.00	Type 2, 249XB.M1750	Type 2, 249XB.M1000
201.90	100	4	1500	60	537.7	21.17	–	–	50.8	2.00	577.9	22.75	Type 2, 249XB.M1750	

*Standard stroke available in increments of 50 mm (2 in). Custom orders can be built up to 2550 mm (100 in) on select models.

**Reduced compression force rating to prevent buckling may exist on long stroke cylinders, contact MTS for more information.

***For static applications, contact MTS for more information.

MTS Series 201 Hydraulic Actuators *(continued)*

Specifications

201 Actuator Specifications


			210 BAR (3000 PSI)						280 BAR (4000 PSI)							
MODEL	ROD DIAMETER		BORE		TENSION		COMPRESSION		BORE		TENSION		COMPRESSION		ROD THREADS	
	mm	in	mm	in	kN	kip	kN	kip	mm	in	kN	kip	kN	kip	Metric	US Customary
201.10	25.4	1.00	34.9	1.38	9.0	2	19	4	31.8	1.25	7.7	2	21	5	M12x1.25 ± 25.4	1/2"-20 ± 1.00
201.11	34.9	1.38	47.6	1.88	16	4	36	8	44.5	1.75	16	4	42	9	M12x1.25 ± 25.4	1/2"-20 ± 1.00
201.17	44.5	1.75	63.5	2.50	32	7	63	14	57.2	2.25	27	6	69	16	M27x2 ± 44.5	1"-14 ± 1.75
201.20	50.8	2.00	82.6	3.25	66	15	107	24	76.2	3.00	68	15	123	28	M27x2 ± 44.5	1"-14 ± 1.75
201.25	63.5	2.50	101.6	4.00	99	22	162	36	95.3	3.75	106	24	192	43	M27x2 ± 44.5	1"-14 ± 1.75
201.30	76.2	3.00	127.0	5.00	162	36	253	57	114.3	4.50	153	34	276	62	M36x2 ± 57.2	1 1/2"-12 ± 2.25
201.35	88.9	3.50	152.4	6.00	241	54	365	82	139.7	5.50	245	55	412	93	M36x2 ± 57.2	1 1/2"-12 ± 2.25
201.40	101.6	4.00	177.8	7.00	334	75	496	112	165.1	6.50	358	80	576	129	M52x2 ± 76.2	2"-12 ± 3.00
201.45	114.3	4.50	203.2	8.00	443	100	648	146	184.2	7.25	440	99	716	161	M52x2 ± 76.2	2"-12 ± 3.00
201.50	127.0	5.00	222.3	8.75	522	117	776	174	203.2	8.00	531	119	872	196	M52x2 ± 76.2	2"-12 ± 3.00
201.60	152.4	6.00	254.0	10.00	648	146	1013	228	228.6	9.00	613	138	1104	248	M76x2 ± 114	3"-12 ± 4.50
201.70	177.8	7.00	304.8	12.00	962	200	1459	328	279.4	11.00	981	221	1649	371	M90x2 ± 133	3 1/2"-12 ± 5.25
201.80	203.2	8.00	355.6	14.00	1337	301	1986	446	323.9	12.75	1343	302	2215	498	M90x2 ± 133	3 1/2"-12 ± 5.25
201.90	228.6	9.00	406.4	16.00	1773	399	2594	583	368.3	14.50	1761	396	2865	644	M125x4 ± 171	4 1/2"-8 ± 6.75



MTS Series 242 Hydraulic Actuators

MTS Series 242 Actuators are used in closed loop control to extend or retract and provide displacement or force into a test article. These fatigue-rated actuators incorporate a coaxially mounted transducer for displacement feedback and internal hydraulic cushions for extended service life. Featuring MTS exclusive bearing technology and low pressure seal system, low-friction Series 242 Actuators are ideal for many applications including high-frequency excitation of a test article. Combined with MTS force transducers, servovalves, and swivels, these actuators generate the precision movements and forces required in test and simulation applications.

- » Double-ended, fatigue-rated design combines balanced dynamic performance and robust side-load tolerance to maximize data fidelity
- » Full-stroke, coaxially-mounted displacement transducer
- » High-capacity, non-metallic bearings bonded directly to the end caps offer high side-load tolerance and resistance to failure from galling and seizure
- » Piston rods machined from a single piece of heat-treated alloy steel and surface coated to reduce the effects of rod banding
- » Cushions protect the actuator from the effects of high-speed and inertial forces, and prevent the actuator from contacting the end caps
- » Proprietary MTS seal and bearing designs and surface-coating processes set the industry standard for durability, longevity and performance.

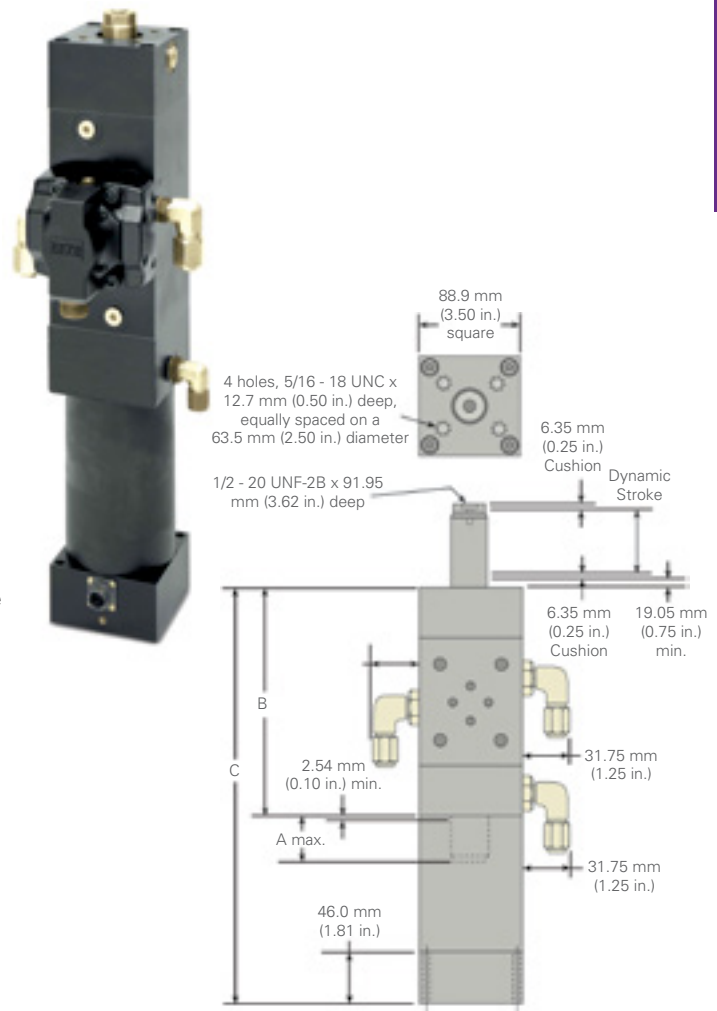

Series 242 Actuator Specifications

MODEL*	FORCE RATING [†]		PISTON AREA		ROD DIAMETER	
	kN	Kip	mm ²	in. ²	mm	in
242	2.7	0.6	137	0.21	28.6	1.12
242.01	4.5	1	270	0.42	28.6	1.12
242.02	9.8	2.2	590	0.91	28.6	1.12
242.03	14.7	3.3	760	1.18	28.6	1.12

* Each model has a standard stroke length of 101.6 or 152.4 mm (4 or 6 in. This value is total stroke minus the total cushion length, 12.7 mm (0.50 in.). Optional stroke lengths are available in 1.00 and 2.00 in. (25.4 and 50.8 mm).

† Nominal force with 17.2 MPa (2500 psi) pressure drop across the piston. Actual force may be up to 20% higher depending on servovalve size and test conditions.

Series 242 Actuator Dimensions*


STROKE (All Models)		A		B		C (Rod-Eye Bearing Swivel)		C (Closed Housing)	
mm	in.	mm	in.	mm	in.	mm	in.	mm	in.
25.4	1	41.4	1.63	215.9	8.5	322.2	12.69	308	12.13
50.8	2	66.8	2.63	215.9	8.5	347.6	13.69	333.4	13.13
101.6	4	117.6	4.63	279.4	11	461.9	18.19	447.7	17.63
152.4	6	168.4	6.63	330.2	13	563.5	22.19	549.3	21.63

* The dimensions listed are shown in drawing.

Dimensions are subject to change without notice. Contact MTS for verification of dimensions critical to your needs

MTS Series 244 Hydraulic Actuators

MTS Series 244 Actuators are double-ended, fatigue-rated and designed with proprietary materials that minimize friction, while maximizing reliability, wear resistance and ease of maintenance. They are hydraulically-powered pistons that can extend or retract (double-acting) to provide displacement or force into a test article or structure. MTS Series 244 Actuators are used worldwide in a variety of demanding applications from vehicle dynamics and structural fatigue to component testing, and are recommended for dynamic, static and fatigue applications.

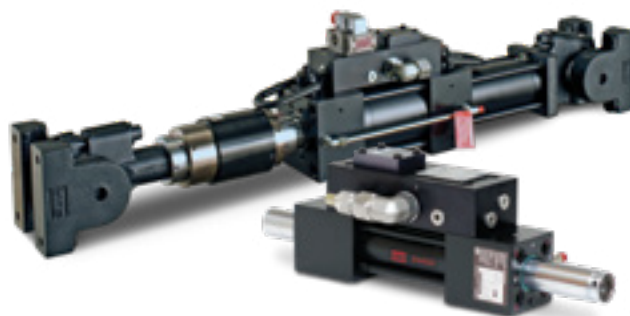

244 Actuator Specifications

MODEL	Force Rating		Rod Threads		HYDRAULIC CONNECTION
	kN	LBS	METRIC	US CUSTOMARY	
244.11	15	3,300	M12x1.25	1/2"-20	Pressure & Return – Pressure and return line connections depend on the servovalve. Configurations with Series 252 Servovalve Manifolds have SAE -16 fittings. Series 256 Servovalve Manifolds have SAE -20 fittings.
244.12	25	5,500	M12x1.25	1/2"-20	
244.21	50	11,000	M27x2	1"-14	
244.20	67	15,000	M27x2	1"-14	
244.22	100	22,000	M27x2	1"-14	
244.23	150	35,000	M27x2	1"-14	
244.31	250	55,000	M36x2	1 1/2"-12	Drain – All 244 Actuators have a SAE -6 drain connection.
244.41	500	110,000	M52x2	2"-12	
244.51	1000	220,000	M76x2	3"-12	

Specifications subject to change without notice. Please contact MTS for specifications critical to your application.

- » Thermal spray-coated piston rods that have 10x greater service life than hard chrome-plated rods
- » Proprietary MTS seal and bearing designs and surface-coating processes set the industry standard for durability, longevity and performance
- » Double-ended, fatigue-rated design (piston rod extends from both sides) combines balanced dynamic performance and robust side-load tolerance to maximize data fidelity
- » Full-stroke, coaxially-mounted displacement transducer
- » High-capacity, non-metallic, ultra-low friction bearings bonded directly to the end caps offer high side-load tolerance and resistance to failure from galling and seizure
- » Cushions protect the actuator from the effects of high-speed and high-mass forces, and prevent the piston from contacting the end caps.


Typical Configurations for Component or Structural Testing Applications

MODEL	SERVOVALVE	LOAD CELL	SWIVEL	SPIRAL WASHERS
244.11	252.2X (qty 1 or 2)	661.19E/F-03	249xx.M25	601.11-20
244.12	252.2X (qty 1 or 2)	661.19E/F-04	249xx.M25	601.11-20
244.21	252.2X (qty 1 or 2)	661.20E/F-02	249xx.M70	601.11-11
244.20	252.2X (qty 1 or 2)	661.20E/F-03	249xx.M70	601.11-11
244.22	252.2X (qty 1 or 2)	661.20E/F-03	249xx.M160	601.11-11
244.23	252.2X (qty 1 or 2)	661.22C/D-01	249xx.M160	601.11-11
244.31	252.2X (qty 1 or 2), 256 (qty 1)	661.22C/D-01	249xx.M340	601.11-19
244.41	252.2X (qty 1 or 2), 256 (qty 1)	661.23E/F-01	249xx.M500	601.11-13
244.51	252.2X (qty 1 or 2), 256 (qty 1)	661.31E/F-01	249xx.M1000	601.11-15

Stroke Length

mm	inch
100	4
150	6
200	8
250	10
300	12
350	14
400	16
450	18
500	20

Each Series 244 Actuator is available in these standard stroke lengths

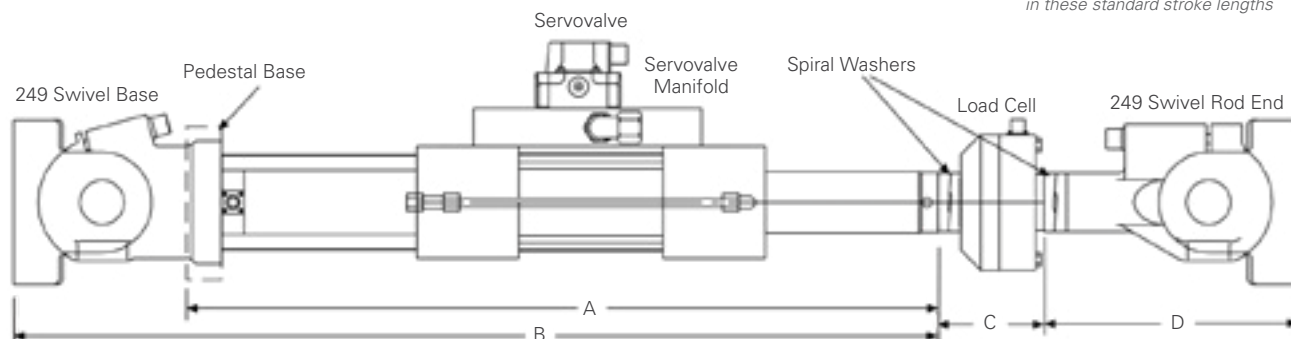

MTS Series 248 Hydraulic Actuators

MTS Series 248 Actuators are heavy-duty, double-ended, linear actuators designed specifically for vibration testing of structures and components. Hydrostatic bearings enable these actuators to withstand high side-load forces without excessive wear on the rod or bearings. The bearings are also coated with a polymer material to protect the actuator if side-load forces exceed the hydrostatic bearing capacity.

- » Recommended for low friction, low distortion, high sideload, vibration testing
- » Double-ended (piston rod extends from both sides of the equal area piston)
- » Double-acting (hydraulic pressure on both sides of the equal area piston)
- » Full stroke coaxially mounted LVDT displacement transducer
- » Hydrostatic pressure-centered bearings for continuous high speed operations
- » Pedestal base allows actuator to be mounted to a reaction mass.



Series 248 Actuator Specifications

MODEL	Force Rating		Piston Area		Rod Diameter	
	kN	KIP	cm ²	in. ²	MM	IN.
248.01	10	2.2	5.23	0.81	80	3.15
248.02	16	3.5	8.32	1.29	80	3.15
248.03	28	6.2	14.52	2.25	80	3.15
248.04	38	8.5	20	3.1	80	3.15
248.05	50	11	25.67	3.98	80	3.15
248.11	100	22	50.77	7.87	100	3.94
248.12	160	35	81.29	12.6	100.1	3.94
248.21	250	55	127.03	19.69	125	4.92

Specifications are subject to change without notice. Contact MTS for verification of specifications critical to your needs.

Series 249 Swivel Base and Rod End

Series 249 Swivels are important load transferring components. Benefits these swivels provide include: minimizing backlash, forgiving offset loads, and providing a pivot to enable planned movements.

TYPICAL APPLICATIONS


Static and dynamic structural or component test systems requiring a backlash-free actuator connection that pivots.

FEATURES

- » Strong, durable, versatile swivels specifically designed for the rigorous demands of closed-loop servohydraulic testing applications.
- » Force ratings available from 25 kN (5,500 lb) to 1550 kN (350,000 lb).
- » Unique positive bearing preload adjustment eliminates backlash that causes load spikes and disturbs loading rates.
- » Design utilizes cast, fatigue resistant, ductile iron for minimum weight and maximum strength.
- » Maintenance free, non-metallic swivel bearings reduce swivel friction, eliminate the need for lubrication, and increase swivel life.

DESCRIPTION

Swivel Base – The swivel base can be used to attach the linear hydraulic actuator or other actuating devices to a reaction mass or frame for structural or component testing. It typically bolts to an actuator end cap, but can be used anywhere a backlash-free swivel connection is required. The mounting holes are designed to fit a variety of MTS actuators and fixtures.

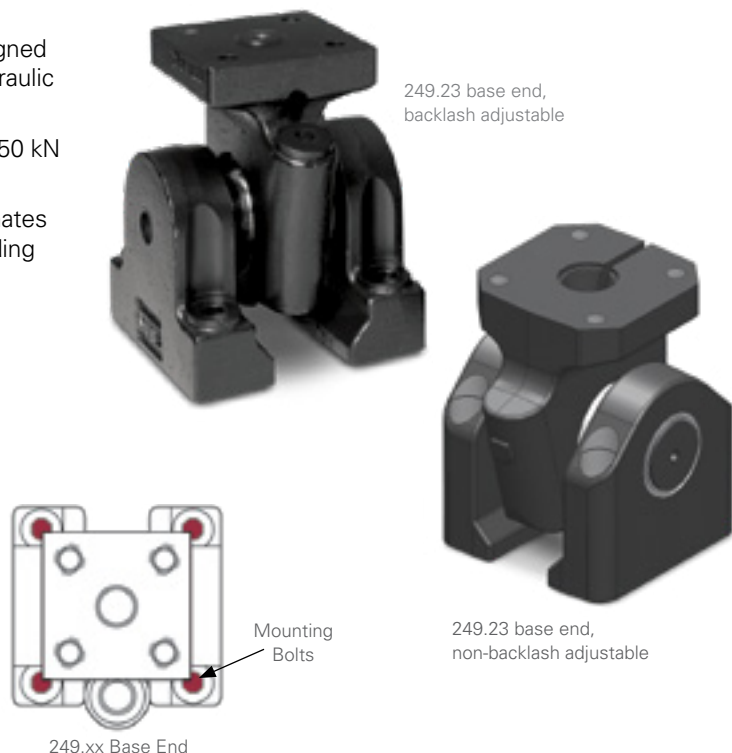

Swivel Base Specifications

MODEL	FORCE RATING	WEIGHT	HEIGHT	ACTUATOR END	MOUNTING BOLT SIZE		PART NUMBER
					SI Metric	US Customary	
249B.M25	25 kN (5.5 kip)	7 kg (15 lb)	127 mm (5 in)	Base	M10x1.50	(3/8 in-16)	100-221-936
249B.M70	70 kN (15 kip)	29 kg (64 lb)	180.0 mm (7.12 in)	Base	M16x2.00	(5/8 in-11)	100-221-932
249B.M160	160 kN (35 kip)	96 kg (211 lb)	238.3 mm (9.38 in)	Base	M16x2.00	(5/8 in-11)	100-221-924
249B.M340	340 kN (77 kip)	222 kg (488 lb)	317.5 mm (12.5 in)	Base	M24x3.50	(1 in-8)	100-221-928
249B.M500	500 kN (110 kip)	334 kg (735 lb)	406.4 mm (16 in)	Base	M30x3.50	(1-1/4 in-7)	100-227-944
249B.M500	730 kN (165 kip)	632 kg (1390 lb)	511.3 mm (20.02 in)	Base	M30x3.50	(1-1/4 in-7)	100-227-948
249B.M1000	1000 kN (220 kip)	995 kg (2189 lb)	546.1 mm (21.5 in)	Base	M40x4.00	(1-1/2 in-6)	100-227-952
249NB.M25	25 kN (5.5 kip)	7 kg (15 lb)	127 mm (5 in)	Non- Adj Base	M10x1.50	(3/8 in-16)	100-221-937
249NB.M70	70 kN (15 kip)	29 kg (64 lb)	180.0 mm (7.12 in)	Non- Adj Base	M16x2.00	(5/8 in-11)	100-221-933
249NB.M160	160 kN (35 kip)	68 kg (150 lb)	238.3 mm (9.38 in)	Non- Adj Base	M16x2.00	(5/8 in-11)	100-221-925
249NB.M340	340 kN (77 kip)	163 kg (359 lb)	317.5 mm (12.5 in)	Non- Adj Base	M24x3.50	(1 in-8)	100-221-929
249NB.M500	500 kN (110 kip)	349 kg (768 lb)	406.4 mm (16 in)	Non- Adj Base	M30x3.50	(1-1/4 in-7)	100-227-945
249NB.M1000	1000 kN (220 kip)	973 kg (2141 lb)	546.1 mm (21.5 in)	Non- Adj Base	M40x4.00	(1-1/2 in-6)	100-227-953

Series 249 Swivel Base and Rod End

DESCRIPTION

Swivel Rod End – The swivel rod end is typically attached to the other components in the force train with a US Customary or Metric thread connector stud. A fatigue-resistant, preloaded connection between the actuator and swivel can be created by using MTS Model 601 Spiral Washers in conjunction with a threaded stud.



Swivel Rod End Specifications

MODEL	FORCE RATING	WEIGHT	HEIGHT	ACTUATOR END	MOUNTING BOLT SIZE		PART NUMBER
					SI Metric	US Customary	
249R.M25	25 kN (5.5 kip)	5 kg (11 lb)	120.7 mm (4.75 in)	Rod	M10x1.50	(3/8 in-16)	100-221-938
249R.M70	70 kN (15 kip)	24 kg (53 lb)	198.2 mm (7.8 in)	Rod	M16x2.00	(5/8 in-11)	100-221-934
249R.M160	160 kN (35 kip)	60 kg (132 lb)	263.7 mm (10.38 in)	Rod	M16x2.00	(5/8 in-11)	100-221-926
249R.M340	340 kN (77 kip)	199 kg (438 lb)	355.6 mm (14 in)	Rod	M24x3.50	(1 in-8)	100-221-930
249R.M500	500 kN (110 kip)	360 kg (792 lb)	406.4 mm (16 in)	Rod	M30x3.50	(1-1/4 in-7)	100-227-946
249R.M500	730 kN (165 kip)	706 kg (1553 lb)	546.1 mm (21.5 in)	Rod	M30x3.50	(1-1/4 in-7)	100-227-950
249R.M1000	1000 kN (220 kip)	1069 kg (2352 lb)	546.1 mm (21.5 in)	Rod	M40x4.00	(1-1/2 in-6)	100-227-954
249NR.M25	25 kN (5.5 kip)	5 kg (11 lb)	120.7 mm (4.75 in)	Non-Adj Rod	M10x1.50	(3/8 in-16)	100-221-939
249NR.M70	70 kN (15 kip)	25 kg (55 lb)	198.2 mm (7.8 in)	Non-Adj Rod	M16x2.00	(5/8 in-11)	100-221-935
249NR.M160	160 kN (35 kip)	60 kg (132 lb)	263.7 mm (10.38 in)	Non-Adj Rod	M16x2.00	(5/8 in-11)	100-221-927
249NR.M340	340 kN (77 kip)	140 kg (308 lb)	355.6 mm (14 in)	Non-Adj Rod	M24x3.50	(1 in-8)	100-221-931
249NR.M500	500 kN (110 kip)	394 kg (867 lb)	406.4 mm (16 in)	Non-Adj Rod	M30x3.50	(1-1/4 in-7)	100-227-947
249NR.M1000	1000 kN (220 kip)	1087 kg (2391 lb)	546.1 mm (21.5 in)	Non-Adj Rod	M40x4.00	(1-1/2 in-6)	100-227-955

Series 252 Servovalves

MTS Series 252 Servovalves are two stage, four-way valves designed for servocontrol of hydraulic systems. These high performance servovalves are available in a range of flow rates from 4 to 227 lpm (1 to 60 gpm).

TYPICAL APPLICATIONS

- » Position, load, or velocity control of a hydraulic actuator in highly dynamic applications
- » Pressure control in a dynamic pressurization system.

STANDARD FEATURES & BENEFITS

- » Rugged design for long operating life in dynamic applications
- » High resolution and low hysteresis for precise flow control
- » 31 MPa (4500 psi) operating pressure for Models 252.2x and 252.4x, 21 MPa (3000 psi) for Model 252.3x
- » Internal filtration to protect nozzle orifices from contamination
- » High spool driving forces reduce the likelihood of a sticking spool due to dirty fluid
- » Standard mounting configuration for Model 252.2x allows interchangeability between valves of different flow ratings.

SERIES 252 SERVOVALVE OPTIONS

External Pilot Pressure – Series 252 Servovalves are available with a fifth port for external pilot pressure. External pilot pressure allows the second stage spool position circuit to be pressurized and controlled without applying system pressure to the actuator. This option is recommended when precise control of the servoactuator must be maintained when system pressure is first started.



Model 252.2x
Servovalve

Series 252 Servovalve Specifications

MODEL	FLOW RATING	PILOT PRESSURE	90° POINT @ 10% COMMAND	PART NUMBER
252.21G-01	4 lpm (1 gpm)	No	240 Hz	056-006-101
252.21G-04	4 lpm (1 gpm)	Yes	240 Hz	056-006-104
252.22G-01	9.5 lpm (2.5 gpm)	No	240 Hz	056-006-201
252.22G-04	9.5 lpm (2.5 gpm)	Yes	240 Hz	056-006-204
252.23G-01	19 lpm (5 gpm)	No	240 Hz	056-006-301
252.23G-04	19 lpm (5 gpm)	Yes	240 Hz	056-006-304
252.24G-01	37 lpm (10 gpm)	No	200 Hz	056-006-401
252.24G-04	37 lpm (10 gpm)	Yes	200 Hz	056-006-404
252.25G-01	56 lpm (15 gpm)	No	170 Hz	056-006-501
252.25G-04	56 lpm (15 gpm)	Yes	170 Hz	056-006-504
252.26G-01	63 lpm (16.5 gpm)	No	120 Hz	056-061-301
252.26G-04	63 lpm (16.5 gpm)	Yes	120 Hz	056-061-304
252.31A-01	93 lpm (25 gpm)	Convertible	80 Hz	004-823-101
252.32A-01	151 lpm (40 gpm)	Convertible	60 Hz	004-823-201
252.33A-01	227 lpm (60 gpm)	Convertible	50 Hz	004-823-301

Weight: 1 kg (2.3 lb)

Electrical Interface: MS3106F14S-2S

Hydraulic Interface: ISO 10372-04-04-0-92

Maximum Operating Pressure: 31 MPa (4500 psi)

Series 215 & 216 Rotary Actuators

MTS Model 215 and 216 Rotary Actuators are heavy duty, torque generating devices designed for precise servo-control applications. These actuators, along with the matching accessories, create a versatile test stand for torsional testing applications.

TYPICAL APPLICATIONS

- » Fatigue and durability testing of drivetrain components such as axles, couplings, and drive shafts
- » Ultimate strength and other material tests on shafts.

STANDARD FEATURES

- » Fatigue-rated design for years of reliable service
- » Heavy duty roller bearings endure high radial loads with minimal friction
- » Hydrostatic thrust bearings to safely withstand inadvertent thrust loads
- » 90 degrees of dynamic rotation and 100 degrees of static rotation
- » Hydraulic cushions protect the vanes during full stroke, high velocity operation
- » An angular displacement transducer (ADT) provides the signal for precise angular position feedback in a closed loop servocontrol application
- » A flange adapter to attach the rotary actuator shaft to the test specimen.



- » The Model 215 includes an actuator mounting bracket with flexures to reduce thrust loads caused by specimen deformation when torsional loads are applied. The Model 216 includes two diaphragm flexures that attach to each end of the test specimen.

HYDRAULIC CONNECTION

Pressure & Return – Pressure and return line connections depend on the servovalve. Configurations with Series 252.2x Servovalve Manifolds have JIC-16 fittings. Series 252.3x and 256.0x Manifolds have JIC-20 fittings.

Drain – All 215 & 216 Actuators have a JIC-6 drain connection.

Rotary Actuator Specifications

MODEL	RATED TORQUE		ANGULAR DISP.		MAX. THRUST LOAD		PART NUMBER
	N-M	LBF-IN	DYNAMIC	STATIC	KN	LB	
215.32	226	2,000	90°	100°	3.3	750	100-191-977
215.35	565	5,000	90°	100°	3.3	750	100-191-978
215.41	1130	10,000	90°	100°	3.3	750	100-191-979
215.42	2260	20,000	90°	100°	3.3	750	100-191-981
215.45	5650	50,000	90°	100°	5.3	1,200	100-191-994
215.51	11,300	100,000	90°	100°	5.3	1,200	100-191-995

Typical Configuration for Component Applications

MODEL	SERVOVALVE	TORQUE TRANSDUCER	REACTION BRACKET	REACTION BASE
215.32	252.2x (qty 1)	663.11A-02	215.35TCRB	215.35RB
215.35	252.2x (qty 1 or 2)	663.12A-01	215.35TCRB	215.35RB
215.41	252.2x (qty 1 or 2)	663.12A-02	215.42TCRB	215.42RB
215.42	252.2x (qty 1 or 2), 252.3x (qty 1)	663.13A-02	215.42TCRB	215.42RB
215.45	252.2x (qty 1 or 2), 252.3x (qty 1)	663.13A-05	215.45TCRB	215.45RB
215.51	252.2x (qty 1 or 2), 252.3x (qty 1), 256.0x (qty 1)	663.14A-01	215.51TCRB	215.51RB
216.10	256.0x (qty 1)	663.15A-02	216.10TCRB	216.10RB
216.20	256.0x (qty 1)	663.16A-01	216.20TCRB	216.20RB
216.30	256.0x (qty 1)	663.18A-01	216.30TCRB	216.30RB

Rotary Actuator Options

Torque Cells

MTS torque cells provide a precise electrical feedback signal that is proportional to the torque applied to the specimen.

Reaction Bracket and Base

A rigid reaction frame is critical for good performance. MTS provides a complete line of reaction fixturing that provide the stiffness necessary for quality testing. When considering any durability testing, remember that MTS equipment is fatigue-rated to ensure your test schedule will not be compromised by your test apparatus.

Diaphragm Flexures

Diaphragm flexures offer the satisfaction of knowing that the test article is not becoming loaded by unwanted externals due to misalignment or an unplanned axial loading component. Be sure the testing you are performing has not become biased by unintended loads.

Series 256 Servovalves

Series 256 Servovalves are high-flow, high-performance valves of three-stage design. They are suited for 3000 psi dynamic application, and provide the user superior position, load, or velocity control resolution when applied to a Model 216 Rotary Actuator. When selecting a Model 256 Servovalve for your Model 216 Actuator, also select manifold part number 052-812-101.

Actuator Configuration

MODEL	CYLINDER	MANIFOLD	ADAPTER FLANGE	ADT	FLEXURE/FOOT MOUNT
215.32		043-932-301		032-610-111	
215.35	047-466-001	043-932-401		032-610-111	
215.41	047-528-601	032-554-801		032-610-111	
215.42	047-535-301	032-554-801		032-610-111	
215.45	049-029-401	032-554-901		032-610-112	036-373-301
215.51	049-029-501	032-554-901		032-610-112	
216.10	050-206-901	034-676-101		032-610-115	NA
216.20		034-676-101	100-016-581	032-610-116	NA
216.30	054-128-401	034-676-101			NA

MODEL	TORQUE CELL	REACTION BRACKET	REACTION BASE	DIAPHRAGM FLEXURE
215.32	006-437-502	043-962-701		032-553-601
215.35	006-437-504	043-962-701	044-527-001	032-553-601
215.41	006-437-505			032-553-701
215.42	006-437-507	043-962-801	044-702-901	032-553-801
215.45	006-437-510		044-703-001	032-553-901
215.51	006-437-511		044-703-001	032-554-001
216.10	006-437-513			032-554-101
216.20	006-437-514			032-554-201
216.30	006-437-516			

SERIES 215 OR 216 ACTUATOR ORDER INFORMATION

When you order a Model 215 or 216 Actuator, you receive a servoactuator assembly including the rotor actuator, a servovalve manifold, an ADT, and an actuator flexure mounting bracket.

See MTS Model 215 Product Specification number 100-007-528a and Model 216 Product Specification number 100-032-134a for additional information.

MTS Series 515 Hydraulic Power Units

- » Variable volume pumps help save energy and reduce operating costs
- » Designed to eliminate ambient heat load, reducing HVAC requirements
- » Compact design fits through a standard doorway
- » Optional Remote-Mounted Water shut-off valve saves water when the unit is shut down
- » Efficient heat transfer technology reduces cooling requirements with standard air and water cooling options
- » Advanced programmable logic controls (PLC) enable options for control and monitoring
- » Engineered for remote monitoring and multi-pump control options
- » High pressure filter available for additional system filtration at the outlet of the pump
- » Output accumulation available as an option for all sizes
- » Designed to relevant EN and UL standards; CE compliant; TÜV certified including TUV CUE Certified to EMC
- » Safety engineered to SIL-3 rating when operated in standalone mode



Model 515.07 – 515.30 Specifications

Operating pressure: 210 bar (3,000 psi)

Pump type: Variable displacement pumps

Filtration: Full flow on the return side

Maximum ambient operating temperature: 40°C (104°F)

Minimum ambient operating temperature: 5°C (40°F)

	MODEL 515.07	MODEL 515.11	MODEL 515.20	MODEL 515.30
Flow rates (for 60 Hz models)	26.5 lpm (7 gpm)	41.6 lpm (11 gpm)	75.7 lpm (20 gpm)	113.5 lpm (30 gpm)
Flow rates (for 50 Hz models)	22.7 lpm (6 gpm)	37.9 lpm (10 gpm)	62.5 lpm (16.5 gpm)	100.7 lpm (26.6 gpm)
Noise level*	58 dB(A)	60 dB(A)	63 dB(A)	63 dB(A)
Reservoir capacity (maximum)	174 L (46 gal)	174 L (46 gal)	341 L (90 gal)	341 L (90 gal)
Unit dimensions				
Width	77.0 cm (30.3 in)	77.0 cm (30.3 in)	89.7 cm (35.3 in)	89.7 cm (35.3 in)
Height	120.1 cm (47.3 in)	120.1 cm (47.3 in)	138.9 cm (54.7 in)	138.9 cm (54.7 in)
Length	114.6 cm (45.1 in)	114.6 cm (45.1 in)	163.1 cm (64.2 in)	163.1 cm (64.2 in)
Weight with maximum oil	522 kg (1,150 lb)	544 kg (1,200 lb)	862 kg (1,900 lb)	907 kg (2,000 lb)
Motor size	11 kW (15 hp)	18.5 kW (25 hp)	30 kW (40 hp)	45 kW (60 hp)

* Sound levels [dB(A)] are expressed as a free field value. Readings may vary with acoustic environment.

Specifications subject to change without notice. Please contact MTS for specifications critical to your application.

Model 515.60 – 515.180 Specifications

Operating pressure: 210 bar (3,000 psi)

Pump type: Variable displacement pumps

Filtration: Full flow on the return side

Maximum ambient operating temperature: 40°C (104°F)

Minimum ambient operating temperature: 5°C (40°F)

	MODEL 515.60	MODEL 515.90	MODEL 515.120	MODEL 515.150	MODEL 515.180
Flow rates (for 60 Hz models)	227 lpm (60 gpm)	340 lpm (90 gpm)	454 lpm (120 gpm)	567 lpm (150 gpm)	681 lpm (180 gpm)
Flow rates (for 50 Hz models)	200 lpm (53.2 gpm)	300 lpm (80 gpm)	400 lpm (106.4 gpm)	500 lpm (133 gpm)	600 lpm (160 gpm)
Noise level*	68 dB (A)	68 dB (A)	70 dB (A)	71 dB (A)	72 dB (A)
Reservoir capacity (maximum)	1211 L (320 gal)	1211 L (320 gal)	2,188 L (578 gal)	2,188 L (578 gal)	2,188 L (578 gal)
Unit dimensions					
Width	103.4 cm (40.7 in)	103.4 cm (40.7 in)	103.4 cm (40.7 in)	103.4 cm (40.7 in)	103.4 cm (40.7 in)
Height	199.4 cm (78.5 in)	199.4 cm (78.5 in)	199.4 cm (78.5 in)	199.4 cm (78.5 in)	199.4 cm (78.5 in)
Length	287.0 cm (113.0 in)	287.0 cm (113.0 in)	430.5 cm (169.5 in)	430.5 cm (169.5 in)	430.5 cm (169.5 in)
Weight with maximum oil	2,835 kg (6,250 lb)	3,289 kg (7,250 lb)	4,876 kg (10,750 lb)	5,330 kg (11,750 lb)	5,783 kg (12,750 lb)
Motor size	45 KW (60 Hp)	45 KW (60 Hp)	45 KW (60 Hp)	45 KW (60 Hp)	45 KW (60 Hp)
Number of Motor/pump units	2, max of 3	3	4, max of 6	5, max of 6	6


* Sound levels [dB(A)] are expressed as a free field value. Readings may vary with acoustic environment.

Specifications subject to change without notice. Please contact MTS for specifications critical to your application.

MTS Multi-Pump Control Manager (MPCM)

Manage the hydraulic power demands of Hydraulic Power Units (HPUs) as a system, rather than unit by unit. The Multi-Pump Control Manager (MPCM) shows all hydraulic power demand in a single touchscreen interface so that it can be automatically or manually managed for the best efficiency. This solution includes the Basic Connection Interface (BCI), Supervisor PC and touchscreen monitor. Control includes all control options from a unit or system level. Similarly, status information can be viewed from a unit or system level.


- » Centralized monitoring and control of up to eight hydraulically-commoned HPUs
- » Input, view or change individual station flow requirements for multiple stations from a single location
- » Monitor and control hydraulic system supply and demand to improve productivity, reduce energy costs and prolong test equipment life
- » Configurations can include up to 8 HPUs, 8 External E-Stop chains and/or Station Flow Managers, E-Stop Out, 32 channels of DIO, and 3 System Status Light Bars
- » Includes cables to connect each HPU to BCI, Supervisor PC and touchscreen monitor

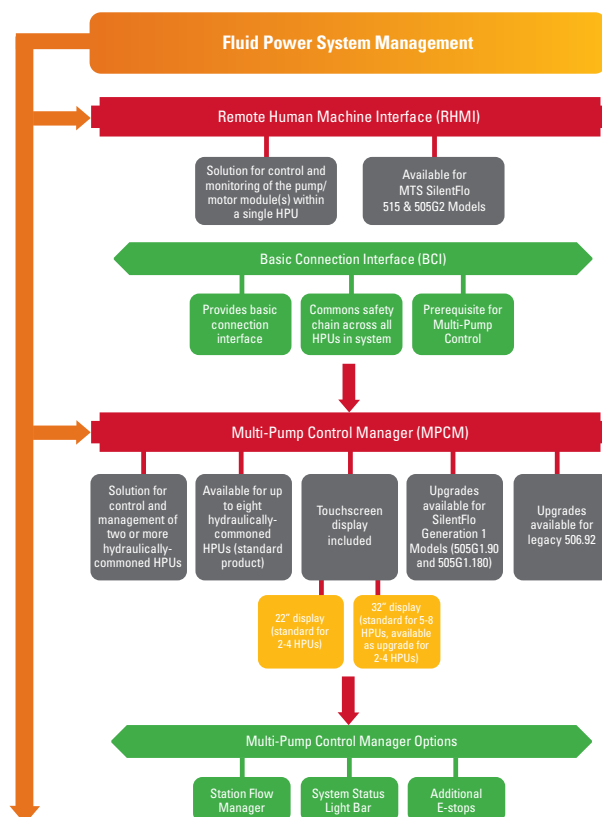

MTS Series 293 Hydraulic Service Manifold

Series 293 Hydraulic Service Manifolds (HSMs) are hydraulic pressure and flow regulation devices that control pressure to a test station independent from the main hydraulic power unit (HPU). Installing the HSM between the HPU and the actuator allows the operator to turn each hydraulic station on and off and set the low pressure level. Smooth, controlled transitions between off, low and high pressure modes enhance safety and permit predictable system control.

- » Minimizes the effects of rapid application and removal of pressure, reducing unexpected actuator movement that could be damaging to the component or assembly being tested
- » Equipped with a hydraulic circuit which provides separately filtered fluid at HPU pressure to operate the pilot stage of multistage servovalves
- » Uses accumulators to reduce pressure fluctuations caused by changing system demands
- » Hydraulic solenoid valves control the hydraulic pressure to the work station.

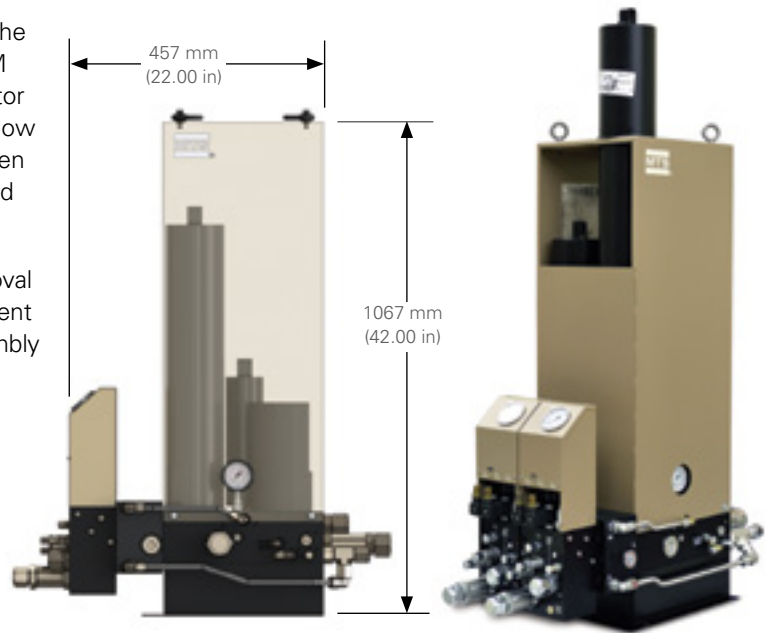

Series 293 Hydraulic Service Manifold Specifications

DESCRIPTION	293.1X	293.22	293.24	293.32	293.34
Operating Pressure	21 MPa (3000 psi)	21 MPa (3000 psi)	28 MPa (4000 psi)	21 MPa (3000 psi)	28 MPa (4000 psi)
Nominal flow	189 lpm (50 gpm)	378 lpm (100 gpm)*	378 lpm (100 gpm)	946 lpm (250 gpm)	946 lpm (250 gpm)
Variable low pressure	1 - 21 MPa (150 - 3000 psi)	1 - 21 MPa (150 - 3000 psi)	1 - 28 MPa (150 - 4000 psi)	1 - 21 MPa (150 - 3000 psi)	1 - 28 MPa (150 - 4000 psi)
Slow on/off ramp time	5.0 - 9.0 seconds	5.0 - 9.0 seconds	5.0 - 9.0 seconds	5.0 - 9.0 seconds	5.0 - 9.0 seconds
Maximum current to solenoid valves	0.3 A at 115 V AC or 1.5 A at 24 V DC	0.3 A at 115 V AC or 1.5 A at 24 V DC	0.3 A at 115 V AC or 1.5 A at 24 V DC	0.3 A at 115 V AC or 1.5 A at 24 V DC	0.3 A at 115 V AC or 1.5 A at 24 V DC

Accumulators

Pressure	min max	0.94 L (0.25 gal) 3.8 L (1 gal)	3.8 L (1 gal) 7.6 L (2 gal)	3.8 L (1 gal) 7.6 L (2 gal)	3.8 L (1 gal) 19.0 L (5 gal)	3.8 L (1 gal) 19.0 L (5 gal)
Return	min max	0.45 l (0.12 gal) 0.94 l (0.25 gal)	1.9 L (0.5 gal) 3.8 L (1 gal)	1.9 L (0.5 gal) 3.8 L (1 gal)	0.94 l (0.25 gal) 7.6 L (2 gal)	0.94 l (0.25 gal) 7.6 L (2 gal)
Pilot pressure	min max	0.45 L (0.12 gal) 0.94 l (0.25 gal)	0.45 L (0.12 gal) 0.94 l (0.25 gal)	0.45 L (0.12 gal) 0.94 l (0.25 gal)	0.45 L (0.12 gal) 0.9 L (0.25 gal)	0.45 L (0.12 gal) 0.94 L (0.25 gal)

Dimensions and Weights

Weight	70 kg (155 lb)	159 kg (350 lb) ¹	182 kg (400 lb) ¹	181 kg (400 lb) ²	227 kg (500 lb) ¹
Height	806 mm (31.75 in)	1067 mm (42 in) max	1067 mm (42 in)	1092 mm (43.5 in) ³	1229 mm (48.38 in)
Width (base)	363 mm (14.3 in)	362 mm (14.25 in)	362 mm (14.25 in)	403 mm (15.9 in) max	387 mm (15.25 in)
Depth (base)	359 mm (14.1 in)	394 mm (15.5 in)	400 mm (15.75 in)	394 mm (15.5 in) max	530 mm (20.9 in)

Footnotes:

* 378 lpm (100 gpm) total, 189 lpm (50 gpm) per station maximum for multi-station flow

† In a multi-station configuration, pressure -20 IN and -16 OUT

‡ In a multi-station configuration, return -16 IN and -20 OUT

** In a multi-station configuration, HSM -06 IN and -12 OUT

*** In a multi-station configuration

¹ Standard size for single-station configuration; 227 kg (500 lb) maximum weight for four-station HSM

² Standard size; with optional accumulator capacity the weight is 454 kg (1000 lb) maximum

³ Standard size; with optional accumulator capacity the height is 2134 mm (84 in) maximum

Specifications are subject to change without notice. Contact MTS for verification of any critical specifications. Note: Accumulator may extend beyond HSM height.

MTS Isolation Service Manifold

The Isolation Service Manifold (ISM) systematically controls hydraulic distribution and is compliant with Machinery Directive 2006/42/EC and sound engineering practice article 3 (3) of pressure equipment directive 97/23/EC. This manifold is placed between the hydraulic power unit (HPU) and the hydraulic service manifold (HSM) to prevent system over-pressurization, provide full system pressure discharge, and allow for full hydraulic power isolation downstream of the device. It is designed to work with MTS Series 293 HSMs, and MTS FlexTest® controllers, and can be adapted to work with other HPUs and HSMs. It is easy to set up in a new system or you may add it to an existing test system using the cabling in the current system.

- » Plug-and-Play service manifold complies with Machinery Directive 2006/42/EC
- » Positive isolation of hydraulic pressure
- » Pressure and valve position feedback in control loop
- » Over-pressurization relief
- » Pressure dump to return line
- » Maintenance valves to drain residual oil
- » Requires Isolation Service Manifold Controller which can be shared with one other ISM.



Isolation Service Manifold Specifications

MODEL	ISM02-210	ISM02-280	ISM04-210	ISM04-280	ISM10-210	ISM10-280
Flow Capacity	200 lpm (50 gpm)	200 lpm (50 gpm)	400 lpm (100 gpm)	400 lpm (100 gpm)	1000 lpm (250 gpm)	1000 lpm (250 gpm)
Pressure Rating	210 bar (3000 psi)	280 bar (4000 psi)	210 bar (3000 psi)	280 bar (4000 psi)	210 bar (3000 psi)	280 bar (4000 psi)
Control Voltage	24 VDC	24 VDC	24 VDC	24 VDC	24 VDC	24 VDC
Required Control Power	70 watts	70 watts	70 watts	70 watts	70 watts	70 watts
DIMENSION	ISM02-210	ISM02-280	ISM04-210	ISM04-280	ISM10-210	ISM10-280
Width	27.30 cm (10.75 in)	27.30 cm (10.75 in)	27.30 cm (10.75 in)	27.30 cm (10.75 in)	36.20 cm (14.25 in)	48.90 cm (19.25 in)
Length	42.24 cm (16.63 in)	42.24 cm (16.63 in)	42.24 cm (16.63 in)	42.24 cm (16.63 in)	60.12 cm (23.67 in)	59.64 cm (23.48 in)
Height	60.60 cm (23.86 in)	60.60 cm (23.86 in)	57.58 cm (22.67 in)	58.85 cm (23.17 in)	61.87 cm (24.36 in)	73.30 cm (28.86 in)
Weight	53.52 kg (118.0 lb)	114.53 kg (252.5 lb)	58.97 kg (130.0 lb)	124.15 kg (273.7 lb)	215.09 kg (474.2 lb)	288.48 kg (636.0 lb)

MTS TestLine™ Modular Components

Hydraulic Hoses


All reinforced hydraulic hoses supplied for the purpose of connecting hydraulic power to MTS Test components, with the exception of low pressure drain hoses, will conform to industry standards: SAE J517–, SAE J343–Test and procedures for SAE 100R, and ISO 11171–Contamination Control.

- » Silicone-free fabrication and testing
- » Oil-resistant labels with part number, date of manufacture, diameter and working pressure.



Reaction Brackets

Reaction brackets attach securely to the reaction base or T-slot table and provide a mounting surface for the torque cell. Each reaction bracket is designed to restrain a specific model torque cell.

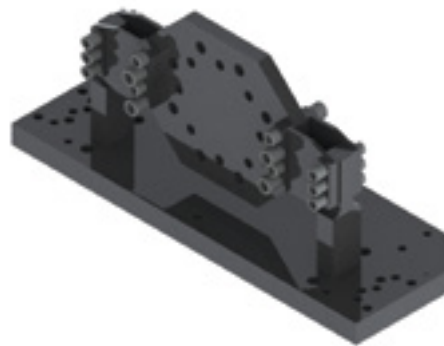

Diaphragm Flexure Assembly

Diaphragm flexures are used at both ends of the specimen if large axial and angular deflections are generated during testing. If reaction forces exceed stated actuator operating limits, diaphragm flexures help reduce the thrust and side loads reacted by the actuator and test specimen.



Load Cells

MTS 661 Series Load Cells are designed for a wide array of static and dynamic testing applications. The cells are designed for cyclic operation in through zero tension/compression modes.



Torque Cells

Torque cells provide a precise electrical feedback signal that is proportional to the torque applied to the specimen. For more information on MTS torque cells, refer to the appropriate MTS product specification.



Additional Components

Model 601 Spiral Washers

MODEL	DESCRIPTION	FORCE	TEMP °C MIN/MAX	TEMP °F MIN/MAX	DIAMETER ID/OD	SPANNER PIN DIAMETER
601.11A/B-20	Spiral Washers	25 kN (5 kip)	-129°/177°	-200°/350°	12 mm (0.50 in)/41 mm (1.62 in)	3.58 mm (0.141 in)
601.11A/B-11	Spiral Washers	100 kN (22 kip)	-129°/177°	-200°/350°	27 mm (1.00 in)/66 mm (2.62 in)	6.73 mm (0.265 in)
601.11A/B-19	Spiral Washers	250 kN (55 kip)	-129°/177°	-200°/350°	36 mm (1.50 in)/92 mm (3.62 in)	6.73 mm (0.265 in)
601.11A/B-13	Spiral Washers	500 kN (110 kip)	-129°/177°	-200°/350°	52 mm (2.00 in)/130 mm (5.12 in)	9.91 mm (0.390 in)
601.11A/B-15	Spiral Washers	1000 kN (220 kip)	-129°/177°	-200°/350°	76 mm (3.00 in)/178 mm (7.00 in)	9.91 mm (0.390 in)

The unique design of patented MTS spiral washers ensures preloading without inducing offsets in the load train. Normally sold for use on load frames, these versatile accessories can be used anywhere a backlash-free, threaded union is required. When setting the spiral washer connection, the stud must be preloaded to a minimum axial load that is 110% of test maximum.

One Set of Two Spiral Washers

MODEL	METRIC PART NUMBER	US CUSTOMARY PIN DIAMETER
601.11A/B-20	040-473-120	040-473-119
601.11A/B-11	040-473-110	040-473-101
601.11A/B-19	040-473-111	040-473-102
601.11A/B-13	040-473-113	040-473-104
601.11A/B-15	040-473-115	040-473-106

Shim Kits

These kits provide multiple shims and thicknesses to allow for 360° of adjustment.

CAPACITY	THREAD SIZE	METRIC PART NUMBER	US CUSTOMARY PIN DIAMETER
25 kN (5.5 kip)	M12 x 1.25 (1/2"-20)	052-105-004	052-105-003
50, 100, 110 kN (11, 22, 25 kip)	M27 x 2 (1"-14)	052-105-006	052-105-005
250 kN (55 kip)	M36 x 2 (1 1/2"-12)	052-105-008	052-105-007
500 kN (110 kip)	M52 x 2 (2"-12)	052-105-010	052-105-009

High-Performance Lubricants for Testing Applications

Room temperature Molykote® lubricant is best suited for smooth operation of metal surfaces that are under pressure or load. It is a good anti-galling agent for grip wedges and other tough lubrication problems.

High-temperature Jet Lube anti-seize compound is best suited for smooth operation of metal surfaces that are under pressure or load at elevated temperatures (up to 2100°F). It is a good anti-galling agent for studs and other tough lubrication problems that include higher temperatures.

DESCRIPTION	SIZE	PART NUMBER
Molykote Lubricant /Grease	17 oz	011-010-207
Anti-Seize Compound	17 oz	011-354-902

MAINTENANCE PARTS

MTS stocks a wide variety of parts to help you maintain your test equipment.

Don't see what you're looking for? Contact your service sales engineer for more information.



Introducing the *MTS Platinum Promise* for Load Cells

Our most popular load cells are in stock and ready to ship.

Look for this symbol to identify Platinum Promise load cells.



MTS Maintenance Parts

Introduction	146
Hydraulic Tools & Sample Bottles	147
Accumulator Parts	148
Actuator Parts	149-153
Load Frame Parts	154-155
Cables	156-159
Replacement Lamps	159
Adapter Cables for 3rd Part Equipment	160
Grip Parts	161
Hydraulic Hoses	162
Hydraulic Power Unit Parts	163-168
Hydraulic Service Manifold Parts	169-170
Filter Elements	171
Flushing & Shut-Off Valves	171
SWIFT® Parts	172-173
FlatTrac® & Rolling Road Systems Parts	174
Series 329 Road Simulator Parts	175
329 System Parts	176
Model 329 4 DOF PC Road Simulator Parts	178-180
Model 329 6 DOF Road Simulator Actuator Parts	181
Model 329 6 DOF PC Road Simulator Parts	182-184
Model 329 4 DOF LT Road Simulator Parts	185-186
Model 329 6 DOF LT Road Simulator Parts	187-189
MAST Parts	190
320 Seal Kits	191

MTS Maintenance Parts

Hydraulic Tools & Sample Bottles

Tools – LVDT

DESCRIPTION	PART NUMBER
Rod-Threaded End-LVDT Adj. (3 kg)	003-116-201
Key-Socket Head-LVDT Adj. (3 kg)	003-115-701
LVDT Plug Wrench (206) (3 kg)	004-059-901

(For 6" Stroke) Note: Consult MTS for additional lengths.

Tools – Spanner Wrench

DESCRIPTION	WASHER SIZE	SPANNER PIN DIAMETER	PART NUMBER
Spanner Wrench – Spiral Washers	19 mm to 50 mm (0.75 in to 2 in)	3.2 mm (0.125 in)	010-099-313
Spanner Wrench – Spiral Washers	32 mm to 76 mm (1.25 in to 3 in)	6.4 mm (0.250 in)	010-099-304
Spanner Wrench – Spiral Washers	50 mm to 120 mm (2 in to 4.75 in)	6.4 mm (0.250 in)	010-099-301
Spanner Wrench – Spiral Washers	115 mm to 158 mm (4.5 in to 6.25 in)	9.5 mm (0.375 in)	010-099-303

Tools – Hydraulic

DESCRIPTION	PART NUMBER
Fluid Sampling Tool Kit – D.I.Y.	055-589-601
Fluid Transfer Pump - Model 590.01	005-149-201

MTS Fluid Analysis Bottle Kit & Results Report – D.I.Y. Customer Performed

DESCRIPTION	PART NUMBER
DIY Bottle - North America	100-030-731
DIY Bottle - Europe	100-188-130
DIY Bottle - South Korea	100-219-318
DIY Bottle - Japan	100-219-319
DIY Bottle - China	100-219-320

MTS Fluid Analysis Bottle Kit & Results Report – D.I.Y. Customer Performed

DESCRIPTION	PART NUMBER
FSE Bottle - North America	100-188-132
FSE Bottle - Europe	100-188-133
FSE Bottle - South Korea	100-219-330
FSE Bottle - Japan	100-219-331
FSE Bottle - China	100-219-332

MTS Maintenance Parts

Accumulator Parts

Accumulator Seal Kits

DESCRIPTION	STYLE	PART NUMBER
2.5 in I.D., MTS Model 111.11A	(old style) (new style)	036-548-301 041-463-501 
4.0 in I.D., MTS Model 111.12A/B	(old style) (new style)	036-555-101 041-463-301 
2.5 in I.D., Made by Parker	(Silver Label)	100-102-494 
4.0 in I.D., Made by Parker	(Silver Label)	100-102-496 

Accumulator Bladders

DESCRIPTION	PART NUMBER
1 Pint Bladder	010-050-509
O-Ring, Plug	101-010-906
1 Quart Bladder	010-050-510
O-Ring, Plug	010-010-908
1 Gallon Bladder	010-050-511
O-Ring, Plug	010-010-911

Bladder Kit contains Washer & O-Ring.

Accumulator Charging Kits (7 kg)

DESCRIPTION	MANUFACTURER	PART NUMBER
Model 590.05 Charging Kit	MTS	037-698-601

**Only available in Europe.*

Miscellaneous Accumulator Parts

DESCRIPTION	PART NUMBER
Gas Valve (Accumulator)	011-012-402
Gas Valve	010-098-205
Hose Assy – Nitrogen	041-922-801
Hose Assy – Extension, 12 Inches Long	041-922-701
Valve-Chuck, High Pressure	010-098-248
Washer-Flat Copper	010-099-901
Valve Core	010-098-202
Valve Stem Extension	010-098-301
Washer-Sealing, Hard Fiber	011-480-202
Pin Removal Tool for Accumulators	011-870-101
Adjustable Spanner Wrench (For 2.5 in & 4 in I.D.)	010-099-317

Estimated weight of each is no more than 2 kg.

MTS Maintenance Parts

Actuator Parts

201.xx Actuator Seal Kits

MODEL	DESCRIPTION	PART NUMBER
201.17	Seal Kit	057-003-401
201.2	Seal Kit	057-003-402
201.25	Seal Kit	057-003-403
201.3	Seal Kit	057-003-404
201.35	Seal Kit	057-003-405
201.40	Seal Kit	057-003-406
201.45	Seal Kit	057-003-407
201.6	Seal Kit	057-003-408
201.7	Seal Kit	057-003-409
201.80	Seal Kit	057-003-410
201.90	Seal Kit	057-003-411

204.xx Seal Kits (Linear)

MODEL	RATING	DESCRIPTION	PART NUMBER
204.08	5 kN/1.1 kip	Old style cream-colored rod seal	008-710-012
204.08	5 kN/1.1 kip	New style blue-colored rod seal	042-451-601
204.09	10 kN/2.2 kip	Old style cream-colored rod seal	008-710-013
204.09	10 kN/2.2 kip	New style blue-colored rod seal	042-451-401
204.11	12 k N/2.5 kip	Seal Kit	008-710-014
204.12	18 kN/4 kip	Seal Kit	008-710-015
204.13	27 kN/6 kip	Seal Kit	008-710-016
204.21	27 kN/6 kip	Seal Kit	008-710-017
204.22	40 kN/9 kip	Seal Kit	008-710-018
204.23	55 kN/12 kip	Seal Kit	008-710-019
204.24	90 kN/20 kip	Seal Kit	008-710-020
204.25	110 kN/24 kip	Seal Kit	008-710-021
204.26	160 kN/35 kip	Seal Kit	008-710-022
204.31	227 kN/50 kip	Seal Kit	008-710-023
204.32	318 kN/70 kip	Seal Kit	008-710-024
204.41*	455 kN/100 kip	Seal Kit	008-710-025
204.42*	682 kN/150 kip	Seal Kit	008-710-026
204.51	15 kN/3.3 kip	Seal Kit	008-710-001
204.52	25 kN/5.5 kip	Seal Kit	008-710-002
204.61	50 kN/11 kip	Seal Kit	008-710-003
204.62	68 kN/15 kip	Seal Kit	008-710-004
204.63	100 kN/22 kip	Seal Kit	008-710-005
204.64	160 kN/35 kip	Seal Kit	008-710-006
204.71	250 kN/55 kip	Seal Kit	008-710-007
204.72	350 kN/77 kip	Seal Kit	008-710-008
204.81	500 kN/110 kip	Seal Kit	008-710-009
204.82	750 kN/165 kip	Seal Kit	008-710-010
204.91	1000 kN/220 kip	Seal Kit	008-710-011

Estimated weight of each is no more than 2 kg.

Each kit contains the parts needed to replace all existing seals.

*Normally not a stocked item. Call Order Services for delivery time.

MTS Maintenance Parts

Actuator Parts

205.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
205.04	2 kip	Seal Kit	008-710-028
205.08	0.5 kip	Seal Kit	008-710-027
205.08B	0.5 kip	Seal Kit	040-226-201
205.09	1.5 kip	Seal Kit	008-710-027
205.09B	1.5 kip	Seal Kit	040-226-201
205.20	3.3 kip	Seal Kit	039-080-101
205.31	6.6 kip	Seal Kit	039-080-201

206.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
206.21	55 kN/12 kip	Seal Kit	008-710-031
206.31*	127 kN/16 kip	Seal Kit	008-710-032
206.41*	327 kN/36 kip	Seal Kit	008-710-034
206.42*	227 kN/50 kip	Seal Kit	008-710-035

*Normally not a stocked item. Call Order Services for delivery time.

231.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
231.02	6.2 kN	Seal Kit	045-895-201

242.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
242.01	All sizes	Seal Kit	042-065-501
		Seal Kit with high pressure seals	042-065-401
242.02	All sizes	Seal Kit	042-065-501
		Seal Kit with high pressure seals	042-065-402
242.03	All sizes	Seal Kit with high pressure seals	042-065-403

Estimated weight of each is no more than 2 kg.

Each kit contains the parts needed to replace all existing seals.

243.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
243.14	30/65 kN	Seal Kit	047-237-715
243.15*	30/65 kN	Seal Kit	047-237-714
243.17	30/65 kN	Seal Kit	047-237-701
243.20	68/109 kN	Seal Kit	047-237-702
243.25	101/166 kN	Seal Kit	047-237-703
243.30	166/259 kN	Seal Kit	047-237-704
243.35	245/373 kN	Seal Kit	047-237-705
243.39*	1523/506 kN	Seal Kit	047-237-712
243.40	298/507 kN	Seal Kit	047-237-706
243.45	455/664 kN	Seal Kit	047-237-707
243.60	663/1035 kN	Seal Kit	047-237-708
243.65*	663/1035 kN	Seal Kit	047-237-713
243.70	982/1491 kN	Seal Kit	047-237-709
243.80*	1364/2027 kN	Seal Kit	047-237-710
243.90*	2250/2727 kN	Seal Kit	047-237-711

MTS Maintenance Parts

Actuator Parts

244.xx Seal Kits (Load Frames)

MODEL	RATING	DESCRIPTION	PART NUMBER
244.11	15 kN/3.3 kip	Seal Kit	008-710-067
244.12	25 kN/5.5 kip	Seal Kit	008-710-068
244.21	50 kN/11 kip	Seal Kit	008-710-069
244.22	100 kN/22 kip	Seal Kit	008-710-070
244.23	160 kN/35 kip	Seal Kit	008-710-071
244.31	250 kN/55 kip	Seal Kit	008-710-072
244.41	500 kN/110 kip	Seal Kit	008-710-073
244.51	1000 kN/220 kip	Seal Kit	008-710-074

Estimated weight of each is no more than 2 kg.

Each kit contains the parts needed to replace all existing seals.

244.xx G1 Seal Kits (Structural Applications)

MODEL	RATING	DESCRIPTION	PART NUMBER
244.11	15 kN/3.3 kip	G1 Seal Kit	008-710-046
244.12	25 kN/5.5 kip	G1 Seal Kit	008-710-047
244.20	68 kN/15 kip	G1 Seal Kit	008-710-075
244.21	50 kN/11 kip	G1 Seal Kit	008-710-048
244.22	100 kN/22 kip	G1 Seal Kit	008-710-049
244.23	160 kN/35 kip	G1 Seal Kit	008-710-050
244.31	250 kN/55 kip	G1 Seal Kit	008-710-051
244.41	500 kN/110 kip	G1 Seal Kit	008-710-052
244.51	1000 kN/220 kip	G1 Seal Kit	008-710-053
O-Ring for Servovalve Manifold			010-010-710

244.xx G2 Seal Kits (Structural Applications)

MODEL	RATING	DESCRIPTION	PART NUMBER
244.11	15 kN/3.3 kip	G2 Seal Kit	056-079-380
244.12	25 kN/5.5 kip	G2 Seal Kit	056-079-381
244.20	68 kN/15 kip	G2 Seal Kit	056-079-382
244.21	50 kN/11 kip	G2 Seal Kit	056-079-383
244.22	100 kN/22 kip	G2 Seal Kit	056-079-384
244.23	160 kN/35 kip	G2 Seal Kit	056-079-385
244.31	250 kN/55 kip	G2 Seal Kit	056-079-386
244.41	500 kN/110 kip	G2 Seal Kit	056-079-387
244.51	1000 kN/220 kip	G2 Seal Kit	056-079-388

244.xx Seal Insertion Tools

MODEL	RATING	DESCRIPTION	PART NUMBER
244.1x	45 mm/1.75 in	Rod Seal Insertion Tool	046-591-901
244.2x	70 mm/2.75 in	Rod Seal Insertion Tool	046-592-001
244.31	95 mm/3.75 in	Rod Seal Insertion Tool	046-592-201
244.41	135 mm/5.25 in	Rod Seal Insertion Tool	046-592-101
244.51	150 mm/6.00 in	Rod Seal Insertion Tool	046-592-301

Estimated weight of seal insertion tool is 5 kg.

MTS Maintenance Parts

Actuator Parts

245.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
245.XX	All sizes	Seal Kit	038-736-601

247.xx Seal Kits (Structural Applications)

MODEL	RATING (TENSION/ COMPRESSION)	DESCRIPTION	PART NUMBER
247.11	15/47 kN	Seal Kit	008-710-054
247.12	25 /59 kN	Seal Kit	008-710-055
247.21	50/130 kN	Seal Kit	008-710-056
247.22*	100/178 kN	Seal Kit	008-710-057
247.23*	159/246 kN	Seal Kit	008-710-058
247.31	250/405 kN	Seal Kit	008-710-059
247.62*	159/373 kN	Seal Kit	008-710-061

*Estimated weight of each seal kit is no more than 2 kg.
Each kit contains the parts needed to replace all existing seals.*

248.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
248.0X	All sizes	Seal Kit	040-262-601
		Seal Kit with static support	040-262-701
248.1X	All sizes	Seal Kit	040-262-801
		Seal Kit with static support	040-262-901
248.2X	All sizes	Seal Kit	040-263-001
		Seal Kit with static support	044-711-801

MTS Maintenance Parts

Actuator Parts

215.xx Seal Kits

MODEL	RATING	ACTUATOR ASSEMBLY #	PART NUMBER
215.35A	565 N.m Rotary Actuator	327308-01	033-283-101
215.41A	1130 N.m Actuator	327309-01	033-305-001
215.42A	2260 N.m Rotary Actuator	327310-01	033-311-901
215.45A	5650 N.m Rotary Actuator	327311-01	033-282-801
215.51A	11300 N.m Rotary Actuator	32731201	033-276-801
215.32B	226 N.m Rotary Actuator	362540-01	036-371-601 
215.35B	565 N.m Rotary Actuator	362541-01	036-371-601 
215.41B	1130 N.m Actuator	362542-01	036-371-601 
215.42B	2260 N.m Rotary Actuator	362543-01	036-371-601 
215.45B	5650 N.m Rotary Actuator	362544-01	036-254-701
215.51B	11300 N.m Rotary Actuator	36254501	036-254-701
215.32B	226 N.m Rotary Actuator	443061-01	036-371-601 
215.35B	565 N.m Rotary Actuator	443062-01	036-371-601 
215.41B	1130 N.m Actuator	443063-01	036-371-601 
215.42B	2260 N.m Rotary Actuator	444064-01	036-371-601 
215.45B	5650 N.m Rotary Actuator	444576-01	044-527-201 
215.51B	11300 N.m Rotary Actuator	44457701	044-527-201 
215.32C	226 N.m Rotary Actuator	474661-01	047-917-101
215.35C	565 N.m Rotary Actuator	474660-01	036-371-601 
215.41C	1130 N.m Actuator	475286-01	036-371-601 
215.42C	2260 N.m Rotary Actuator	475353-01	036-371-601 

Swivel Head/Base Spares

MODEL	RATING	DESCRIPTION	PART NUMBER
249.12	25 kN/5.5 kip	Spares Package	051-852-401
249.20	66 kN/15 kip	Spares Package	051-852-402
249.23	155 kN/35 kip	Spares Package	051-852-403
249.31	250 kN/55 kip	Spares Package	051-852-404
249.41	500 kN/100 kip	Spares Package	051-852-405
249.42	750 kN/165 kip	Spares Package	051-852-406
249.42S	750 kN/165 kip	Spares Package	051-852-407
249.51	1000 kN/220 kip	Spares Package	051-852-408

Estimated weight of each is no more than 2 kg.
Each kit contains the parts needed to replace all existing seals.

MTS Maintenance Parts

Load Frame Parts

LVDT

DESCRIPTION	PART NUMBER
Transducer LVDT/2 in stroke	039-075-101
Transducer LVDT/4 in stroke	039-075-102
Transducer LVDT/6 in stroke	039-075-103
Transducer LVDT/8 in stroke	039-075-111
Transducer LVDT/10 in stroke	039-075-104
Transducer LVDT/14 in stroke	039-075-105
Transducer LVDT/20 in stroke	039-075-106
Transducer LVDT/2 in stroke Transformer length 3.75 in	006-839-901
Transducer LVDT/2 in stroke Transformer length 4.63 in	006-839-902
Transducer LVDT/4 in stroke	006-839-903
Transducer LVDT/6 in stroke	006-839-904
Transducer LVDT/8 in stroke	Contact MTS
Transducer LVDT/10 in stroke	Contact MTS
252/253/256 Servovalve-LVDT	043-229-601

Estimated weight of each is no more than 2 kg.

Each kit contains the parts needed to replace all existing seals.

Spares Kits Description

MTS INSIGHT	PART NUMBER
Spares Kit - Insight 1K/2K	100-187-792
Spares Kit - Insight 5K/10K	100-187-793
Spares Kit - Insight 30K/50K	100-187-794
Spares Kit - Insight 100K	100-187-795
Spares Kit - Insight 150K	100-187-796
Spares Kit - Insight 200K	100-187-797
Spares Kit - Insight 300K	100-187-798
Spares Kit - Insight 50K Wide	100-187-799
Spares Kit - XLT Extensometer 100 R Std	100-187-800
Spares Kit - XLT Extensometer 100 R Ext	100-187-801

311.xx Seal Kits

MTS INSIGHT	PART NUMBER	PART NUMBER
311.11 & .21	Lift Seal Kit	032-853-001
	Lock Seal Kit	032-853-301
311.31 & .41	Lift Seal Kit	032-853-101
	Lock Seal Kit (new style)	032-853-701

312.xx Seal Kits

MTS INSIGHT	PART NUMBER	PART NUMBER
312.21	Lift Seal Kit	032-853-201
	Lock Seal Kit	032-853-301
312.31	Lift Seal Kit	032-853-401
	Lock Seal Kit	032-853-501
312.41	Lift Seal Kit	032-853-201
	Lock Seal Kit	032-853-701

Estimated weight of each is no more than 3 kg.

Each kit contains the parts needed to replace the existing seals.

MTS Maintenance Parts

Load Frame Parts

318.xx Seal Kits (810 Series)

MODEL	DESCRIPTION	PART NUMBER
318.10	Lift Seal Kit	040-602-101
	Lock Seal Kit	040-601-801
	25 kN/5.5 kip Actuator Seal Kit	038-171-201
	50 kN/11 kip Actuator Seal Kit	038-171-101
	100 kN/22 kip Actuator Seal Kit	038-171-001
318.25 / 319.25	Lift Seal Kit	040-602-101
	Lock Seal Kit	040-601-901
	100 kN/22 kip Actuator Seal Kit	038-171-401
	250 kN/55 kip Actuator Seal Kit	038-171-301
318.50	Lift Seal Kit	040-602-201
	Lock Seal Kit	040-602-001
	250 kN/55 kip Actuator Seal Kit	038-171-601
	500 kN/110 kip Actuator Seal Kit	038-171-501

Estimated weight of each is no more than 10 kg.

Each kit contains the parts needed to replace the existing seals.

322.xx Seal Kits

MODEL	DESCRIPTION	PART NUMBER
322.21	Lift Seal Kit	032-853-801
	Lock Seal Kit	032-853-901
322.31	Lift Seal Kit	032-853-401

370.xx Seal Kits (MTS Landmark® Systems)

MODEL	DESCRIPTION	PART NUMBER
370.02	Lift Seal Kit (load frames manufactured after January 2010)	100-222-825
	Lift Seal Kit (load frames manufactured before January 2010)	Contact MTS
	Lock Seal Kit	100-216-447
	15 kN / 25 kN Linear Actuator Seal Kit	100-215-462
370.10	Lift Seal Kit (load frames manufactured after January 2010)	100-222-825
	Lift Seal Kit (load frames manufactured before January 2010)	Contact MTS
	Lock Seal Kit	100-216-447
	15 kN standard and hydrostatic bearing	100-195-267
	25 kN standard and hydrostatic bearing	100-195-268
	50 kN standard and hydrostatic bearing	100-195-271
	67 kN standard and hydrostatic bearing (piston hub seal N/A)	100-195-273
	100 kN standard and hydrostatic bearing	100-195-274
370.25	Lift Seal Kit (load frames manufactured after January 2010)	100-222-825
	Lift Seal Kit (load frames manufactured before January 2010)	Contact MTS
	Lock Seal Kit	100-216-448
	100 kN standard and hydrostatic bearing	100-195-274
	250 kN standard and hydrostatic bearing	100-195-275
370.50	Lift Seal Kit (load frames manufactured after January 2010)	100-222-825
	Lift Seal Kit (load frames manufactured before January 2010)	Contact MTS
	Lock Seal Kit	100-216-449
	250 kN standard and hydrostatic bearing	100-195-275
	500 kN standard and hydrostatic bearing	100-195-276

370.02 Seal Kits (Bionix Servohydraulic System)

MODEL	DESCRIPTION	PART NUMBER
370.02	Lift Seal Kit (load frames manufactured after January 2010)	100-222-825
	Lift Seal Kit (load frames manufactured before January 2010)	Contact MTS
	Lock Seal Kit	100-216-447
	15 kN / 25 kN Linear Actuator Seal Kit	100-215-462
	Seal work on rotary actuator not recommended in the field	

MTS Maintenance Parts

Cables

Cables for FlexTest 40, 60, 100, & 200 Controllers

DESCRIPTION	CONNECTOR	LENGTH	CONTROLLER	PART NUMBER
Adapter, RJ50-D15, 4-Wire		0.3 m (1 ft)	FlexTest 40	100-187-439
Adapter, RJ50-D15, 8-Wire		0.3 m (1 ft)	FlexTest 40	057-303-001
Adapter, RJ50-D15, 4-Wire		1.5 m (5 ft)	FlexTest 60; 100; 200	057-294-001
Adapter, RJ50-D15, 8-Wire		1.5 m (5 ft)	FlexTest 60; 100; 200	057-303-002
ADT 120 Degrees		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-309-104
ADT 120 Degrees		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-309-109
ADT 300 Degrees		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-309-204
ADT 300 Degrees		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-309-209
Delta P Cell		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-144-504
Delta P Cell		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-144-509
Extensometer	PT	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-230-704
Extensometer	PT	15 m (50 ft)	FlexTest 40; 60; 100; 200	057-230-709
Load Cell	PT	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-099-204
Load Cell	PT	15 m (50 ft)	FlexTest 40; 60; 100; 200	057-099-209
Load Cell	MS	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-221-604
Load Cell	MS	15 m (50 ft)	FlexTest 40; 60; 100; 200	057-221-609
HPU, Adapter Cable (Dual Control, Single Pump) 24 V DC		5 m (1.5 ft)	FlexTest 40; 60; 100; 200	005-407-801
HSM Adapter, D9P-CPC4S		0.2 m (10 in)	FlexTest 40; 60; 100; 200	100-175-734
HSM Service Manifold 290/3/4, On/Off, 115 V AC		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	050-053-904
HSM Service Manifold 290/3/4, On/Off, 115 V AC		15 m (50 ft)	FlexTest 40; 60; 100; 200	050-053-909
HSM Service Manifold 290/3/4, Hi/Low, 115 V AC		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	039-718-304
HSM Service Manifold 290/3/4, Hi/Low, 115 V AC		15 m (50 ft)	FlexTest 40; 60; 100; 200	039-718-309
HSM Service Manifold 290/3/4, Hi/Low, 24 V DC		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-119-204
HSM Service Manifold 290/3/4, Hi/Low, 24 V DC		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-119-209
HSM Service Manifold 298, On/Off, 24 V DC		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-172-904
HSM Service Manifold 298, On/Off, 24 V DC		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-172-909
HSM, Proportional Valve 298.12		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	056-126-404
HSM, Proportional Valve 298.12		15 m (50 ft)	FlexTest 40; 60; 100; 200	056-126-409
LVDT		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-099-104
LVDT		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-099-109
Pressure - Sensotec		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-238-804
Pressure - Sensotec		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-238-809
RJ50-RJ50 for 494.16/25/26/47		30 m (100 ft)	FlexTest 40; 60; 100; 200	057-267-511
RJ50-RJ50 for 494.45/46/49		30 m (100 ft)	FlexTest 40; 60; 100; 200	057-144-111
RVDT		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-255-604
RVDT		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-255-609
Servo Valve Adapter, RJ50-D15		0.3 m (1 ft)	FlexTest 40; 60; 100; 200	100-108-342
Servo Valve LVDT, 3-Stage		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-204-104
Servo Valve LVDT, 3-Stage		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-204-109
Servo Valve, 252, 2-Stage Dual, in phase		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-255-504
Servo Valve, 252, 2-Stage Dual, in phase		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-255-509
Servo Valve, 252, 2-Stage Dual, out of phase		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-255-404
Servo Valve, 252, 2-Stage Dual, out of phase		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-255-409
Servo Valve, 252,/6 2-Stage Single		7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-193-704
Servo Valve, 252/6, 2-Stage Single		15 m (50 ft)	FlexTest 40; 60; 100; 200	057-193-709
Torque Cell	MS	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-230-504
Torque Cell	MS	15 m (50 ft)	FlexTest 40; 60; 100; 200	057-230-509
Torque Cell	PT	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-255-304
Torque Cell	PT	15 m (50 ft)	FlexTest 40; 60; 100; 200	057-255-309

MTS Maintenance Parts

Cables

TEDS Cables for FlexTest 40, 60, 100, 200 Controllers

DESCRIPTION	CONNECTOR	LENGTH	CONTROLLER	PART NUMBER
ADT 120 Degrees, TEDS Module	MS	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-702
ADT 300 Degrees, TEDS Module	MS	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-902
Delta P Cell, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-102
Extensometer, 632.11x-90 (3 units, 1 output), 634.28, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-610
Extensometer, 632.11x-90 (2 units, 1 output), 632.12, 632.23, 632.24, 632.92H-05, 634.12, 634.25, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-608
Extensometer, 632.03, 632.13, 632.18, 632.20, 632.65, 632.85E Axial (001), 632.90, 632.92H-03, 632.92H-04, 634.11, 634.33, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-606
Extensometer, 632.17. (001), 632.85F Axial (001), 634.31, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-605
Extensometer, 632.02, 632.05, 632.06H-2x, 632.06H-3x (002) thru (006), 632.11x-90 (single unit), 632.19, 632.26-2x, 632.26-4x, 632.27-2x, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-613
Extensometer, 632.85 (002), TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-611
Extensometer, 632.06H-3x (001), 632.17 (002), 632.27-3x, 632.29, 632.53, 632.54, 632.85 Transverse, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-604
Extensometer, 632.26-3x TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-603
Extensometer, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-602
Extensometer, TEDS Module	Amphenol	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-302
Load Cell, 661.09A/B-2X, 661.19E/F-02/04, 661.20E/F-01/02/03, 661.36C/D-03, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-513
Load Cell, 661.11A/B-0X, 661.18E/F-01/02, 661.19E/F-01, 661.22C/D, 661.23E/F, 661.31E/F, 661.34E/F, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-511
Load Cell, A/T, 662.10A/B, Axial, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-513
Load Cell, A/T, 662.10A/B, Torsional, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-813
Load Cell, A/T, 662.20C/D-01, Axial, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-514
Load Cell, A/T, 662.20C/D-01, Torsional, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-814
Load Cell, A/T, 662.20C/D-03/04/05/06, Axial, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-512
Load Cell, 662.20C/D-03/06, Torsional, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-804
Load Cell, A/T, 662.20C/D-04/05, Torsional, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-803
Load Cell, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-502
Jumper Plug Kit		1.5 m (5 ft)	FlexTest 40; 60; 100; 200	100-182-496
LVDT, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-202
Pressure Sensor, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-274-002
RVDT, TEDS Module	MS	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-272-502
System Cable, ADT 120 degrees with TEDS	JT to D9/D15	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-316-204
System Cable, ADT 300 degrees with TEDS	JT to D9/D15	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-316-304
System Cable, Transducer with TEDS	JT to RJ50	7.5 m (25 ft)	FlexTest 40; 60; 100; 200	057-241-404
Tempo G, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-274-202
Torque, TEDS Module	MS	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-902
Torque, TEDS Module	PT	1.5 m (5 ft)	FlexTest 40; 60; 100; 200	057-273-802

MTS Maintenance Parts

Cables

Cables for FlexTest SE & GT Controllers

DESCRIPTION	CONNECTOR	LENGTH	CONTROLLER	PART NUMBER
ADT 120 Degrees		7.5 m (25 ft)	FlexTest SE; FlexTest GT	056-234-804
ADT 120 Degrees		15 m (50 ft)	FlexTest SE; FlexTest GT	056-234-809
ADT 300 Degrees		7.5 m (25 ft)	FlexTest SE; FlexTest GT	056-253-104
ADT 300 Degrees		15 m (50 ft)	FlexTest SE; FlexTest GT	056-253-109
Delta P Cell – Stability		7.5 m (25 ft)	FlexTest SE; FlexTest GT	047-927-604
Delta P Cell – Stability		15 m (50 ft)	FlexTest SE; FlexTest GT	047-927-609
Extensometer	PT	7.5 m (25 ft)	FlexTest SE; FlexTest GT	050-120-004
Extensometer	PT	15 m (50 ft)	FlexTest SE; FlexTest GT	050-120-009
Load Cell	PT	7.5 m (25 ft)	FlexTest SE; FlexTest GT	046-440-204
Load Cell	PT	15 m (50 ft)	FlexTest SE; FlexTest GT	046-440-209
Load Cell	MS	7.5 m (25 ft)	FlexTest SE; FlexTest GT	046-440-604
Load Cell	MS	15 m (50 ft)	FlexTest SE; FlexTest GT	046-440-609
HSM Service Manifold; 290/3/4 On/Off, 115 V AC		7.5 m (25 ft)	FlexTest SE; FlexTest GT	050-053-904
HSM Service Manifold; 290/3/4 On/Off, 115 V AC		15 m (50 ft)	FlexTest SE; FlexTest GT	050-053-909
HSM Service Manifold 290/3/4, Hi/Low, 115 V AC		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-718-304
HSM Service Manifold 290/3/4, Hi/Low, 115 V AC		15 m (50 ft)	FlexTest SE; FlexTest GT	039-718-309
HSM Service Manifold 290/3/4, Hi/Low, 24 V DC		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-701-404
HSM Service Manifold 290/3/4, Hi/Low, 24 V DC		15 m (50 ft)	FlexTest SE; FlexTest GT	039-701-409
HSM Service Manifold 298, On/Off, 24 V DC		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-701-504
HSM Service Manifold 298, On/Off, 24 V DC		15 m (50 ft)	FlexTest SE; FlexTest GT	039-701-509
HSM, Proportional Valve 298.12		7.5 m (25 ft)	FlexTest SE; FlexTest GT	056-126-404
HSM, Proportional Valve 298.12		15 m (50 ft)	FlexTest SE; FlexTest GT	056-126-409
LVDT		7.5 m (25 ft)	FlexTest SE; FlexTest GT	046-440-304
LVDT		15 m (50 ft)	FlexTest SE; FlexTest GT	046-440-309
Pressure - Sensotec		7.5 m (25 ft)	FlexTest SE; FlexTest GT	056-253-004
Pressure - Sensotec		15 m (50 ft)	FlexTest SE; FlexTest GT	056-253-009
Servovalve & Valve LVDT, Y Adapter		5 m (1.5 ft)	FlexTest SE; FlexTest GT	039-710-501
Servovalve LVDT, 3-Stage		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-708-604
Servovalve LVDT, 3-Stage		15 m (50 ft)	FlexTest SE; FlexTest GT	039-708-609
Servovalve, 252, 2-Stage Dual, in phase		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-708-504
Servovalve, 252, 2-Stage Dual, in phase		15 m (50 ft)	FlexTest SE; FlexTest GT	039-708-509
Servovalve, 252, 2-Stage Dual, out of phase		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-708-404
Servovalve, 252, 2-Stage Dual, out of phase		15 m (50 ft)	FlexTest SE; FlexTest GT	039-708-409
Servovalve, 252, 2-Stage Single		7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-708-304
Servovalve, 252, 2-Stage Single		15 m (50 ft)	FlexTest SE; FlexTest GT	039-708-309
Servovalve, 256, 3-Stage		7.5 m (25 ft)	FlexTest SE; FlexTest GT	055-439-604
Servovalve, 256, 3-Stage		15 m (50 ft)	FlexTest SE; FlexTest GT	055-439-609
Torque Cell	PT	7.5 m (25 ft)	FlexTest SE; FlexTest GT	039-704-204
Torque Cell	PT	15 m (50 ft)	FlexTest SE; FlexTest GT	039-704-209
Torque Cell	MS	7.5 m (25 ft)	FlexTest SE; FlexTest GT	051-649-504
Torque Cell	MS	15 m (50 ft)	FlexTest SE; FlexTest GT	051-649-509

MTS Maintenance Parts

Cables & Replacement Lamps

Cables for Model 407 Controllers

DESCRIPTION	CONNECTOR	LENGTH	CONTROLLER	PART NUMBER
E-STOP to HPU Control		7.5 m (25 ft)	Model 407	049-474-704
E-STOP to HPU Control		9.1 m (30 ft)	Model 407	049-474-705
E-STOP, Remote		7.5 m (25 ft)	Model 407	049-474-604
Extensometer	PT	7.5 m (25 ft)	Model 407	050-120-004
Load Cell	MS	7.5 m (25 ft)	Model 407	046-440-604
Load Cell	PT	7.5 m (25 ft)	Model 407	046-440-204
HPU Control 24 V DC		7.5 m (25 ft)	Model 407	039-708-704
HPU Control 24 V DC		9.1 m (30 ft)	Model 407	039-708-705
HSM Hi/Low, 115 V AC		7.5 m (25 ft)	Model 407	049-474-504
HSM Hi/Low, 24 V DC		7.5 m (25 ft)	Model 407	049-474-404
HSM Hi/Low, 24 V DC		9.1 m (30 ft)	Model 407	049-474-405
HSM On/Off Valve, 298.xx		7.5 m (25 ft)	Model 407	049-740-304
HSM, Proportional Valve 298.12		7.5 m (25 ft)	Model 407	049-474-304
LVDT		7.5 m (25 ft)	Model 407	046-440-304
LVDT		9.1 m (30 ft)	Model 407	046-440-305
Servovalve, 252 Dual		7.5 m (25 ft)	Model 407	046-440-404
Servovalve, 252 Dual		9.1 m (30 ft)	Model 407	046-440-405
Servovalve, 252 Single		7.5 m (25 ft)	Model 407	046-440-104
Servovalve, 252 Single		9.1 m (30 ft)	Model 407	046-440-105
System Ground Cable		7.5 m (25 ft)	Model 407	039-709-204
System Ground Cable		9.1 m (30 ft)	Model 407	039-709-205

Cables for FlexTest IIm Controllers

DESCRIPTION	CONNECTOR	LENGTH	CONTROLLER	PART NUMBER
Load Cell	MS	7.5 m (25 ft)	FlexTest IIm	046-440-604
Load Cell	PT	7.5 m (25 ft)	FlexTest IIm	046-440-204
LVDT		7.5 m (25 ft)	FlexTest IIm	046-440-304
Servovalve, 252 Dual		7.5 m (25 ft)	FlexTest IIm	046-440-404
Servovalve, 252 Single		7.5 m (25 ft)	FlexTest IIm	046-440-104

Replacement Lamps

DESCRIPTION	PART NUMBER
421.21	010-076-901
436.11 (Type 334)	010-076-901
Oven Lamp	011-322-312

Estimated weight of each is no more than 2 kg.

MTS Maintenance Parts

Adapter Cables

Adapter Cables for 3rd Party Equipment

DESCRIPTION	CONNECTOR	LENGTH	CONTROLLER	PART NUMBER
Instron Extensometer to Insight Controller Cable	MS to RJ50	3m (10 ft)	Insight	100-220-824
Instron 133X Xhead Interlock to 493.01 Cable	MS to CPC	15m (50 ft)	FlexTest SE; FlexTest GT	052-750-709
Instron 3340 Act Off/Low/High to 493.08 Cable	MS to CPC	7.5m (25 ft)	FlexTest SE; FlexTest GT	052-752-304
Instron 3340 Act Off/Low/High to 493.08 Cable	MS to CPC	15m (50 ft)	FlexTest SE; FlexTest GT	052-752-309
Instron 3340 Act Off/High to 493.08 Cable	MS to CPC	7.5m (25 ft)	FlexTest SE; FlexTest GT	052-751-204
Instron 3340 Act Off/High to 493.08 Cable	MS to CPC	15m (50 ft)	FlexTest SE; FlexTest GT	052-751-209
Instron 210 5/10/20 HPS to 493.08 Cable (Replacing 2150)	MS	7.5m (25 ft)	FlexTest SE; FlexTest GT	054-145-504
Instron 210 5/10/20 HPS to 493.08 Cable (Replacing 8500)	MS	7.5m (25 ft)	FlexTest SE; FlexTest GT	054-145-704
Instron Extensometer to MTS PT Cable Adapter	D to PT	.3m (1 ft)		100-140-617
MTS Extensometer to Instron 8500 Cable	D to PT	7.5m (25 ft)		049-751-304
Instron 133X Xhead Interlock to FT40/60/100/200 Cable	D to MS	7.5m (25 ft)	FlexTest SE; FlexTest GT	057-048-504
Instron 133X Xhead Interlock to FT40/60/100/200 Cable	D to MS	15m (50 ft)	FlexTest SE; FlexTest GT	057-048-509
Instron 133X & 851X SelfID Extensometer to FT40/60/100/200 Cable	D to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-335-604
Instron 133X & 851X SelfID Extensometer to FT40/60/100/200 Cable	D to RJ50	15m (50 ft)	FlexTest 40; 60; 100; 200	057-335-609
Instron 133X & 851X Non-SelfID Extensometer to FT40/60/100/200 Cable	MS to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-335-704
Instron 133X & 851X SelfID L/C to FT40/60/100/200 Cable	D to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-335-404
Instron 133X & 851X SelfID L/C to FT40/60/100/200 Cable	D to RJ50	15m (50 ft)	FlexTest 40; 60; 100; 200	057-335-409
Instron 133X Lebow L/C to FT40/60/100/200 Cable	MS to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-221-604
Instron 133X Lebow L/C to FT40/60/100/200 Cable	MS to RJ50	15m (50 ft)	FlexTest 40; 60; 100; 200	057-221-609
Instron 133X S/V (Single) to FT40/60/100/200 Cable	MS to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-193-704
Instron 133X S/V (Single) to FT40/60/100/200 Cable	MS to RJ50	15m (50 ft)	FlexTest 40; 60; 100; 200	057-193-709
Instron 133X S/V (Dual) to FT40/60/100/200 Cable	(2) MS to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-255-404
Instron 133X S/V (Dual) to FT40/60/100/200 Cable	(2) MS to RJ50	15m (50 ft)	FlexTest 40; 60; 100; 200	057-255-409
Instron 133X LVDT (2150 & 8500 w/ MS Connector) to FT40/60/100/200 Cable	MS to RJ50	7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-335-504
Instron 133X LVDT (2150 & 8500 w/ MS Connector) to FT40/60/100/200 Cable	MS to RJ50	15m (50 ft)	FlexTest 40; 60; 100; 200	057-335-509
Instron 133X & 851X System Ground to FT40/60/100/200 Cable		7.5m (25 ft)	FlexTest 40; 60; 100; 200	039-709-204
Instron 133X & 851X System Ground to FT40/60/100/200 Cable		15m (50 ft)	FlexTest 40; 60; 100; 200	039-709-209
Instron 2GPM 1500PSI HPS to FT40/60/100/200		7.5m (25 ft)	FlexTest 40; 60; 100; 200	054-872-204
Instron 2GPM 1500PSI HPS to FT40/60/100/200		15m (50 ft)	FlexTest 40; 60; 100; 200	054-872-209
Instron 347X Mod 718 HPS to FT40/60/100/200		7.5m (25 ft)	FlexTest 40; 60; 100; 200	054-872-004
Instron 347X Mod 718 HPS to FT40/60/100/200		15m (50 ft)	FlexTest 40; 60; 100; 200	054-872-009
Instron 851X LF Load Rel Pres Switch to MTS Load Release Interface Box Cable	CPC to DIN	7.5m (25 ft)		057-140-404
Instron 851X LF Load Rel Pres Switch to MTS Load Release Interface Box Cable	CPC to DIN	15m (50 ft)		057-140-409
Instron 851X LF Pres Red Man Opt to MTS Load Release Interface Box Cable	CPC to DIN	7.5m (25 ft)		055-859-804
Instron 851X LF Pres Red Man Opt to MTS Load Release Interface Box Cable	CPC to DIN	15m (50 ft)		055-859-809
Instron 851X LF Load Rel Sol to MTS Load Release Interface Box Cable	CPC to DIN	7.5m (25 ft)		055-859-704
Instron 851X LF Load Rel Sol to MTS Load Release Interface Box Cable	CPC to DIN	15m (50 ft)		055-859-709
Load Release Interface Box for Instron 851X LF - Floor Standing w/ 10' Cable		3m (10 ft)		053-045-301
Load Release Interface Box for Instron 851X LF - Rack Mounted w/ 15' Cable		4.5m (15 ft)		053-045-302
Instron 851X LF Load Unit Interlock Adapter Cable	CPC to DIN	1.5m (5 ft)		053-047-201
Instron 851X LF Load Unit Interlock Adapter Cable to FT40/60/100/200		7.5m (25 ft)	FlexTest 40; 60; 100; 200	056-126-504
Instron 851X LF Load Unit Interlock Adapter Cable to FT40/60/100/200		15m (50 ft)	FlexTest 40; 60; 100; 200	056-126-509
Instron 851X S/V (Single) to FT40/60/100/200 Cable		7.5m (25 ft)	FlexTest 40; 60; 100; 200	052-752-204
Instron 851X S/V (Single) to FT40/60/100/200 Cable		15m (50 ft)	FlexTest 40; 60; 100; 200	052-752-209
Instron 851X S/V (Dual) to FT40/60/100/200 Cable		7.5m (25 ft)	FlexTest 40; 60; 100; 200	039-708-404
Instron 851X S/V (Dual) to FT40/60/100/200 Cable		15m (50 ft)	FlexTest 40; 60; 100; 200	039-708-409
Instron 851X LVDT w/ Burndy Connector to FT40/60/100/200 Cable		7.5m (25 ft)	FlexTest 40; 60; 100; 200	057-335-804
Instron 851X LVDT w/ Burndy Connector to FT40/60/100/200 Cable		15m (50 ft)	FlexTest 40; 60; 100; 200	057-335-809

MTS Maintenance Parts

Grip Parts

641.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
641.35	100 kN	Seal Kit/Metric (Two required)	020-718-301
641.37	500 kN	Seal Kit/Metric (Two required)	020-718-701
641.35	100 kN	Seal Kit	036-792-301
641.36	250 kN	Seal Kit	038-753-801
641.37	500 kN	Seal Kit	036-809-201
641.38	1000 kN	Seal Kit	038-753-901

Metric grip sets have a knurled ring on the main body.

646.xx Seal Kits

MODEL	RATING	DESCRIPTION	PART NUMBER
646.10A	100 kN/22 kip	Seal Kit	038-400-501
646.10B	100 kN/22 kip	Seal Kit	042-739-101
646.25A	250 kN/55 kip	Seal Kit	038-400-701
646.25B	250 kN/55 kip	Seal Kit	042-739-201

647.xx Seal Kits (seals for upper and lower grip)

MODEL	RATING	DESCRIPTION	TEMPERATURE RANGE	PART NUMBER
647.02	25 kN/5.5 kip	Seal Kit	-17/+65°C (0/+150°F)	047-810-701
647.02	25 kN/5.5 kip	Replacement Springs		010-086-754
6647.02	25 kN/5.5 kip	Spring Anchors		038-529-001
647.02	25 kN/5.5 kip	Seal Kit (high temp)	-40/+175°C (-40/+350°F)	047-810-702
647.02B	25 kN/5.5 kip	Seal Kit	-40/+120°C (-40/+250°F)	056-143-501
647.02B	25 kN/5.5 kip	Seal Kit (high temp)	-40/+175°C (-40/+350°F)	056-143-502
647.10B	100 kN/22 kip	Seal Kit	-17/+65°C (0/+150°F)	047-810-901
647.10	100 kN/22 kip	Replacement Springs		010-086-752
647.10	100 kN/22 kip	Spring Anchors		038-310-501
647.10B	100 kN/22 kip	Seal Kit (high temp)	-40/+175°C (-40/+350°F)	047-810-902
647.25B	250 kN/55 kip	Seal Kit	-17/+65°C (0/+150°F)	047-810-801
647.25	250 kN/55 kip	Replacement Springs		010-086-753
647.25	250 kN/55 kip	Spring Anchors		038-310-501
647.25B	250 kN/55 kip	Seal Kit (high temp)	-40/+175°C (-40/+350°F)	047-810-802
647.50	500 kN/110 kip	Seal Kit	-17/+65°C (0/+150°F)	047-811-001
647.50	500 kN/110 kip	Replacement Springs		010-086-753
647.50	500 kN/110 kip	Spring Anchors		038-310-501
647.50	500 kN/110 kip	Seal Kit (high temp)	-40/+175°C (-40/+350°F)	047-811-002

680.xx Seal Kits

MODEL	DESCRIPTION	PART NUMBER
680.01	Seal Kit	040-985-901

Estimated weight of each is no more than 2 kg. Each kit contains parts needed to replace all existing seals.

Electromechanical Grip Parts

ADVANTAGE PNEUMATIC GRIP SPARES KIT	PART NUMBER
Spares Kit – Advantage Pneumatic 100 N/200 N	100-044-616
Spares Kit – Advantage Pneumatic Grip 2 kN	100-044-617
Spares Kit – Advantage Pneumatic 10 N	100-044-618
Spares Kit – Advantage Pneumatic 10 kN	100-044-619
ADVANTAGE WEDGE ACTION REPLACEMENT KIT	PART NUMBER
Replacement Kit Wedge Action 30 kN/50 kN	056-011-601
Replacement Kit Wedge Action 100 kN/150 kN	056-011-602

Grip Supplies and Environmental Simulation Parts

GRIP SUPPLY SPARES KIT	PART NUMBER
Kit-Spares, Grip Supply 685.10/685.22	100-109-989
653.XX FURNACE SERVICE KIT	PART NUMBER
Heater Element Service Kit – 653.01 Furnace	056-543-201
Heater Element Service Kit – 653.02 &.03 Furnace	056-543-202
Heater Element Service Kit – 653.04 Furnace	056-543-203

MTS Maintenance Parts

Hydraulic Hoses

Pressure and Return Hoses

FITTING SIZE	DIAMETER IN INCHES	LENGTH (FT/M)	PRESSURE RATING	PRESSURE HOSE PART NUMBER	RETURN RATING	RETURN HOSE PART NUMBER
-6	0.375	5/1.5	3000 psi	006-890-305	2000 psi	006-890-405
-6	0.375	10/3.048	3000 psi	006-890-310	2000 psi	006-890-410
-6	0.375	15/4.5	3000 psi	006-890-311	2000 psi	006-890-411
-6	0.375	20/6.096	3000 psi	006-890-312	2000 psi	006-890-412
-6	0.375	25/7.6	3000 psi	006-890-313	2000 psi	006-890-413
-6	0.375	30/9.144	3000 psi	006-890-314	2000 psi	006-890-414
-8	0.5	5/1.5	3000 psi	006-890-605	2000 psi	006-890-705
-8	0.5	10/3.048	3000 psi	006-890-610	2000 psi	006-890-710
-8	0.5	15/4.5	3000 psi	006-890-611	2000 psi	006-890-711
-8	0.5	20/6.096	3000 psi	006-890-612	2000 psi	006-890-712
-8	0.5	25/7.6	3000 psi	006-890-613	2000 psi	006-890-713
-8	0.5	30/9.144	3000 psi	006-890-614	2000 psi	006-890-714
-12	0.75	5/1.5	4000 psi	006-890-905	2000 psi	006-891-005
-12	0.75	10/3.048	4000 psi	006-890-910	2000 psi	006-891-010
-12	0.75	15/4.5	4000 psi	006-890-911	2000 psi	006-891-011
-12	0.75	20/6.096	4000 psi	006-890-912	2000 psi	006-891-012
-12	0.75	25/7.6	4000 psi	006-890-913	2000 psi	006-891-013
-12	0.75	30/9.144	4000 psi	006-890-914	2000 psi	006-891-014
-16	1	5/1.5	4000 psi	006-891-205	2000 psi	006-891-305
-16	1	10/3.048	4000 psi	006-891-210	2000 psi	006-891-310
-16	1	15/4.5	4000 psi	006-891-211	2000 psi	006-891-311
-16	1	20/6.096	4000 psi	006-891-212	2000 psi	006-891-312
-16	1	25/7.6	4000 psi	006-891-213	2000 psi	006-891-313
-16	1	30/9.144	4000 psi	006-891-214	2000 psi	006-891-314
-20	1.25	5/1.5	3000 psi	006-891-505	2000 psi	006-891-605
-20	1.25	10/3.048	3000 psi	006-891-510	2000 psi	006-891-610
-20	1.25	15/4.5	3000 psi	006-891-511	2000 psi	006-891-611
-20	1.25	20/6.096	3000 psi	006-891-512	2000 psi	006-891-612
-20	1.25	25/7.6	3000 psi	006-891-513	2000 psi	006-891-613
-20	1.25	30/9.144	3000 psi	006-891-514	2000 psi	006-891-614

Drain Hoses

FITTING SIZE	DIAMETER IN INCHES	LENGTH (FT/M)	PRESSURE RATING	PART NUMBER
-6	0.375	5/1.5	100 psi	006-890-505
-6	0.375	10/3.048	100 psi	006-890-510
-6	0.375	15/4.5	100 psi	006-890-511
-6	0.375	20/6.096	100 psi	006-890-512
-6	0.375	25/7.5	100 psi	006-890-513
-6	0.375	30/9.144	100 psi	006-890-514
-8	0.5	5/1.5	100 psi	006-890-805
-8	0.5	10/3.048	100 psi	006-890-810
-8	0.5	15/4.5	100 psi	006-890-811
-8	0.5	20/6.096	100 psi	006-890-812
-8	0.5	25/7.5	100 psi	006-890-813
-8	0.5	30/9.144	100 psi	006-890-814

Actual hoses used depend on specific equipment in the installation.

To determine desired hose size: 1) Determine servovalve and actuator manifold; 2) determine HSM;

3) determine HPU; and 4) then refer to those sections for proper fitting sizes to complete your hose selection.

MTS Maintenance Parts

Hydraulic Power Unit Parts

505.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
505.07	High Pressure Filter Kit (Adds High Pressure Filtration)	100-008-737
	High Pressure Element & Seal Kit	100-030-194
	High Pressure Element	100-030-007
	Return Element & Seal Kit	100-030-195
	Return Element	100-009-495
505.11	High Pressure Filter Kit (Adds High Pressure Filtration)	100-008-737
	High Pressure Element & Seal Kit	100-030-194
	High Pressure Element	100-030-007
	Return Element & Seal Kit	100-030-195
	Return Element	100-009-495
505.20/.30	High Pressure Filter Kit for 505.20 & .30	054-728-801
	Return Element	010-053-349
	High Pressure Element for 505.20 & .30	011-970-652
505.60	Return Element**	**MFG prior to 2007 100-015-519
	Return Element w/ O-ring Seal**	**MFG prior to 2007 100-030-199
	Return Element w/ O-ring Seal***	***MFG in 2007 & later 100-029-989
	High Pressure Element (two required)	010-053-305
505.90	Return Element**	**MFG prior to 2007 100-015-519
	Return Element w/ O-ring Seal**	**MFG prior to 2007 100-030-199
	Return Element w/ O-ring Seal***	***MFG in 2007 & later 100-029-989
	High Pressure Element (two required)	010-053-305
505.120	Return Element**	**MFG prior to 2007 100-015-519
	Return Element w/ O-ring Seal**	**MFG prior to 2007 100-030-199
	Return Element w/ O-ring Seal***	***MFG in 2007 & later 100-029-989
	High Pressure Element (four required)	010-053-305
505.120	Return Element**	**MFG prior to 2007 100-015-519
	Return Element w/ O-ring Seal**	**MFG prior to 2007 100-030-199
	Return Element w/ O-ring Seal***	***MFG in 2007 & later 100-029-989
	High Pressure Element (four required)	010-053-305
505.180	Return Element**	**MFG prior to 2007 100-015-519
	Return Element w/ O-ring Seal**	**MFG prior to 2007 100-030-199
	Return Element w/ O-ring Seal***	***MFG in 2007 & later 100-029-989
	High Pressure Element (four required)	010-053-305

**Machines before January 2007: Clamp style filter cover. Filter part number: 100-015-519 (approx. 24.25" long).

- Kit (filter & seals): 100-030-199

***Machines from January 2007 on: Bolt on filter cover. Filter part number: 100-029-989 (includes seals) (approx. 18-25" long).

506.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
506.00-01	Filter Element-Main Pressure	010-088-309
506.00-02	Filter Element-Inlet Pressure **	011-068-418
All serial numbers	Filter Kit-Level 1	038-752-401
	Filter Kit-Level 2	038-094-801
	Filter Kit-Level 4	038-750-901
506.01 Serial numbers 100-470, 472-474, 481 and 482	Filter Element-Main Pressure	010-088-334
	Filter Element-Inlet Pressure **	011-129-901
	Filter Kit-Level 1	035-843-601
	Filter Kit-Level 2	031-875-301
	Filter Kit-Level 4	031-875-401
506.01 Serial numbers 471, 475-480, 484 and up	Filter Element-Main Pressure	010-088-328
	Filter Element-Inlet Pressure **	011-129-901
	Filter Kit-Level 1	038-750-301
	Filter Kit-Level 2	038-094-901
	Filter Kit-Level 4	038-751-001

*** Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 10 kg on kits.

Level 1 kits contain filters and seals to perform one maintenance requirement.

Level 2 kits contain filters and seals to perform two maintenance requirements, one cleanable filter.

Level 4 kits contain filters and seals to perform four maintenance requirements, one cleanable filter.

MTS Maintenance Parts

Hydraulic Power Unit Parts

506.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
506.02 Serial numbers 100-545 4 in Element, Silver Bowl	Filter Element-Main Pressure	010-088-334
	Filter Element-Inlet Pressure **	011-129-901
	Filter Kit-Level 1	035-843-601
	Filter Kit-Level 2	031-875-301
	Filter Kit-Level 4	031-875-401
506.02 Serial numbers 100-545 8 in Element, Silver Bowl	Filter Element-Main Pressure	010-088-354
	Filter Element-Inlet Pressure **	011-129-901
	Filter Kit-Level 1	035-843-701
	Filter Kit-Level 2	032-845-001
	Filter Kit-Level 4	032-845-101
506.02 Serial numbers 546 and up 4 in Element, Black Bowl	Filter Element-Main Pressure	010-088-328 
	Filter Element-Inlet Pressure **	011-129-901
	Filter Kit-Level 1	038-750-301 
	Filter Kit-Level 2	038-094-901
	Filter Kit-Level 4	038-751-001
506.03 All serial numbers	Filter Element-Main Pressure	010-088-354
	Filter Element-Inlet Pressure **	011-129-901
	Filter Kit-Level 1	035-843-701
	Filter Kit-Level 2	032-845-001
	Filter Kit-Level 4	032-845-101
506.10/.20 All serial numbers 4 in Element	Filter Element-Main Pressure	010-088-328 
	Filter Element-Inlet Pressure **	011-068-406
	Filter Kit-Level 1	035-843-801
	Filter Kit-Level 2	031-875-701
	Filter Kit-Level 4	036-911-001
506.10/.20 All serial numbers 8 in Element	Filter Element-Main Pressure	010-088-323 
	Filter Element-Inlet Pressure **	011-068-406
	Filter Kit-Level 1	038-514-001
	Filter Kit-Level 2	036-911-101
	Filter Kit-Level 4	031-875-801
506.22 All serial numbers	Filter Element-Inlet Pressure **	011-068-406
	Filter Kit-Level 1	038-514-001
	Filter Kit-Level 2	036-911-101
	Filter Kit-Level 4	031-875-801

** Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 10 kg on kits.

Level 1 kits contain filters and seals to perform one maintenance requirement.

Level 2 kits contain filters and seals to perform two maintenance requirements, one cleanable filter.

Level 4 kits contain filters and seals to perform four maintenance requirements, one cleanable filter.

MTS Maintenance Parts

Hydraulic Power Unit Parts

506.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
506.41 Serial numbers 100-186 and 193	Filter Element-Main Pressure (two required)	010-053-305
	Filter Element-Inlet Pressure	011-068-403
	Filter Element-Bypass	010-088-323
	Filter Kit-Level 1	035-843-901
	Filter Kit-Level 2	031-875-901
	Filter Kit-Level 4	031-876-001
506.41 Serial numbers 187-192, 194 and up	Filter Element-Main Pressure (two required)	010-053-305
	Filter Element-Inlet Pressure	011-402-819
	Filter Kit-Level 1	038-750-401
	Filter Kit-Level 2	038-105-201
	Filter Kit-Level 4	038-751-101
506.51 Serial numbers 100-143	Filter Element-Main Pressure (two required)	010-053-305
	Filter Element-Inlet Pressure	011-068-403
	Filter Element-Bypass	010-088-323
	Filter Kit-Level 1	035-843-901
	Filter Kit-Level 2	031-875-901
	Filter Kit-Level 4	031-876-001
506.51 Serial numbers 144 and up	Filter Element-Main Pressure (two required)	010-053-305
	Filter Element-Inlet Pressure	011-402-819
	Filter Kit-Level 1	038-750-401
	Filter Kit-Level 2	038-105-201
	Filter Kit-Level 4	038-751-101
Filter Element-Main Pressure	Filter Element-Main Pressure	011-951-524
	Filter Element-Inlet Pressure	011-178-619
506.61 Serial numbers 100-189, and 205	Filter Element-Main Pressure (two required)	010-053-305
	Filter Element-Inlet Pressure	011-068-403
	Filter Element-Bypass	010-088-323
	Filter Kit-Level 1	035-843-901
	Filter Kit-Level 2	031-875-901
	Filter Kit-Level 4	031-876-001
506.61 Serial numbers 190-204, 206 and up	Filter Element-Main Pressure (two required)	010-053-305
	Filter Element-Inlet Pressure	011-402-819
	Filter Kit-Level 1	038-750-401
	Filter Kit-Level 2	038-105-201
	Filter Kit-Level 4	038-751-101

** Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 10 kg on kits.

Level 1 kits contain filters and seals to perform one maintenance requirement.

Level 2 kits contain filters and seals to perform two maintenance requirements, one cleanable filter.

Level 4 kits contain filters and seals to perform four maintenance requirements, one cleanable filter.

MTS Maintenance Parts

Hydraulic Power Unit Parts

506.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
506.62 All serial numbers	Filter Element-Main Pressure (two required)	010-053-305 
	Filter Element-Inlet Pressure	011-402-819 
	Filter Kit-Level 1	038-750-401
	Filter Kit-Level 2	038-105-201
	Filter Kit-Level 4	038-751-101
506.71 Serial numbers 100-165	Filter Element-Main Pressure (four required)	010-053-305 
	Filter Element-Inlet Pressure	011-068-403
	Filter Element-Inlet Pressure	011-068-412
	Filter Element-Bypass	010-088-323
	Filter Kit-Level 1	038-514-101
	Filter Kit-Level 2	032-845-201
	Filter Kit-Level 4	032-845-301
506.71 Serial numbers 166 and up	Filter Element-Main Pressure (four required)	010-053-305 
	Filter Element-Inlet Pressure	011-402-817 
	Filter Kit-Level 1	038-750-501
	Filter Kit-Level 2	038-750-801
	Filter Kit-Level 4	038-751-201
506.81 Serial numbers 100-141	Filter Element-Main Pressure (four required)	010-053-305 
	Filter Element-Inlet Pressure	011-068-403
	Filter Element-Inlet Pressure	011-068-412
	Filter Element-Bypass	010-088-323
	Filter Kit-Level 1	038-514-101
	Filter Kit-Level 2	032-845-201
	Filter Kit-Level 4	032-845-301
506.81 Serial numbers 142 and up	Filter Element-Main Pressure (four required)	010-053-305 
	Filter Element-Inlet Pressure	011-402-817 
	Filter Kit-Level 1	038-750-501
	Filter Kit-Level 2	038-750-801
	Filter Kit-Level 4	038-751-201
506.82/92 All serial numbers	Filter Element-Main Pressure (four required)	010-053-305 
	Filter Element-Inlet Pressure	011-402-817 
	Filter Kit-Level 1	038-750-501
	Filter Kit-Level 2	038-750-801
	Filter Kit-Level 4	038-751-201

** Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 10 kg on kits.

Level 1 kits contain filters and seals to perform one maintenance requirement.

Level 2 kits contain filters and seals to perform two maintenance requirements, one cleanable filter.

Level 4 kits contain filters and seals to perform four maintenance requirements, one cleanable filter.

MTS Maintenance Parts

Hydraulic Power Unit Parts

510.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
510.10 All serial numbers	Filter Element-Main Pressure	010-088-323
	Filter Element-Inlet Pressure **	011-068-416
	Filter Kit-Level 1	038-750-601
	Filter Kit-Level 2	038-105-001
	Filter Kit-Level 4	038-751-301
510.21/.23 – Rev A & B Built through April, 1991 Top Level Assy #'s – 380821-xx; 380823-xx	Filter Element-Main Pressure	011-402-827
	Filter Element-Inlet Pressure **	011-068-417
	Filter Kit-Level 1	038-750-701
	Filter Kit-Level 2	038-105-101
	Filter Kit-Level 4	038-751-401
510.21/.23 - Rev C Built from May, 1991 – present Top Level Assy #'s – 480821-xx; 480823-xx	Filter Element-Main Pressure	011-402-827
	Filter Element-Inlet Pressure **	011-068-406
	Filter Kit-Level 1	038-750-701
	Filter Kit-Level 2	048-603-101
	Filter Kit-Level 4	048-603-201
510.30 - Rev A Built from May, 1991 – present Top Level Assy #'s – 480830-xx	Filter Element-Main Pressure	011-402-827
	Filter Element-Inlet Pressure **	011-068-406
	Filter Kit-Level 1	038-750-701
	Filter Kit-Level 2	048-603-101
	Filter Kit-Level 4	048-603-201

** Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 10 kg on kits.

Level 1 kits contain filters and seals to perform one maintenance requirement.

Level 2 kits contain filters and seals to perform two maintenance requirements, one cleanable filter.

Level 4 kits contain filters and seals to perform four maintenance requirements, one cleanable filter.

515.xx Filter Kits

MODEL	DESCRIPTION	PART NUMBER
515.07	High Pressure Filter Kit (3000 psi)	100-322-895
	High Pressure Filter Kit (4000 psi)	100-322-896
	High Pressure Seal Kit (3000 psi)	100-030-031
	High Pressure Element (3000 psi)	100-030-007
	High Pressure Seal Kit (4000 psi)	010-010-771 & 100-351-864
	High Pressure Element (4000 psi)	100-337-062
	Return Seal Kit	100-030-009
	Return Element	100-009-495
515.11	High Pressure Filter Kit (3000 psi)	100-322-895
	High Pressure Filter Kit (4000 psi)	100-322-896
	High Pressure Seal Kit (3000 psi)	100-030-031
	High Pressure Element (3000 psi)	100-030-007
	High Pressure Seal Kit (4000 psi)	010-010-771 & 100-351-864
	High Pressure Element (4000 psi)	100-337-062
	Return Seal Kit	100-030-009
	Return Element	100-009-495
515.20/.30	High Pressure Filter Kit for 505.20 & .30	100-322-897
	Return Element (two required)	010-053-349
	High Pressure Element for 505.20 & .30	010-053-305
515.60/90/90S-1	Return Element	100-029-989
	Return Element w/ O-ring Seal	100-322-894
	High Pressure Element w/ O-ring Seal	100-322-898
	High Pressure Element (two required)	010-053-305
515.120-180/180S-1/2/3	Return Element	100-029-989
	Return Element w/ O-ring Seal	100-322-894
	High Pressure Element w/ O-ring Seal	100-322-899
	High Pressure Element (four required)	010-053-305

MTS Maintenance Parts

Hydraulic Power Unit Parts

512 and 530 Filter Elements

MODEL	DESCRIPTION	PART NUMBER
512.01	Filter Element 5-in spin on	011-951-521 
512.04	Filter Element 5-in spin on	011-951-521 
530.07	Filter Element 5-in spin on	011-951-522
530.10	Filter Element 5-in spin on	011-951-521 
530.20	Filter Element 8-in	010-053-349 

Seal Kit for Radial Piston Pump - RKP 63

MODEL	DESCRIPTION	PART NUMBER
504.11/.12/.21	Seal Kit RKP 63	Contact MTS

Please refer to the HPS manual for detailed information.
Contact MTS Technical Support in case of difficulties.

Filter Elements / European Standards

MODEL/TYPE	DESCRIPTION	PART NUMBER
HC 9021 FUT 4H	Filter Element	010-088-302 
HC 9020 FDP 4H	Filter Element	010-088-309 
HC 9021 FDP 4H	Filter Element	011-402-821
HC 9600 FDP 4H	Filter Element	010-088-328 
HC 9601 FDP 4H	Filter Element	010-088-361
HC 9600 FDP 8H	Filter Element	010-088-323 
HC 9600 FDP 13H	Filter Element	011-402-827 
HC 9600 FDP 16H	Filter Element	011-402-880
HC 9601 FDP 16H	Filter Element	011-402-883
HC 9800 FDP 8H	Filter Element	011-402-824
927 965 10BXLV 10-Micron	Filter Element	011-395-936
927 964 03BXCW 3-Micron	Filter Element	011-395-937 

** Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 10 kg on kits.

Level 1 kits contain filters and seals to perform one maintenance requirement.

Level 2 kits contain filters and seals to perform two maintenance requirements, one cleanable filter.

Level 4 kits contain filters and seals to perform four maintenance requirements, one cleanable filter.

Filter Bowl Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
SB 9020 SKH*	Seal Kit	027-001-052
SB 9600 SKH*	Seal Kit	027-001-049
SB 9800 SKH*	Seal Kit	027-001-051

Heat Exchangers

MODEL/TYPE	DESCRIPTION	PART NUMBER
504.21	Heat Exchanger	027-000-353
504.22	Heat Exchanger	027-001-973

Upgrade of heat exchanger required. Contact MTS.

Each Package contains the parts needed to replace existing seals, filter elements, and accumulators.

MTS Maintenance Parts

Hydraulic Service Manifold Parts

290.xx Filter & Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
290.11	Filter Element-Main Pressure **	011-130-208
	Filter Kit – Level 2	038-926-201
	Seal Kit	038-926-301
290.12/.14	Filter Element-Main Pressure **	011-130-208
	Filter Element-Pilot Pressure	010-088-309
	Filter Kit – Level 1	038-802-901
	Filter Kit – Level 2	038-803-601
	Seal Kit	038-806-001
290.13	Filter Element-Main Pressure **	011-130-208
	Filter Kit – Level 2	038-926-201
	Seal Kit	038-806-001

**Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 5 kg on kits.

Level 1 kit contains the paper filters and filter seals required to perform one complete filter change.

Level 2 kit contains one cleanable element, paper filters & filter seals required for two filter changes.

290.xx Filter & Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
290.21/.22	Filter Element-Main Pressure **	011-130-205
	Filter Element-Pilot Pressure	010-088-309
	Filter Kit – Level 1	038-802-901
	Filter Kit – Level 2	038-803-701
	Seal Kit	038-806-101
290.31	Filter Element-Main Pressure **	011-130-205
	Filter Element-Pilot Pressure	010-088-309
	Filter Kit – Level 1	038-802-901
	Filter Kit – Level 2	038-803-701
	Seal Kit	038-806-201
290.32	Filter Element-Main Pressure **	011-130-205
	Filter Element-Pilot Pressure	010-088-309
	Filter Kit – Level 1	038-802-901
	Filter Kit – Level 2	038-803-701
	Seal Kit	038-806-301
290.4X	Seal Kit – 7-in Accumulator	011-484-307

**Cleanable element. Estimated weights: No more than 2 kg on elements; no more than 5 kg on kits.

Level 1 kit contains the paper filters and filter seals required to perform one complete filter change.

Level 2 kit contains one cleanable element, paper filters & filter seals required for two filter changes.

293.xx Filter & Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
293.1X	Seal Kit, 1-station	052-816-601
	Seal Kit, 2-station	052-816-602
	Seal Kit, 3-station	052-816-603
	Seal Kit, 4-station	052-816-604
	Filter Kit, Main Element, 3-Micron	052-957-601
	Filter Kit, Main Element, 10-Micron	052-957-501
	Filter Kit, Pilot Pressure	052-957-701
293.22	Bowl O-ring	010-010-947
	Seal Kit, 1-station	052-816-701
	Seal Kit, 2-station	052-816-801
	Seal Kit, 3-station	052-816-802
	Seal Kit, 4-station	052-816-803
	Filter Kit, Main Element, cleanable	052-957-801
	Filter Kit, Pilot Pressure	052-957-701
293.32	Seal Kit	052-816-901
	Filter Kit, Main Element, cleanable	052-957-901
	Filter Kit, Pilot Pressure	052-957-701

Note: All Series 293 Elements are disposable, except part numbers 052-957-801 & 052-957-901.

MTS Maintenance Parts

Hydraulic Service Manifold Parts

294.xx Filter & Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
294.11	Filter Element-Main Pressure	010-088-302
	Filter Kit – Level 1	038-803-001
	Filter Kit – Level 2	038-803-801
	Seal Kit	038-806-401
294.12	Filter Element-Main Pressure	010-088-302
	Filter Kit – Level 1	038-803-001
	Filter Kit – Level 2	038-803-801
	Seal Kit	038-806-501

Level 1 kit contains the paper filters and filter seals required to perform one complete filter change.

Level 2 kit contains one cleanable element, paper filters, & filter seals required for two filter changes.

295.1x & 295.2x Filter & Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
295.1x HSM	Spares Kit, HPF	Filter Kit, Main Elements, 10 Micron
	Spares Kit, PPF	Filter Kit, Pilot Pressure Elements, 3-Micron
	Spares Kit, Seals, Accumulator	Seal Kit, Pressure Accumulator
	Spares Kit, Seals, Control	Seal Kit, Pressure Accumulator
	Spares Kit, Seals, Control, Multi	Seal Kit, Multi-Channel Control Manifold
	Spares Kit, Seals, Pilot, PLD/PLE	Seal Kit, Pilot Pressure Manifold PLD/PLE
295.2x HSM	Spares Kit, HPF	Filter Kit, Main Elements, 25-Micron
	Spares Kit, PPF	Filter Kit, Pilot Pressure Elements, 3-Micron
	Spares Kit, Seals, ACC, MFLD	Seals Kit, Accumulators and Manifolds

298.xx Filter & Seal Kits

MODEL/TYPE	DESCRIPTION	PART NUMBER
298.10/11/12	Filter Element-Main Pressure, 10-Micron	011-395-936
	Filter Kit, Level 1, 10-Micron	044-205-201
	Filter Kit, Level 1, 3-Micron	044-205-301
298.12	Seal Kit	047-651-201
298.10/11	Seal Kit	047-651-301

Model 298 Notes: Filter Assy is an option on Series 298 Manifolds. All Series 298 Elements are disposable.

Level 1 kit contains the paper filters and filter seals required to perform one complete filter change.

Level 2 kit contains one cleanable element, paper filters & filter seals required for two filter changes.

Marotta Valve Replacement Kit

MODEL/TYPE	PART NUMBER
24 Volt	047-792-401
110 Volt	047-792-402

Whenever possible, we recommend the purchase of replacement kits instead of the original Marotta valves.

Replacement kits provide all the components required to replace the existing Marotta valves.

Replacement kits are less expensive to purchase, less expensive to maintain, and readily available from MTS.

Installation is not included in this price. Not intended for gas use.

MTS Maintenance Parts

Filter Elements and Flushing & Shut-Off Valves

Hydraulic Fluid Filter Elements by P/N

MATERIAL DESCRIPTION	PART NUMBER
Filter Element – 10 Micron*	008-725-501
Filter Element – Cotton, 10 in 5 Micron, Tin P	010-052-802
Filter Element – 10 Micron, for 506 Pumps	010-053-305
Filter Element – 10 Micron	010-053-308
Filter Element – 3 Micron 8 in Length	010-053-349
Filter Element – 4 in Large B15=75,3000 psid, 10 gp	010-088-302
Filter Element – 4 in Large B3=75,150 psid, 12 gp	010-088-309
Filter Element – 8 in Large B3=75,150 psid, 40 gp	010-088-323
Filter Element – HC 9600 FKP 4H	010-088-328
Filter Element – Dispos, 4 in B6=75,150 psid	010-088-334
Filter Element – 8 in Large B6=75,150 psid, 20 gp	010-088-354
Filter Element – Paper 40 Micron	011-068-403
Filter Element – Strainer 100 Mesh, 50 gpm*	011-068-406
Filter Element – 40 Micron	011-068-412
Filter Element – Strainer 60 Mesh, 10 gpm	011-068-416
Filter Element – Strainer 200 Mesh, 20 gpm*	011-068-417
Filter Element – Strainer 60 Mesh, 5 gpm*	011-068-418
Filter Element – 6 in Large 74 Micron, 10 gpm*	011-129-901
Filter Element – 8 in Large 74 Micron, 2000 psid*	011-130-205
Filter Element – 4 in Large 74 Micron, 2000 psid*	011-130-208
Filter Element – 19 in Large 25 Micron, 2000 psid	011-130-224
Filter Element – Strainer 60 Mesh, 100 gpm	011-178-619
Filter Element – B10=75,3000 psid, 30 gpm	011-395-936
Filter Element – B3=75,3000 psid, 30 gpm	011-395-937
Filter Element – 39 in Large B3=75,150 psid, 175	011-402-817
Filter Element – 16 in Large B3=75,150 psid, 150	011-402-819
Filter Element – 13 in Large B3=75,150 psid,	011-402-827
Filter Element – 5 in Spin On B3=75, Microgl	011-951-521
Filter Element – 5 in Spin On B11=75, Cellulo	011-951-522
Filter Element – 18 in Large B3=200, 150 psid	011-951-524
Filter Element – Criterion C64 IOP	100-248-238
Filter Element – 8 in Large B3=200, 150 psid	011-951-528

Controller Chassis Air Filter Elements by P/N

MATERIAL DESCRIPTION	PART NUMBER
Guard & Filter – TestStar IIs – Model 493.01	010-087-541
Filter – FlexTest SE & FlexTest 40 – Models 493.02 & 494.04	100-057-809
Filter – FlexTest GT & TestStar IIIm & FlexTest 100 – Models 493.10 [10-Slot] & 494.10	100-020-714
Filter – FlexTest GT & Aero ST – Model 493.20 [20-Slot]	100-169-943
Filter – FlexTest 60 – Model 494.06 [6-Slot]	100-168-328
Filter – FlexTest 200 – Model 494.20 [20-Slot]	100-169-890

*Cleanable element.

Estimated weights: No more than 2 kg on elements; no more than 5 kg on kits.

System Flushing Valves

MODEL	FOR SERVOVALVE	PART NUMBER
291.01B-01	252.2x Series	031-622-001
291.02B-01	252.3x Series	032-484-401
291.08B-01	256.0x Series	033-874-001

Port-Shut-off-Valves

MODEL	FOR SERVOVALVE	PART NUMBER
254.02A-01	252.2x Series	031-101-001
291.01B-01	Seal kit	049-160-701
254.02B-01	252.3x Series	006-540-801

Estimated weight of each is no more than 2 kg.

MTS Maintenance Parts

SWIFT Parts

SWIFT 20 Ultra Parts

DESCRIPTION	QTY NEEDED	PART NUMBER
Storage Case	1	100-025-776
Beam Hole Plug Assembly 20T	4	100-026-667
Beam Hole Plug Assembly 20A	4	100-019-585
Outer Label SWIFT 20T	1	100-027-373
Outer Label SWIFT 20A	1	100-019-296
Leveling Assembly – with inclinometer*	1	100-014-450
Slip Ring Assembly (with bracket and slip ring encoder)	1	100-019-727
Slip Ring Assembly (with bracket and ratcheting slip ring encoder)	1	100-052-351
Slip Ring Assembly (with bracket and slip ring encoder)	1	100-062-128

SWIFT 30 PC Parts

DESCRIPTION	QTY NEEDED	PART NUMBER
Storage Case	1	054-924-501
Flexure Covers	8	054-405-901
O-Ring Type Beam Hole Cover like SWIFT LT Assembly	4	Contact MTS
Slip Ring Bracket Assembly – with slip ring/encoder	1	100-052-353
Inner Label	1	054-286-401
Outer Label	1	054-286-201
Leveling Assembly – with inclinometer*	1	100-010-124

SWIFT 40 LT Parts

MODEL	QTY NEEDED	PART NUMBER
Storage Case	1	055-975-401
Flexure Covers	8	054-947-301
Beam Hole Cover Assembly	4	Contact MTS
Slip Ring Bracket Assembly – with slip ring/encoder	1	100-027-993
Inner Label	1	055-235-701
Outer Label	1	055-235-801
Leveling Assembly – with inclinometer*	1	100-010-509

SWIFT 45 MT Parts

MODEL	QTY NEEDED	PART NUMBER
Storage Case	1	100-211-561
Beam Hole Covers	16	100-203-509
Slip Ring Bracket	1	100-205-155
Slip Ring/Encoder	1	100-172-139
Outer Label	1	100-205-328
Level Bracket Ref 700-004-754	1	100-211-570
Level Bracket	1	100-211-585
Extension Assembly	1	100-219-464

MTS Maintenance Parts

SWIFT Parts

SWIFT 50 GLP Parts

DESCRIPTION	QTY NEEDED	PART NUMBER
Storage Case	1	100-162-342
Beam Hole Plug	8	100-155-050
Modified Set Screws needed for Beam Hole Plug	8	100-157-769
Spinning Output Cable with D-style Connector-25ft	1	056-997-604
SWIFT 50GLP has special low profile feedback cables	1	056-997-604
Outer Label	1	100-147-491
Leveling Assembly - with inclinometer	1	100-148-483
Encoder/Slip Ring with D-style Connector	1	100-143-325
Flat Spider Bracket	1	100-163-204
Dished-in-Spider Bracket	1	100-154-409
Anti-Rotate Hinge/Tube Assembly	1	100-149-775
Non-Spinning Connector Adapter for Slip Ring Replacement (also requires special cable 056-493-704)	1	100-166-648
Shunt Cable with D-style Connector-25ft	2	056-997-504
Modified Lug Nut-M22	10	056-984-401
Lug Nut Shim Washer	10	057-090-201

Shared Parts between ULTRA(20), PC(30) & SWIFT LT(40), & SWIFT 50GLP SWIFT 30 PC Parts

DESCRIPTION	QTY NEEDED	PART NUMBER
‡Small Oval Connector Cover on top of slip ring brkt. assy.	2/spider assembly	053-323-701
‡Encoder/Slip Ring – connector closer to spinning tire - non ratcheting**	1	055-554-001
‡Encoder/Slip Ring – ratcheting connector**	1	100-043-522
‡Shunt Cal Cover Assembly	2/spider assembly	Contact MTS
‡Bottom Extruded Aluminum Cover for Slip Ring assembly	1/spider assembly	053-324-001
‡Anti-Rotate Hinge/Tube Assembly	1	055-862-401
‡Anti-Rotate Steel Sleeve*	1	054-405-201
‡Anti-Rotate Modified Delrin Ball*	1	054-405-301
‡Spinning Output Cable – non ratcheting - 25 ft***	1	054-145-004
‡Spinning Output Cable – non ratcheting – Low Profile – 25 ft***	1	056-493-604
‡Spinning Output Cable – ratcheting - 25 ft***	1	056-221-904
‡Spinning Output Cable – ratcheting – Low Profile 25 ft***	1	056-753-804
‡Non-Spinning Ratcheting Output Cable – 50 ft***	1	054-145-106
‡Non-Spinning Ratcheting Output Cable – Low Profile – 50 ft***	1	056-493-709
Comm Cable to PC	1	051-074-101
Power to TI Box	1	054-870-906
Power to TI Box – Low Profile	1	056-494-606
Fast Blow Fuse*	1	010-052-509
‡Shunt Cal Cables – non ratcheting – 25 ft	1	054-871-004
Shunt Cal Cables – non ratcheting – Low Profile – 25 ft	1	056-493-504
‡Shunt Cal Cable – non-spinning ratcheting for use with connector housing, 50 ft	1	054-871-105
Non-Spinning Connector Adaptor for slip ring replacement	1	100-192-580
AC/DC Lab Power Supply (for up to 4 SWIFTS):	1	054-864-903
‡Lug Nut Washer for M14 or 9/16 studs	depends on vehicle	055-879-401
‡Lug Nut Washer for M12x1.5 or 1/2-20 studs	depends on vehicle	055-522-801
‡Lug Nut - M12x1.5, 1/2-20, M14x2.0 or 9/16-18 thread	depends on vehicle	545-521-1xx

Mini TI Section

DESCRIPTION	QTY NEEDED	PART NUMBER
Mini TI	1	100-182-635
Cable SWIFT 10 or 45	1	052-702-9xx
Mini Cable Adapter to SWIFT 20, 30, 40, and 50	1	100-224-052

*) Recommended spare part.

**) Typically only one of the three encoders/Slip rings are required.

***) Typically only one of the output cables are required per transducer.

‡) Recommended spare part for the Ultra(20), PC(30), & SWIFT LT(40) only. Not for use with the SWIFT50 GLP.

MTS Maintenance Parts

FlatTrac & Rolling Road Systems Parts

Flat Trac System Endless Belts

DESCRIPTION	PART NUMBER
Belt-Endless, Painted, FT LTR	100-047-038
Belt-Endless, Painted, FT CT Plus	100-225-198
Belt-Endless, Painted, FT 2	049-050-801
Belt-Endless, Painted, FT 1	100-064-608

Note 1: Always verify correct belt part number to specific system application.

Rolling Road System Endless Belts - Uncoated

DESCRIPTION	PART NUMBER
Belt-Endless, 0.9M, RRS	100-247-104
Belt-Endless, 1M, RRS	100-136-088
Belt-Endless, 1M, RRS	100-153-637
Belt-Endless, 1M, RRS	100-210-255
Belt-Endless, 1.1M, RRS	100-153-639
Belt-Endless, 1.1M, RRS	100-237-053
Belt-Endless, 1.1M, RRS	100-247-105
Belt-Endless, 1.1M, RRS	056-624-301
Belt-Endless, 1.1M, RRS	056-624-302
Belt-Endless, 1.1M, RRS	056-624-303
Belt-Endless, 1.1M, RRS	056-624-304
Belt-Endless, 2M, RRS	057-374-405
Belt-Endless, 2M, RRS	057-374-406
Belt-Endless, 2.4M, RRS	057-374-401
Belt-Endless, 2.4M, RRS	057-374-402
Belt-Endless, 2.4M, RRS	057-374-403
Belt-Endless, 2.4M, RRS	057-374-404
Belt-Endless, 2.4M, RRS	100-127-230
Belt-Endless, 3.2M, RRS	100-048-656
Belt-Endless, 3.2M HS, RRS	100-236-024
Belt-Endless, 6.0M, RRS	100-223-993

Note 1: Always verify correct belt part number to specific system application.

Note 2: Contact MTS for special belt coating requests.

Rolling Road System Air Bearings - Standard

DESCRIPTION	PART NUMBER
AirBearing 3.0 LG Strip	100-031-076
AirBearing 6.0 LG Strip	100-014-067
AirBearing 9.0 LG Strip	100-014-068
AirBearing 12.0 LG Strip	100-014-069

Note 1: For custom system air bearings contact MTS for pricing and availability.

MTS Maintenance Parts

Series 329 Road Simulator Parts and Durability Systems

Overview

Series 329 Road Simulator maintenance parts kits are designed for several levels of service, from routine maintenance to repair. Replace just some of the pieces or the entire component. All kits are manufactured using the latest design features and highest quality parts for maximum performance. Most of these kits are normally stocked. Contact MTS for availability.

Component Kit

A component kit enables you to replace the wear parts of a fixture, commonly known as bushing/bearing kits. These are bushings and bearings that keep the fixture operating smoothly. When used as part of a routine maintenance program, these kits help to prolong the repeatability and life of your test rig.

Component Assembly Kit

A component assembly is designed to completely replace a given component in your test rig. This gives you a fast and convenient way to perform routine maintenance and repairs. Each component assembly comes with the appropriate bushings and bearings already installed (for attachment pins and shafts, see connecting kit). Replacing a component is simple: remove the old component and put the new one in its place.

Connecting Kit

The connecting kit is designed to replace the pieces that connect two adjacent components, commonly referred to as pin kits. These kits consist of pins, shafts, and grease fittings. These parts are not included in a component assembly or kit because they are not always needed for a repair.

Material Handling Trays

The wear parts shown in the trays below are included in the Series 329 rebuild kit. These original OEM replacement parts are used to rebuild the 11 unique structural fixture assemblies used within one MTS Series 329 4DOF corner.

Factory Maintenance

Not every Series 329 Road Simulator component is listed on the following pages. There are some components that can be serviced with the use of custom-made tools and special processes available only at an MTS facility. Please contact your MTS service specialist for parts needs not found in this catalog.


Series 329 Rebuild Program

It's almost like getting a brand new system. The Series 329 rebuild program was developed to make it efficient and economical to accomplish extensive routine maintenance during a single maintenance period. The Series 329 rebuild program will disassemble your corners at our factory and rebuild them to like-new condition. This program can be tailored to meet the maintenance requirements of your system. Please contact us for an evaluation and a discussion of these options.

MTS offers a planned approach for fast cycle times in completing complex 329 Rig overhauls either onsite or in the factory with our efficient 329 rebuild program. The two trays shown below contain all pin, bushings, bearings, and shims to completely rebuild all moving parts on one corner of 329 LT or PC. The trays are shipped stacked together in a box. The kits are compact and easily stored at the customer's site for future use if desired.



Small parts trays

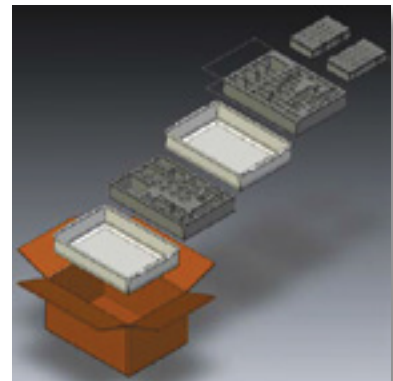

Series 329 rebuild kits help organize the many small parts needed for a successful rebuild.

MTS Maintenance Parts

329 System Parts

329 6DOF Hardware Kits for Level 2 Rebuild

Pre 2007 329 PC 6DOF REF 562040-XX

DESCRIPTION	PART NUMBER
Hardware Kits-Need to be created	

Pre 2007 329 LT 6DOF REF 562045-XX

2007-09 329 PC/LT REF 571210-XX

DESCRIPTION	PART NUMBER
Spindle Service ASSY 329-6D (HDWR)	056-736-651
VERT STRUT SVC ASSY-329 6DOF TUBE (HDWR)	057-460-351
SERVICE ASSY-VERT HDWR KIT, 329-6DOF LT	056-736-452
INPUT ASSY-LONG (HDWR), 329LT-6DOF SVC	056-736-751
SERVICE ASSY,329-6D PC/LT,LAT/STR (HDWR)	056-736-851

2010-Current 329 PC/LT REF 575000-XX

DESCRIPTION	PART NUMBER
SPINDLE SVC ASSY 329-6D (HDWR)	056-736-651
VERT STRUT SVC ASSY-329 6DOF TUBE (HDWR)	057-460-351
SERVICE ASSY-VERT HDWR KIT, 329-6DOF LT	056-736-452
INPUT ASSY-LONG (HDWR), 329LT-6DOF SVC	056-736-751
SERVICE ASSY-LAT/STR (HDWR)329-6D PC/LT	057-460-551

MTS Maintenance Parts

Model 329 4 DOF Passenger Car Road Simulator Parts

329 Bearings

DESCRIPTION	PART NUMBER
Rod End Bearing, 3/4-16, R.H. Thread, 0.750 Bore	042-798-601
Rod End Bearing, 3/4-16, L.H. Thread, 0.750 Bore	042-798-602
Rod End Bearing, 3/4-16, R.H. Thread, 0.625 Bore	042-795-401
Rod End Bearing, 1/2-20, R.H. Thread, 0.500 Bore	042-798-603

329 Load Cells

DESCRIPTION	PART NUMBER
661.57x-xx Load Cell for 329 & 329LT Longitudinal (Includes attachment screws)	048-104-273
661.57x-xx Load Cell for 329 & 329LT Lateral (Includes attachment screws)	048-102-974
661.19x-xx Load Cell-Modified for 329 Brake Actuator	440396-XX

329 Actuator Seal Kits

DESCRIPTION	PART NUMBER
Vertical Actuator Seal Kit for 244.20 Spcl (Static Support) This Seal Kit provides seals for both the Model 244.20 Actuator and the Static Support Option.	048-705-606
Vertical Actuator Seal Kit for 244.22 Spcl (Static Support) This Seal Kit provides seals for both the Model 244.22 Actuator and the Static Support Option.	047-785-401
Lateral Actuator Seal Kit for 244.12	008-710-047 
Longitudinal Actuator Seal Kit for 244.21	008-710-048 
Brake Actuator Seal Kit for 242.0x	042-065-404

Miscellaneous 329 Parts

DESCRIPTION	PART NUMBER
Water Dummy, Torso ONLY	011-860-904
Water Dummy Assembly (w/valves)	045-962-001
Brake-Strut Shaft Service Kit	050-668-961
Brake-Strut Bushing/Bearing Service Kit	050-668-981
Seat Mounted Steer Robot (cables & amplifier enclosure not included)	100-107-417
Dynamic Brake Intensifier (brake fluid not included, EPR seals compatible with most brake fluids)	045-233-306
Replacement Seal Kit for Dynamic Brake Intensifier	100-110-498
Dynamic Brake Intensifier Air Inlet Filter	100-130-988
HSM Pilot Pressure Filter	052-957-701
AEROSHELL14 GREASE (tube for grease gun)	011-010-225

MTS Maintenance Parts

Model 329 4 DOF Passenger Car Road Simulator Parts

329 Structural Parts

DESCRIPTION	DRAWING NUMBER	CLEVIS P/N	LORD P/N
Spindle Adapter (Component Assembly)	1	048-102-691	048-103-091
Spindle Adapter @ Vertical Strut Bearings and Bushings (Component Kit)		048-102-681	048-103-081
Spindle Adapter @ Vertical Strut Clamp and Shaft Kit (Connecting Kit)		048-102-662	048-103-062
Spindle Adapter @ Lateral Strut Clamp and Shaft Kit (Connecting Kit)		048-102-661	048-103-061
Vertical Strut Assembly (Component Assembly)	2	048-103-191	048-103-191
Vertical Strut-Bearings (Component Kit)		048-103-181	048-103-181
Vertical Strut-Pin Kit @ Center (Connecting Kit)		048-103-161	048-103-161
Vertical Strut-Special/For Limited Clearance Applications (Component Assembly)		044-038-101	044-038-101
Longitudinal Attach Beam Assembly (Component Assembly)	3	049-071-891	049-071-691
Longitudinal Attach –Bushings & Bearings (Component Kit)		049-071-881	049-071-681
Longitudinal Attach – Triangular Cast Adapters (Connecting Kit)		049-071-871	049-071-671
Longitudinal Attach – Pin Kit @ Center (Connecting Kit)		049-071-661	049-071-661
Delta Crank Assembly (Component Assembly)	4	048-102-891	048-103-391
Delta Crank @ Vertical Strut Bushing w/Bearing (Single Bushing Only)		048-102-882	048-103-382
Delta Crank @ Vertical Struts Bushing w/Bearing (Component Kit)		048-102-881	048-103-381
Delta Crank @ Vertical Bellcrank/Bearings (Component Kit)		048-102-871	048-103-371
Delta Crank @ Vertical Bellcrank/Pin Kit (Connecting Kit)		048-102-862	048-103-362
Delta Crank @ Vertical Strut/ Pin Kit (Single Pin Only)		048-102-861	048-103-361
Brake Lever (Component Assembly)	5	048-103-491	048-103-491
Brake Lever/Bearings (Component Kit)		048-103-481	048-103-481
Brake Lever/Pin Kit @ Vertical Bellcrank (Connecting Kit)		048-103-461	048-103-461
Brake Input Strut Assembly (Component Assembly)	6	048-103-591	048-103-591
Brake Input Strut/Rod End Bearing/ 3/4"-16 R.H.		042-798-601	042-798-601
Brake Input Strut/Rod End Bearing/ 3/4"-16 L.H.		042-798-602	042-798-602
Vertical Bellcrank Assembly (Component Assembly)	7	048-103-691	048-103-691
Vertical Bellcrank & Flexure Bearings (Component Kit)		048-103-681	048-103-681
Vertical Flexure Housing (Component Assembly)		048-103-671	048-103-671
Vertical Flexure Housing/Pin Kit (Connecting Kit)		048-103-662	048-103-662
Vertical Bellcrank/Pin Kit @ Pivot (Connecting Kit)		048-103-661	048-103-661
Longitudinal Bellcrank (Component Assembly)	8	048-103-791	048-103-791
Longitudinal Strut (Component Assembly)	9	048-104-291	048-103-891
Longitudinal Strut Bearings & Bushings (Component Kit)		048-104-282	048-103-882
Longitudinal Strut Bearings (Component Kit)		048-104-281	048-103-881
Longitudinal Strut Load Cells (Connecting Kit)		048-104-273	048-103-873
Longitudinal Strut Standard Housing Assy (Component Assembly)		048-104-272	048-103-872
Longitudinal Strut Swivel Housing Assy (Component Assembly)		048-104-271	048-103-871
Longitudinal Strut/Pin Kit @ Bellcrank (Connecting Kit)		048-104-261	048-103-861

MTS Maintenance Parts

Model 329 4 DOF Passenger Car Road Simulator Parts

329 Structural Parts (continued)

DESCRIPTION	DRAWING NUMBER	CLEVIT P/N	LORD P/N
Lateral Strut (Component Assembly)	10	048-102-991	048-103-991
Lateral Strut Wrist Joint (Component Kit)		048-102-982	
Lateral Strut Bearings & Bushings (Component Kit)		048-102-981	048-103-981
Lateral Strut Bearing Housing Assy @ Spindle (Component Assembly)		048-102-975	048-103-975
Lateral Strut Load Cell (Connecting Kit)		048-102-974	048-103-974
Lateral Strut Bearings/Bushings/Wrist (Component Kit)		048-102-973	048-103-973
Lateral Strut Bushing Housing Assy @ Bellcrank (Component Assembly)		048-102-972	048-103-972
Lateral Strut Tube & Tie Rod (Connecting Kit)		048-102-971	048-103-971
Lateral Strut Pin Kit @ Spindle (Connecting Kit)		048-102-962	048-103-962
Lateral Strut Pin Kit @ Bellcrank (Connecting Kit)		048-102-961	048-103-961
Lateral Bellcrank (Component Assembly)	11	048-104-091	048-104-091

(329 4 DOF PC drawing on next page)

Component Assembly:

- Is usually the -x91 part numbers
- Involves a preassembly of parts for fast “plug & play” and to reduce customer downtime.

Connecting Kit


- Usually involves a smaller group of parts than a Component Assembly, which may or may not be preassembled.
- Focus of a connecting kit is on repair or replacement of a failed part.
- Customer uptime on a system is usually not as critical; involves system teardown to a single component level.

Component Kit

- Involves replacement of parts with a predetermined life span. Components of this nature are all bearings, bushings, and so forth.

MTS Maintenance Parts


Model 329 4 DOF Passenger Car Road Simulator Parts

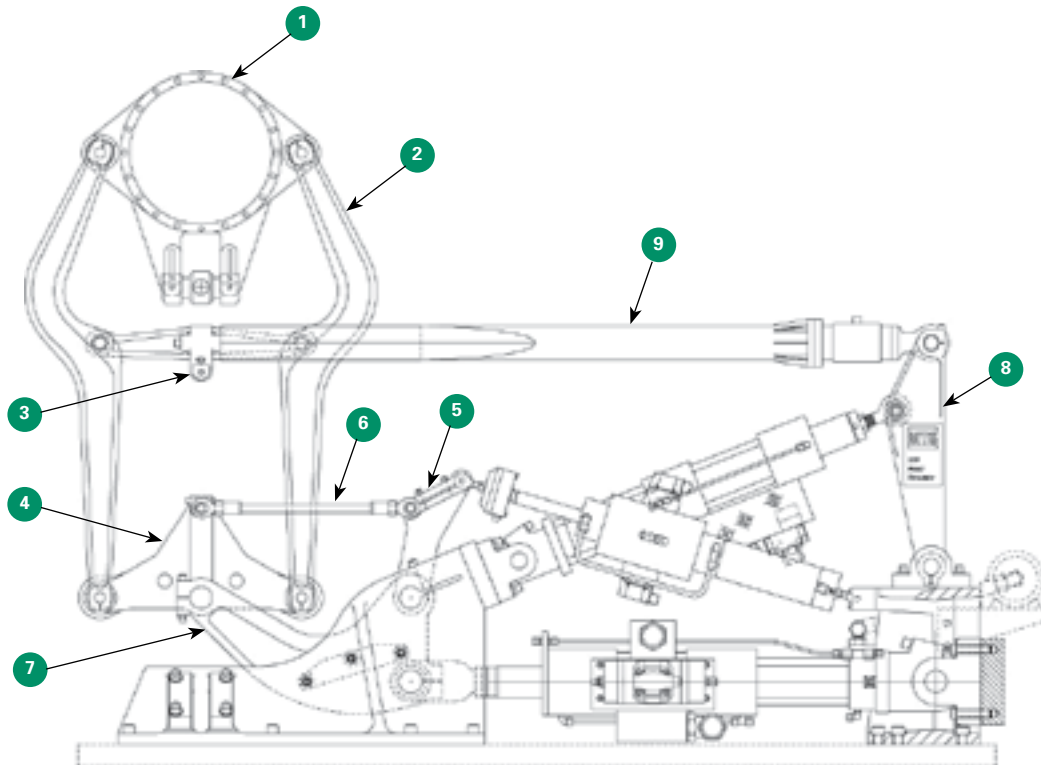

To determine which part number is correct for your Series 329 Road Simulator, you must determine if you need assemblies with Lord bushings or Clevite bushings. MTS keeps detailed records of your system to help determine which bushings you have. Make sure you discuss which bushing is correct for your system when placing your order.

Guidelines for choosing the correct bushing:

- » If you are replacing just one assembly, you should replace it with the same type of bushing that was in the original assembly.
- » If you are replacing multiple assemblies as part of a major overhaul of your system, you should consider assemblies with Clevite bushings.
- » If you decide to replace your existing Lord bushing assembly or assemblies with an assembly with Clevite bushings, you may need to adjust (reiterate) your drive files to compensate for the slightly different performance characteristics of the Clevite bushings.

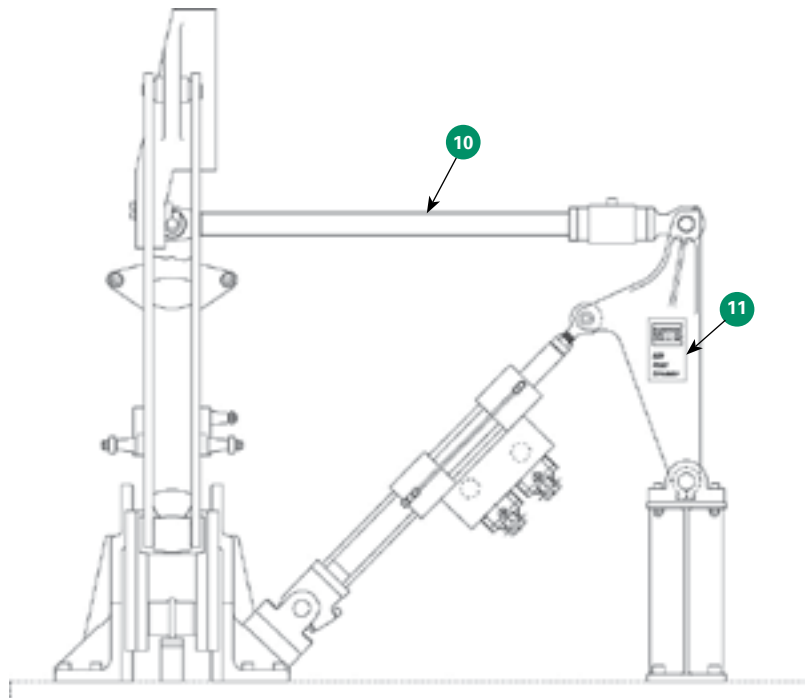

Note: Some maintenance and part replacement requires special tools to complete. MTS field service can also complete these tasks on site with the appropriate tools and expertise.

MTS Maintenance Parts

Model 329 6 DOF Road Simulator Actuator Parts

329 6 DOF Road Simulator Actuator Parts

**To determine the correct service parts locate the actuator part number on cylinder assembly.
This applies to all 329 6 DOF Passenger Car and Light Truck. Contact MTS if your actuator is not shown here.
The last two digits indicate left or right version and will typically be a -01 or -02.*

VERTICAL ACTUATOR PART NUMBER	DESCRIPTION	PART NUMBER
571599-XX	LVDT	039-075-103
	3 Stage Servo Valve	033-561-002
	2 Stage Servo Valve	056-006-201
	Seal Kit	008-710-050
564342-XX	LVDT	048-821-201
	3 Stage Servo Valve	033-561-002
	2 Stage Servo Valve	056-006-201
	Seal Kit	056-079-353
571599-XX	LVDT	039-075-103
	3 Stage Servo Valve	033-561-002
	2 Stage Servo Valve	056-006-201
	Seal Kit	008-710-050
561982-XX	LVDT	048-821-201
	3 Stage Servo Valve	033-561-002
	2 Stage Servo Valve	056-006-201
	Seal Kit	048-978-504
569889-XX	LVDT	048-821-201
	3 Stage Servo Valve	033-561-002
	2 Stage Servo Valve	056-006-201
	Seal Kit	056-079-348
LONGITUDINAL ACTUATOR PART NUMBERS	DESCRIPTION	PART NUMBER
564247-XX	LVDT	039-075-111
	2 Stage Servo Valve	056-006-501
	Seal Kit	008-710-048
560805-XX	LVDT	039-075-111
	3 Stage Servo Valve	033-366-501
	2 Stage Servo Valve	056-006-101
	Seal Kit	008-710-075
LATERAL ACTUATOR PART NUMBERS	DESCRIPTION	PART NUMBER
100-063-299	LVDT	039-075-111
	2 Stage Servo Valve	056-006-501
	Seal Kit	056-079-305
100-064-331	LVDT	039-075-111
	3 Stage Servo Valve	033-366-508
	2 Stage Servo Valve	056-348-201
	Seal Kit	008-710-048
STEER ACTUATOR PART NUMBERS	DESCRIPTION	PART NUMBER
564234-XX	LVDT	039-075-104
	2 Stage Servo Valve	056-006-501
	Seal Kit	008-710-046
564248-XX	LVDT	039-075-104
	2 Stage Servo Valve	056-006-501
	Seal Kit	008-710-047
BRAKE ACTUATOR PART NUMBERS	DESCRIPTION	PART NUMBER
564235-XX	LVDT	039-075-104
	2 Stage Servo Valve	056-006-501
	Seal Kit	056-079-304
564324-XX	LVDT	039-075-104
	2 Stage Servo Valve	056-006-501
	Seal Kit	048-705-609

MTS Maintenance Parts

Model 329 6 DOF Passenger Car Road Simulator Parts

329 6 DOF Passenger Car Road Simulator Parts

**The following part numbers apply to 329 6DOF Passenger Car Road Simulators purchased pre December 1, 2006.
See 329 6DOF Light Truck section for purchases after this date.*

SPINDLE HOUSING ASSEMBLY*	DRAWING NUMBER	DESCRIPTION	PART NUMBER
329-6DOF SVC Assembly Spindle (HSG/BRG)	1	Spindle Housing, Shafts & Bearings	056-736-692
VERTICAL STRUT ASSEMBLY*	2	DESCRIPTION	PART NUMBER
Vertical Strut SVC Assembly-329-6D (Strut)		Vertical Struts, Bearings, and Bushings	056-737-092
Vertical Strut SVC Assembly-329-6D (BRG)		Vertical Strut Assembly Bearings	056-737-082
Vertical Strut SVC Assembly-329-6D (BSHG)		Vertical Strut Assembly Bushings	056-737-084
Vertical Strut SVC Assembly-329-6D (PIN)		Vertical Strut Assembly Pin Kit	056-737-072
Spherical Ball & Pin Upgrade Assembly		Updates ball and pin to include a grease zerk for lubrication of the ball joint	056-736-972
VERTICAL STRUT ASSEMBLY*	3	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly (J Crank Component Kit)		J-Cranks only	056-736-492
329-6 DOF SVC Assembly (J Crank SVC Assembly)		J-Crank Replacement Set; (one corner) no pins	056-736-491
329-6 DOF SVC Assembly (J Crank Bearing)		J-Crank Bearings and Wear Components	056-736-481
329-6 DOF SVC Assembly (Pin Kit) x		Vertical Input Assembly Pin Kit	056-736-471
329-6 DOF SVC Assembly (Flexure)		Vertical Input Flexure and Bearings/Wear Components	056-736-392
329-6 DOF SVC Assembly (Flexure Bearing)		Vertical Input Flexure Bearings	056-736-381
329-6 DOF SVC Assembly (Delta Crank)		Delta Crank and Bearings/Wear Components	056-736-291
329-6 DOF SVC Assembly (Delta Crank Bearing)		Delta Crank Bearings	056-736-281
329-6 DOF SVC Assembly (Brake)		Brake Bellcrank and Bearings/Wear Components	056-736-191
329-6 DOF SVC Assembly (Brake Bearing)		Brake Bearings	056-736-181
329-6 DOF SVC Assembly (Brake Bushing)		Brake Input Bushing	056-736-172
329-6 DOF SVC Assembly (Brake Lever)		Brake Input Links and Bearings/Wear Components	056-736-171
LONGITUDINAL INPUT ASSEMBLY*	4	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly Longitudinal Input (Bellcrank and Bearing)		Longitudinal Bellcrank and Bearings/Wear Components	056-736-791
329-6 DOF SVC Assembly Longitudinal Input (Bearing)		Longitudinal Bearings and Wear Components	056-736-781
329-6 DOF SVC Assembly Longitudinal Input (Pin Kit)		Longitudinal Pin Kit	056-736-771
LATERAL INPUT ASSEMBLY*	5	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly Lateral 1 (Bellcrank and Bearing)		Lateral 1 Bellcrank and Bearings/Wear Components	056-736-891
329-6 DOF SVC Assembly Lateral 2 (Bellcrank and Bearing)		Lateral 2 Bellcrank and Bearings/Wear Components	056-736-892
329-6 DOF SVC Assembly Lateral 3 (Bellcrank and Bearing)		Lateral 3 Bellcrank and Bearings/Wear Components	056-736-893
329-6 DOF SVC Assembly Lateral/Steer (Bearing)		Lateral 1, Lateral 2, Lateral 3 Bearings/Wear Components	056-736-881
329-6 DOF SVC Assembly Lateral/Steer (Pin Kit)		Lateral 1, Lateral 2, Lateral 3 Pin Kit	056-736-771
LONGITUDINAL STRUT ASSEMBLY*	6	DESCRIPTION	PART NUMBER
Strut Assembly Longitudinal SVC Assembly (Strut)		Longitudinal Strut and End caps, Assembled	056-737-391
Strut Assembly Longitudinal SVC Assembly (Bushing)		Longitudinal Strut Bearings/Wear Components	056-737-381
Strut Assembly Longitudinal SVC Assembly (Clevis)		Longitudinal Strut Clevis and Bearings	056-737-382
Strut Assembly Longitudinal SVC Assembly (U-Joint)		Longitudinal U-Joint and Bearings	056-737-383

MTS Maintenance Parts

Model 329 6 DOF Passenger Car Road Simulator Parts


329 6 DOF Passenger Car Road Simulator Parts (continued)

CENTER LATERAL STRUT ASSEMBLY*	DRAWING NUMBER 7	DESCRIPTION	PART NUMBER
Strut Assembly Lateral LT SVC Assembly (Strut)		Lateral Strut and End Caps, Assembled	056-737-191
Strut Assembly Lateral LT SVC Assembly (Bearing)		Lateral Strut Rod-Eye and Bearings/Wear Components	056-737-181
Strut Assembly Lateral LT SVC Assembly (Bearing)		Upgrade for Lateral Input Spherical Bearing	100-178-317
STEER STRUT ASSEMBLY*	8	DESCRIPTION	PART NUMBER
Strut Assembly Steer SVC Assembly (Strut)		Steer Strut and End caps, Assembled	056-737-291
Strut Assembly Steer SVC Assembly (Bearing)		Steer Strut Bearings/Wear Components	056-737-281
Strut Assembly Steer SVC Assembly (HSG/Bushing)		Wrist Joint Assembly for all Horizontal Struts	056-735-761
ACTUATOR ASSEMBLIES		DESCRIPTION	PART NUMBER
Vertical			
LVDT Assembly-Closed Housing, Static Support, 7 in, 244.3		Vertical LVDT Assembly	044-965-901
Servovalve-2.5 gpm 252.22G-01 4 Port		Two-Stage Vertical Servovalve Assembly	056-006-201
244.3x Seal Kit 45 Kip, 1 in Cushion, STCSUPR		Vertical Actuator Seal Kit	008-710-085
Longitudinal			
LVDT Assembly Closed Housing		Longitudinal LVDT Assembly	036-738-931
Servovalve Assembly-50 gpm 256.05A-01, 1 gpm PL		Three-Stage Longitudinal Servovalve Assembly (includes two-stage)	Contact MTS
Servovalve-252.xx		Two-Stage Longitudinal Servovalve Assembly	056-006-501
Shipping Assembly-Seal Kit, 244.20, Structural		Longitudinal Actuator Seal Kit	008-710-048
Lateral			
LVDT Assembly Closed Housing		Lateral LVDT Assembly	036-738-931
Servovalve Assembly-256.05A-005, 1/2 gpm PLT		Three-stage Lateral Servovalve Assembly (includes two-stage)	Contact MTS
Servovalve- 252.xx		Two-stage Lateral Servovalve Assembly	056-006-501
Seal Kit 244.21, Structural		Lateral Actuator Seal Kit	056-079-305
Steer			
LVDT Assembly Closed Housing 244.11/12 8 in SP		Steer LVDT Assembly	036-738-821
Servovalve-15 gpm 252.25G-01 4 Port		Steer Servovalve Assembly	056-006-501
Seal Kit, 244.12 Structural		Steer Actuator Seal Kit	008-710-046
Brake			
LVDT Assembly Closed Housing 244.11/12 8 in SP		Brake LVDT Assembly	036-738-821
Servovalve-15 gpm 252.25G-01 4 Port		Brake Servovalve Assembly	056-006-501
Seal Kit, 244.10		Brake Actuator Seal Kit	056-079-304
HSM FILTER KIT		DESCRIPTION	PART NUMBER
Filter Element (Pilot Pressure)		8 in lg b3 = 200, 150 psid	100-216-412
LUBRICANT			PART NUMBER
AEROSHELL14 GREASE (tube for grease gun)			011-010-225
KRYTOX GREASE (tube for grease gun)			100-123-981
DYNAMIC BRAKE INTENSIFIER			PART NUMBER
Dynamic Brake Intensifier Air Inlet Filter			100-130-988
Replacement Seal Kit for Dynamic Brake Intensifier			100-110-498

(329 6 DOF PC drawing on next page)

MTS Maintenance Parts

Model 329 6 DOF Passenger Car Road Simulator Parts

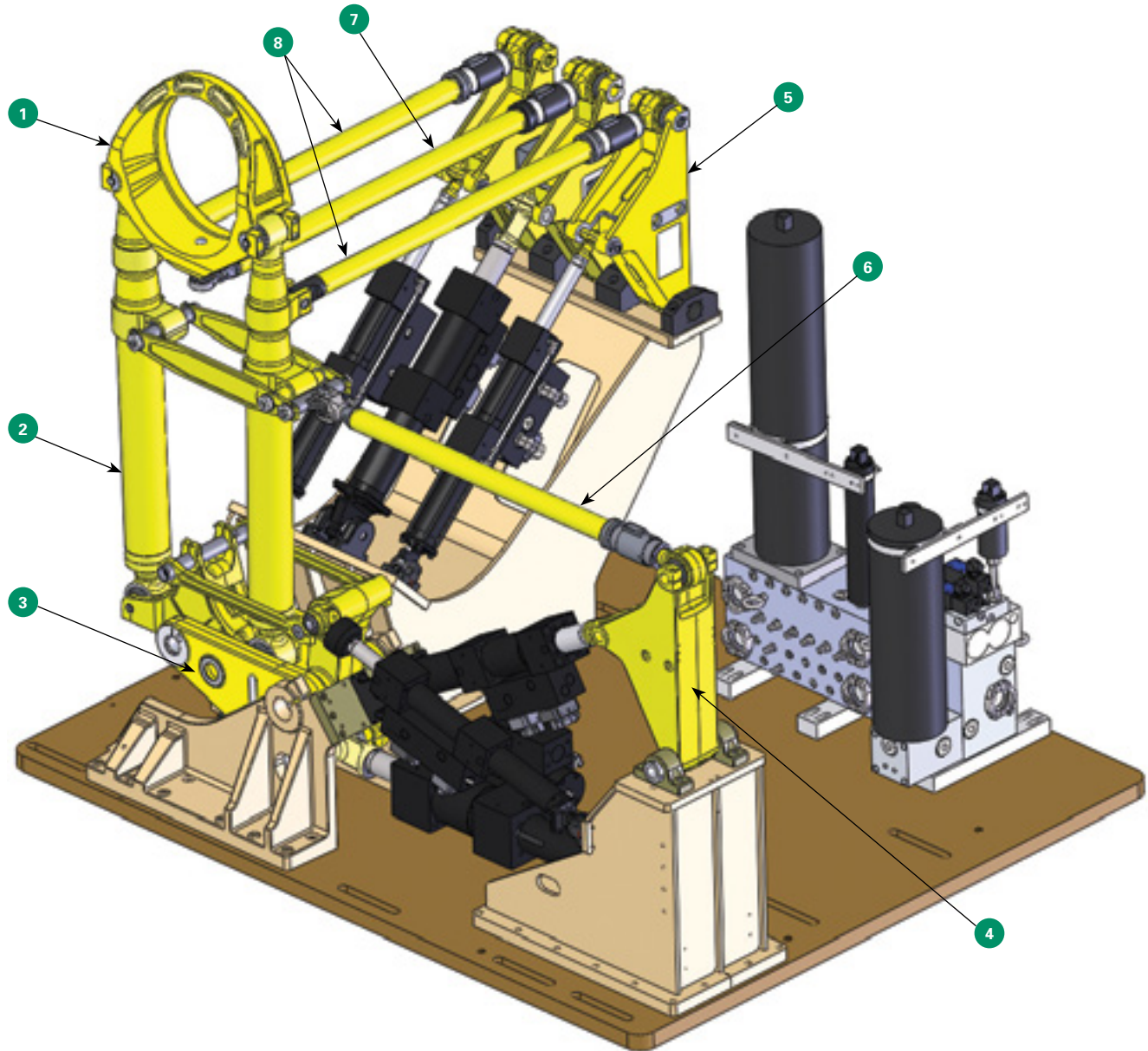

MTS Maintenance Parts

Model 329 4 DOF Light Truck Road Simulator Parts

329 4 DOF Light Truck Structural Parts

DESCRIPTION	DRAWING NUMBER	QUANTITY PER CORNER	PART NUMBER
LT Spindle Adapter (Component Assembly)	1	1	045-538-791
LT Tire Patch Input (Component Assembly)		1	044-453-391
LT Vertical Strut Assembly (Component Assembly)	2	2	047-419-791
LT Longitudinal Strut Assembly (Component Assembly)	3	1	045-844-491
LT Longitudinal Wrist Joint Component Assembly	4	1	045-516-791
LT Brake Strut (No-Brake option) Connecting Kit		1	044-452-391
LT Brake Housing (w/Brake option) Component Assembly	5	1	045-836-191
LT Brake Lever (Component Assembly)	6	1	044-451-591
LT Brake Lever (Component Kit)		1	044-451-581
LT Brake Link Component Assembly	7	2	044-452-491
LT Brake Link Component Kit		2	044-452-481
LT Delta Crank Component Assembly	8	1	045-324-491
LT Delta Crank Component Kit		1	045-324-481
LT Vertical J-Crank Component Assembly	9	1	045-538-891
LT Vertical J-Crank Component Kit		1	045-538-881
LT Flexure-Vert-45 Kip Connecting Kit	10	1	045-324-791
LT Flexure-Vert-45 Kip Component Kit		1	045-324-781
LT Flexure-Vert-35 Kip Connecting Kit		1	045-324-792
LT Flexure-Vert-35 Kip Component Kit		1	045-324-782
LT Flexure-Lat/Long Connecting Kit	11	2	045-324-991
LT Flexure-Lat/Long Component Kit		2	045-324-981
LT Bellcrank-Lat/Long Component Assembly	12	2	045-537-691
LT Bellcrank-Lat/Long Component Kit		2	045-537-681
LT Lateral Strut Component Assembly	13	1	045-844-291
LT Lateral Wrist Joint Service Assembly	14	1	045-516-691
LT Pillow Block Bearing-Lat/Long Connecting Kit	15	4	043-683-991

MISCELLANEOUS	PART NUMBER
AEROSHELL14 GREASE (tube for grease gun)	011-010-225
HSM Pilot Pressure Filter	100-216-412
Replacement Seal Kit for Dynamic Brake Intensifier	100-110-498
Dynamic Brake Intensifier Air Inlet Filter	100-130-988

(329 4 DOF LT drawing on next page)


Series 329 4 DOF LT Strut and Wrist Joint Repair

Most Series 329 4 DOF LT struts (drawing item 2,3,5 and 13) and wrist joints (drawing item 4 and 14) can be more economically repaired than buying whole new assemblies. Lead times vary.

Contact MTS for an evaluation and a quote.


MTS Maintenance Parts

Model 329 4 DOF Light Truck Road Simulator Parts

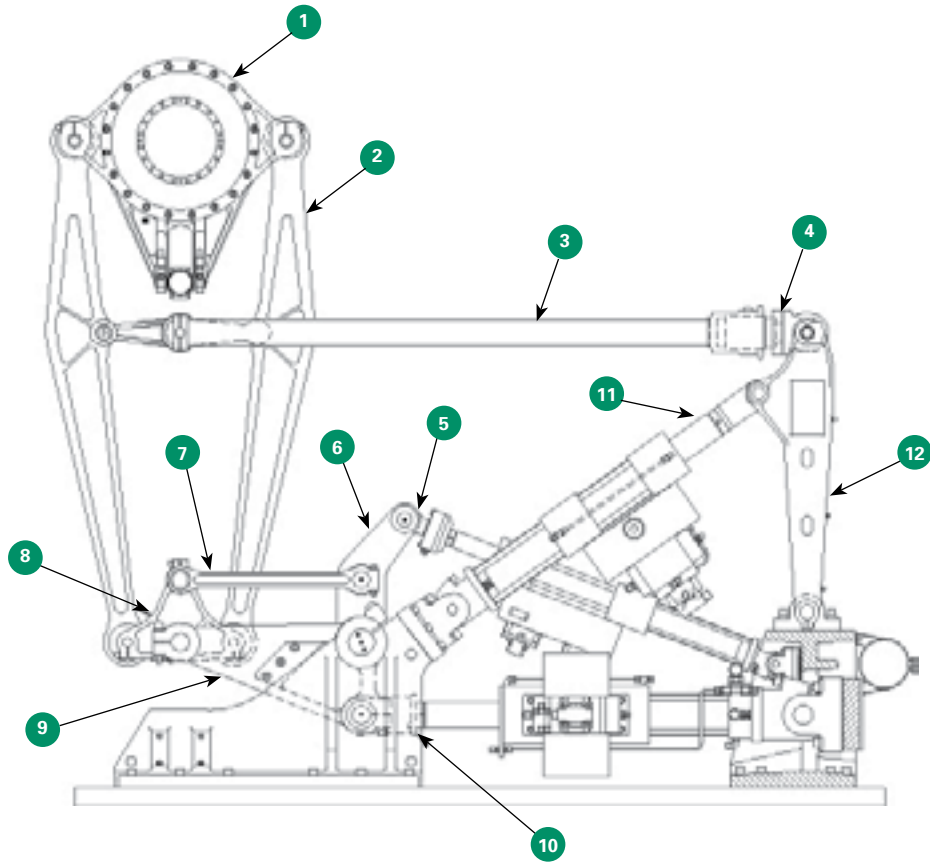

Note: It is recommended that any 329 4DOF LT bushing that is replaced to be done so at MTS. This can be part of a total system overhaul or with select components/housings being sent back for bushing replacement.

Note: The longitudinal strut (item 3) should only be disassembled if absolutely necessary. The tensioning rods surrounding each strut tube must always be replaced if taken apart. Also, it is likely that disassembly will disrupt the bond between the strut tubes and either housing end. This requires re-machining the housings, replacing the strut tubes and completing the bonding process. This must also be completed at MTS.

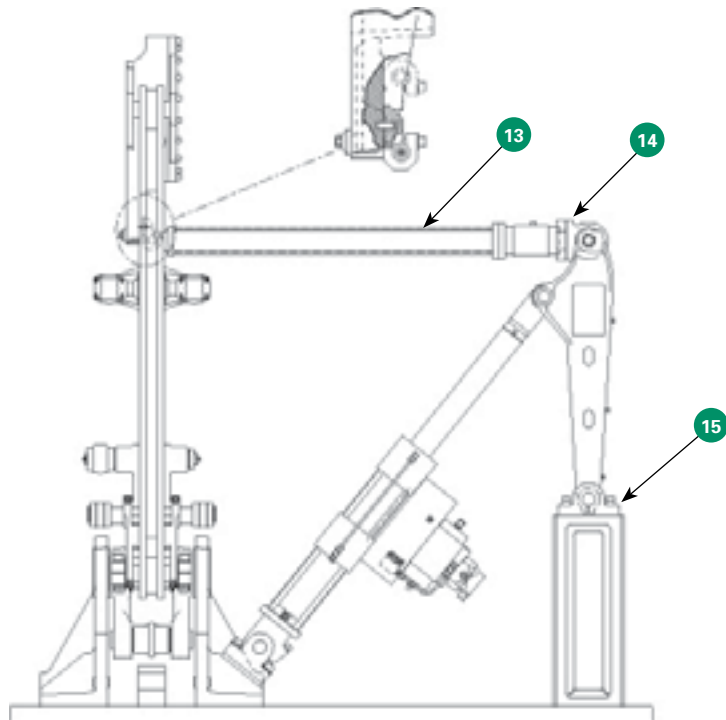

MTS Maintenance Parts

Model 329 6 DOF Light Truck Road Simulator Parts

329 6 DOF Light Truck Road Simulator Parts

The following part numbers apply to 329 6DOF Light Truck Road Simulators purchased pre September 1, 2009 and 329 6DOF Passenger Car Road Simulators (excluding) purchased between December 1, 2006 and September 1, 2009.*

SPINDLE HOUSING ASSEMBLY*	DRAWING NUMBER	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly Spindle (Housing/Bearing)	1	Spindle Housing, Shafts & Bearings	056-736-691
VERTICAL STRUT ASSEMBLY*	2	DESCRIPTION	PART NUMBER
Vertical Strut SVC Assembly-329-6D (Pin)		Vertical Strut Assembly Pin Kit	057-460-372
Vertical Strut SVC Assembly-329-6D (Bearing)		Vertical Strut Assembly Bearings	057-460-382
Vertical Strut SVC Assembly-329-6D (Bushing)		Vertical Strut Assembly Bushings	057-460-384
Vertical Strut SVC Assembly-329-6D (RF/LR Strut)		Vertical Struts, Bearings, and Bushings	057-460-391
Vertical Strut SVC Assembly-329-6D (LF/RR Strut)		Vertical Struts, Bearings, and Bushings	057-460-392
VERTICAL INPUT ASSEMBLY*	3	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly (J Crank Component Kit)		J-Crank only	056-736-492
329-6 DOF SVC Assembly (J Crank SVC Assembly)		J-Crank Replacement Set; (one corner) no pins	056-736-491
329-6 DOF SVC Assembly (J Crank Bearing)		J-Crank Bearings and Wear Components	056-736-481
329-6 DOF SVC Assembly (Pin Kit)		Vertical Input Assembly Pin Kit	056-736-471
329-6 DOF SVC Assembly (Flexure)**		Vertical Input Flexure and Bearings/Wear Components	056-736-391
329-6 DOF SVC Assembly (Flexure Bearing)		Vertical Input Flexure Bearings	056-736-381
329-6 DOF SVC Assembly (Delta Crank)		Delta Crank and Bearings/Wear Components	056-736-291
329-6 DOF SVC Assembly (Delta Crank Bearing)		Delta Crank Bearings	056-736-281
329-6 DOF SVC Assembly (Brake)		Brake Bellcrank and Bearings/Wear Components	056-736-191
329-6 DOF SVC Assembly (Brake)		Brake Bearings	056-736-181
329-6 DOF SVC Assembly (Brake Bushing)		Brake Input Bushing	056-736-172
329-6 DOF SVC Assembly (Brake Lever)		Brake Input Links and Bearings/Wear Components	056-736-171
LONGITUDINAL INPUT ASSEMBLY*	4	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly Longitudinal Input (Bellcrank and Bearing)		Longitudinal Bellcrank and Bearings/Wear Components	056-736-791
329-6 DOF SVC Assembly Longitudinal Input (Bearing)		Longitudinal Bearings and Wear Components	056-736-781
329-6 DOF SVC Assembly Longitudinal Input (Pin Kit)		Longitudinal Pin Kit	056-736-771
LATERAL INPUT ASSEMBLY*	5	DESCRIPTION	PART NUMBER
329-6 DOF SVC Assembly Lateral 1 (Bellcrank and Bearing)		Lateral 1 Bellcrank and Bearings/Wear Components	056-736-891
329-6 DOF SVC Assembly Lateral 2 (Bellcrank and Bearing)		Lateral 2 Bellcrank and Bearings/Wear Components	056-736-892
329-6 DOF SVC Assembly Lateral 3 (Bellcrank and Bearing)		Lateral 3 Bellcrank and Bearings/Wear Components	056-736-893
329-6 DOF SVC Assembly Lateral/Steer (Bearing)		Lateral 1, Lateral 2, Lateral 3 Bearings/Wear Components	056-736-881
329-6 DOF SVC Assembly Lateral/Steer (Pin Kit)		Lateral 1, Lateral 2, Lateral 3 Pin Kit	056-736-871
LONGITUDINAL STRUT ASSEMBLY*	6	DESCRIPTION	PART NUMBER
Strut Assembly Longitudinal SVC Assembly (Strut)		Longitudinal Strut and End caps, Assembled	056-735-791
Strut Assembly Longitudinal SVC Assembly (U-Joint)		Longitudinal Joint Assembly with bearings/pins	100-219-013
Strut Assembly Longitudinal SVC Assembly (Pin)		Longitudinal Strut Pin Kit	056-735-772
Strut Assembly Longitudinal SVC Assembly (Bearing)			056-735-781
Strut Assembly Longitudinal SVC Assembly (Clevis)			056-735-783
Strut Assembly Longitudinal SVC Assembly (U-Joint)			056-735-783

MTS Maintenance Parts

Model 329 6 DOF Light Truck Road Simulator Parts

329 6 DOF Light Truck Road Simulator Parts (continued)

CENTER LATERAL STRUT ASSEMBLY*	7	DESCRIPTION	PART NUMBER
Strut Assembly Lateral LT SVC Assembly (Strut)		Lateral Strut and End caps, Assembled	056-735-991
Strut Assembly Lateral LT SVC Assembly (Bearing)		Lateral Strut Rod-Eye and Bearings/Wear Components	056-735-981
Rod End Bearing Upgrade, 329 6DOF LT		Upgrade for Lateral Input Spherical Bearing	100-178-334

STEER STRUT ASSEMBLY*	8	DESCRIPTION	PART NUMBER
Strut Assembly Steer SVC Assembly (Strut)		Steer Strut and End caps, Assembled	056-737-291
Strut Assembly Steer SVC Assembly (Bearing)		Steer Strut Bearings/Wear Components	056-737-281
Strut Assembly Steer SVC Assembly (Housing/Bushing)		Wrist Joint Assembly for all Horizontal Struts	056-735-761

ACTUATOR ASSEMBLIES	DESCRIPTION	PART NUMBER
Vertical**		
LVDT Assembly-Closed Housing, Static Support, 7 in, 244.3	Vertical LVDT Assembly	045-362-001
Servo Valve Assembly-250 gpm 256.25A-02, 2.5 gpm	Three-Stage Vertical Servo Valve Assembly (includes two-stage)	033-561-002
Servo Valve-2.5 gpm 252.22G-01 4 Port	Two-Stage Vertical Servo Valve Assembly	056-006-201
244.3x Seal Kit 45 Kip, 1 in Cushion, Static Support	Vertical Actuator Seal Kit	048-978-504
Longitudinal**		
LVDT Assembly Closed Housing	Longitudinal LVDT Assembly	036-738-931
Servo Valve Assembly-50 gpm 256.05A-01, 1 gpm PL	Three-stage Longitudinal Servo Valve Assembly (includes two-stage)	033-366-501
Servo Valve-252.xx	Two-stage Longitudinal Servo Valve Assembly	056-006-101
Shipping Assembly-Seal Kit, 244.20, Structural	Longitudinal Actuator Seal Kit	008-710-075
Lateral**		
LVDT Assembly Closed Housing	Lateral LVDT Assembly	036-738-931
Servo Valve Assembly-256.05A-005, 1/2 gpm PLT	Three-Stage Lateral Servo Valve Assembly (includes two-stage)	033-366-508
Servo Valve- 252.xx	Two-stage Lateral Servo Valve Assembly	056-348-201
Seal Kit 244.21, Structural	Lateral Actuator Seal Kit	008-710-048
Steer**		
LVDT Assembly-Closed Housing 244.11/12 8 in SP	Steer LVDT Assembly	036-738-821
Servo Valve-15 gpm 252.25G-01 4 Port	Steer Servo Valve Assembly	056-006-501
Seal Kit, 244.12 Structural	Steer Actuator Seal Kit	008-710-047
Brake**		
LVDT Assembly-Closed Housing 244.11/12 8 in SP	Brake LVDT Assembly	036-738-821
Servo Valve-15 gpm 252.25G-01 4 Port	Brake Servo Valve Assembly	056-006-501
Seal Kit, 244.10	Brake Actuator Seal Kit	048-705-609


HSM FILTER KIT	DESCRIPTION	PART NUMBER
Filter Element (Pilot Pressure)	8 in lg b3 = 200, 150 psid	100-216-412

LUBRICANT	PART NUMBER
AEROSHELL 14 GREASE (tube for grease gun)	011-010-225
KRYTOX GREASE (tube for grease gun)	100-123-981

DYNAMIC BRAKE INTENSIFIER	PART NUMBER
Dynamic Brake Intensifier Air Inlet Filter	100-130-988
Replacement Seal Kit for Dynamic Intensifier	100-110-498

(329 6 DOF LT drawing on next page)

Model 329 6 DOF Light Truck Road Simulator Parts

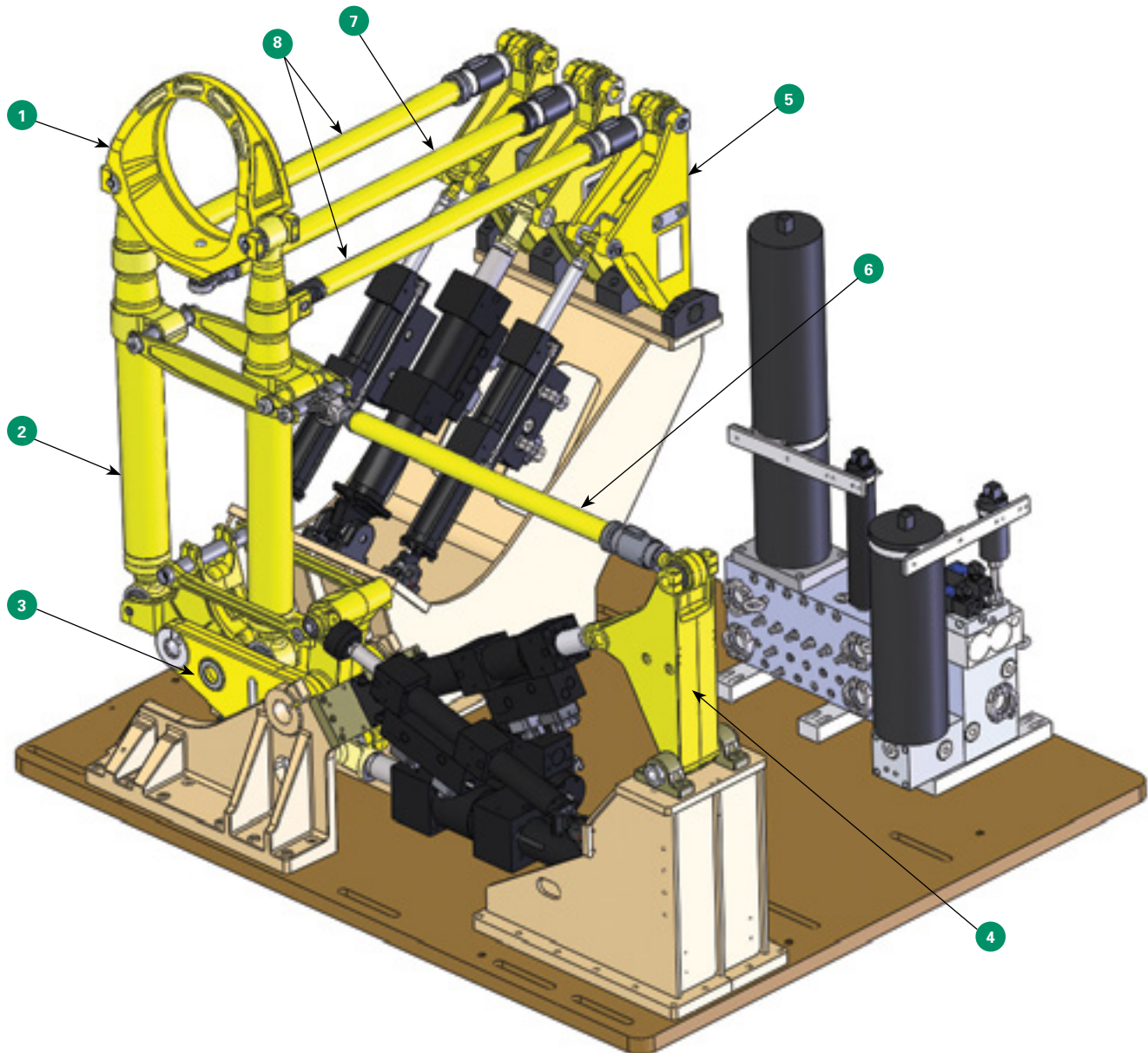

MTS Maintenance Parts

MAST Systems Parts

MAST

MODEL	DESCRIPTION	PART NUMBER
323.10	Vert table swivels	051-107-001
323.10	Lat/long actuator rod swivel	050-043-101
323.10	Actuator base swivels	049-577-401
323.10	Strut swivel- table end	050-043-102
323.10	Strut swivel- bellcrank end	050-043-101
323.10	Pillow block bearings	048-894-802

MODEL	DESCRIPTION	PART NUMBER
323.20	Vert table swivels	052-173-001
323.20	Lat/long actuator rod swivel	050-043-101
323.20	Actuator base swivels	049-577-401
323.20	Strut swivel- table end	051-107-002
323.20	Strut swivel- bellcrank end	050-043-101
323.20	Pillow block bearings	048-894-802

MODEL	DESCRIPTION	PART NUMBER
323.30	Vert table swivels	051-107-001
323.30	Lat/long actuator rod swivel	050-043-101
323.30	Actuator base swivels	049-577-401
323.30	Strut swivel- table end	050-043-102
323.30	Strut swivel- bellcrank end	050-043-101
323.30	Pillow block bearings	048-894-802

MODEL	DESCRIPTION	PART NUMBER
323.40	Vert table swivels	052-173-001
323.40	Lat/long actuator rod swivel	050-043-101
323.40	Actuator base swivels	049-577-401
323.40	Strut swivel- table end	051-107-002
323.40	Strut swivel- bellcrank end	050-043-101
323.40	Pillow block bearings	048-894-802

MODEL	DESCRIPTION	PART NUMBER
353.10	Bearing Assy- U-Joint (STD)	100-137-817
353.10	Boot - Rolling, 3" Stroke, 353.10	100-220-297

MODEL	DESCRIPTION	PART NUMBER
353.20	Swivel Assy-353.20-S	100-174-868
353.20	Boot - Convolute, 8.25" Stroke, 353.20	100-291-907
353.20	Boot - Upper, Convolute, 353.20	100-214-399
353.20	Upgrade Kit - Bellows, 353.20 Strut	100-301-438

MTS Maintenance Parts

320 Seal Kits

320.025 First Road and Small Car Systems

DESCRIPTION	PART NUMBER
Actuator Seal Kit	040-262-601 
Actuator Seal Kit for Environmental Systems	056-079-329
Servovalve	042-795-401

320.035 Large Car Systems

DESCRIPTION	PART NUMBER
Actuator Seal Kit	040-262-601 
Actuator Seal Kit for Environmental Systems	056-079-329
Servovalve	033-366-601

320.050 Light Truck Standard Performance Systems


DESCRIPTION	PART NUMBER
Actuator Seal Kit	040-262-601 
Actuator Seal Kit for Environmental Systems	056-079-329
Servovalve	033-366-601

320.050 Light Truck High Performance Systems

DESCRIPTION	PART NUMBER
Actuator Seal Kit	040-262-601 
Actuator Seal Kit for Environmental Systems	056-079-329
Servovalve	033-206-921

EXCHANGE & REPAIR

The exchange programs provide a cost-effective turnaround option to extend the life of your test equipment. You'll receive a fully-warranted, remanufactured component to replace the worn component. A service contract and installation is required. The field service engineer will return the worn component to MTS as a core. This offer will be sold through a service contract.

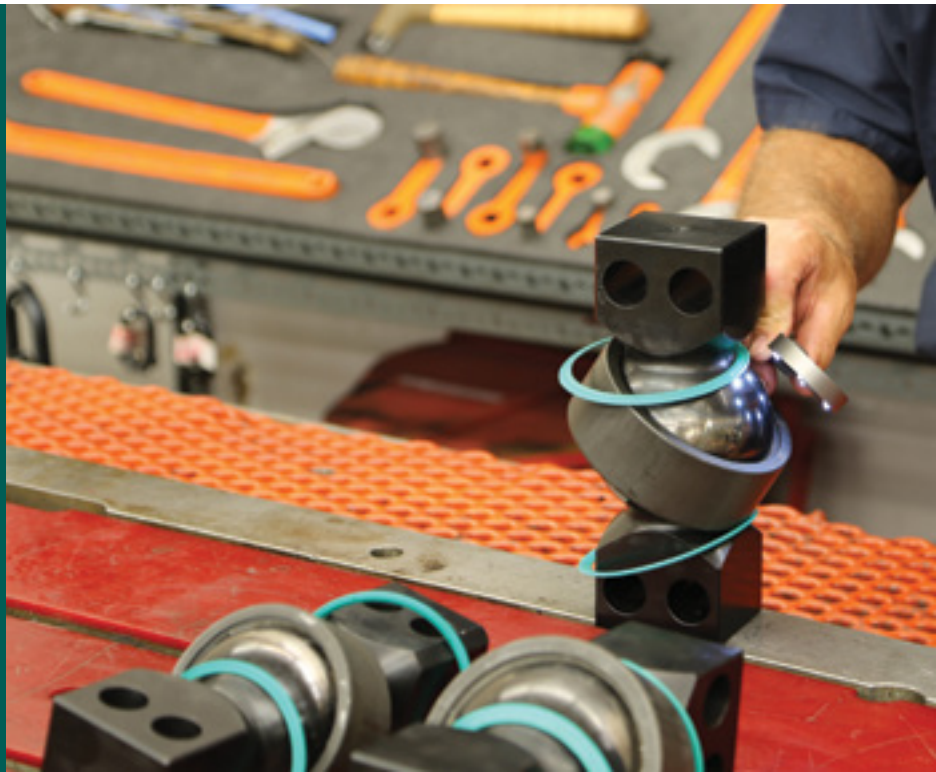

Exchange & Repair Programs

Introduction	192
Just In Case™	193
Electronics	194-201
Onsite Exchange Programs	202-204
MTS Series 252 & 256 Servovalves	205-206
Remanufacturer Programs - Actuators	207

Exchange & Repair Programs

FlexTest Controllers/Insight Electromechanical Load Frames

PRODUCTIVITY INSURANCE FOR YOUR MTS DIGITAL CONTROLLERS



Controller Spare Kits – Just in Case

FLEXTEST 40, 60, 100, & 200	PART NUMBER
Spares Kit - Just in Case, FlexTest 40 -NV350/60 w/o Handset, Basic	100-246-613
Spares Kit - Just in Case, FlexTest 40 -NV350/60 w/Handset & Accel Cards, Advanced	100-246-634
Spares Kit - Just in Case, FlexTest 40 -NV350/60—Multiple Content w/Handset, Extended	100-246-635
Spares Kit - Just in Case, FlexTest 40 -NV350/ 60—Multiple Content w/o Handset, Extended	100-246-636
Spares Kit - Just in Case, FlexTest 100 & 200—Multiple Content w/o Handset, Basic	100-200-812
Spares Kit - Just in Case, FlexTest 200 for Aero Structures	100-195-413

We make custom spares kits for the FlexTest IIIm, FlexTest 40, 60, 100, & 200 Controllers and other systems. Contact your sales engineer for details.

Insight Spares Kits

NEW MTS DESCRIPTION	PART NUMBER
Spares Kit - Insight 1K / 2K	100-187-792
Spares Kit - Insight 5K / 10K	100-187-793
Spares Kit - Insight 30K / 50K	100-187-794
Spares Kit - Insight 100K	100-187-795
Spares Kit - Insight 150K	100-187-796
Spares Kit - Insight 200K	100-187-797
Spares Kit - Insight 300K	100-187-798
Spares Kit - Insight 50K Wide	100-187-799
Spares Kit - XLT Extensometer 100R STD	100-187-800
Spares Kit - XLT Extensometer 100R Ext	100-187-801

All Insight spares kits contain both consumable and essential components that are prone to wear including: fuses, bolts, grease, switches, bellows and brushes.

Exchange & Repair Programs

Electronics

Electronic Exchange Program

The MTS electronic exchange program provides a cost-effective, quick turnaround, replacement option for electronic boards. You can purchase a fully warranted and tested product in exchange for your failed unit.

FEATURES

- » Good-as-new performance for less than replacement price
- » Shipment within 24 hours (Monday–Friday), if item is in stock
- » A fixed price for your repair
- » A 12 month warranty
- » The latest product revision levels
- » Each product is cosmetically remanufactured
- » MTS provides a convenient return shipping carton and label to make it easy for you to return the failed module to MTS.

HOW THE EXCHANGE PROGRAM WORKS

1. Call MTS and ask for Order Services.
2. You will be connected with an order specialist who will help you select the proper product. You will be required to provide a secured payment method equal to the core charge by one of the following methods:
 - » A purchase order
 - » A credit card
 - » An MTS Value Payment Plan contract number or Pay-As-You-Go Plan.
3. MTS will ship your exchange part, if it is in stock, within 24 hours.
4. Once you receive the exchange product, you must return your defective item to MTS within 30 days of initial order placement.
5. Ship the defective item in the same package in which you received the exchanged item. Make sure the Return Authorization Number (RAN) is clearly marked on the shipping labels and enclosed within the package.
6. If MTS receives your defective item within 30 days of initial order placement and if it can be remanufactured, MTS will invoice the exchange price.
7. In addition to the exchange price, if MTS does not receive your defective item within 30 days of order placement, or if it does not qualify for remanufacture, MTS will invoice the core charge. (In a warranty exchange only, a core charge would apply if the item is unreturned.)

Electronic Repair Program

An alternative to the exchange program is to send your failed unit to our factory for repair.

FEATURES 3 OR 10 DAY TURNAROUND

The electronic repair program offers a 3-day or 10-day repair service. Troubleshooting and repair are performed at an MTS facility. Under our standard offering, your product will be shipped back to you in 10 business days. We also have an expedited service available that provides a 3-day turnaround after your product is received at the MTS facility. The expedite fee is 20% of the normal repair price. Items repaired have a 90 day warranty. Some items may only be available as part of the exchange program. Contact your sales engineer for details.

HOW THE REPAIR PROGRAM WORKS

1. Call MTS and ask for order services. MTS will issue a Return Authorization Number.
2. Mark your shipping package with the RAN.
3. Return your item to MTS for repair.

REPAIR PROGRAM INSPECTION FEES

When a product to be repaired is received by MTS, the first step we take is to inspect it to determine the scope of the repair. If the repair is not possible, or if the estimated repair cost significantly exceeds the scope originally quoted, you will be contacted with that information and offered the opportunity to reconsider the original repair request.

MTS reserves the right to charge a fee for Pre-Repair Inspections (fees may vary). If it is decided that the repair should not be completed, you will be asked to decide whether the product should be returned or if MTS should arrange for appropriate disposal. In some cases MTS may charge for the disposal of unrepared or unrepairable products.

Exchange & Repair Programs

Electronics

Ordering

To place an order for either an exchange or repair product, please call MTS and ask for Order Services.

BE PREPARED WITH THE FOLLOWING INFORMATION WHEN YOU CALL MTS

- » Your MTS site number
- » The MTS model number of the product you want to exchange or repair
- » The MTS assembly number of the product you want to repair
- » The MTS serial number of the product you want to exchange or repair
- » The MTS system in which the product is used
- » The symptoms of the problem you are experiencing
- » Your PO number or method of payment.

COMPATIBILITY

Inclusion in this catalog does not imply that a product is electrically or physically compatible with your existing equipment. Many MTS product configurations have been customized to meet specific system testing needs. If you are unsure if your part can be exchanged, contact your service engineer or the MTS Customer Care Center.

Established Systems

We attempt to repair or exchange all possible systems. However, we may not be able to repair some of our more established systems for various reasons. For these situations we will make a best effort at the time of the call. Please contact your application engineer for pricing and availability.

Exchange & Repair Programs

Electronics

Model 407

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
407	Base Unit, Pwr Sup, Proc, PIDF, CE	407BASE-EX
407	Base Unit, Pwr Sup, Proc, PIDF, AMP, CE	407BASEAMP-EX
407.12	DC Conditioner	049-334-101-EX
407.14	AC Conditioner (repair only)	
407.14B	AC Conditioner	100-012-737-EX
407.15	Three-Stage Valve Driver	049-334-501-EX
407.16	Valve Driver (repair only)	
407.16B	Valve Driver	053-306-801-EX
407.36	Communications Module	049-333-901-EX

Models 408-413 are no longer supported.

Motorola Power PC CPUs

Limited Supply (Contact your MTS application engineer)

MODEL	CONFIG'D MATERIAL #	TOP LEVEL #	EXCHANGE PART NUMBER
498.93-1	052-244-101	052-244-101	011-964-709-EX
498.93-1	052-748-001	052-748-001	011-964-709-EX
498.94-1	054-396-401	054-396-401 (kit only)	011-970-862-EX
498.94-1	054-396-401	054-396-402 (kit only)	011-970-862-EX
498.96-1	056-027-601	100-040-837	100-003-121-EX
498.96-1	056-027-602	100-040-837	100-003-121-EX
498.96-1	100-009-996	100-009-996	100-003-121-EX
498.96-2	056-023-701	100-040-838	011-973-553-EX
498.96-2	056-023-702	100-040-838	011-973-553-EX
498.96-2	056-023-703	100-040-838	011-973-553-EX
498.96-1B			100-100-908-EX
498.96-3	056-309-601	100-053-993	100-044-858-EX
498.96-3	056-309-602	100-053-993	100-044-858-EX
498.96-3	056-309-603	100-053-993	100-044-858-EX
498.96-4			100-053-122-EX
498.96-5B			100-100-909-EX
498.96-5D			057-174-301-EX
498.96-5			100-074-939-EX
498.96-6			100-074-940-EX
498.96-7			100-102-106-EX
498.96-7B			057-159-301-EX

Model 497 & 498 Modules are shipped unconfigured.

Exchange & Repair Programs

Electronics

Flat-Trac II, III

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
497.13B	AC Conditioner	100-005-333-EX
497.15B	Three-Stage Valve Driver	100-019-313-EX
497.22C	Dual DC Conditioner	055-943-401-EX
497.26B	Dual Channel Valve Driver	054-102-401-EX
498.05	Serial Port (Used with the 498.32)	011-487-360-EX
498.06	Serial Port with Ethernet	046-778-901-EX
498.47	PAMUX Bus Controller	011-487-343-EX
498.55	MIOP II	046-555-501-EX
498.65	ADDA	048-077-001-EX
498.70	16 Channel Digital I/O	044-127-701-EX
498.00	16 Channel Digital I/O Transition	049-808-401-EX
498.00	Serial Transition	054-170-601-EX
498.00	Eight-Channel D/A Daughter Card	047-849-602-EX
498.00	Eight-Channel A/D Daughter Card	100-042-515-EX

Note: 497 modules are backward compatible for letter releases.

Flat-Trac III

Limited Supply (Contact your MTS application engineer)

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
497.15	Three-Stage Valve Driver	047-178-601-EX
497.22A	Dual DC Conditioner	046-359-501-EX
497.22B	Dual DC Conditioner	048-757-601-EX
498.06	Serial Port	046-779-201-EX
498.32	Processor Module	011-487-365-EX
498.36	Processor Module	011-487-392-EX
498.71	GRES II	049-808-501-EX
498.00	System Digital I/O Transition	053-112-101-EX
498.00	Eight-Channel Filtered Analog Input	051-064-610-EX
498.00	Eight-Channel Filtered Analog Input	051-064-630-EX
498.00	Eight-Channel Filtered Analog Input	051-064-670-EX
498.00	Eight-Channel A/D Daughter Card	049-676-901-EX

RPC III

Limited Supply (Contact your MTS application engineer)

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
498.60	DCON	043-669-201-EX
498.60B	DCON	052-736-901-EX
498.36	Processor Module	011-487-392-EX
498.37	MJME 147SA-1 Processor	011-487-372-EX
498.37	MJME 147-012 Processo	011-963-434-EX
498.61	DSPAD	043-669-101-EX
498.61B	DSPAD	050-106-101-EX
498.70	16 Channel Digital I/O	044-127-701-EX

Exchange & Repair Programs

Electronics

FlexTest IIm, IIs, CTM & CTC

Limited Supply (Contact your MTS application engineer)

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
497.13B	AC Conditioner	100-005-333-EX
497.15B	Three-Stage Valve Driver	100-019-313-EX
497.22C	Dual DC Conditioner	055-943-401-EX
497.26B	Dual Channel Valve Driver	054-102-401-EX
498.34B	Station Controller	100-048-996-EX
497.36	FlexTest Communications	046-791-101-EX
493.50	ADDA II	100-029-330-EX
498.65	ADDA	048-077-001-EX
498.70	16 Channel Digital I/O	044-127-701-EX
498.71B	GRES III	100-111-372-EX
498.00	16 Channel Digital I/O Transition	049-808-401-EX
498.00	System Digital I/O Transition	053-112-101-EX
498.00	Serial Transition	054-170-601-EX

Note: 497 modules are backward compatible for letter releases.

FlexTest IIm, IIs, CTM & CTC

Limited Supply (Contact your MTS application engineer)

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
497.12	DC Conditioner	045-771-801-EX
497.14	AC Conditioner	046-481-501-EX
497.15A	Three-Stage Valve Driver	047-178-601-EX
497.22A	Dual DC Conditioner	046-359-501-EX
497.22B	Dual DC Conditioner	048-757-601-EX
497.26	Dual Channel Valve Driver	045-771-701-EX
497.34A	Station Controller	045-772-001-EX
498.71	GRES II Transition Module	049-808-301-EX
498.71	GRES II	049-808-501-EX
498.71B	GRES III	054-170-401-EX
498.71C	GRES III	100-111-372-EX
498.00	Eight-Channel Filtered Analog Input	051-064-610-EX
498.00	Eight-Channel Filtered Analog Input	051-064-630-EX
498.00	Eight-Channel Filtered Analog Input	051-064-670-EX
498.00	Eight-Channel Filtered Analog Input (repair only)	051-064-690-EX
498.00	Eight-Channel A/D Daughter Card	049-676-901-EX
498.00	Eight-Channel A/D Daughter Card	100-042-515-EX
498.00	Eight-Channel D/A Daughter Card	047-849-602-EX

Note: 497 modules are backward compatible for letter releases (A, B, C...).

Exchange & Repair Programs

Electronics

FlexTest 40, 60, 100, & 200

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
493.80	493.80, Encoder Interface Board	056-085-301-EX
494.04	Power Supply, 494.04 Chassis (depending on what the customer has for an assembly)	100-173-146-EX or 100-231-674-EX
494.06	494.06 Power Supply	100-178-428-EX
494.10	494.10 Power Supply (1 round fan)	100-155-696-EX
494.10	494.10 Power Supply (4 fans)	100-205-894-EX
494.16	494.16 DUC/Valve Driver	100-174-701-EX
494.20	494.20 Power Supply (2 round fans)	100-167-738-EX
494.20	494.20 Power Supply (7 fans)	100-206-233-EX
494.21	494.21 EL DUC	100-176-671-EX
494.25	494.25 Single Conditioner	100-181-455-EX
494.26	494.26, Dual DUC	100-173-510-EX
494.40	494.40 VME I/O Car, IEEE+/-15	100-169-815-EX
495.40B	494.40 VME I/O Car, IEEE+/-15	100-206-704-EX
494.41	EEE 1011 Model 494.41 System I/O Mod	100-179-562-EX
494.44	494.44 Module	100-184-696-EX
494.45	494.45 Eight-Channel A/D	100-210-991-EX
494.46	494.46 Eight-Channel D/A	100-173-095-EX
494.47	494.47 Two-Channel UART/Encoder	100-161-906-EX
494.49	494.49 Encoder	100-170-372-EX
493.72	493.72 Digital I/O Transition SHLD	100-003-853-EX
493.73	493.73 HPU Transition Board	100-003-855-EX
493.74	493.74 HSM Transition, w/o Label	100-039-452-EX
494.74	494.74 Dual HSM w/o Label	100-169-789-EX
494.75	494.75, Analog In, No Label	057-203-101-EX
494.76	494.76, Analog Out, No Label	057-203-201-EX
494.76B	494.76, Analog Out, No Label	057-442-101-EX
494.79	494.79, 8 Channel, Univ Driver	100-153-566-EX
494.96-1B	Processor – 494.96-1B, 5110 IEEE VX5.5 CPU0	057-204-601-EX
494.96-2A	Processor – 494.96-2A, 5500 IEEE, CPU 0	057-204-701-EX

Note: 497 modules are backward compatible for letter releases (A, B, C...).

Exchange & Repair Programs

Electronics

FlexTest GT

Limited Supply (Contact your MTS application engineer)

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
493.10*	Chassis Assembly – (repair only)	
493.10	Power Supply Assembly – Lambda	100-031-143-EX
493.14B	Valve Driver for 252.x or 5x Valves	100-109-855-EX
493.15B	Three-Stage Valve Driver	100-109-974-EX
493.20	493.20 Power Supply	100-031-183-EX
493.21BC	Digital Universal Conditioner	100-109-976-EX
493.25B	Digital Universal Conditioner	100-109-659-EX
493.25C	Digital Universal Conditioner	100-207-810-EX
493.40	I/O Carrier Board	100-003-854-EX
493.45	Analog Input Module – 6 Channel A/D	051-904-401-EX
493.46	Analog Output Module – 6 Channel D/A	052-907-701-EX
493.47	High Speed Encoder Module	055-351-501-EX
493.48	Accelerometer Conditioner Module	100-011-578-EX
493.71	RS485	100-003-832-EX
493.72	Digital I/O Transition Module	100-003-853-EX
493.73	HPU Transition Board	100-003-855-EX
493.74	493.74 HSM Transition Station 1 & 2	100-039-452-EX
493.75	Analog Input BNC Transition Module, 6 Channel	100-008-984-EX
493.76	Analog Output BNC Transition Module, 6 Channel	100-008-986-EX
493.77	Filter Input Transition Board, 6 Channel 300 HZ	056-013-705-EX
493.78	Accelerated Input	100-013-651-EX
498.71C	Global resource III Module	100-111-372-EX

FlexTest GT

Limited Supply (Contact your MTS application engineer)

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
493.02	Chassis Assembly – Repair only	100-061-333-EX
493.07	Hydraulic Interface Box 24 V	049-969-401-EX
493.07	Hydraulic Interface Box 24 V (407 Controller Version)	049-969-402-EX
493.07	Hydraulic Interface Box 115 V	049-969-403-EX
493.07	Hydraulic Interface Box 115 V (407 Controller Version)	049-969-404-EX
493.14	Two-Stage Valve Driver	048-934-001-EX
493.14B	Valve Driver for 252.x or 5x Valves	100-109-855-EX
493.15	Three-Stage Valve Driver	053-782-701-EX
493.15B	Three-Stage Valve Driver	100-109-974-EX
493.21B	Digital Universal Conditioner	053-476-101-EX
493.21BC	Digital Universal Conditioner	100-109-976-EX
493.25	Digital Universal Conditioner	100-033-468-EX
493.25B	Digital Universal Conditioner	100-109-659-EX
493.25C	Digital Universal Conditioner	100-207-810-EX
493.40	I/O Carrier Board	100-003-854-EX
493.42	I/O Module Assembly	100-056-617-EX
493.45	Analog Input Module – 6 Channel A/D	051-904-401-EX
493.46	Analog Output Module – 6 Channel D/A	052-907-701-EX

Note: 493 modules are backward compatible for letter releases (A, B, C...), with the exception of the 493.21A.

Exchange & Repair Programs

Electronics

T/RAC & FlexTest DSC

Limited Supply (Contact your MTS application engineer)

MODEL NUMBER	DESCRIPTION	EXCHANGE PART NUMBER
497.12	DC Conditioner	045-771-801-EX
497.14	AC Conditioner	046-481-501-EX
497.15A	Three-Stage Valve Driver	047-178-601-EX
497.15B	Three-Stage Valve Driver	100-019-313-EX
497.22A	Dual DC Conditioner	046-359-501-EX
497.22B	Dual DC Conditioner	048-757-601-EX
497.22C	Dual DC Conditioner	055-943-401-EX
497.26	Dual Channel Valve Driver	045-771-701-EX
497.34	Station Controller	045-772-001-EX
497.34B	Station Controller	100-048-996-EX
497.35	Analog Interface	045-223-401-EX
497.36	FlexTest Communications	046-791-101-EX

T/RAC & FlexTest DSC

Limited Supply (Contact your MTS application engineer)

MODEL NUMBER	DESCRIPTION	EXCHANGE PART NUMBER
498.04	Global Resource Board	046-480-301-EX
498.04	Global Resource Board (FlexTest)	047-025-401-EX
498.33	Processor Module	011-487-337-EX
498.34	Processor Module	011-487-359-EX
498.42	1 Mb EPROM Memory	011-487-350-EX
498.54	MIOP I	043-774-501-EX
498.55	MIOP II	046-555-501-EX
498.70	16 Channel Digital I/O	044-127-701-EX

Note: 497 modules are backward compatible for letter releases.

Exchange & Repair Programs

Onsite Exchange Programs

The onsite exchange programs provide a cost-effective quick turnaround option to extend the life of your test equipment. You'll receive a fully-warranted, remanufactured component to replace the worn component. A service contract and installation by MTS-trained personnel is required. The field service engineer will return the worn component to MTS as a core. This offer will be sold through a service contract. No core charge is needed or charged.

ADVANTAGES OF AN ONSITE EXCHANGE

- » Good-as-new performance for less than new replacement cost
- » Latest product revision levels
- » Exchange components are remanufactured and tested to stringent standards
- » Minimizes downtime as this service may be scheduled with other onsite service activities
- » Extends the overall life of your test equipment
- » Saves labor and downtime costs by exchanging components instead of replacing equipment
- » Provides a perfect opportunity to add MTS Routine Maintenance and Fluid Assessment Services to monitor system performance moving forward. These services extend equipment life and help prevent untimely failures
- » Field delivery through MTS personnel only, eliminating your staff's need to package the core for return.



HPU pump/
motor module

Series 329 System Actuator Onsite Exchange

Available for 4DOF and 6DOF systems. Please consult your service sales representatives for more information.

CAPABILITY

Inclusion of a product in this catalog does not imply that it is mechanically or physically compatible with your existing equipment. Many MTS product configurations have been customized to meet specific system testing needs. Let us help you determine the compatibility required for an exchange.



329 system

Exchange & Repair Programs

Onsite Exchange Programs

505.90 & 505.180 G1 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.90 G1	380-480V	T.H.	3000	057-492-601EX
505.180 G1	380-480V	T.H.	3000	057-492-602EX
505.90 G1	380-480V	S.S.	3000	057-492-603EX
505.180 G1	380-480V	S.S.	3000	057-492-604EX
505.90 G1	575V	T.H.	3000	057-492-605EX
505.180 G1	575V	T.H.	3000	057-492-606EX
505.90 G1	575V	S.S.	3000	057-492-607EX
505.180 G1	575V	S.S.	3000	057-492-608EX

505.07 G1 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.07 G1	190/220-240V	T.H.	3000	057-492-701EX
505.07 G1	200-208V	T.H.	3000	057-492-702EX
505.07 G1	380-480V	T.H.	3000	057-492-703EX
505.07 G1	575V	T.H.	3000	057-492-704EX
505.07 G1	190/220-240V	S.S./T.H.	3000	057-492-705EX
505.07 G1	200-208V	S.S./T.H.	3000	057-492-706EX
505.07 G1	380-480V	S.S./T.H.	3000	057-492-707EX
505.07 G1	575V	S.S./T.H.	3000	057-492-708EX

505.11 G1 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.11 G1	190/220-240V	T.H.	3000	057-492-801EX
505.11 G1	200-208V	T.H.	3000	057-492-802EX
505.11 G1	380-480V	T.H.	3000	057-492-803EX
505.11 G1	575V	T.H.	3000	057-492-804EX
505.11 G1	190/220-240V	S.S./T.H.	3000	057-492-805EX
505.11 G1	200-208V	S.S./T.H.	3000	057-492-806EX
505.11 G1	380-480V	S.S./T.H.	3000	057-492-807EX
505.11 G1	575V	S.S./T.H.	3000	057-492-808EX

505.20 G1 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.20 G1	190/220-240V	T.H.	3000	057-492-901EX
505.20 G1	200-208V	T.H.	3000	057-492-902EX
505.20 G1	380-480V	T.H.	3000	057-492-903EX
505.20 G1	575V	T.H.	3000	057-492-904EX

505.30 G1 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.30 G1	190/220-240V	S.S.	3000	057-493-001EX
505.30 G1	200-208V	S.S.	3000	057-493-002EX
505.30 G1	380-480V	S.S.	3000	057-493-003EX
505.30 G1	575V	S.S.	3000	057-493-004EX
505.30 G1	200-208V	T.H.	3000	057-493-006EX

Exchange & Repair Programs

Onsite Exchange Programs

505.90 & 505.180 G2 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.90 & 180 G2	380-480V	T.H.	3000	057-492-601EX
505.90 & 180 G2	380-480V	S.S.	3000	057-492-603EX
505.90 & 180 G2	575V	T.H.	3000	057-492-605EX
505.90 & 180 G2	575V	S.S.	3000	057-492-607EX

505.07 G2 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.07 G2	190/220-240V	T.H.	3000	057-492-701EX
505.07 G2	200-208V	T.H.	3000	057-492-702EX
505.07 G2	380-480V	T.H.	3000	057-492-703EX
505.07 G2	575V	T.H.	3000	057-492-704EX
505.07 G2	190/220-240V	S.S./T.H.	3000	057-492-705EX
505.07 G2	200-208V	S.S./T.H.	3000	057-492-706EX
505.07 G2	380-480V	S.S./T.H.	3000	057-492-707EX
505.07 G2	575V	S.S./T.H.	3000	057-492-708EX

505.11 G2 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.11 G2	190/220-240V	T.H.	3000	057-492-801EX
505.11 G2	200-208V	T.H.	3000	057-492-802EX
505.11 G2	380-480V	T.H.	3000	057-492-803EX
505.11 G2	575V	T.H.	3000	057-492-804EX
505.11 G2	190/220-240V	S.S./T.H.	3000	057-492-805EX
505.11 G2	200-208V	S.S./T.H.	3000	057-492-806EX
505.11 G2	380-480V	S.S./T.H.	3000	057-492-807EX
505.11 G2	575V	S.S./T.H.	3000	057-492-808EX

505.20 G2 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.20 G2	190/220-240V	T.H.	3000	057-492-901EX
505.20 G2	200-208V	T.H.	3000	057-492-902EX
505.20 G2	380-480V	T.H.	3000	057-492-903EX
505.20 G2	575V	T.H.	3000	057-492-904EX

505.30 G2 Series

HPU SERIES	VOLTAGE	T.H./S.S.	PRESSURE (PSI)	EXCHANGE PART NUMBER
505.30 G2	190/220-240V	S.S.	3000	057-493-001EX
505.30 G2	200-208V	S.S.	3000	057-493-002EX
505.30 G12	380-480V	S.S.	3000	057-493-003EX
505.30 G2	575V	S.S.	3000	057-493-004EX

Exchange & Repair Programs

MTS Series 252 & 256 Servovalves

The MTS servovalve exchange program provides a cost effective, rapid turnaround replacement option for non-functioning servovalves. You can purchase a remanufactured, tested and fully-warranted product in exchange for your used unit to ensure top quality performance of your complete system.

FEATURES

- » Good-as-new performance for less than new product replacement price
- » Shipment within 24-hours (Monday-Friday)
- » 12 Month Warranty
- » Remanufactured to stringent MTS specifications using the latest product revision levels where applicable
- » Convenient return shipping carton and label provided for easy return of the used servovalve core.

COMPATIBILITY

Inclusion of a product in this catalog does not imply that it is mechanically or physically compatible with your existing equipment. Many MTS product configurations have been customized to meet specific system testing needs. Let us help you determine the compatibility required for exchange or if functional replacement/upgrade is the best solution. MTS provides fit-form-function interchange exchange using the G-Type Model servovalves of the same 4 or 5 port configuration and flow size. The G-Type Model servovalves are completely interchangeable with the legacy C-Type Model servovalves.



EXCLUSIONS

The following types of exchanges are not supported/allowed:

- » Exchange of 1970's vintage Model A series servovalves
- » Exchange of vintage 1980's A-Series copper-body servovalves
- » Change in 4 or 5 port configuration servovalves
- » Change in servovalve flow size
- » Exchange of Series 252.5x servovalves
- » Exchange of high response (HR) servovalves
- » Repair of A or C-Type Model servovalves. (C-Type servovalves can be Exchanged for G-Type).



How the valve exchange program works

1. Contact the MTS Order Services Team to help you select the proper remanufactured exchange valve. You will be asked to provide a secured payment method to cover the value of both the remanufactured valve and the value of your used core to be returned. Payment can be secured by one of the following methods.
 - » A purchase order
 - » A credit card
 - » An MTS service plan contract number
2. MTS will ship the remanufactured exchange valve that includes the used valve core return packaging.
3. Upon receiving the exchange valve, you must return your used valve core to MTS within 30 days of order placement.
4. Ship your used valve core back to MTS using the white inner carton provided. Ensure that the Return Authorization Number (RAN) is clearly indicated on the shipping labels and enclosed within the package.
5. When MTS receives your used valve core within 30 days of order placement, your final invoice will have only the exchange valve price.
6. When MTS does not receive your used valve core within 30 days of order placement or the used core is found to not qualify for remanufacture, your final invoice will have the exchange valve price plus the core value price as set in the initial purchase order.

NOTE FOR ALL EXCHANGES

All MTS remanufactured product exchanges require a core value to be associated with the purchase. Your purchase order must explicitly include both the exchange product price plus the core value. Please consult with your order services representative to answer any questions about exchange product and related core pricing. Upon returning your good condition and qualifying used core to MTS within 30 days of issuing the purchase order, MTS will not levy the core charge and you will pay for only the remanufactured exchange product.

FACTORY REPAIR SERVICE 256.XX

1. Pilot
2. Basic cleaning/tuning/validation
3. Full rebuild

Exchange & Repair Programs

MTS Series 252 Servovalves

Servovalves Models 252.2x - Exchange

MODEL	DESCRIPTION	ALSO APPLIES	EXCHANGE PART NUMBER
252.21G-01	1 gpm, 4 Port Servovalve	252.21C-01	056-006-101-EX
252.21G-04	1 gpm, 5 Port Servovalve	252.21C-04	056-006-104-EX
252.22G-01	2.5 gpm, 4 Port Servovalve	252.22C-01	056-006-201-EX
252.22G-04	2.5 gpm, 5 Port Servovalve	252.22C-04	056-006-204-EX
252.23G-01	5 gpm, 4 Port Servovalve	252.23C-01	056-006-301-EX
252.23G-04	5 gpm, 5 Port Servovalve	252.23C-04	056-006-304-EX
252.24G-01	10 gpm, 4 Port Servovalve	252.24C-01	056-006-401-EX
252.24G-04	10 gpm, 5 Port Servovalve	252.24C-04	056-006-404-EX
252.25G-01	15 gpm, 4 Port Servovalve	252.25C-01	056-006-501-EX
252.25G-04	15 gpm, 5 Port Servovalve	252.25C-04	056-006-504-EX

Servovalves Models 256.xx - Repair Services

MODEL	DESCRIPTION	EXCHANGE PART NUMBER
256.09	Pilot	Contact MTS
256.09	Clean/Tunning/Validate	Contact MTS
256.09	New Spool/Full Rebuild	Contact MTS
256.18	Pilot	Contact MTS
256.18	Clean/Tunning/Validate	Contact MTS
256.18	New Spool/Full Rebuild	Contact MTS
256.25	Pilot	Contact MTS
256.25	Clean/Tunning/Validate	Contact MTS
256.25	New Spool/Full Rebuild	Contact MTS

Remanufacture Programs

Actuators

MTS offers you an alternative to buying new actuators or constantly repairing your existing ones. The MTS Actuator Remanufacture Program follows a strict protocol to restore the performance and reliability of your old actuators. We are so confident in the performance and quality of our remanufactured actuators that we provide the same warranty that you would receive when purchasing a new actuator.

Here's What You Get When You Choose to Remanufacture:

Your actuator will be updated to meet the newest design standards, as required. Improvements that were incorporated after your unit was built are retrofitted, bringing performance to "better than original" in many instances.

All remanufacturing work is performed in our factory to our specifications, using MTS parts. This is why we can assure that your actuator will perform as intended after the remanufacture. Components subjected to wear are resurfaced or replaced. That's how we can promise you that you'll receive the same or better service life from your remanufactured actuator.

Service Contracts Eliminate Procurement Hassle


You can include actuator remanufacturing, as well as other MTS services, on an MTS service contract. This solution saves money and time, and avoids the hassle of issuing separate purchase orders.

The Actuator Remanufacture Program covers the basic cylinder assembly.

HERE ARE THE SPECIFICS:

- » End cap bearings are replaced or resurfaced
- » Piston rods are stripped and thermal spray coated or rechromed
- » Seals are replaced with the latest design
- » Other damaged or worn components are repaired or replaced
- » Clearances are inspected and brought up to new MTS actuator specifications
- » Remanufactured units are pressure checked and performance tested

Other components such as servovalves, LVDTs, manifolds, accumulators, and swivels may be repaired or replaced at an additional charge.

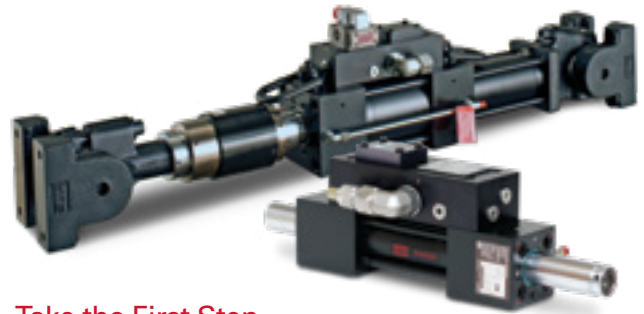

Take the First Step

Restore the performance of your worn hydraulic actuators by contacting MTS about actuator remanufactures today. Call us at 1-800-328-2255 or 1-952-937-4000 or email info@mts.com for more information.

BENEFITS

- » Economical solution
- » Extended service life
- » Like-new, or better than original performance
- » MTS quality and reliability
- » Quick turnaround service for some models
- » Expedited service upon request


QUICK TURN SERVICE

At our U.S. facility, we can process your remanufacture within approximately 10 business days of actuator receipt, and then send it back to you. The following models and stroke lengths currently eligible for this service:

- » Stroke length 150 mm (6 in)
- » Stroke length 250 mm (10 in)
- » 244.11 – 15 kN (3.3 kip)
- » 244.12 – 25 kN (5.5 kip)
- » 244.20 – 67 kN (15 kip)
- » 244.21 – 50 kN (11 kip)
- » 244.22 – 100 kN (22 kip)
- » 244.23 – 150 kN (35 kip)
- » 244.31 – 250 kN (55 kip)

EXPEDITE SERVICE

MTS may be able to offer expedited service at an additional charge. Our goal is to provide you with a quality product at a reliable lead time. Please discuss this option with your Service Sales Representative.

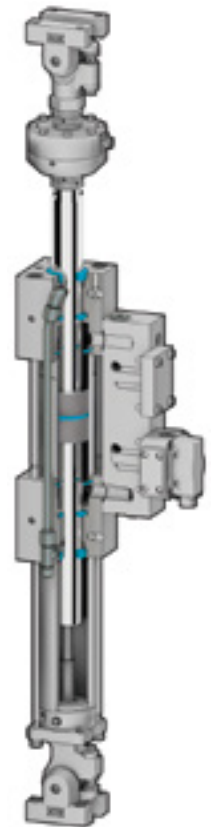

Company Locations

Regional Business Centers

THE AMERICAS

MTS Systems Corporation

14000 Technology Drive
Eden Prairie, MN 55344-2290
USA
Telephone: 952-937-4000
Toll Free: 800-328-2255
E-mail: info@mts.com
Internet: www.mts.com

EUROPE

MTS Systems France

BAT EXA 16
16/18 rue Eugène Dupuis
94046 Créteil Cedex
France
Telephone: +33-(0)1-58 43 90 00
E-mail: contact.france@mts.com

MTS Systems GmbH

Hohentwielsteig 3
14163 Berlin
Germany
Telephone: +49-(0)30 81002-0
E-mail: euinfo@mts.com

MTS Systems S.R.L. socio unico

Strada Pianezza 289
10151 Torino
Italy
Telephone: +39-(0)11 45175 11 sel. pass.
E-mail: mtstorino@mts.com

MTS Systems Norden AB

Datavägen 37b
SE-436 32 Askim
Sweden
Telephone: +46-(0)31-68 69 99
E-mail: norden@mts.com

MTS Systems Limited. UK

40 Alan Turing Road
Surrey Research Park
Guildford
Surrey
GU2 7YF
United Kingdom
Telephone: +44 (0) 1483 446500
E-mail: mtsuksales@mts.com

ASIA/PACIFIC

MTS Japan Ltd.

ArcaCentral Bldg. 8F
1-2-1 Kinshi, Sumida-ku
Tokyo 130-0013
Japan
Telephone: 81-3-6658-0901
E-mail: mts-jinfo@mts.com

MTS Korea, Inc.

4th F, ATEC Tower, 289,
Pankyo-ro, Bundang-gu
Seongnam-si
Gyeonggi-do 13488,
Korea
Telephone: 82-31-728-1600
E-mail: mtsk-info@mts.com

MTS Systems (China) C., Ltd.

Floor 34, Building B,
New Caohejing International,
Business Center,
No.391, Guiping Road,
Xuhui District,
Shanghai 200233,
P.R.China
Telephone: +021-24151000
E-mail: MTSC-Info@mts.com

MTS Testing Solutions Pvt Ltd.

Unit No. 201 & 202, Second Floor
Donata Radiance,
Krishna Nagar Industrial Layout,
Koramangala, Bangalore - 560029
Karnataka, India
Telephone: +91 80 46254100
Email: mts.india@mts.com

Training Centers

Americas Training Center

MTS Systems Corporation

14000 Technology Drive
Eden Prairie, MN 55344-2290
USA
Telephone: 952-937-4000
Toll Free: 800-328-2255
E-mail: info@mts.com
Internet: www.mts.com

Asia Training Center

MTS Korea, Inc.

Sungkyunkwan Univ. Natural Sciences Campus:
#83158 Research Complex 2, 300 Cheoncheon-dong
Jangan-gu, Suwon-si,
Gyeonggi-do 440-746 Korea
Office Phone: +82-31-714-7151
E-mail: mtsk-info@mts.com

Europe Training Center

MTS Systems GmbH

Hohentwielsteig 3
14163 Berlin
Germany
Telephone: +49-(0)30 81002-222
E-mail: training.europe@mts.com

Company Locations

Sales & Service Organizations

NORTH AMERICA

Tech Support:

Toll Free: 1-800-328-2255
Outside US: 1-952-937-4000
E-mail: techsupport@mts.com

Canada

Telephone: 800-328-2255

Mexico

Telephone: 52-449-922-9200

United States

Telephone: 800-328-2255
1-952-937-4000

CENTRAL AND SOUTH AMERICA

Argentina

Telephone: (54) (11) 4331-5766
or (54) (11) 4331 4503
or (54) (11) 4331 4504

Brazil

Telephone: 55-11-5501-5300

Chile

Telephone: 55-11-5501-5300

Colombia

Telephone: 55-11-5501-5300

Costa Rica

Telephone: 2256-178

Ecuador

Telephone: 55-11-5501-5300

Peru

Telephone: 55-11-5501-5300

Venezuela

Telephone: 58-212-573-0270
or 58-212-573-8687

EUROPE

Tech Support:

Toll Free: +800-81002-222
Germany: +49-(0)30-81002-8222
France: +33-(0)1-58 43 90 05
E-mail: techsupport.europe@mts.com

Belarus

DL Kinesis, LLC
37/1 Surganova St, Suite 107
220013, Minsk
Republic of Belarus
Telephone: +375 (17) 235-39-57,
+375 (33) 322-41-83 or +375 (29) 761-76-27

Belgium (Sales – Electromechanical)

MTS Systems SAS
BAT EXA 16
16/18 rue Eugène Dupuis
94046 Créteil Cedex
France
Telephone: +33-(0)1-58 43 90 00

Belgium (Sales – Servohydraulic)

Technisch Advies- en
Handelsburo J.J. Bos B.V.
Marconistraat 11
NL-2809 PH Gouda
The Netherlands
Telephone: +31-(0)182-619 333

Czech Republic

JD DVORAK, s.r.o.
V Holešovičkách 1448/14
180 00 Praha 8
Czech Republic
Telephone: +420 284 688 050

France (Sales – Service)

MTS Systems SAS
BAT EXA 16
16/18 rue Eugène Dupuis
94046 Créteil Cedex
France
Telephone: +33-(0)1-58 43 90 00

Germany, Austria, Benelux, Southeast Europe (Sales – Service)

MTS Systems GmbH
Hohentwielsteig 3
14163 Berlin
Germany
Telephone: +49-(0)30 81002-0

Greece (Sales)

Neotek OE
P. Xystris & Co.
E. Venizelou 105
N. Smyrni 17123
Athens, Greece
Telephone: 30-210-9341533 or
30-210-9359142

Italy (Sales – Service)

MTS Systems S.r.l. socio unico
Strada Pianezza 289
10151 Torino, Italy
Telephone: +39-(0)11 45175 11 sel. pass.

Netherlands (Sales)

Technisch Advies- en
Handelsburo J.J. Bos B.V.
Marconistraat 11
NL-2809 PH Gouda
The Netherlands
Telephone: 31-(0)182-619 333

Northern Europe (Sales – Service)

MTS Systems Norden AB
Datavägen 37b
SE-436 32 Askim
Sweden
Telephone: +46-(0)31-68 69 99

Poland (Sales – Service)

Elhys Sp. Z o.o.
ul. Naukowa 45
PL 02-463
Warsaw, Poland
Telephone: 00-48-22-8633049 or
00-48-22-8633052

Russia GRV

LLC "BLM Synergie"
Electrozavodskaya st. 24
Moscow, 107023
Telephone: +7 (495) 781-3939

Russia MAT

SYNERCON, LLC
Varshavskoe shosse, 32
Moscow, 117105
Telephone: +7 (499) 611-15-09,
+7 (499) 611-52-89
or +7 (495) 741-59-04

Spain/Portugal (Sales – Service)

SEM, SA (Sistemas de Ensayo de
Materiales SA)
Moscou, 30, 1-1
E-08005 Barcelona
Spain
Telephone: +34 (0)93 221 0065

Turkey (Sales – Service)

MEGA Danışmanlık, Temsilcilik ve
Dis Ticaret Ltd.
Uzuncayir Cad. No: 31, A Blok D: 18-19
34722 - Hasanpasa
Istanbul, Turkey
Telephone: +90-216-326-4535

United Kingdom (Sales – Service)

MTS Systems Ltd.
40 Alan Turing Road
Surrey Research Park
Guildford
Surrey
GU2 7YF
United Kingdom
Telephone: +44 (0) 1483 446500

AFRICA AND THE MIDDLE EAST

Egypt

Cairo
Telephone: 20-22-642-3604

Israel

Lod
Telephone: 972-73-2570000

U.A.E., Oman, Qatar, Kuwait, Pakistan, Jordan

RAK
Telephone: 971 50 6559539

ASIA/PACIFIC

India

Karnataka
Telephone: 91-80-46254100

Indonesia

Jakarta
Telephone: +62-21-6251690

Japan

Nagoya
Telephone: 81-52-709-6328
Fax: 81-52-709-6434
Tokyo
Telephone: 81-3-6658-0901

People's Republic of China

Beijing
Telephone: +86-10-65900980
Fax: +86-10-65900942
Shanghai
Telephone: 86-21-5495 1122

Republic of Korea

Seongnam City
Telephone: 82-31-728-1600

Saudi Arabia, Bahrain

Al Khobar
Telephone: +966-54-882-4125

Singapore

Singapore
Telephone: 65-6778-3655

Taiwan

Kaohsiung
Telephone: 886-7-554-1800

Taichung

Telephone: 886-4-2245 4700

Taipei

Telephone: 886-2-2792 2440

AUSTRALIA AND NEW ZEALAND

Australia

Chatswood
Telephone: +61-(02) 9496-4200

New Zealand

Auckland
Telephone: +64-9-366-3999

Services & Accessories

	Introduction
	Services
	Load Cells
	Extensometers
	Grips & Fixtures
	Furnace & Chambers
	Testing System Components
	MTS Maintenance Parts
	Exchange & Repair Programs
	Company Locations



MTS Systems Corporation
14000 Technology Drive
Eden Prairie, Minnesota 55344-2290 USA
Toll Free: 800-328-2255
Telephone: 952-937-4000
Fax: 952-937-4515
E-mail: info@mts.com
<http://www.mts.com>

ISO 9001 CERTIFIED QMS

©2020 MTS Systems Corporation
100-522-590 English • Printed in U.S.A. • 1/2020

MTS, Bionix, Flat-Trac, FlexTest, MTS Acumen, MTS Landmark, MTS EM Extend, MTS Insight, MTS Echo, RPC, SWIFT, TestWare, and TestWorks are registered trademarks; Advantage, Aero-90, Elastomer Express, Just In Case, MAST, MicroProfilers, MTS Fundamental, SilentFlo, TestStar, and T/RAC are trademarks of MTS Systems Corporation within the United States. These trademarks may be protected in other countries.

Android is a trademark of Google Inc.

Apple is a registered trademark of Apple, Inc., registered in the U.S. and other countries. iOS is a registered trademark of Cisco Systems, Inc.

Products or company names other than trademarks of MTS Systems Corporation are used for identification purposes only and may be trademarks of their respective owners.

All specifications are subject to change without notice. Please contact MTS for specifications critical to your application. MTS reserves the right to change products, services and prices without notice.